



Sigma 6-16S

from serial no. 147414



Laboratory Centrifuge

Operating Manual

Please retain for later use!

© Copyright by
Sigma Laborzentrifugen GmbH
An der Unteren Söse 50
37520 Osterode am Harz
Germany

Tel.: +49 (0) 5522 / 5007-0
Fax: +49 (0) 5522 / 5007-12
Web: www.sigma-zentrifugen.de
E-mail: info@sigma-zentrifugen.de

Table of contents

1	General information	9
1.1	Importance of the operating manual	9
1.2	Intended use	9
1.2.1	Centrifuge	9
1.2.2	Rotors, buckets, carriers and accessories	9
1.2.3	Starting products	10
1.3	Warranty and liability	10
1.4	Copyright	10
1.5	Standards and regulations	11
1.6	Scope of Supply	11
2	Layout and mode of operation	12
2.1	Layout of the centrifuge	12
2.1.1	Functional and operating elements	12
2.1.2	Name plate	13
2.2	Mode of operation	14
2.2.1	Principle of centrifugation	14
2.2.1.1	Speed, radius, and relative centrifugal force	14
2.2.1.2	Density	14
3	Safety	15
3.1	Marking of the unit	15
3.2	Explanation of the symbols and notes	16
3.3	Responsibility of the operator	17
3.4	Requirements concerning the personnel	18
3.5	Informal safety instructions	19
3.6	Safety instructions	20
3.6.1	Electrical safety	20
3.6.2	Mechanical safety	20
3.6.3	Fire prevention	21
3.6.4	Chemical and biological safety	21
3.6.5	Safety instructions for centrifugation	22
3.6.6	Resistance of plastics	22
3.6.7	Safety of rotors and accessories	23
3.6.7.1	Marking of rotors and accessories	23
3.6.7.2	Service life of rotors and accessories	24
3.6.7.3	Service life of the Sigma "Comfort" rotor coating	25
3.7	Safety devices	27
3.7.1	Lid lock device	27
3.7.2	Standstill monitoring system	27
3.7.3	System check	27
3.7.4	Earth conductor check	27
3.7.5	Imbalance monitoring system	27
3.7.6	Rotor monitoring system	27
3.8	Measures in the event of hazards and accidents	28

Table of contents

3.9	Remaining hazards	28
4	Storage and transport	29
4.1	Dimensions and weight	29
4.2	Storage conditions	29
4.3	Notes on transport	29
4.4	Packaging	30
4.5	Transport safety device	30
5	Set-up and connection	32
5.1	Installation site	32
5.2	Power supply	32
5.2.1	Type of connection	32
5.2.2	Customer-provided fuses	33
6	Using the centrifuge	34
6.1	Initial start-up	34
6.2	Switching the centrifuge on	34
6.2.1	Opening and closing the lid	34
6.2.2	Installation of rotors and accessories	34
6.2.2.1	Installation of a rotor	34
6.2.2.2	Installation of an angle rotor with a hermetic lid	35
6.2.2.3	Installation of accessories	37
6.2.2.4	Adapters	39
6.2.2.5	Vessels	39
6.2.2.6	Blood bag systems	40
6.2.2.7	Centrifugation with tubes for mineral oil	41
6.3	Control system "Spincontrol S"	42
6.3.1	User interface	42
6.3.2	Manual mode	43
6.3.2.1	Starting a centrifugation run	43
6.3.2.2	Interrupting a centrifugation run	43
6.3.2.3	Interrupting a deceleration process	43
6.3.2.4	Selection, display, and modification of data	43
6.3.2.5	Standard menu	44
6.3.2.6	Process library menu	48
6.3.2.7	Parameters menu	49
6.3.2.8	Setup menu	52
6.3.2.9	Curve menu	55
6.3.2.10	Help menu	56
6.3.2.11	Changing the contrast	57
6.3.3	Program mode	57
6.3.3.1	Saving a program	58
6.3.3.2	Loading a program	58
6.3.3.3	Executing a program	59
6.3.3.4	Deleting a program	59
6.3.3.5	Automatic program rotation	60

Table of contents

6.3.4	Options for data input and output.....	61
6.4	Switching the centrifuge off.....	61
7	Malfunctions and error correction.....	62
7.1	General malfunctions.....	62
7.1.1	Emergency lid release	63
7.2	Table of error codes.....	64
7.3	Service contact	65
8	Maintenance and service.....	66
8.1	General notes	66
8.1.1	Cleaning	67
8.1.1.1	Cleaning of the centrifuge	68
8.1.1.2	Cleaning of rotors and buckets.....	68
8.1.1.3	Cleaning of adapters	68
8.1.2	Disinfection	69
8.1.2.1	Disinfection of the centrifuge	69
8.1.2.2	Disinfection of rotors and buckets	69
8.1.2.3	Disinfection of adapters.....	70
8.1.3	Autoclaving	71
8.2	Maintenance tasks.....	72
8.2.1	Maintenance of the centrifuge.....	72
8.2.2	Maintenance of rotors, buckets and adapters	72
8.2.2.1	Plastic accessories	73
8.2.2.2	Load-bearing bolts	73
8.2.3	Broken glass.....	75
8.3	Maintenance and service	76
8.4	Return of defective centrifuges or parts	77
9	Disposal.....	79
9.1	Disposal of the centrifuge	79
9.2	Disposal of the packaging.....	79
10	Technical data.....	80
10.1	Ambient conditions	80
10.2	Technical documentation	81
11	Appendix	82
11.1	Range of accessories	82
11.1.1	Rotor radii.....	82
11.2	Acceleration and deceleration curves	83
11.3	Table of the service life of rotors and accessories.....	85
11.4	Resistance data	86
11.5	EC declaration of conformity	91
11.6	Declaration of conformity – China RoHS 2.....	93
12	Index	95

Table of contents

1 General information

1.1 Importance of the operating manual

A fundamental requirement for the safe and trouble-free operation of the centrifuge is to be familiar with the fundamental safety instructions and all possible hazards.

The operating manual includes important information concerning the safe operation of the centrifuge.

This operating manual and, in particular, the notes on safety and hazards must be observed by all persons operating the centrifuge.

In addition, the local rules and regulations for the prevention of accidents must be complied with.

1.2 Intended use

This chapter and its subchapters describe the centrifuge's intended area of application. Any use other than the one described here or any use going beyond this description is regarded as improper use.

Sigma Laborzentrifugen GmbH cannot be held liable for any damage resulting from this.

The intended use also includes

- Observation of all of the notes and instructions that are included in the operating manual,
- Compliance with the care, cleaning and repair instructions.

1.2.1 Centrifuge

Centrifuges are power-driven machines that use centrifugal force to separate mixtures with a maximum density of 1.2 g/cm³ into their components. They are intended for this specific purpose only.

The centrifuge may only be operated by trained and authorised specialist personnel within closed spaces if the following conditions are fulfilled:

- The centrifuge has been properly set up and connected.
- The centrifuge is in a perfect technical state.
- Maintenance and service intervals are adhered to.
- The centrifuge is not located in an area with an explosion hazard.

1.2.2 Rotors, buckets, carriers and accessories

- The rotors, buckets, carriers and accessories have been produced by Sigma Laborzentrifugen GmbH and are authorised for use in this centrifuge.
- The maintenance intervals are adhered to.
- The service life specifications and, in particular, the specifications concerning the cycle limitations are complied with (see chapter 8.2.2 - "Maintenance of rotors, buckets and adapters").

1 General information

1.2.3 Starting products

The centrifuge is intended for the centrifugation of solid substances as well as mixtures of liquids or solids. The substances to be centrifuged must fulfil the following requirements:

- The substances must not damage the material of the centrifuge, rotors or inserts nor impair its mechanical strength.
- Biological material may only be centrifuged in special, certified containment systems with a biological seal to prevent the release of the material.

1.3 Warranty and liability

The warranty and liability are subject to our "General Conditions" that were distributed to the operator upon the conclusion of the contract.

Warranty and liability claims are excluded if they are due to:

- improper use.
- non-compliance with the safety instructions and hazard warnings in the operating manual.
- improper installation, start-up, operation, or maintenance of the centrifuge.

1.4 Copyright

The copyright concerning the operating manual remains with Sigma Laborzentrifugen GmbH.

The operating manual is solely intended for the operator and their personnel. It includes instructions and information that must not be

- duplicated,
- distributed, or
- communicated in any other way.

Non-compliance may be prosecuted under criminal law.

1.5 Standards and regulations

These operating instructions have been created in accordance with the relevant European standards and regulations (see chapter 11.5 - "EC declaration of conformity").

1.6 Scope of Supply

The centrifuge comprises:

- | | |
|--|------------------|
| • 1 rotor wrench, size 13 | Part no. 930 102 |
| • 1 hexagon socket wrench (size 5) | Part no. 930 051 |
| • 1 tube (30 g) heavy-duty grease for load-bearing bolts | Part no. 71 401 |

Documentation

Operating manual incl. EC declaration of conformity
(see chapter 11.5 - "EC declaration of conformity")

Accessories

according to your order, our order confirmation, and our delivery note.

2 Layout and mode of operation

2 Layout and mode of operation

2.1 Layout of the centrifuge

2.1.1 Functional and operating elements

- 1 Lid
- 2 User interface (see chapter 6.3.1 - "User interface")
- 3 Mains switch

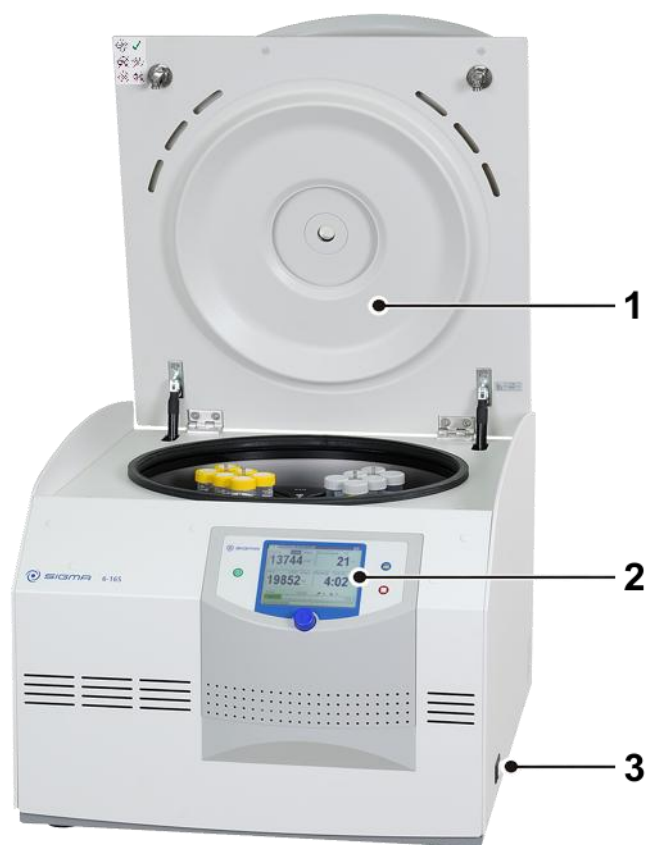


Fig. 1: Total view of the centrifuge

- 4 Name plate (see chapter 2.1.2 - "Name plate")
- 5 Fuse holder
- 6 Mains cable

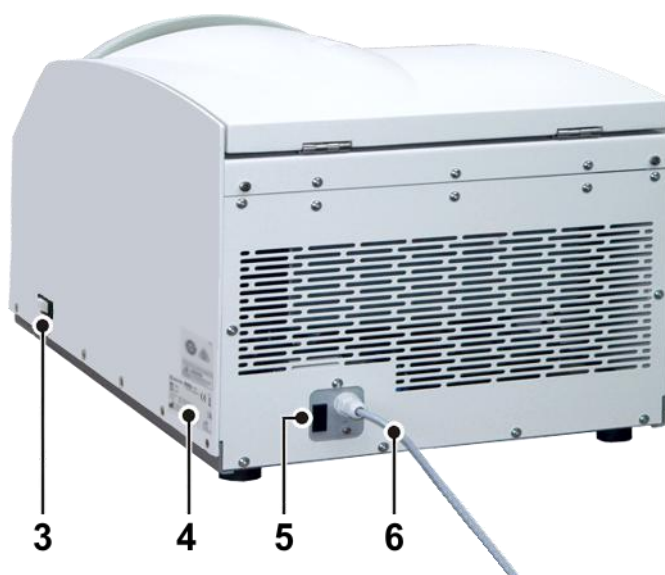


Fig. 2: Rear view of the centrifuge

2 Layout and mode of operation

2.1.2 Name plate

- 1 Manufacturer
- 2 Power consumption
- 3 Max. speed
- 4 Max. kinetic energy
- 5 Product name
- 6 Part number
- 7 Serial number
- 8 Nominal voltage
- 9 Product designation
- 10 CE mark in compliance with the directive 2006/42/EC
- 11 Symbol for special disposal (see chapter 9 - "Disposal")
- 12 Date of manufacture
- 13 Consult operating manual
- 14 Max. permissible density

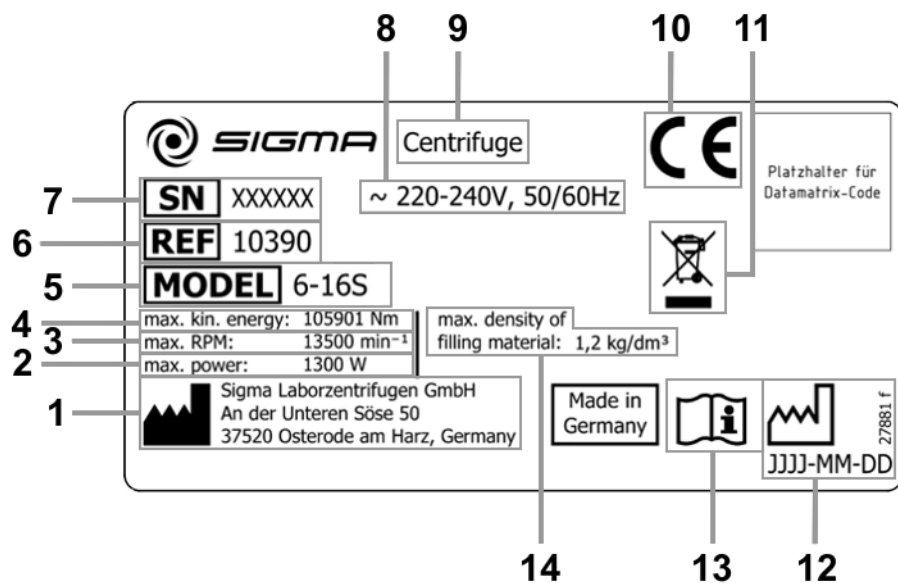


Fig. 3: Example of a name plate

2 Layout and mode of operation

2.2 Mode of operation

2.2.1 Principle of centrifugation

Centrifugation is a process for separating heterogeneous mixtures of substances (suspensions, emulsions or gas mixtures) into their components. Rotating on a circular trajectory, the mixture of substances is subjected to centrifugal acceleration many times greater than the acceleration due to gravity.

Centrifuges harness inertia in the rotor chamber to separate substances. Due to their higher inertia, particles or media with a higher density migrate outwards, thereby displacing any components with a lower density which are forced towards the centre.

The centrifugal acceleration of an object in a centrifuge depends on the distance of the object from the axis of rotation and on the angular velocity. The larger the radius of the rotor chamber and the higher the speed are, the greater the centrifugal acceleration will be.

For detailed information about the various areas of application, see www.sigma-zentrifugen.de → [Applications].

2.2.1.1 Speed, radius, and relative centrifugal force

The acceleration g , which the samples are subject to, can be increased by increasing the radius in the rotor chamber and by increasing the speed. These three parameters are interdependent and linked with each other via the following formula:

$$\text{Relative centrifugal force RCF} = 11.18 \times 10^{-6} \times r \times n^2$$

r = radius in cm

n = speed in rpm

RCF without any dimension

If two values are entered, the third value is determined by way of the stated formula. If, afterwards, the speed or the radius is changed, the resulting relative centrifugal force will be recalculated automatically by the control unit. If the RCF is changed, the speed will be adapted while the specified radius is maintained.

The speed-gravitational-field-diagram provides an overview of the relationship between speed, radius, and RCF.

2.2.1.2 Density

The laboratory centrifuge is suitable for the separation of constituents of different densities in mixtures with a maximum density of 1.2 g/cm^3 . All information concerning the speed of rotors and accessories refers to liquids with a density corresponding to this specification. If the density is above this value, the maximum permissible speed of the centrifuge must be reduced based on the following formula:


$$n = n_{\max} \times \sqrt{(1.2/\rho)}$$

ρ = density in g/cm^3

3 Safety

3.1 Marking of the unit

The following symbols are used on this centrifuge:

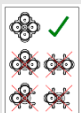	On (Power)		Arrow indicating the direction of rotation
	Off (Power)		Rotor loading information (see chapter 6.2.2.3 - "Installation of accessories")
	Name plate (see chapter 2.1.2 - "Name plate")		CE mark in compliance with the directive 2006/42/EC
	Do not dispose as part of domestic waste		Consult operating manual
	RCM mark (only for Australia)		China RoHS 2 mark (only for China)
	California Proposition 65 mark (only for the USA)		



Safety indications on the centrifuge must be kept readable at all times. If necessary, they must be replaced.



The marking varies depending on the version and country of destination of the centrifuge.

3 Safety

3.2 Explanation of the symbols and notes

In this operating manual, the following names and symbols to indicate hazards are used:



DANGER

This symbol stands for a **direct** hazard to the life and health of persons.

Non-observance of these symbols **causes** serious health problems up to life-endangering injuries.



DANGER

This symbol stands for a **direct** hazard to the life and health of persons due to electrical voltage.

Non-observance of these symbols **causes** serious health problems up to life-endangering injuries.



WARNING

This symbol stands for a **potential** hazard to the life and health of persons.

Non-observance of these symbols **can** cause serious health problems up to life-endangering injuries.



CAUTION

This symbol indicates a potentially hazardous situation.

Non-observance of these notes can cause minor injuries or damage to property.



NOTE

This symbol indicates important information.

3.3 Responsibility of the operator

The operator is the person who uses the centrifuge for commercial or economic purposes or lets a third party use it and who bears the legal product responsibility for the protection of personnel or third parties during operation.

Personnel

The operator is obliged:

- to comply with the specifications concerning the personnel (see chapter 3.4 - "Requirements concerning the personnel")
- to clearly define the areas of responsibility of the personnel concerning the operation, cleaning and care.
- to ensure that any personnel undergoing training as well as persons with physical, sensory or mental disabilities only work on the centrifuge under the supervision of an experienced person.
- to provide training for the tasks on the centrifuge, safety at the workplace and accident prevention.

Work environment

The operator is obliged:

- to perform a risk assessment concerning potential accidents in connection with the centrifuge and take design-related countermeasures, if necessary.
- to perform a risk assessment concerning the specific hazards with regard to the processing of the various products in the centrifuge (e.g. leakage of harmful substances, contaminated surfaces, fire hazard) and to take appropriate protective measures, if necessary.
- to perform compatibility tests of the substances to be processed in the centrifuge. Substances that damage the materials of the centrifuge, rotors or inserts or that weaken their mechanical strength must not be used.
- to provide standard operating procedures (SOPs) for the processing of the various products in the centrifuge.
- to check the safety-conscious work of the personnel in compliance with the operating instructions, other applicable documents and SOPs as well as with the relevant EU directives and the national laws concerning health and safety and the prevention of accidents at regular intervals (e.g. every month).
- to ensure that a sufficient quantity of suitable fire-extinguishing equipment is provided and ready for use.

3 Safety

Maintenance

The operator is obliged:

- to have the centrifuge checked by authorised specialist personnel at regular intervals to ensure its correct operation.
- to immediately stop using a damaged device and to have defective parts replaced without delay.
- to keep all safety and hazard notes on the centrifuge readable at all times and to replace them, if necessary.
- to check – at their own accord and at regular intervals – the number of cycles performed and the condition of the rotors and accessories.
- to comply with the inspection obligation as per the information FBRCI-025 of the German Berufsgenossenschaft Rohstoffe und chemische Industrie for centrifuges operated in Germany (see chapter 8.3 - "Maintenance and service").

3.4 Requirements concerning the personnel



Risk of injury if the personnel are not sufficiently qualified

If unqualified personnel perform work on the centrifuge or are present in the danger zone of the centrifuge, hazards result that can cause serious injuries and considerable damage to property.

- Ensure that all the tasks are performed by personnel with the corresponding qualifications.
- Ensure that unqualified personnel stay clear of the danger zones.



Risk of fatal injury to unauthorised persons due to hazards in the danger zone or work area

Unauthorised persons who do not fulfil the requirements described herein are not aware of the hazards in the work area. This is why there is a risk of serious or even fatal injuries for unauthorised persons.

- Ensure that unauthorised persons stay clear of the danger zone and work area.
- If in doubt, address these persons and instruct them to leave the danger zone and work area.
- Interrupt any running work if unauthorised persons are present in the danger zone or work area.

This manual uses the following personnel qualifications for various areas of activity:

Qualified electrician

Due to their special training, knowledge, experience and familiarity with the relevant standards and regulations, qualified electricians are in the position to perform work on electrical systems and to autonomously identify and prevent possible hazards.

Qualified electricians have been specifically trained for the environment in which they work and they are familiar with all the relevant standards and regulations.

Qualified electricians must fulfil the requirements as set out in the applicable legal provisions concerning the prevention of accidents.

Specialised personnel

Due to their special training, knowledge, experience and familiarity with the relevant regulations, specialised personnel are in the position to perform any tasks assigned to them and to autonomously identify and prevent possible hazards.

Operating personnel

Only trained, specialised personnel are authorised to operate the unit. The persons operating the unit must

- be familiar with the fundamental health, safety, and accident prevention regulations,
- have read and understood this operating manual, in particular the safety sections and warning notes, and confirmed this with their signature,
- have been instructed in the operation and maintenance of this centrifuge.

3.5 Informal safety instructions

- This operating manual is a part of the product.
- The operating manual must be kept at the location of use of the centrifuge. Ensure that it is accessible at all times.
- The operating manual must be handed over to any subsequent owner or operator of the centrifuge.
- Any changes, additions or updates received must be added to the operating manual.
- In addition to the operating manual, the general and local rules and regulations concerning the prevention of accidents and the protection of the environment must also be supplied.
- Safety and danger indications on the centrifuge must be kept readable at all times. If necessary, they must be replaced.

3 Safety

3.6 Safety instructions

3.6.1 Electrical safety

As protection against electric shock, the centrifuge is equipped with an earthed mains power cable and connector. To ensure the effectiveness of this safety feature, the following must be ensured:



DANGER

- Ensure that the wall socket is properly wired and grounded.
- Check that the mains voltage agrees with the nominal voltage listed on the name plate.
- Ensure that the mains power cable is intact prior to using the centrifuge. Damaged or faulty mains power cables must be replaced immediately.
- Do not place vessels containing liquid on the centrifuge lid or within the safety distance of 30 cm around the centrifuge. Spilled liquids may get into the centrifuge and damage electrical or mechanical components.
- Service tasks or repairs of the electrical system for which the housing needs to be removed must only be carried out by authorised specialist personnel.
- Inspect the electrical equipment of the unit regularly. Defects such as loose or burnt cables must be eliminated immediately.
- Following the completion of any type of repair or service, the authorised specialist personnel must perform final inspection and testing in compliance with the relevant standards (see chapter 8.3 - "Maintenance and service").

3.6.2 Mechanical safety

In order to ensure the safe operation of the centrifuge, observe the following:



WARNING

- Do not open the lid when the rotor is in motion!
- Do not reach into the rotor chamber when the rotor is in motion!
- Do not use the centrifuge if it was installed incorrectly.
- Do not use the centrifuge without panels.
- Do not use the centrifuge if the rotors and inserts show signs of corrosion or other defects.
- Only use the centrifuge with rotors and accessories that have been approved by the manufacturer. In case of doubt, contact the manufacturer (see chapter 7.3 - "Service contact").
- Do not hold your fingers between the lid and the housing when closing the lid. Risk of crushing!
- Defective lid relieving devices could cause the centrifuge lid to fall (contact the service department, if necessary). Risk of crushing!
- Do not hit or move the centrifuge during its operation.
- Do not lean against or rest on the centrifuge during its operation.

3 Safety



WARNING

- Do not spin any substances that could damage the material of the rotors and buckets of the centrifuge in any way. Highly corrosive substances, for example, damage the material and affect the mechanical strength of the rotors and buckets.
- Stop the centrifuge immediately in the event of a malfunction. Eliminate the malfunction (see chapter 7 - "Malfunctions and error correction") or inform the service department of the manufacturer (see chapter 7.3 - "Service contact").
- Ensure that all repairs are performed only by authorised and specialised personnel.
- Prior to any start-up, check the centrifuge, rotor, and accessories for signs of damage that can be discerned from the outside. Special attention must be paid to all of the rubber parts (e.g. motor cover, lid seal, and adapters) in terms of visible structural changes. Defective parts must be replaced immediately.
- Open the centrifuge when it is not in use so that moisture can evaporate.

3.6.3 Fire prevention



DANGER

- Do not spin explosive or inflammable substances.
- Do not use the centrifuge within hazardous locations.

3.6.4 Chemical and biological safety

If pathogenic, toxic, or radioactive samples are intended to be used in the centrifuge, it is in the responsibility of the user to ensure that all necessary safety regulations, guidelines, precautions, and practices are adhered to accordingly.



DANGER

- Infectious, toxic, pathogenic, and radioactive substances may only be used in special, certified containment systems with a bio-seal in order to prevent the material from being released.
- Take suitable precautions for your own safety if there is a risk of toxic, radioactive, or pathogenic contamination
- Materials that chemically react with each other with a high level of energy are prohibited.



WARNING

- Keep informed about local measures to avoid harmful emissions (depending on the substances to be centrifuged).
- Protective clothing is not required for the operation of the centrifuge. The materials to be centrifuged may, however, require special safety measures (e.g. centrifugation of infectious, toxic, radioactive, or pathogenic substances).

3 Safety

3.6.5 Safety instructions for centrifugation

For safe operation, observe the following before starting the centrifuge:



WARNING

- Ensure that the centrifuge was set up properly (see chapter 5 - "Set-up and connection").
- Keep a safety range of at least 30 cm free around the centrifuge as well as with regard to walls and other devices.
- Do not store any dangerous goods in the centrifuge area.
- Do not stay in the safety area longer than what is absolutely necessary for the operation of the centrifuge.
- Only use the centrifuge with rotors and accessories that have been approved by the manufacturer. We explicitly warn against the use of equipment of poor quality. Breaking glass or bursting vessels can cause dangerous imbalances at high speeds.
- Ensure that the rotor and buckets are correctly fitted (see chapter 6.2.2.1 - "Installation of a rotor").
- Observe the instructions on the installation of accessories (see chapter 6.2.2.3 - "Installation of accessories").
- The rotor must be loaded in a rotationally symmetrical manner at equal weights.
- If liquids with a density $> 1.2 \text{ g/cm}^3$ are used, reduce the speed (see chapter 2.2.1.2 - "Density").
- Do not use the centrifuge if the rotor is loaded asymmetrically.
- Do not use the centrifuge with tubes that are excessively long.

3.6.6 Resistance of plastics

Chemical influences have a strong effect on the polymeric chains of plastics, and, therefore, on their physical properties. Plastic parts can be damaged if solvents, acids, or alkaline solutions are used.



NOTE

- Refer to the resistance data (see chapter 11.4 - "Resistance data")!

3.6.7 Safety of rotors and accessories

3.6.7.1 Marking of rotors and accessories

Batch number and serial number

During production, every rotor and bucket receives a batch number enabling conclusions to be drawn concerning the production process and the subsequent quality inspection.

Some rotors also have an additional serial number providing further detailed information.

The batch number and serial number is engraved on the rotor as follows:

- 1 Batch number
- 2 Serial number

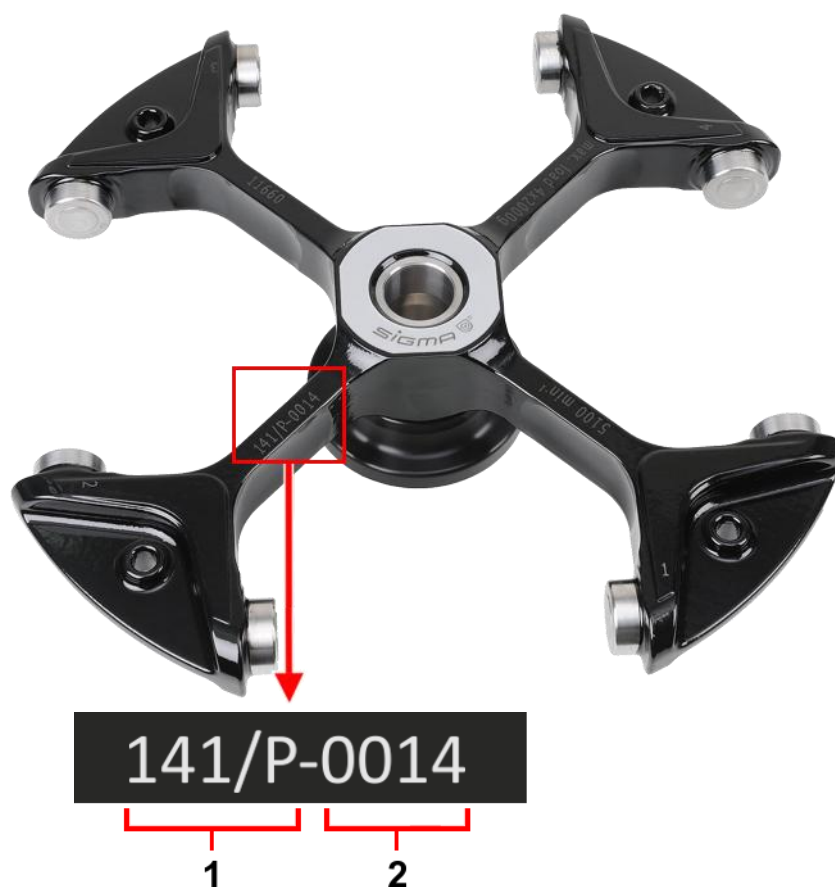


Fig. 4: Rotor with an engraved batch number and serial number (example)



If there are any enquiries concerning the rotor, please state the batch number and serial number!

Marking of the service life of rotors and accessories

see chapter 3.6.7.2 - "Service life of rotors and accessories"

Marking of rotors with Sigma "Comfort" rotor coating

see chapter 3.6.7.3 - "Service life of the Sigma "Comfort" rotor coating"

3 Safety

3.6.7.2 Service life of rotors and accessories

The rotors and accessories have a limited service life.



WARNING

- Perform regular checks (at least once per month) for safety reasons!
- Pay special attention to changes, such as corrosion, cracks, material abrasion, etc.

- After 10 years, they must be inspected by the manufacturer.
- After 50,000 cycles, the rotor must be scrapped for reasons of safety.
- If other data concerning the service life are engraved on the rotor or bucket, these data shall apply accordingly. For example, a bucket with the engraving "max. cycles = 15000" has a service life of 15,000 cycles, and a rotor with the engraving "Exp. date 01/27" must be scrapped in January 2027 at the latest (see figure).
- If a specification concerning the maximum number of cycles **and** a specification concerning the service life (i.e. a date) are provided, the specification that occurs first shall apply.

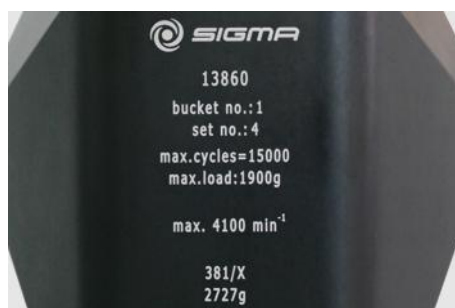


Fig. 5: Different service life – engraving on the bucket/rotor

i
NOTE

The information engraved on the rotor or bucket regarding the service life is binding. The information provided by the control system is for informational purposes only and, thereby, non-binding.

i
NOTE

- Refer to the table of the service life of rotors and accessories (see chapter 11.3 - "Table of the service life of rotors and accessories")!

3.6.7.3 Service life of the Sigma "Comfort" rotor coating

Marking

All "Comfort" coated swing-out rotors have a "C" shown after their number. These rotors do not need to be greased during the service life of the coating.

- 1 Marking indicating the "Sigma Comfort" rotor coating
- 2 Load-bearing bolt



Fig. 6: Marking of a rotor with "Comfort" coating (example)



A list of the available rotors with a "Comfort" rotor coating as well as information about the service life of the coating can be found in the "Sigma Comfort rotor coating" document that comes supplied with every rotor having a "Comfort" rotor coating.



It is not permissible to combine "Comfort" coated rotors with buckets or carriers with a cleanroom coating.

Service life

The service life of the coating varies and depends on the degree of utilisation of the rotor. Tests have shown a service life of 7,000 to 40,000 cycles. The service life of the "Comfort" rotor coating can be increased by way of specific measures:

The following factors have an influence on the service life:

- Speed/load:
The service life of the coating can be increased by a factor of 3.5 if the speed or load is decreased by approximately 15%. In case of a load decrease by 20%, the service life increases approximately by a factor of 5.5.
- Temperature:
The service life of the coating decreases by approximately 30% at lower temperatures (approximately 4°C) compared to room temperature (approximately 20°C).
- Use of buckets:
The buckets and rotor form a joint unit. For the coating to reach the maximum possible service life, ensure to always use the same set of buckets and to always install the individual buckets in their specific position within the rotor.

3 Safety

- **Cleaning frequency:**
Regular cleaning of the accessories is indispensable. The more frequently the rotor needs to be cleaned, the shorter the service life of the coating will be.



Traces of use on the black coating do not affect the anti-friction properties (see the following illustration).



Fig. 7: Load-bearing bolt with traces of use – no greasing required

End of the service life of the "Comfort" rotor coating

When the coating has worn off nearly completely and the metallic surface of the load-bearing bolt becomes visible (see the following illustration), the anti-friction effect decreases. As a result, the buckets will swing out irregularly, thereby potentially resulting in unwanted imbalance.

- From this moment on, the rotor must be used with greased load-bearing bolts (heavy-duty grease for load-bearing bolts, ref. no. 71401) until it reaches the end of its service life.
- Mark the rotor accordingly to prevent it from being used without grease on the load-bearing bolts.



Fig. 8: Load-bearing bolt with completely worn-off anti-friction coating – greasing required

3.7 Safety devices

3.7.1 Lid lock device

The centrifuge can only be started when the lid is properly closed. The electrical lock must be locked. The lid can only be opened when the rotor has stopped. If the lid is opened by way of the emergency release system during operation, the centrifuge will immediately switch off and decelerate brakeless. If the lid is open, the drive is completely separated from the mains power supply, i.e. the centrifuge cannot be started (see chapter 7.1.1 - "Emergency lid release").

3.7.2 Standstill monitoring system

Opening of the centrifuge lid is only possible if the rotor is at a standstill. This standstill is checked by the microprocessor.

3.7.3 System check

An internal system check monitors the data transfer and sensor signals with regard to plausibility. Errors are detected with extreme sensitivity and displayed as error messages in a dialog box (see chapter 7.2 - "Table of error codes").

3.7.4 Earth conductor check

An earth conductor check can be carried out by authorised and specialised personnel using a suitable measuring instrument. Please contact the Sigma service department (see chapter 7.3 - "Service contact").

3.7.5 Imbalance monitoring system

A dialog box may pop up or emit a sound signal in order to indicate that the centrifuge is in the inadmissible imbalance range. If the rotor is loaded unevenly, the drive will be switched off in the acceleration phase or during the run.

3.7.6 Rotor monitoring system

When a rotor number and, if applicable, a bucket number are selected, the computer will automatically check whether the entered speed or the entered gravitational field are permissible for the selected rotor.

3 Safety

3.8 Measures in the event of hazards and accidents



DANGER

- If an emergency arises, switch off the centrifuge immediately!
- If in doubt, call the emergency doctor!

3.9 Remaining hazards

The centrifuge was built in accordance with the state of the art and in compliance with the generally recognized safety rules. However, danger to life and limb of the operator, or of third parties, or impairments of the unit or other material assets cannot be completely excluded when the unit is being used.

- Use the unit only for the purpose that it was originally intended for (see chapter 1.2 - "Intended use").
- Use the unit only if it is in a perfect running state.
- Immediately eliminate any problems that can affect safety.

4 Storage and transport

4.1 Dimensions and weight

	Sigma 6-16S
Height:	483 mm
Height with open lid:	990 mm
Width:	581 mm
Depth:	711 mm
Weight:	116 kg

4.2 Storage conditions

The centrifuge can be stored in its original packaging for up to a year.

- Store the centrifuge only in dry rooms.
- The permissible storage temperature is between -20°C and +60°C.
- If you would like to store it for more than one year, or if you intend to ship it overseas, please contact the manufacturer.

4.3 Notes on transport

- Install the transport safety device (see chapter 4.5 - "Transport safety device")
- Always lift the centrifuge with a lifting device.
- When lifting the centrifuge, always reach under the centrifuge from the side.


CAUTION

The centrifuge weighs approx. 116 kg!

- For transport, use suitable packaging and, if at all possible, the original packaging (see chapter 4.4 - "Packaging").

4 Storage and transport

4.4 Packaging

The centrifuge is packaged in a wooden crate.

- After taking off the lid, remove the side panels.
- Remove the packaging material.
- Lift the centrifuge upwards with a lifting device to lift it safely. When lifting the centrifuge, always reach under the centrifuge from the side.



CAUTION

The centrifuge weighs approx. 116 kg!

- Retain the packaging for any possible future transport of the centrifuge.

4.5 Transport safety device

The transport safety device consists of two plastic screws. They are located on the base plate and can be accessed from below.



CAUTION

The transport safety device must be removed prior to start-up because the screws lock the motor bearings!

Removal

- Lift the centrifuge upwards at the front side. Always reach under from the side.
- Position a suitable object, e.g. a wooden block, between the tabletop and centrifuge. The two plastic screws are now visible on the base plate.
- Unscrew the screws by hand anti-clockwise.

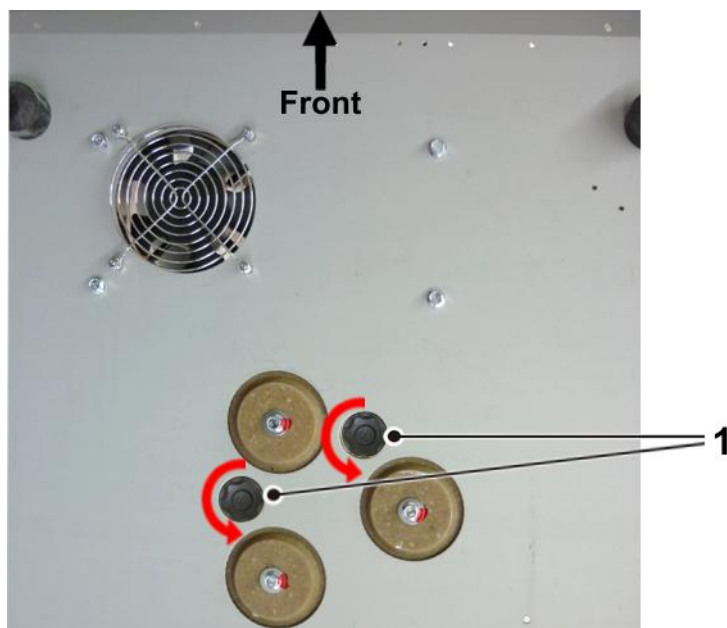
4 Storage and transport**1 Locking screws**

Fig. 9: Unscrew the locking screws by turning them anti-clockwise

- Retain the transport safety device for the possibility of the return of the centrifuge.

5 Set-up and connection

5 Set-up and connection

5.1 Installation site

Operate the centrifuge only in closed and dry rooms.

All the energy supplied to the centrifuge is converted into heat and emitted to the ambient air.

- Ensure sufficient ventilation.
- Keep a safety range of at least 30 cm free around the centrifuge as well as with regard to walls or other devices so that the vents in the machine remain unobstructed and fully effective.
- Do not subject the centrifuge to thermal stress, e.g. by positioning it near heat generators.
- Avoid direct sunlight (UV radiation).
- The table must be stable and have a solid, even surface.
- Attention: During transport from cold to warmer places, condensational water will collect inside the centrifuge. It is important to allow sufficient time for drying (min. 24 h) before the centrifuge can be used again.

5.2 Power supply

5.2.1 Type of connection



DANGER

The operating voltage on the name plate must correspond to the local supply voltage!



CAUTION

The mains power plug is an isolating device which is why it must be accessible at all times.

Sigma centrifuges are units of protection class I. The centrifuges of this model series have a three-wire power cord with a fixed cable. They are equipped with a mains power switch with an integrated thermal circuit breaker.

- Switch the unit off by actuating the mains power switch.
- If it has tripped, let the circuit breaker cool for approximately 2 minutes.
- Switch the unit on.

The centrifuge is now ready for operation.

5 Set-up and connection

5.2.2 Customer-provided fuses

Typically, the centrifuge must be protected with 16 Amp B fuses that are to be provided by the customer.



To ensure safe disconnection in the event of a fault, an AC/DC-sensitive RCD (residual current device) must be integrated in the wiring system of the building.

6 Using the centrifuge

6 Using the centrifuge

6.1 Initial start-up



DANGER

- Before the initial start-up, please ensure that your centrifuge is properly set up and installed (see chapter 5 - "Set-up and connection").

6.2 Switching the centrifuge on

- Press the mains power switch.

The display then illuminates. The centrifuge is ready for operation.

6.2.1 Opening and closing the lid

The lid can be opened if the centrifuge is at a standstill and if the lid key is illuminated.

- Press the lid key in order to open the lid.

The centrifuge cannot be started if the lid is opened.

- To close, press with both hands slightly on the lid until the lock is audibly locked.



WARNING

Do not place your fingers between the lid and the housing when closing the lid. Risk of crushing!

6.2.2 Installation of rotors and accessories

6.2.2.1 Installation of a rotor

- Open the centrifuge lid by pressing the lid key.
- Unscrew the rotor tie-down screw from the motor shaft (counter-clockwise).
- Lower the rotor with its central bore straight down onto the motor shaft.
- Tighten the rotor tie-down screw clockwise with the supplied rotor wrench with 10 Nm (Option: torque wrench 10 Nm, part no. 17060). In doing so, hold the rotor at its outer rim.
- Follow the safety instructions and hazard warnings (see chapter 3 - "Safety")!



WARNING

Loosen and tighten the rotor tie-down screw at regular intervals (depending on the frequency of use) to ensure a proper connection between the rotor holder and the motor shaft.

6 Using the centrifuge



CAUTION

When using rotors for microtiter plate formats:

Ensure that the plate holders are inserted together with the plates into the buckets.



NOTE

The lid screw serves for the fastening of the lid onto the rotor only, not for the fastening of the rotor onto the motor shaft.

Removing a rotor

- Loosen the rotor tie-down screw by turning it anti-clockwise and remove the rotor.

6.2.2.2 Installation of an angle rotor with a hermetic lid

- Open the centrifuge lid by pressing the lid key.
- Unscrew the rotor tie-down screw from the motor shaft (counter-clockwise).
- Screw the hermetic lid onto the rotor and tighten it clockwise by hand.
- Lower the rotor with the lid with its central bore straight down onto the motor shaft.
- Insert the rotor tie-down screw into the motor shaft. Tighten the rotor tie-down screw clockwise with the supplied rotor wrench with 10 Nm (Option: torque wrench 10 Nm, part no. 17060). In doing so, hold the rotor at its outer rim.
- Follow the safety instructions and hazard warnings (see chapter 3 - "Safety")!



NOTE

The rotor can also be used without a lid.



NOTE

The lid screw serves for the fastening of the lid onto the rotor only, not for the fastening of the rotor onto the motor shaft.

6 Using the centrifuge

Removing a rotor

- Loosen the rotor tie-down screw anti-clockwise and remove the rotor together with the lid.
- Loosen the lid by turning it anti-clockwise and then remove it.

- 1 Rotor tie down screw
- 2 Lid
- 3 Rotor

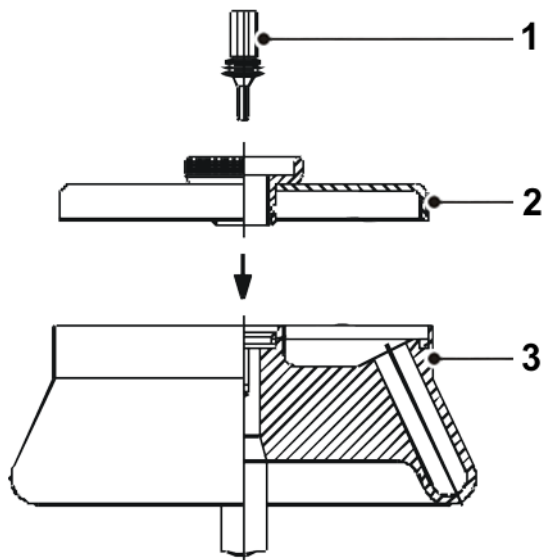


Fig. 10: Angle rotor with a hermetic lid

6 Using the centrifuge


6.2.2.3 Installation of accessories

- Only use vessels that are suitable for the rotor.
- In the case of swing-out rotors, all of the compartments must be equipped with buckets.
- Always load rotationally symmetrical compartments of the rotors with the same accessories and fill to avoid imbalance.

Centrifugation with vessels of various sizes

Working with vessel of various sizes is possible. In this case, however, it is very important that the rotationally symmetrical inserts are identical (see figure).

correct

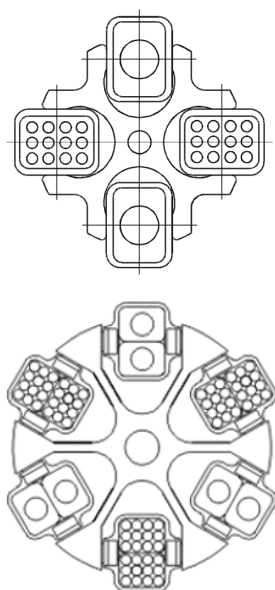

incorrect

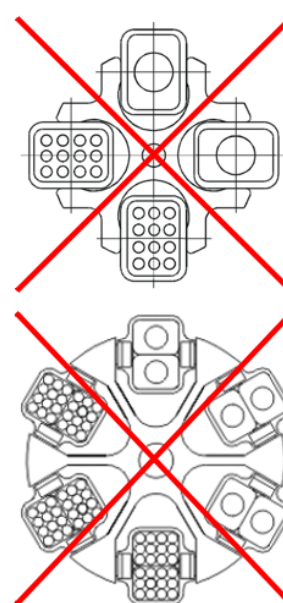



Fig. 11: Permissible and impermissible loading of a swing-out rotor with vessels of various sizes (example illustration)

6 Using the centrifuge

Centrifugation with low capacity

- Install the sample vessels in a rotationally symmetrical manner so that the buckets and their suspensions are loaded evenly.
- It is not permissible to load angle rotors on only one axis.

correct

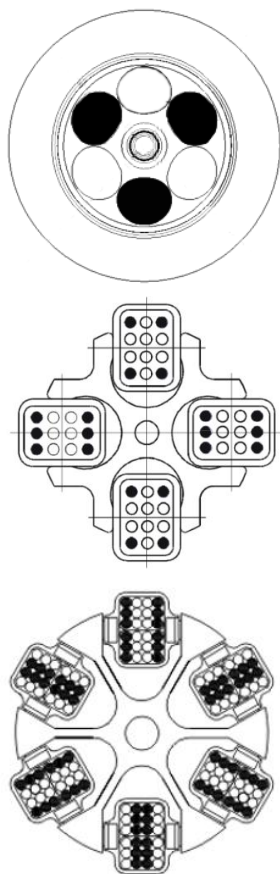

incorrect

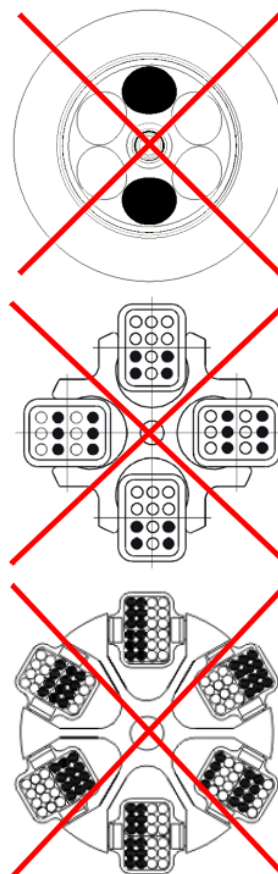


Fig. 12: Permissible and impermissible loading of an angle rotor and a swing-out rotor (example illustrations)



Pay attention to the marking of the centrifuge (see the illustration below)! Safety indications on the centrifuge must be kept readable at all times. If necessary, they must be replaced.

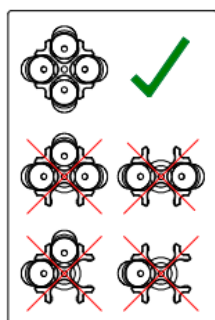


Fig. 13: Safety indication on the centrifuge: Loading of a swing-out rotor

6 Using the centrifuge

Double rectangular bucket 13420 and bucket 13621

The double rectangular bucket 13420 and the bucket for microtiter plates 13621 are dynamically balanced at the factory. Therefore, the buckets must be installed in sets of 4 and in the prescribed order (see the numbers on the buckets and on the rotor).



CAUTION

Mixed use with only two buckets on opposite sides will cause increased imbalance!

6.2.2.4 Adapters

In order to ensure easy handling, even if vessels of various sizes are used, carrier systems were developed.

- Load the opposite adapters with the same number of vessels and with the same weights in order to avoid imbalance.
- If all of the compartments of a carrier are not used, the buckets must be loaded evenly. Loading the edges of a bucket only is not permissible.

6.2.2.5 Vessels

- Load the vessels outside of the centrifuge. Liquids in the buckets or multiple carriers cause corrosion.
- Fill the vessels carefully and arrange them according to their weight. Imbalances result in the excessive wear of the bearings.
- Always fill the tubes up to their useful volume (= the volume that is stated for the tube).
- After the centrifugation, remove the vessels carefully in order to prevent the samples from mixing.
- Follow the safety instructions and hazard warnings (see chapter 3 - "Safety")!

Maximum speed for tubes

Some tubes, such as centrifuge glass tubes, microtubes, culture tubes, PTFE tubes and especially high-volume tubes can be used in our rotors, buckets, and adapters at higher speeds than their breaking limit.



NOTE

When using glass vessels, the maximum value of 4,000 x g must not be exceeded (except special high-strength glass tubes; please refer to the information provided by the manufacturer).



NOTE

When installing the 500 ml bottles use the supplied supporting rings

6 Using the centrifuge

6.2.2.6 Blood bag systems

The centrifuge tolerates an imbalance of approx. 20 g. Therefore, it is important to balance the buckets carefully.

- All four places of a rotor must be loaded with buckets. It is permissible to fill just two opposite buckets with blood bag systems.
- The opposite buckets including the filled blood bag systems must have an equal weight.
- The blood bags have to be put into the opposite buckets mirror-inverted.
- If in one bucket the main blood bag is on the left side, then the main blood bag inserted in the opposite bucket has to be on the right side.

- 1 Filter
2 Blood bag

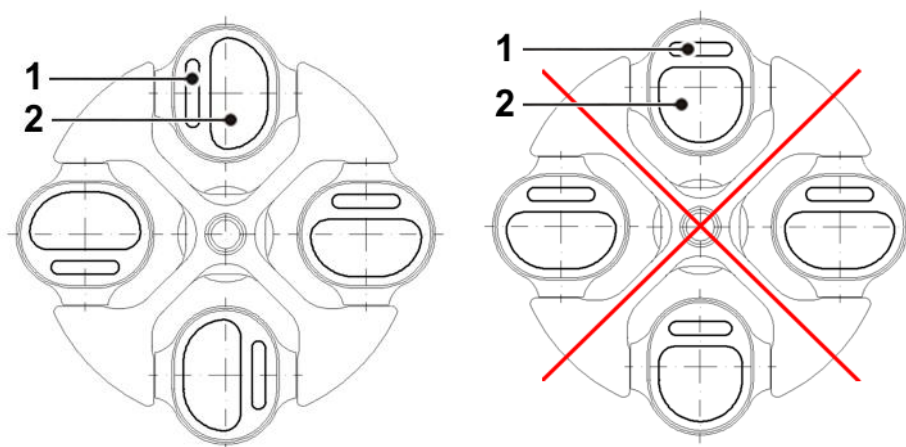


Fig. 14: Loading of blood bag systems

- When using smaller bag systems or in case of incompletely filled blood bags adapters (e.g. part no. 17750) have to be inserted together with the blood bag systems. This will avoid slipping of the bags which could result in an imbalance.
- For the centrifugation of human blood suitable acceleration and deceleration curves must be used. Curve no. 12 is the most suitable for this purpose.



If the centrifuge does not reach the mentioned RCF values, the maximum gravitational field must be set and the run time must be extended.

6 Using the centrifuge

6.2.2.7 Centrifugation with tubes for mineral oil

In the previous ASTM standard, a maximum gravitational field of 800 x g is required for the centrifugation with tubes for mineral oil. Practically, it is possible to reach higher RCF-values with the rotor 11165 and the buckets 13113, 13114 or 13116.

**CAUTION**

The required gravitational field must not be exceeded.
Values above 800 x g increase the risk of glass breakage!

- Use sealing caps (part no. 17107).

6 Using the centrifuge

6.3 Control system "Spincontrol S"

6.3.1 User interface

The centrifuge is operated via three buttons with integrated light-emitting diodes and one function knob. The display is divided into several different fields. The various functions of the system can be called up by pressing and turning the function knob.

- 1 Start key
- 2 Display
- 3 Function knob
- 4 Stop key
- 5 Lid key

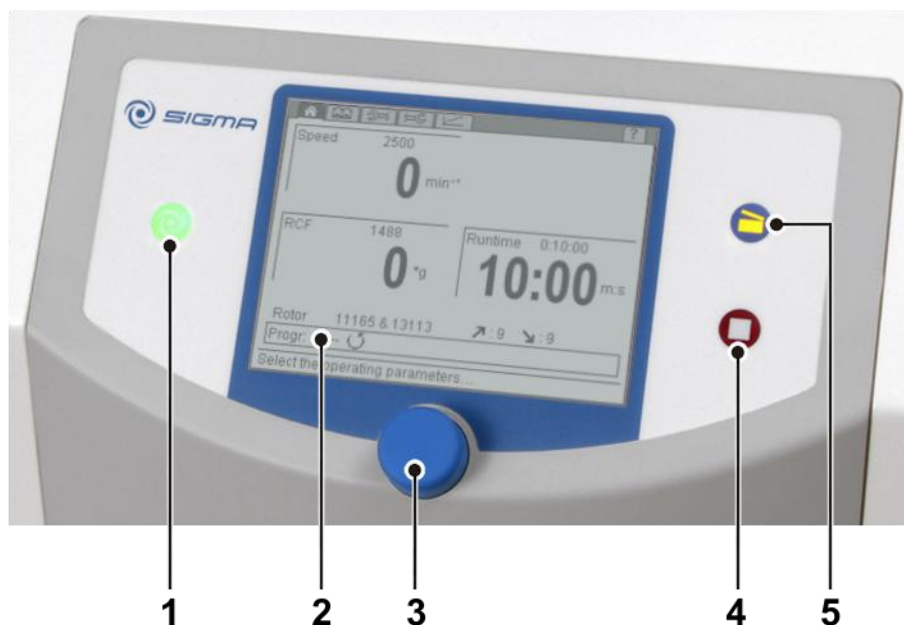


Fig. 15: User interface of the Spincontrol S control system

Display

The centrifuge display has the following fields:

- 1 Menu bar
- 2 Speed field
- 3 RCF field
- 4 Rotor field
- 5 Program field
- 6 Status bar
- 7 Runtime field
- 8 Acceleration curve
- 9 Deceleration curve

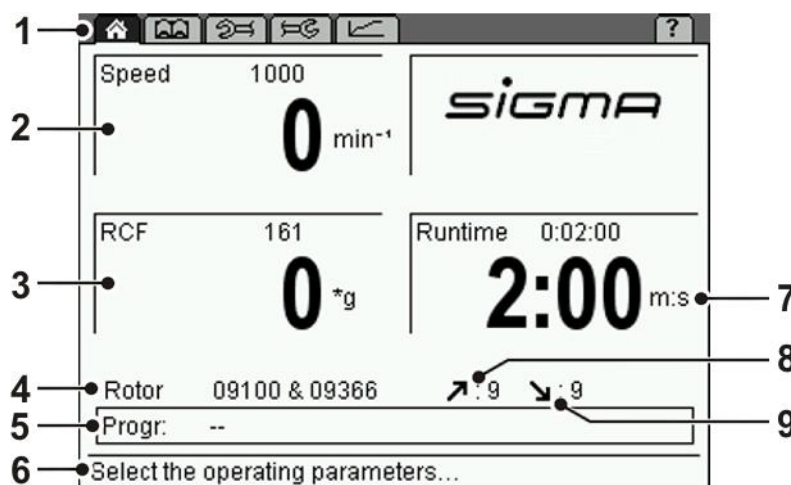


Fig. 16: Display of the Spincontrol S control system

6 Using the centrifuge

6.3.2 Manual mode

6.3.2.1 Starting a centrifugation run

The centrifuge is ready for operation when the start key is illuminated.

- Press the start key in order to start a centrifugation run.

6.3.2.2 Interrupting a centrifugation run

- Press the stop key in order to interrupt a centrifugation run. The centrifugation run will be terminated prematurely.

Quick stop

- Press the stop key for more than three seconds.

The centrifuge decelerates with the maximum deceleration curve.

After a quick stop, the centrifuge lid must be opened before a new centrifugation run can be started.

A quick stop can also be triggered during a normal deceleration, e.g. in order to speed up the deceleration.

When a quick stop is triggered, "Quick stop" will be displayed in the speed field.



A quick stop can be performed even if the centrifuge is blocked against unauthorised use.

6.3.2.3 Interrupting a deceleration process

- Press the start key during a deceleration process in order to interrupt it and to restart the centrifuge.

6.3.2.4 Selection, display, and modification of data

The "Standard" menu is displayed.

- Turn the function knob in order to select a field. The selected field is inverted.
- Press the function knob. The display starts to flash and the modification mode is active.
- Turn the function knob in order to modify the set value of the selected field.
- Press the function knob again to confirm the entry and to quit the modification mode.

6 Using the centrifuge

6.3.2.5 Standard menu

The "Standard" menu is symbolised by the icon "⌂" on the menu bar and it is displayed a few seconds after the centrifuge has been switched on. In this menu, the parameters of a centrifugation run can be displayed and modified.

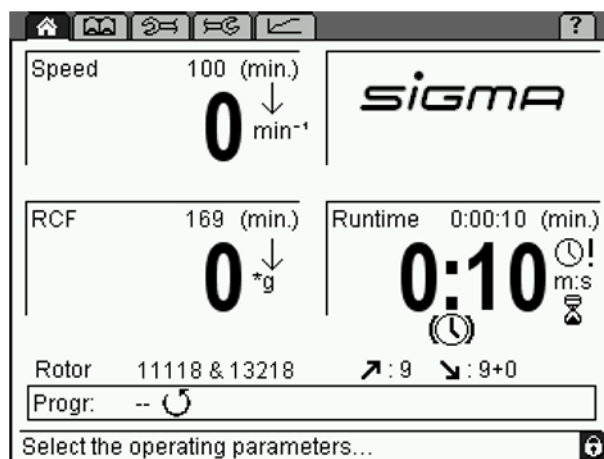


Fig. 17: Standard menu; here shown with all of the possible symbols

Speed

In the upper section of the field, the set speed of the centrifuge is displayed. The actual speed is displayed below this value. The values are stated in revolutions per minute ($\text{min}^{-1} = \text{rpm}$) and depend on the RCF values (see chapter 2.2.1.1 - "Speed, radius, and relative centrifugal force"). The maximum speed values depend on the rotor that is used.

Relative centrifugal force (RCF)

The relative centrifugal force is the acceleration that the sample is subjected to during the centrifugation run. The set value of this parameter is displayed in the upper section of this field, with the actual value shown below. The values are stated in g (gravitational acceleration) and they depend on the speed values (see chapter 2.2.1.1 - "Speed, radius, and relative centrifugal force"). The maximum RCF values depend on the rotor that is used.

Runtime

The set runtime is displayed in the upper section of this field, with the remaining runtime shown below. The runtime is defined as the period from the start of the centrifuge to the beginning of the deceleration phase. The maximum value is 99 h 59 min 59 sec.

In the "Setup" menu , it can be specified that the runtime is not to be started until the set speed is reached (see chapter 6.3.2.8 - "Setup menu"). In this case, the symbol "⌚!" appears in the runtime field.

6 Using the centrifuge

Continuous run

During the continuous run, the runtime of the centrifuge is unlimited and must be stopped manually. The centrifuge accelerates during the continuous run until the set speed is reached.

- Select the field "Runtime" and press the function knob. The display flashes when it is activated.
- Turn the function knob from the time 0:00:10 anti-clockwise or from the time 99:59:59 clockwise. "Infinite" will be displayed. After the start of the centrifuge, the elapsed time will be displayed.
- Deactivate the continuous run by pressing the stop button or by entering a specific runtime.

Short run

A short run can be started if no run is active.

- Keep the start button pressed during the short run.

During the short run, the centrifuge accelerates with acceleration curve 9 (maximum acceleration) until the maximum speed of the rotor is reached. The runtime is counted and in the speed field the message "Short run" flashes.

When the start button is released, the centrifuge decelerates with the maximum deceleration curve to a standstill.



The parameters speed, RCF, and runtime can be changed during the centrifugation.



If the centrifuge is locked with level 02 or higher, it is not possible to start a short run.

Rotor: rotor selection list

This field shows the rotor that is currently being used.

- Select the field "Rotor" and confirm the selection. A list with all of the possible rotors without buckets is displayed.
- Select the desired rotor.
 - If an angle rotor is selected, additional information concerning this rotor will be displayed.
 - In the case of swing-out rotors, a list with all of the possible rotor/bucket combinations will be displayed. Select an item from the list so that the additional information concerning the combination will be displayed.
- Press the function knob in order to accept the data.

6 Using the centrifuge

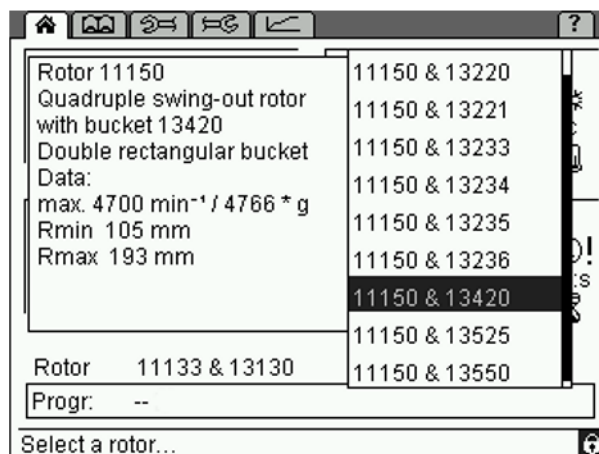


Fig. 18: Rotor selection list; here for a swing-out rotor

Automatic rotor identification

The centrifuge automatically identifies the rotor that is being currently used.

- If the system identifies a different rotor than the one that is set and if there are no different buckets for this rotor, the rotor input will be adapted automatically. The system will not display a message.
- If the system identifies a different rotor than the one that is set, and if there are different rotor/bucket combinations for this rotor, the system will automatically identify the correct rotor and select the rotor/bucket combination with the lowest speed. The system will display a corresponding message so that the combination can be adapted manually.
- If the system cannot identify the rotor, a message will be displayed. The rotor cannot be used in the centrifuge.

This prevents the maximum permissible speed from being exceeded.

Acceleration ↗

This function is used to select an acceleration curve. One can select a linear rise (curves 0-9) or a quadratic rise (curves 10-19). The acceleration curves 20-29 can be programmed as desired (see chapter 11.2 - "Acceleration and deceleration curves").

Deceleration (brake) ↘

This function is used to select a curve that decelerates the centrifuge to a standstill. Deceleration curves are inverted images of the acceleration curves and are labelled with identical numbers. Deceleration curve no. 0 represents a brakeless deceleration.

Progr.: program list

This field in the "Standard" menu shows the program that is currently loaded. When the field is selected, the program list is displayed (for information on how to work with the programs, please see chapter 6.3.3 - "Program mode").

6 Using the centrifuge

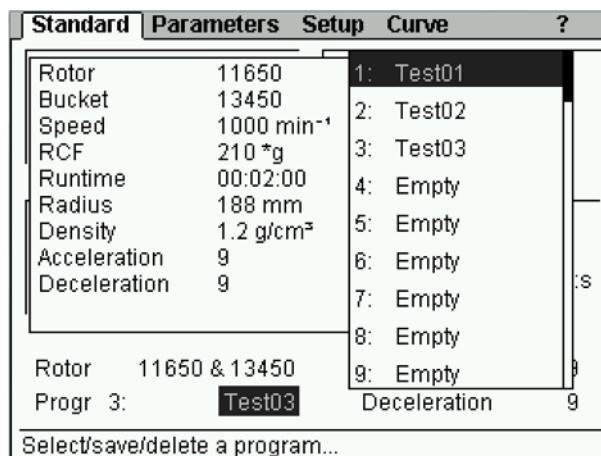


Fig. 19: Program list

Progress indicator

The progress indicator provides a quick overview of the remaining runtime of the running centrifugation run. For this purpose, a green progress bar and percentage value are displayed in the program field.

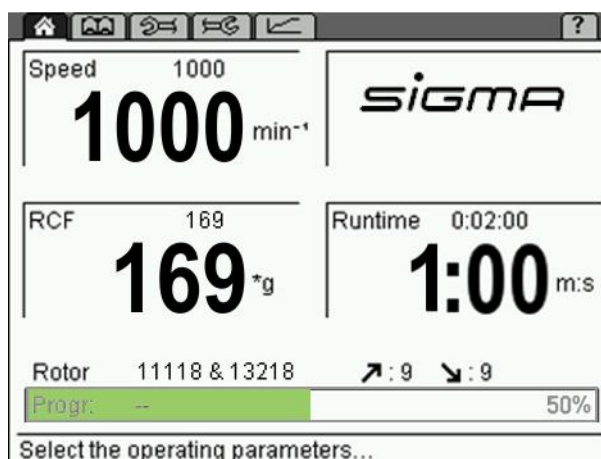


Fig. 20: Progress indicator during a centrifugation run

After the completion of the centrifugation run, the progress indicator remains at 100% until

- the lid is opened,
- a parameter is changed,
- a program is loaded, saved, or deleted, or
- a new centrifugation run is started.


6 Using the centrifuge

6.3.2.6 Process library menu

The process library can be opened via the "📖" symbol on the menu bar. It provides the user with an overview of all of the stored programs as processes.

The processes are listed with their name, RCF, and runtime.

The order of the processes corresponds to their storage locations on the program list. Empty storage locations will not be displayed. If more than 11 programs have been stored, the user can scroll through the list.

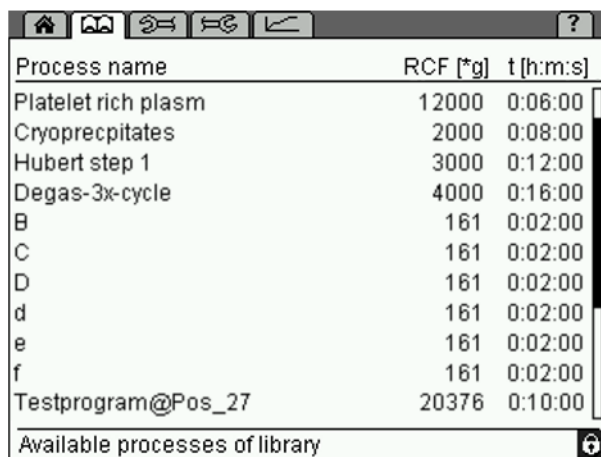

Process name	RCF [*g]	t [h:m:s]
Platelet rich plasm	12000	0:06:00
Cryoprecipitates	2000	0:08:00
Hubert step 1	3000	0:12:00
Degas-3x-cycle	4000	0:16:00
B	161	0:02:00
C	161	0:02:00
D	161	0:02:00
d	161	0:02:00
e	161	0:02:00
f	161	0:02:00
Testprogram@Pos_27	20376	0:10:00

Available processes of library

Fig. 21: Process library menu

Loading a process

- In the process library, select the desired process by turning the function knob. Press the function knob in order to confirm the selection.

The process will be loaded and the "Standard" menu 🏠 will be displayed.

Starting a process

- In the process library, select the desired process by turning the function knob. Then, press the start button.

The process will be loaded and started. The "Standard" menu 🏠 will be displayed.

6 Using the centrifuge

6.3.2.7 Parameters menu

The "Parameters" menu is symbolised by the "≡" symbol on the menu bar. It is used to specify various conditions for the centrifugation. These conditions are used to monitor the process and to control access to the centrifuge.

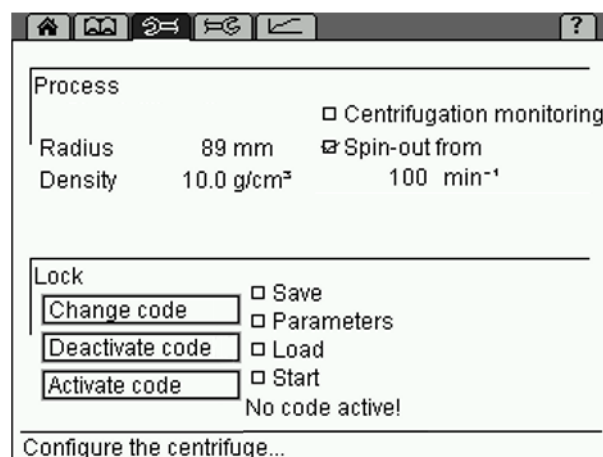


Fig. 22: Menu "Parameters"

Process

Radius

The radius determines the values of the relative centrifugal force (RCF) that the sample is subjected to. Normally, the maximum RCF value is displayed. If the value is reduced manually, a downward facing arrow "↓" will be displayed in the RCF field.

Density

This setting is useful for glass vessels. If the density of the liquid to be centrifuged is higher than 1.2 g/cm³, the value must be adapted manually in order to prevent the glass vessel from breaking. This will reduce the maximum possible final speed (see chapter 2.2.1.2 - "Density"). The reduction will be represented by a downward facing arrow "↓" in the speed field. Values between 1.2 and 10.0 g/cm³ are possible.

Centrifugation monitoring

The centrifugation monitoring function enables the continuous monitoring of the speed and runtime parameters during the centrifugation.

- Activate the centrifugation monitoring function by clicking.



If the function is activated during a centrifugation run, the monitoring process will not be started until the start of the next centrifugation run.

The centrifugation monitoring function compares the speed values of the current run with the reference values that are stored in the control unit. After every run, it issues a corresponding message.

6 Using the centrifuge

The runtime is considered faulty if the centrifugation run had to be stopped prematurely.

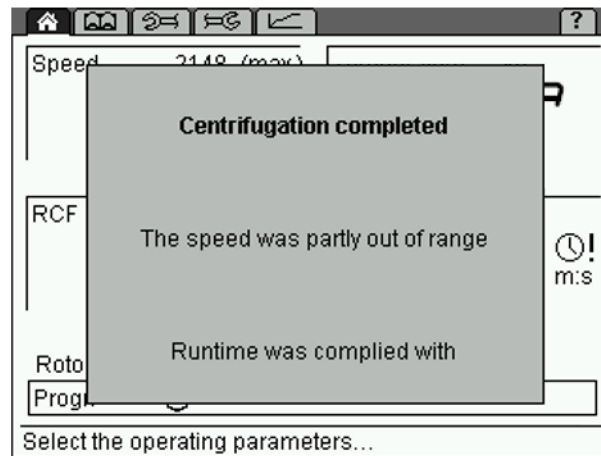


Fig. 23: Example of a centrifugation monitoring message



If the centrifugation monitoring system is used in combination with free acceleration or deceleration curves, unjustified error messages may result in certain cases.

Spin-out from...

If this function is active, the drive will be disconnected if the actual speed is below the set speed. As a result, the rotor decelerates and stops in a brakeless manner.



A spin-out, in particular with heavy rotors and at high speeds, can take a lot of time! (Depending on the rotor and load, the speed will be reduced by approximately 0.5 to 1 rpm per second.)

If the spin-out is active, "+0" is displayed next to the deceleration curve.

- The spin-out can be interrupted by a quick stop or by restarting the centrifuge.

Lock

In order to prevent any unauthorised use of the centrifuge, the following functions can be blocked:

- Saving of programs (level 01)
- Changing of parameters, short run (level 02)
- Loading of programs (level 03)
- Start button (level 04)

6 Using the centrifuge

Blocking a function

- Select the function that is to be blocked. The lower levels will also be automatically selected (if, for example, the "Parameters" function is selected, the "Save" function will also be selected).
- Select the button "Activate code".
- Enter a four-digit code and confirm the entry.

The lock is now active. The symbol "🔒" will be displayed in the status line and the lock level will be indicated.

If changes are made after a function has been blocked, the system will ask for the code prior to executing the change.

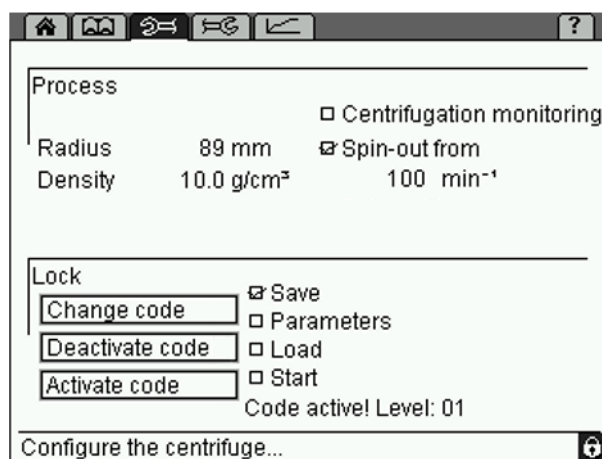


Fig. 24: "Parameters" menu with an active lock (level 01)

Unblocking a function

- Select the button "Deactivate code".
- Enter the code and confirm the entry.

The lock is now deactivated.

Changing the code

- Select the button "Change code".
- Enter the old code and confirm the entry.
- Enter the new code.
- For safety reasons, the code must be entered a second time.

The code is now changed.

6 Using the centrifuge

6.3.2.8 Setup menu

The "Setup" menu is symbolised by the "⚙️" symbol on the menu bar. It is used to perform basic settings concerning the control system of the centrifuge. It enables the optimum adaptation of the centrifuge to its specific area of application.

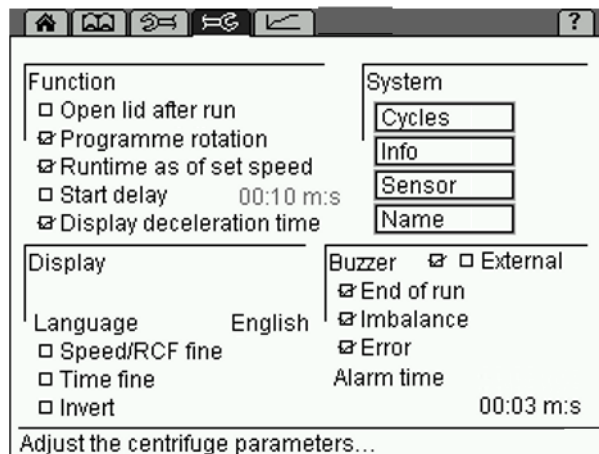


Fig. 25: Menu "Setup"

Function

Open lid after run

The automatic lid opening function ensures that the lid opens when the rotor has stopped.

Program rotation

See chapter 6.3.3.5 - "Automatic program rotation"

Runtime as of set speed

If this function is active, the runtime will not be measured until the set speed is reached. In the "Standard" menu, the symbol "⌚!" will be displayed in the runtime field.



In the program mode, this function can be saved separately for every program. In this case, the symbol "⌚!" will be displayed under the runtime of the detailed program description.

Start delay

If the start delay function is active, the centrifuge will not start until the pre-set time has elapsed. The symbol "⌚" will be displayed in the runtime field.

Display deceleration time

In the activated mode, the deceleration time will be displayed instead of the runtime during and after the deceleration process. Below the time display, the symbol "⌚" will be displayed. During a deceleration process, the symbol flashes. Once the deceleration is complete, it is displayed in a permanent manner.

6 Using the centrifuge

Display

Language

The control system can be used in various language versions.

If a language is selected by mistake, it can be changed on any screen as follows:

- Press and hold the stop button.
- Turn the function knob one notch to the left and then one notch to the right.
- Release the stop button. The "Language" window will be displayed.
- Select the desired language.

Speed/RCF fine and Time fine

This menu item can be used to preselect the set speed in steps of 1 rpm (instead of 100 rpm), the RCF value in steps of 1 x g (instead of 10 x g), and the set time in steps of 1 min or 1 sec (instead of 10 min or 10 sec).



Regardless of the fine adjustment, the step size increases when the function knob is turned quickly.

NOTE

Invert

If this function is activated, the display switches from the standard setting with a bright background and dark writing to a dark background with bright writing.

System

Cycles

This field shows the number of cycles as well as the runtime of the rotor and buckets that are used.



The information engraved on the rotor or bucket regarding the service life is binding. The information provided by the control system is for informational purposes only and, thereby, non-binding.

NOTE

Info

This menu provides information on the software versions that are used in this centrifuge.

Sensor

The sensor mode is reserved for service personnel.

6 Using the centrifuge



Values can neither be entered nor changed in the menus "Cycles", "Info", and "Sensor".

Name

In this field, an identification will be assigned to the centrifuge.

- The letters and characters can be entered when the cursor flashes in the text field. Turn the function knob in order to select a character and press it to confirm the selection. Then, press the knob again in order to enter the next character. Pressing the arrow button ← will delete the last character. The maximum number of characters is 19.
- When the name is complete, select the option "Accept" and confirm it.

Buzzer (signal)

With this function, a warning sound signal can be selected for

- the end of a centrifugation run,
- an imbalance message,
- an error message.

The duration of the warning signal can be specified.

External

This function is only available if the centrifuge is equipped with the option for the input and output of data (external signal, floating switch) (see chapter 6.3.4 - "Options for data input and output").

6 Using the centrifuge

6.3.2.9 Curve menu

This menu is used to create and edit customised acceleration and deceleration curves (see chapter 11.2 - "Acceleration and deceleration curves"). It is symbolised by the "⏮" symbol on the menu bar.

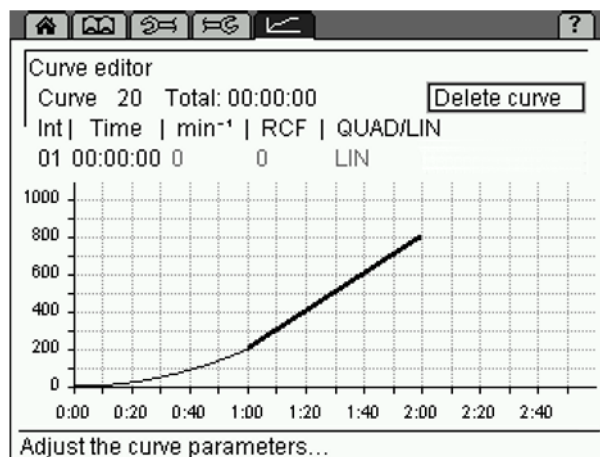


Fig. 26: Menu "Curve"

Creating or changing an acceleration curve



During a centrifugation run, curves can only be displayed. They cannot be changed or deleted.

- Select the "Curve" menu. The curve editor will be displayed.
- Select a curve number between 20 and 29 in the input field "Curve". If the curve number is already used, the stored curve will be displayed.
- The input field "Int" is used to specify the interval number of the process. Up to ten intervals can be entered for a curve.
- Enter the interval time of the current interval into the input field "Time". While doing so, certain restrictions must be taken into consideration (see below).
- Enter the desired acceleration under "min⁻¹" (rpm) or "RCF". While doing so, certain restrictions must be taken into consideration (see below). The values are interdependent.
- In the first interval, "QUAD/LIN" can be used to select a linear or quadratic rise. All of the other intervals are linear.

The field "Total" shows the total runtime of the process. The maximum total runtime of a curve depends on the slope of the curve and on the final speed of the rotor.



Only the last curve interval can be changed retroactively.

6 Using the centrifuge

Restrictions

- Acceleration and deceleration curves can include intervals with a positive slope as well as intervals with a negative slope and also intervals with 0 slope.
- The slope of the curve intervals can be $1 \text{ min}^{-1}/\text{sec}$ (rpm/sec) minimum and $1000 \text{ min}^{-1}/\text{sec}$ (rpm/sec) maximum.
- Quadratic curve intervals are only possible between 0 and 1000 min^{-1} (rpm) maximum. If a final speed $> 1000 \text{ min}^{-1}$ (rpm) is selected, this interval will automatically become linear above 1000 min^{-1} (rpm).
- The possible runtime results from the maximum possible speed (depending on the rotor) and from the slope limitation.

Example 1: Start speed 0 min^{-1} (rpm), final speed 100 min^{-1} (rpm), runtime 1 hour not possible, since the necessary slope $< 0.03 \text{ min}^{-1}$ (rpm), which means that it is beyond the defined range.

Example 2: Start speed 0 min^{-1} (rpm), final speed 15000 min^{-1} (rpm), runtime 10 sec. not possible, since the necessary slope is 1500 min^{-1} (rpm), which means that it is beyond the defined range.

6.3.2.10 Help menu

The help function is symbolised by the "?" symbol on the menu bar. It provides a short description of the control elements of the selected option.

Activating and deactivating the help function

- Select the question mark on the menu bar and press the function knob.
- Quit the help function by selecting the question mark and by pressing the function knob again.

Parameters can still be changed when the help function is activated.

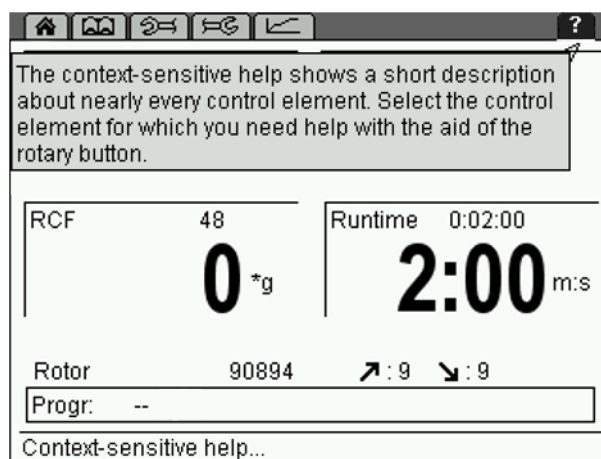


Fig. 27: Menu "Help"

6 Using the centrifuge

6.3.2.11 Changing the contrast

To change the contrast:

- Press and hold the stop button and turn the function knob one notch to the left. A dialog box will be displayed once the stop button is released.
- Adjust the contrast of the centrifuge display and confirm the change.

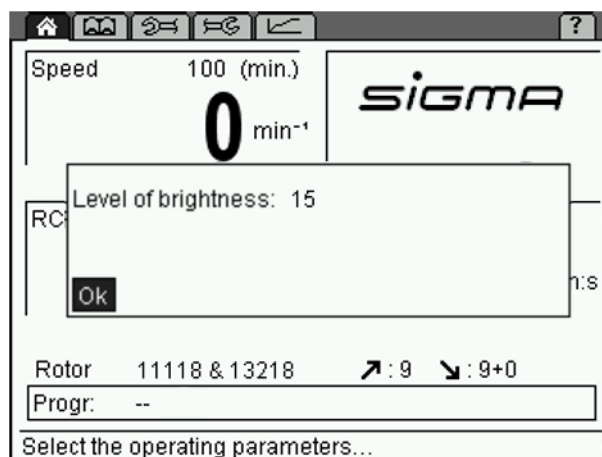


Fig. 28: Dialog box for changing the contrast

6.3.3 Program mode

A program contains all of the data that are required for a centrifuge run. Certain sedimentation results can be repeated under identical conditions. Programs can be loaded, executed, edited, and deleted when the centrifuge is at a standstill.

A maximum of 60 programs can be stored under the numbers 1 - 60.

"--" means that the values that are currently set are not a stored program.

The programs can be protected against unauthorised use, modification, or deletion with the aid of a code (see chapter 6.3.2.7 - "Parameters menu").

Stored programs are listed in the "Process library" menu .

6 Using the centrifuge

6.3.3.1 Saving a program

- Enter the parameters that are to be included in the program.
- Select the option "Progr" in the "Standard" menu  and confirm the selection. The program list will be displayed.
- Select a storage location from the program list.
- Save the program under the desired name. The letters and characters can be entered when the cursor flashes in the text field.
 - Turn the function knob in order to select a letter and press it to confirm the selection. Then, the next character can be selected.
 - Pressing the arrow button ← will delete the last character.
- When the program name is complete, select "OK" and confirm it.

The program will be saved and the "Standard" menu will be displayed.

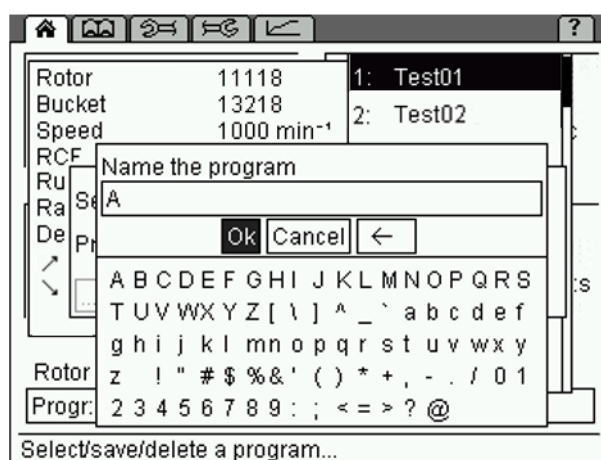


Fig. 29: Assignment of a program name prior to saving the program

6.3.3.2 Loading a program

- Select the option "Progr" from the "Standard" menu  and confirm the selection by pressing the function knob. The program list will be displayed.
- Select the desired program from the list and confirm the selection by pressing the function knob.

Or:

- Open the process library , select the desired program, and confirm the selection by pressing the function knob.

The program will be loaded and the "Standard" menu will be displayed.

6 Using the centrifuge

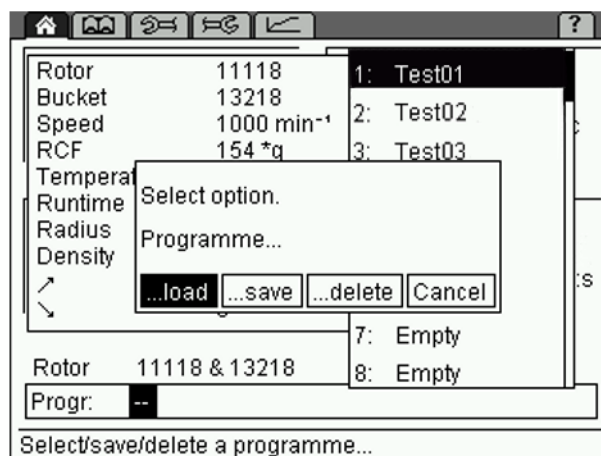


Fig. 30: Loading a program

6.3.3.3 Executing a program

- Select the option "Progr" from the "Standard" menu  and confirm the selection by pressing the function knob. The program list will be displayed.
- Select the desired program from the list and confirm the selection by pressing the function knob.
- Press the start button.

Or:

- Open the process library , select the desired program, and press the start button.

The program will be executed and the "Standard" menu will be displayed.

6.3.3.4 Deleting a program

- Select the option "Progr" in the "Standard" menu  and confirm the selection. The program list will be displayed.
- Select the program that is to be deleted.
- Select the option "Delete" and confirm it.

The program will be deleted and the "Standard" menu will be displayed.

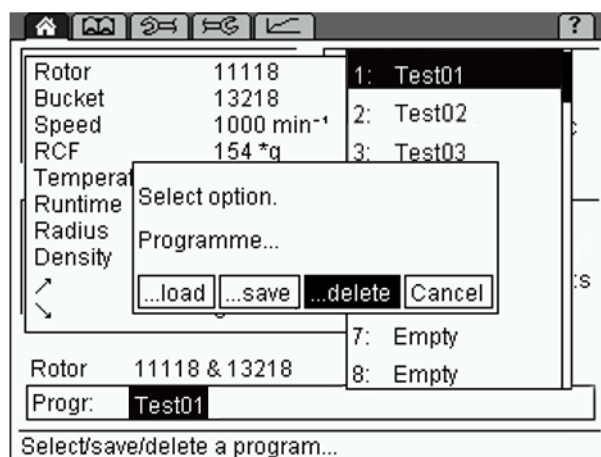


Fig. 31: Deleting a program

6 Using the centrifuge

6.3.3.5 Automatic program rotation

With the automatic program rotation, several programs can be executed directly one after the other.

- Activate the "Program rotation" function in the "Setup" menu .

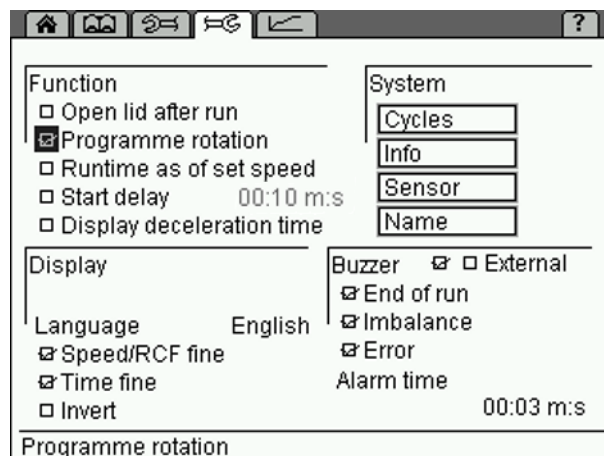


Fig. 32: Program rotation function

When a program is loaded while the program rotation function is active, this program will be used as the start program for the rotation. After the completion of the program, the next program on the program list will be loaded automatically. The rotation continues up to the next empty storage location and then restarts from the beginning (see the following illustration).

Example 1: Loading of Test04

Rotation: Test04, Test05, Test06, Test04,...

Example 2: Loading of Test05

Rotation: Test05, Test06, Test05,...

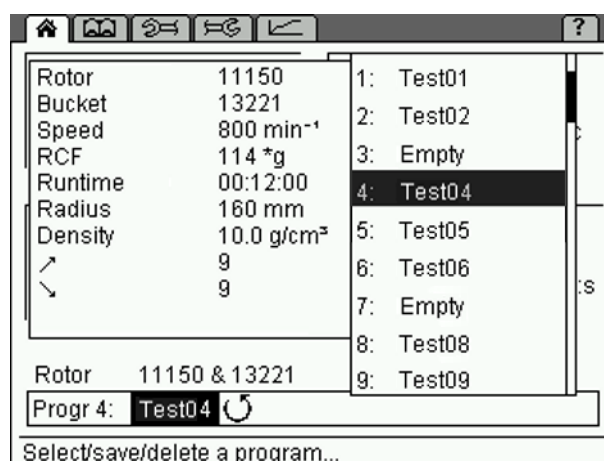


Fig. 33: Automatic program rotation

While the program rotation function is active, the arrow "↺" is displayed in the program line in the "Standard" menu.

6.3.4 Options for data input and output

- Connection for a serial interface (depending on the centrifuge type, partly standard).
- External signal active DC 24V, 0.5 A max. (part no. 17701)
- Floating switch AC 250V max., 6 A (part no. 17702)

6.4 Switching the centrifuge off

- Open the centrifuge when it is not in use so moisture can evaporate.
- Switch the centrifuge off by pressing the mains power switch.

7 Malfunctions and error correction

7 Malfunctions and error correction

7.1 General malfunctions

Malfunctions are indicated by a dialog box. If the acoustic signal is activated, it sounds when the error message is displayed.

- Eliminate the source of the problem (see table below).
- Acknowledge the error messages by pressing the lid key.



Error messages can be eliminated by pressing the lid key. The error itself will not be eliminated, but the centrifuge can be operated again.

Type of error	Possible reason	Correction
No indication on the display	No power in the mains supply	Check fuse in the mains supply
	Power cord is not plugged in	Plug in power cord correctly
	Fuses have tripped	Reactivate temperature fuse (see chapter 5.2.1 - "Type of connection")
	Mains power switch off	Switch mains power switch on
Centrifuge cannot be started: start key LED is not illuminated	Several possible causes	Power off/on. If the error occurs again, contact service
Centrifuge cannot be started: lid key LED flashes	The lid lock is not closed correctly	Open and close the lid. If the error occurs again, contact service
Centrifuge decelerates during operation	Brief mains power failure	Press start key in order to restart the centrifuge
	System error	Power off/on. If the error occurs again, contact service
Centrifuge decelerates during operation, imbalance dialog box is displayed	– Improper loading – Centrifuge is inclined – Drive problem – Centrifuge was moved during run	Balance load and restart the centrifuge. If the error occurs again, contact service
	– Ungreased load- bearing bolts	Clean and grease load- bearing bolts
Lid cannot be opened	Lid lock has not released	Unlock the lid manually (see chapter 7.1.1 - "Emergency lid release") and contact service
	Lid seal sticks	Clean the lid seal and apply talcum powder
Temperature value cannot be reached (only for refrigerated centrifuges)	Condenser dirty (only air-cooled units)	Clean the condenser. If the error occurs again, contact service
Hard running noise during the centrifugation	Screws of the transport safety device are not removed	Remove screws of the transport safety device (see chapter 4.3 - "Notes on transport")

7 Malfunctions and error correction

7.1.1 Emergency lid release

In the event of a power failure, it is possible to manually open the centrifuge lid.

- Switch off the mains power switch and disconnect the power cord from the socket.
- Remove the plugs (see figure, item 1) from the opening on the left and on the right side of the control panel, e.g. with a screwdriver.



Fig. 34: Position of the openings for the emergency lid release

- Insert the supplied hexagon socket key horizontally into the hole. The key will be guided through a funnel-shaped tube to the shaft of the lid lock motor.

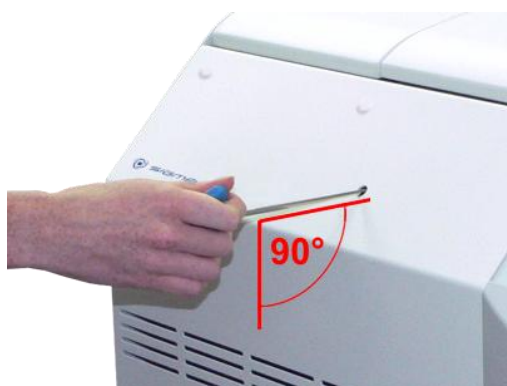


Fig. 35: The emergency lid release key must be inserted horizontally

- Unlock the motorised lid locks as follows:
 - Turn the left lid lock anti-clockwise.
 - Turn the right lid lock clockwise.
- Then, reinsert the plugs.



WARNING

Do not unlock or open the lid unless the rotor is at a standstill.

If the lid is opened via the emergency lid release system during a centrifuge run, the centrifuge will be switched off immediately and decelerate in an unbraked manner.

7 Malfunctions and error correction

7.2 Table of error codes

Error no.	Kind of error	Measures	Note
1-9	System error	- Allow to slow down - Power off/on	All these errors stop the centrifuge or cause it to decelerate brakeless
10-19	Speedometer error	- Allow to slow down - Power off/on	
20-29	Motor error	- Power off - Ensure ventilation	
30-39	EEPROM error	- Allow to slow down - Power off/on	With error 34, 35, and 36, the centrifuge will stop; with error 37 and 38 only an error message will be given
40-45	Temperature error (only for refrigerated centrifuges)	- Allow to slow down - Power off - Allow to cool down - Provide better ventilation (only air-cooled centrifuges) - Provide sufficient water throughput (only water-cooled centrifuges)	
46-49	Imbalance error (only for centrifuges with imbalance monitoring system)	- Allow to slow down - Power off - Eliminate the imbalance	
50-59	Lid error	- Press lid key - Close lid - Remove foreign matter from the opening of the lid lock device	With error 50 and 51, the centrifuge will stop
60-69	Process error	- Allow to slow down - Power off/on	With error 60, the message "power failure during run" will be displayed, with error 61, the message "stop after power on" will be displayed
70-79	Communication error	- Allow to slow down - Power off/on	
80-89	Parameter error	- Power off - Allow to cool down - Provide for better ventilation	With error 83, error message only
90-99	Other errors	- Check connections - Provide sufficient water throughput (only water-cooled centrifuges)	



If it is impossible to eliminate the errors, contact the service!

7.3 Service contact

In the event of queries, malfunctions, or spare part enquiries:

From Germany:

Contact

Sigma Laborzentrifugen GmbH
An der Unteren Söse 50
37520 Osterode (Germany)
Tel. +49 (0) 55 22 / 50 07-44 44
E-mail: support@sigma-zentrifugen.de

Outside Germany:

Contact our agency in your country. All agencies are listed at
www.sigma-zentrifugen.de → [Sales Partners]



- If you would like to utilise our service, please state the type of your centrifuge and its serial number.

8 Maintenance and service

8 Maintenance and service

The centrifuge, rotor, and accessories are subject to high mechanical stress. Thorough maintenance performed by the user extends the service life and prevents premature failure.



CAUTION

If corrosion or other damage occurs due to improper care, the manufacturer cannot be held liable or subject to any warranty claims.

8.1 General notes

The centrifuge and its accessories are subject to high mechanical stress. In order to withstand the continuous, extreme levels of stress within their area of application with high centrifugal force, the products are manufactured with the highest levels of precision. The service life of the centrifuge depends on its cleaning, maintenance, correct storage, frequency of use and other specific conditions of use. Diligent care and maintenance by the user will extend the service life of the products and prevent premature wear.

The centrifuge and accessories are made from different materials. Please check for any possible incompatibilities.

Prior to using cleaning agents and disinfectants other than those recommended, the user must consult with the manufacturer to ensure that they will not damage the centrifuge or its accessories.

For autoclaving, the long-term temperature resistance of the individual materials must be observed (see chapter 8.1.3 - "Autoclaving").

Protect the centrifuge and its accessories from intensive UV radiation (e.g. direct exposure to sunlight) and thermal stress (e.g. caused by heat generators).



DANGER

When using hazardous substances (e.g. infectious or pathogenic substances), it is mandatory to disinfect the centrifuge and its accessories.



WARNING

For your own protection, always take appropriate precautions if there is a risk of toxic, radioactive or pathogenic contamination. Wear personal protective equipment.

8 Maintenance and service

Frequency of the required activities

Chapter	Activity	Frequency
8.1.1	Cleaning	
8.1.1.1	Cleaning of the centrifuge	As required
8.1.1.2	Cleaning of rotors and buckets	As required
8.1.1.3	Cleaning of adapters	As required
8.1.2	Disinfection	
8.1.2.1	Disinfection of the centrifuge	As required
8.1.2.2	Disinfection of rotors and buckets	As required
8.1.2.3	Disinfection of adapters	As required
8.2	Maintenance tasks	
8.2.1	Inspection of the centrifuge for signs of damage	Monthly
8.2.1	Greasing of the motor shaft	After cleaning
8.2.2	Check of the rotors, buckets and adapters for signs of damage	Monthly
8.2.2	Greasing of the rotor tie-down screw	After cleaning
8.2.2.2	Greasing of the load-bearing bolts	As required/ after cleaning



If corrosion or other damage occurs due to lack of care, the manufacturer cannot be held liable and no warranty claims can be made.

8.1.1 Cleaning

- Use soapy water or other water-soluble, mild cleaning agents with a pH value between 6 and 8 for cleaning the centrifuge and its accessories.
- Clean the centrifuge by hand only.
- The water temperature must be < 40°C.
- Avoid corrosive and aggressive cleaning agents.
- Do not use any solvents.
- Do not use any agents with abrasive substances.



Ensure that no liquid penetrates the inside of the centrifuge. Otherwise, electrical or mechanical components can be damaged.

8 Maintenance and service

8.1.1.1 Cleaning of the centrifuge

1. Open the lid.
2. Disconnect the centrifuge from the power supply.
3. Remove the rotors and buckets (see chapter 6.2.2 - "Installation of rotors and accessories").
4. Clean the housing and the rotor chamber with the specified cleaning agents and a soft, lint-free cloth.
5. Remove any cleaning agent residues completely with water and a soft, lint-free cloth.
6. After cleaning, dry the surfaces with a soft, lint-free and absorbent cloth.

8.1.1.2 Cleaning of rotors and buckets



CAUTION

Do not clean the accessories in a dishwasher!

The dishwasher will remove the anodised coating of the adapters and carriers, thereby leading to cracks in areas under stress.

1. Clean the rotor and buckets with the specified cleaning agents and a soft, lint-free cloth.
2. To clean narrow holes, use a test tube brush with a wool tip and then rinse the holes with water.
3. Remove any cleaning agent residues completely with water and a soft, lint-free cloth.
4. Dry the rotor and buckets with a soft, lint-free and absorbent cloth or in a drying cabinet at 50°C maximum.
5. To dry the holes, let the rotor dry lying upside down on a surface.

8.1.1.3 Cleaning of adapters

1. Clean the adapters with the specified cleaning agents and a soft, lint-free cloth.
2. To clean narrow holes, use a test tube brush with a wool tip and then rinse the holes with water.
3. Remove any cleaning agent residues completely with water and a soft, lint-free cloth.
4. Dry the adapters with a soft, lint-free and absorbent cloth.
5. To dry the holes, let the adapter dry lying upside down on a surface.

8.1.2 Disinfection

Prior to a disinfection, the cleaning steps described in chapter 8.1.1 - "Cleaning" must have been performed.

8.1.2.1 Disinfection of the centrifuge

- To disinfect the centrifuge (control panel, contact surfaces and rotor chamber), use an alcohol-based surface disinfectant with a concentration between 30% and 75% (e.g. Bacillol® AF).
 - Disinfect the lid window and the display (acrylic glass), motor cover (silicone) and (screen-printed) labelling with a disinfectant that is suitable for the material (e.g. Bacillol® 30 sensitive tissues). Never use a surface disinfectant with a high alcohol concentration!
 - Follow the manufacturer's instructions regarding the application time and concentration of the disinfectant.
1. Open the lid.
 2. Disconnect the centrifuge from the power supply.
 3. Remove the rotors and buckets (see chapter 6.2.2 - "Installation of rotors and accessories").
 4. Wipe down the housing and rotor chamber with the above-mentioned disinfectant and a soft, lint-free cloth.
 5. Comply with the specified application time.
 6. Let the disinfected surfaces dry completely.

8.1.2.2 Disinfection of rotors and buckets

- To disinfect the rotors and buckets, use an alcohol-based surface disinfectant with a concentration between 30% and 75% (e.g. Bacillol® AF).
 - Follow the manufacturer's instructions regarding the application time and concentration of the disinfectant.
1. Wipe down the rotor and buckets with the above-mentioned disinfectant and a soft, lint-free cloth.
 2. To disinfect narrow holes, fill them with the disinfectant without entrapping any air.
 3. Ensure complete wetting.
 4. Comply with the specified application time.
 5. Let the disinfected rotor and buckets dry completely.

8 Maintenance and service

8.1.2.3 Disinfection of adapters

- To disinfect the adapters, use an alcohol-based surface disinfectant with a concentration between 30% and 75% (e.g. Bacillol® AF).
 - Follow the manufacturer's instructions regarding the application time and concentration of the disinfectant.
1. Wipe down the adapters with the above-mentioned disinfectant and a soft, lint-free cloth.
 2. To disinfect narrow holes, fill them with the disinfectant without entrapping any air.
 3. Ensure complete wetting.
 4. Comply with the specified application time.
 5. Let the disinfected adapters dry completely.

8 Maintenance and service

8.1.3 Autoclaving

The service life of the accessories primarily depends on the autoclaving frequency and their general frequency of use.

- Replace the accessories immediately if you notice any leaks or any changes in the colour or structure, etc.
- When autoclaving, in order to prevent the vessels from being deformed, ensure that the caps are not screwed onto the vessels.



NOTE

The possibility of plastic parts, such as covers or racks, becoming deformed during autoclaving cannot be ruled out.



NOTE

The specified autoclaving parameters exclusively refer to the material resistance and not to the degree of sterility.

Category	Type of accessory	Material abbreviation	121°C 20 min	Notes
Rotors and covers	Aluminium rotors	AL	yes	
	Polypropylene rotors	PP	no	
	Polycarbonate covers for angle rotors	PC	no	
	Polyallomer covers for angle rotors	PA	no	
	Polysulfone covers for angle rotors	PSU	yes	100 cycles max.
Buckets and caps	Aluminium buckets	AL	yes	
	Polyamide buckets	PA	no	13035, 13296, 13299
	Polyphenylsulfone caps	PPSU	yes	100 cycles max.
	Polysulfone caps	PSU	yes	100 cycles max.
Adapters	Polyallomer racks	PA	no	
	Polyallomer racks	PC	no	
	Polyallomer racks	PP	no	
Vessels	Stainless steel vessels and bottles	--	yes	
	Glass vessels	--	yes	
	Polyethylene vessels	PE	no	
	Polyflor vessels	PF	yes	100 cycles max.
	Polycarbonate vessels	PC	no	
	Polypropylene copolymer vessels	PPCO	yes	20 cycles max.
	Polystyrene vessels	PS	no	
Additional equipment	Stainless steel taring weights for blood bag systems	--	yes	

8 Maintenance and service

8.2 Maintenance tasks

8.2.1 Maintenance of the centrifuge

- Liquids, such as water, solvents, acids and lyes, must be removed from the rotor chamber immediately using a cloth. This prevents the motor bearings from being damaged.
- In the event of a contamination with toxic, radioactive or pathogenic substances, clean the inside of the centrifuge immediately with a suitable decontamination agent depending on the type of contamination (see chapter 8.1.1 - "Cleaning" and chapter 8.1.2 - "Disinfection").
- After every cleaning process, grease the motor shaft slightly with a small amount of heavy-duty grease for load-bearing bolts and distribute the grease with a cloth so that it forms a thin layer.
- Perform a visual inspection of the centrifuge at least monthly to detect the following signs of damage:
 - Cracks
 - Visible surface damage
 - Pressure marks
 - Signs of corrosion
 - Other changes

8.2.2 Maintenance of rotors, buckets and adapters



CAUTION

The special precautions for the care of the accessories must be strictly observed.

These are measures to ensure operational safety!

- Liquids that may cause corrosion must be immediately removed completely from the rotor, buckets and accessories by way of a cloth.
- In the event of a toxic, radioactive or pathogenic contamination of the rotors and accessories, clean them immediately with a suitable decontamination agent depending on the type of contamination. For your own protection, always take appropriate precautions if there is a risk of toxic, radioactive or pathogenic contamination (see chapter 8.1.1 - "Cleaning" and chapter 8.1.2 - "Disinfection").

Rotors, buckets and multiple carriers are manufactured with the highest levels of precision in order to withstand the continuous, extreme levels of stress with their area of application with high centrifugal force.

Chemical reactions and pressure corrosion (a combination of alternating pressure and chemical reactions) can attack or destroy the structure of metals. Barely detectable cracks on the surface expand and weaken the material without leaving any clearly visible signs.

8 Maintenance and service

- This is why the material must be checked regularly, at least monthly, for the following:
 - Cracks
 - Visible surface damage
 - Pressure marks
 - Signs of corrosion
 - Other changes
- Check the holes of rotors and multiple carriers.
- Damaged parts must be replaced immediately for your own safety.
- After every cleaning process, grease the rotor tie-down screw slightly with a small amount of heavy-duty grease for load-bearing bolts and distribute the grease with a cloth so that it forms a thin layer.
- Comply with the special requirements for the use of plastic accessories (see chapter 8.2.2.1 - "Plastic accessories").

8.2.2.1 Plastic accessories

The chemical resistance of plastic decreases with rising temperatures (e.g. during drying) (see chapter 11.4 - "Resistance data").

- If solvents, acids, or alkaline solutions have been used, clean the plastic accessories thoroughly.



WARNING

Plastic accessories must not be greased!

8.2.2.2 Load-bearing bolts

Rotors with Sigma "Comfort" rotor coating



NOTE

The load-bearing bolts of some types of swing-out rotors have an anti-friction coating. This coating prevents friction between the buckets and bolts.

The bolts do not require any greasing during the service life of the coating (see chapter 3.6.7.3 - "Service life of the Sigma "Comfort" rotor coating")!

All "Comfort" coated swing-out rotors have a "C" shown after their number (see chapter 3.6.7.1 - "Marking of rotors and accessories").



NOTE

A list of the available rotors with a "Comfort" rotor coating as well as information about the service life of the coating can be found in the "Sigma Comfort rotor coating" document that comes supplied with every rotor that has anti-friction coating.



NOTE

It is not permissible to combine "Comfort" coated rotors with buckets or carriers with a cleanroom coating.

8 Maintenance and service

Rotors without "Sigma Comfort" rotor coating

The load-bearing bolts of rotors that do not (or no longer) have anti-friction coating must be greased. This is the only way to ensure a uniform swing-out motion of the buckets and, therefore, the smooth operation of the centrifuge. Load-bearing bolts that are insufficiently greased may cause the centrifuge to stop due to an imbalance.



Fig. 36: Sufficient quantity of grease for one bolt

Greasing the load-bearing bolts

1. Clean the load-bearing bolts and bucket groove in order to remove the old grease.
2. Apply a small amount of heavy-duty grease for load-bearing bolts (see the picture above) to both load-bearing bolts of a bucket.
3. Install the bucket and swing it manually back and forth once in order to distribute the grease.
4. Repeat this process with all the other buckets.

8.2.3 Broken glass

In the event of glass breakage, all shards must be removed immediately and completely (e.g. with a vacuum cleaner). Rubber cushions must be replaced since even thorough cleaning cannot remove all the glass particles.



There is a risk of injury due to sharp shards when removing the broken glass from the motor chamber.

- Wear personal protective equipment (e.g. gloves).

Potential consequences of broken glass

- Broken glass damages the surface coating (e.g. anodised aluminium) of the buckets, thereby causing corrosion.
- Broken glass on the rubber cushions of the buckets causes more glass breakage.
- Broken glass in the swivel bearings of the load-bearing bolts prevents the buckets and multiple carriers from swinging out evenly, thereby causing imbalance.
- Broken glass in the rotor chamber causes metallic abrasion due to the strong air circulation. This fine metal dust not only strongly contaminates the rotor chamber, rotor and samples, but it also damages the surface of the accessories, rotors and rotor chamber.

Removing broken glass and metal dust from the rotor chamber

1. Apply a thick layer of petroleum jelly or similar to the upper third of the rotor chamber.
2. Then, let the rotor rotate for a few minutes at medium speed (approx. 2,000 rpm). During this process, dust and glass particles will be bound in the grease.
3. Then, remove the grease together with the dust and glass particles by way of a cloth.
4. If necessary, repeat this process.

8 Maintenance and service

8.3 Maintenance and service



DANGER

In the case of maintenance or service requiring the removal of the housing panels, there is a risk of electric shock or mechanical injury.

- Maintenance and service tasks must only be carried out by authorised and specialised personnel.
- Following the completion of any type of maintenance or service, the authorised and specialised personnel must perform final inspection and testing in compliance with the relevant standards (see below).

Centrifuges used in Germany are subject to mandatory testing as per the information sheet of the German Berufsgenossenschaft Rohstoffe und chemische Industrie FBRCI-025 (November 2023). The test intervals must ensure that the centrifuge can be used safely between tests. The operator is responsible for specifying the test intervals.

When specifying the test intervals, the following time frames should not be exceeded:

- Every 1-2 years in operating condition
- Every 3-4 years in dismantled condition

We recommend the following:

- The motor damping elements (rubber parts of the motor suspension) should be replaced every 3 to 4 years or every 10,000 cycles.

Sigma Laborzentrifugen GmbH performs the inspection and maintenance of centrifuges as per FBRCI-025 and in accordance with the German DGUV regulation 3 (DGUV = German statutory accident insurance) (electrical systems and equipment).

These services can be agreed as part of a maintenance contract.

www.sigma-zentrifugen.de → [Service] → [Overhaul and repair]

Scan QR-Code

If you have any questions to this subject, please contact our service department:



8.4 Return of defective centrifuges or parts

Although we exercise great care during the production of our products, it may be necessary to return a unit or accessory to the manufacturer. In order to ensure the quick and economical processing of returns of centrifuges, spare parts, or accessories, we require complete and extensive information concerning the process. Please fill in the following forms completely, sign them, enclose them with the return package, and send them together with the product to:

Sigma Laborzentrifugen GmbH
An der Unteren Söse 50
37520 Osterode (Germany)

1. Declaration of decontamination

As a certified company and due to the legal regulations for the protection of our employees and of the environment, we are obliged to certify the harmlessness of all incoming goods. For this purpose, we require a declaration of decontamination.

- The form must be filled in completely and signed by authorised and specialised personnel only.
- Affix the original form in a clearly visible manner to the outside of the packaging.



We will return the part/unit if no declaration of decontamination is provided!

2. Form for the return of defective parts

This form is for the product-related data. They facilitate the assignment, and they enable the quick processing of the return. If several parts are returned together in one packaging, please enclose a separate problem description for every defective part.

- A detailed problem description is necessary in order to perform the repair quickly and economically.



If the form does not include a description of the malfunction, neither a refund nor a credit note can be issued. In this case, we reserve the right to return the part/unit to you at your expense.

- Upon request, we will prepare and submit to you a cost estimate prior to performing the repair. Please confirm such cost estimate within 14 days. If the cost estimate has still not been confirmed after 4 weeks, we will return the defective part/unit. Please note that you must bear the incurred costs.

8 Maintenance and service

3. Pick-up order

Upon request, we will commission a forwarding agent to collect the unit. In this case, complete the pick-up order and return it to us by e-mail.



The defective part/unit must be packaged in a transport-safe manner. Please use the original packaging for the unit, if at all possible.

If the product is dispatched to us in unsuitable packaging, you will be charged the cost for returning it to you in new packaging.

The forms can be downloaded online from www.sigma-zentrifugen.de → [Service] → [Overhaul and repair].

9 Disposal

9.1 Disposal of the centrifuge



In accordance with the directive 2012/19/EU, SIGMA centrifuges are marked with the symbol shown to the left. This symbol means that it is not permissible to dispose of the unit among household waste.

- You can return these centrifuges free of cost to Sigma Laborzentrifugen GmbH.
- Ensure that the unit is decontaminated. Fill in a declaration of decontamination (see chapter 8.4 - "Return of defective centrifuges or parts").
- Comply with any other applicable local rules and regulations.

9.2 Disposal of the packaging

- Use the packaging to return the centrifuge for disposal or
- dispose of the packaging, after having separated the individual materials.
- Comply with all local rules and regulations.

10 Technical data

10 Technical data

Manufacturer	Sigma Laborzentrifugen GmbH An der Unteren Söse 50 37520 Osterode (Germany)
Type:	6-16S
<u>Connection requirements</u> Electr. connection: Protection class: IP code: Power consumption (W): Input fuse (AT):	see name plate I 20 1 300 16.0 (at 220-240 V / 50 Hz), temperature fuses
<u>Performance data</u> Max. speed (rpm): Max. capacity (ml): Max. gravitational field (x g): Max. kin. energy (Nm):	13 500 3 200 20 376 105 901
<u>Other parameters</u> Time range: Storage locations:	10 sec – 99 h 59 min, short run, continuous run 60
<u>Physical data</u> Height (mm): Height with open lid (mm): Width (mm): Depth (mm): Weight (kg): Noise level (dB(A)):	483 990 581 711 116 < 73 (at maximum speed)

10.1 Ambient conditions

- The figures are valid for an ambient temperature of +23°C and a nominal voltage $\pm 10\%$.



At a nominal voltage of 100V or 200V, a tolerance of +10% / -5% applies.

- For indoor use only.
- Maximum altitude 2,000 m above sea level.
- Allowable ambient temperature +5°C to +40°C.
- Max. allowable relative humidity of air 80% from 5°C up to 31°C with a linear decrease to 50% relative humidity of air at 40°C.
- Pollution degree 2.

10.2 Technical documentation

For environmental reasons, the comprehensive technical documentation of the centrifuge (e.g. circuit diagrams) and the safety data sheets of the manufacturers of refrigerants and lubricants are not attached to this documentation.

You can order these documents from our service department.

11 Appendix

11 Appendix

11.1 Range of accessories

The complete list of accessories can be downloaded from www.sigma-zentrifugen.de.



Some accessories come supplied together with a data sheet that includes important information and notes on safety. This data sheet must be added to the operating manual.

11.1.1 Rotor radii

The information in the accessories table concerning the radius refers to the values of the respective rotor as shown below. The radius calculation is described in chapter 2.2.1.1 - "Speed, radius, and relative centrifugal force".

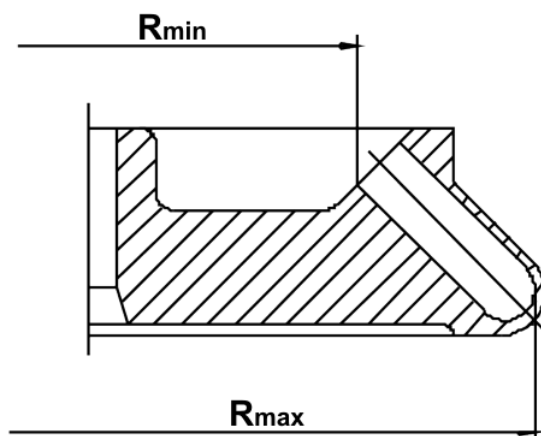


Fig. 37: Minimum and maximum radius of an angle rotor

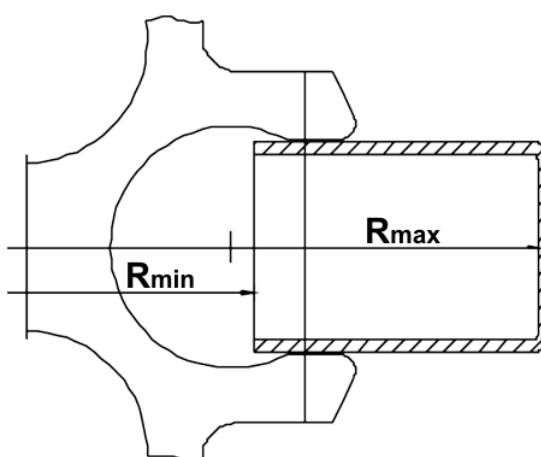


Fig. 38: Minimum and maximum radius of a swing-out rotor

11.2 Acceleration and deceleration curves

Linear as well as quadratic curves are numbered in the direction of increasing acceleration (from right to left).

The deceleration curves are inverted images of the acceleration curves and are assigned the same numbers. An exception is curve 0. It decelerates brakeless (spin-out).

In general, the runtime, until the set speed is reached, depends on the moment of inertia of the rotor.

Linear curves

The slope of the fixed acceleration curves defines the time that is required to accelerate the rotor by 1,000 rpm.

Curve 9 is a special case compared to the other curves. The centrifuge accelerates with maximum power. The runtime, until the set speed is reached, depends solely on the moment of inertia of the rotor.

Linear curve no.	Slope
0	4 [rpm/sec]
1	6 [rpm/sec]
2	8 [rpm/sec]
3	17 [rpm/sec]
4	25 [rpm/sec]
5	33 [rpm/sec]
6	50 [rpm/sec]
7	100 [rpm/sec]
8	200 [rpm/sec]
9	1.000 [rpm/sec]

Fig. 39: Slope of linear curves

Quadratic curves

Curve 19 is a special case compared to the other curves. The centrifuge accelerates with maximum power. The runtime depends solely on the moment of inertia of the rotor.

Quadratic curve no.	Time until 1,000 rpm	Slope as of 1,000 rpm
10	500 sec	4 [rpm/sec]
11	333 sec	6 [rpm/sec]
12	250 sec	8 [rpm/sec]
13	118 sec	17 [rpm/sec]
14	80 sec	25 [rpm/sec]
15	60 sec	33 [rpm/sec]
16	40 sec	50 [rpm/sec]
17	20 sec	100 [rpm/sec]
18	10 sec	200 [rpm/sec]
19	2 sec	1.000 [rpm/sec]

Fig. 40: Slope of quadratic curves

11 Appendix

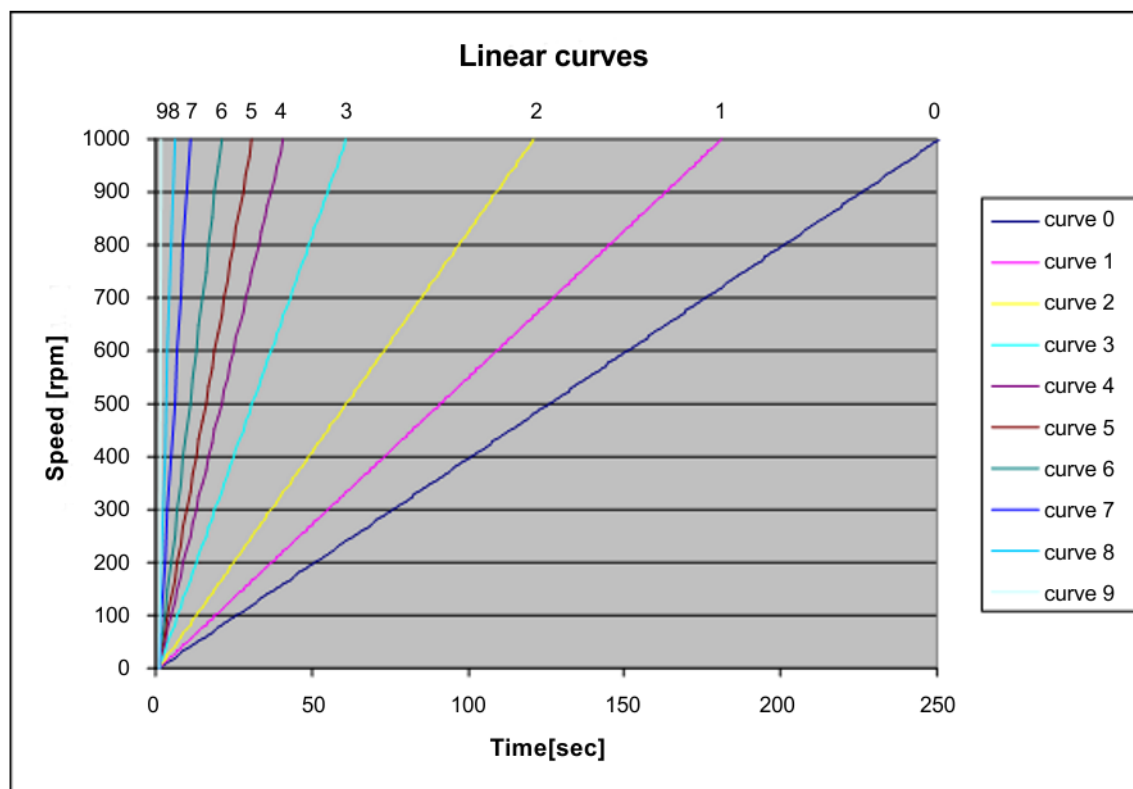


Fig. 41: Diagram of linear curves

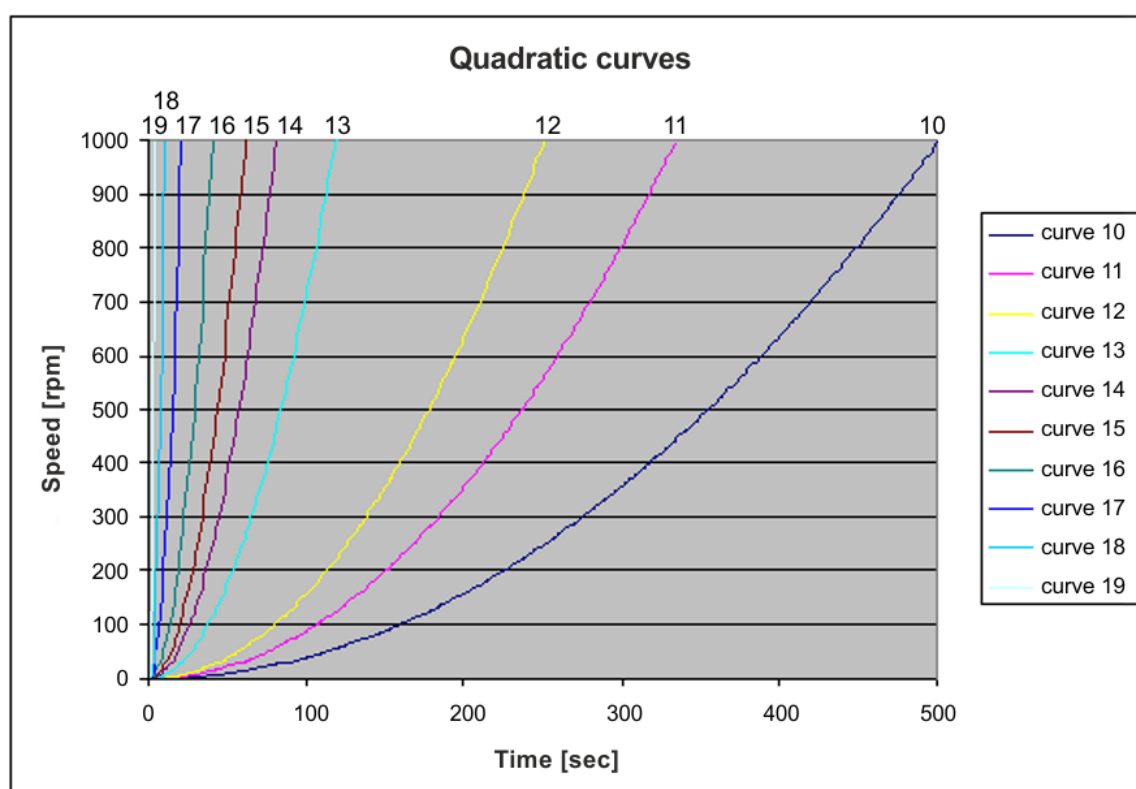


Fig. 42: Diagram of quadratic curves

11.3 Table of the service life of rotors and accessories

- The rotors and accessories must be put out of service after 10 years. Any use after this period may be permissible in individual cases after an inspection performed by the manufacturer.
- If a specification concerning the maximum number of cycles and a specification concerning the service life (i.e. a date) are provided, the specification that occurs first shall apply.
- After 50,000 cycles, rotors must be scrapped for safety reasons.

Rotor / bucket	Cycles	Service life ("Exp.Date")	Suitable for centrifuge	Remarks
12082		7 years	1-14, 1-14K	
12084		7 years	1-14, 1-14K	
12092		5 years	1-14, 1-14K	
12094		5 years	1-14, 1-14K	
12096		5 years	1-14, 1-14K	
12097		5 years	1-14, 1-14K	
12134		5 years	1-16, 1-16K	
12135		5 years	1-16, 1-16K	
12137		5 years	1-16, 1-16K	
91429		7 years	1-7	
13035			2-7	Do not grease the load-bearing bolts of the rotor
13296	35,000	5 years	2-7, 2-16KL, 2-16KHL	Do not grease the load-bearing bolts of the rotor
13299		5 years	2-7, 2-7 Cyto, 2-16KL, 2-16KHL, 3-30KS, 3-30KHS	Do not grease the load-bearing bolts of the rotor
13221	10,000		4-16S, 4-16KS, 4-16KHS, 6-16S, 6-16HS, 6-16KS, 6-16KHS	
9100	35,000		4-5KL, 4-16S, 4-16KS, 4-16KHS, 6-16S, 6-16HS, 6-16KS, 6-16KHS	
9366	15,000		4-5KL, 4-16S, 4-16KS, 4-16KHS, 6-16S, 6-16HS, 6-16KS, 6-16KHS	
13218	20,000		4-5KL, 4-16S, 4-16KS, 4-16KHS, 6-16S, 6-16HS, 6-16KS, 6-16KHS	
12600		7 years	6-16S, 6-16HS, 6-16KS, 6-16KHS	
13625	20,000		6-16S, 6-16HS, 6-16KS, 6-16KHS	
13635	25,000		6-16S, 6-16HS, 6-16KS, 6-16KHS	
13650	20,000		6-16S, 6-16HS, 6-16KS, 6-16KHS	
13860	15,000	10 years	8KBS	
11805	15,000	10 years	8KBS, 8KS	
11806		10 years	8KBS, 8KS	
12505	30,000		8KS	
13845	20,000		8KS	
13850	10,000	10 years	8KS	
13890	25,000		8KS Clinic	

11 Appendix

11.4 Resistance data



The data refer to resistance at 20°C.

Medium	Formula	Concentration [%]	High Density Polyethylene HDPE	Polyamide PA	Polycarbonate PC	Polyoxymethylene POM	Polypropylene PP	Polysulfone PSU	Polyvinyl chloride, hard PVC	Polyvinyl chloride, soft PVC	Polytetrafluorethylene PTFE	Acrylonitrile-butadiene- caoutchouc NBR	Aluminium AL
- no data 1 resistant 2 practically resistant 3 partially resistant 4 not resistant													
Acetaldehyde	C ₂ H ₄ O	40	3	2	4	2	3	4	4	-	1	4	1
Acetamide	C ₂ H ₅ NO	saturated	1	1	4	1	1	4	4	-	1	-	1
Acetone	C ₃ H ₆ O	100	1	1	4	1	1	4	4	-	1	4	1
Acrylonitrile	C ₃ H ₃ N	100	1	1	4	3	3	4	4	4	1	4	1
Allyl alcohol	C ₃ H ₆ O	96	1	3	3	2	2	2	2	4	1	1	1
Aluminium chloride	AlCl ₃	saturated	1	3	2	4	1	-	1	-	1	1	4
Aluminium sulfate	Al ₂ (SO ₄) ₃	10	1	1	1	3	1	1	1	1	1	1	1
Ammonium chloride	(NH ₄)Cl	aqueous	1	1	1	2	1	1	1	1	1	1	3
Ammonium hydroxide	NH ₃ + H ₂ O	30	1	3	4	1	1	2	1	-	1	-	1
Aniline	C ₆ H ₇ N	100	1	3	4	1	2	4	4	4	1	4	1
Anisole	C ₇ H ₈ O	100	3	4	4	1	4	4	2	-	1	4	1
Antimony trichloride	SbCl ₃	90	1	4	1	4	1	-	1	-	1	-	4
Benzaldehyde	C ₇ H ₆ O	100	1	3	4	1	1	3	4	4	1	4	1
Benzene	C ₆ H ₆	100	3	2	4	1	3	4	4	-	1	4	1
Boric acid	H ₃ BO ₃	aqueous	1	3	1	2	1	-	-	-	1	1	1
Butyl acrylate	C ₇ H ₁₂ O ₂	100	1	2	4	2	3	4	4	4	1	-	1
Butyl alcohol, normal	C ₄ H ₁₀ O	100	1	1	2	1	1	2	2	4	1	1	1
Calcium chloride	CaCl ₂	alcoholic	1	4	2	3	1	-	-	4	1	1	3
Carbon disulfide	CS ₂	100	4	3	4	2	4	4	4	4	1	3	1
Carbon tetrachloride (TETRA)	CCl ₄	100	4	4	4	2	4	4	4	4	1	3	1
Chlorine	Cl ₂	100	4	4	4	4	4	4	4	4	1	-	3
Chlorine water	Cl ₂ x H ₂ O		3	4	4	4	3	-	3	3	1	-	4
Chlorobenzene	C ₆ H ₅ Cl	100	3	4	4	1	3	4	4	4	1	4	1
Chloroform	CHCl ₃	100	3	3	4	4	3	4	4	4	1	4	3

11 Appendix

Medium	Formula	Concentration [%]	High Density Polyethylene HDPE	Polyamide PA	Polycarbonate PC	Polyoxymethylene POM	Polypropylene PP	Polysulfone PSU	Polyvinyl chloride, hard PVC	Polyvinyl chloride, soft PVC	Polytetrafluorethylene PTFE	Acrylonitrile-butadiene- caoutchouc NBR	Aluminium AL
- no data 1 resistant 2 practically resistant 3 partially resistant 4 not resistant													
Chromic acid	CrO ₃	10	1	4	2	4	1	4	1	-	1	4	1
Chromic potassium sulphate	KCr(SO ₄) ₂ x 12H ₂ O	saturated	1	2	1	3	1	-	1	-	1	-	3
Citric acid	C ₆ H ₈ O ₇	10	1	1	1	2	1	1	1	1	1	1	1
Citric acid	C ₆ H ₈ O ₇	50	1	3	1	2	1	-	-	-	1	1	1
Copper sulphate	CuSO ₄ x 5H ₂ O	10	1	1	1	1	1	1	1	1	1	1	4
Cyclohexanol	C ₆ H ₁₂ O	100	1	1	3	1	1	1	1	4	1	2	1
Decane	C ₁₀ H ₂₂	100	-	1	2	1	3	-	-	-	1	2	1
Diaminoethane	C ₂ H ₈ N ₂	100	1	1	3	1	1	-	3	4	1	1	1
Diesel fuel	—	100	1	1	3	1	1	-	1	3	1	1	1
Dimethyl formamide (DMF)	C ₃ D ₇ NO	100	1	1	4	1	1	4	3	-	1	3	1
Dimethyl sulfoxide (DMSO)	C ₂ H ₆ SO	100	1	2	4	1	1	4	4	-	1	-	1
Dimethylaniline	C ₈ H ₁₁ N	100	-	3	4	2	4	-	-	-	1	-	1
Dioxane	C ₄ H ₈ O ₂	100	2	1	4	1	3	2	3	4	1	3	1
Dipropylene glycol (mono)methyl ether	C ₄ H ₁₀ O	100	3	1	4	1	4	4	4	4	1	-	1
Ethyl acetate	C ₄ H ₈ O ₂	100	1	1	4	1	1	4	4	4	1	4	1
Ethylene chloride	C ₂ H ₄ Cl ₂	100	3	3	4	1	3	4	4	4	1	-	1
Ferrous chloride	FeCl ₂	saturated	1	3	1	3	1	1	1	1	1	-	4
Formaldehyde solution	CH ₂ O	30	1	3	1	1	1	-	-	-	1	2	1
Formic acid	CH ₂ O ₂	100	1	4	3	4	1	3	3	1	1	2	1
Furfural	C ₅ H ₄ O ₂	100	1	3	3	2	4	-	-	-	1	4	1
Gasoline	C ₅ H ₁₂ - C ₁₂ H ₂₆	100	2	1	3	1	3	3	2	-	1	1	1
Glycerol	C ₃ H ₈ O ₃	100	1	1	3	1	1	1	1	2	1	1	1
Heptane, normal	C ₇ H ₁₆	100	2	1	1	1	2	1	2	4	1	1	1
Hexane, n-	C ₆ H ₁₄	100	2	1	2	1	2	1	2	4	1	1	1
Hydrogen chloride	HCl	5	1	4	1	4	1	1	1	-	1	2	4
Hydrogen chloride	HCl	concentrated	1	4	4	4	1	1	2	3	1	4	4
Hydrogen peroxide	H ₂ O ₂	3	1	3	1	1	1	1	1	-	1	3	3
Hydrogen peroxide	H ₂ O ₂	30	1	4	1	4	1	1	1	-	1	3	3
Hydrogen sulphide	H ₂ S	10	1	1	1	1	1	1	1	3	1	3	1
Iodine, tincture of	I ₂		1	4	3	1	1	-	4	4	1	1	1

11 Appendix

- no data 1 resistant 2 practically resistant 3 partially resistant 4 not resistant			Concentration	High Density Polyethylene	Polyamide	Polycarbonate	Polyoxymethylene	Polypropylene	Polysulfone	Polyvinyl chloride, hard	Polyvinyl chloride, soft	Polytetrafluorethylene	Acrylonitrile-butadiene-caoutchouc	Aluminium
Medium	Formula	[%]	HDPE	PA	PC	POM	PP	PSU	PVC	PVC	PTFE	NBR	AL	
Isopropyl alcohol	C3H8O	100	1	1	1	1	1	1	1	1	4	1	-	2
Lactic acid	C3H6O3	3	1	3	1	2	1	1	2	-	1	1	1	1
Magnesium chloride	MgCl2	10	1	1	1	1	1	1	1	1	1	1	1	1
Mercuric chloride	HgCl2	10	1	4	1	3	1	1	1	1	1	1	1	4
Mercury	Hg	100	1	1	1	1	1	1	1	3	1	1	1	3
Methyl acetate	C3H6O2	100	1	1	4	2	1	-	4	4	1	-	1	1
Methyl alcohol	CH4O	100	1	2	4	1	1	3	1	3	1	2	1	1
Methyl benzene	C7H8	100	3	1	4	1	3	4	4	4	1	4	1	1
Methyl ethyl ketone (MEK)	C4H8O	100	1	1	4	1	1	4	4	4	1	4	1	1
Methylene chloride	CH2Cl2	100	4	3	4	3	3	4	4	4	1	-	1	1
Mineral oil	—	100	1	1	1	1	1	1	1	-	1	1	1	1
Nitric acid	HNO3	10	1	4	1	4	1	1	1	-	1	4	3	1
Nitric acid	HNO3	100	4	4	4	4	4	-	4	-	1	4	1	1
Nitrobenzene	C6H5NO2	100	3	4	4	3	2	4	4	4	1	4	1	1
Oleic acid	C18H34O2	100	1	1	1	2	1	-	1	-	1	3	1	1
Oxalic acid	C2H2O4 x 2H2O	100	1	3	1	4	1	1	1	1	1	2	1	1
Ozone	O3	100	3	4	1	4	3	1	1	-	1	4	2	1
Petroleum	—	100	1	1	3	1	1	1	1	3	1	1	1	1
Phenol	C6H6O	10	1	4	4	4	1	4	1	3	1	3	1	1
Phenol	C6H6O	100	2	4	4	4	1	3	4	3	1	3	1	1
Phosphoric acid	H3PO4	20	1	4	2	4	1	-	-	-	1	2	4	1
Phosphorus pentachloride	PCl5	100	-	4	4	4	1	-	4	4	1	-	1	1
Potassium hydrogen carbonate	CHKO3	saturated	1	1	2	1	1	-	-	-	1	-	4	1
Potassium hydroxide	KOH	30	1	1	4	3	1	1	1	1	1	-	4	1
Potassium hydroxide	KOH	50	1	1	4	3	1	1	1	1	1	-	4	1
Potassium nitrate	KNO3	10	1	1	1	1	1	-	-	-	1	1	1	1
Potassium permanganate	KMnO4	100	1	4	1	1	1	-	1	-	1	3	1	1
Pyridine	C5H5N	100	1	1	4	1	3	4	4	4	1	4	1	1
Resorcinol	C6H6O2	5	1	4	2	3	1	4	2	-	1	-	2	1
Silver nitrate	AgNO3	100	1	1	1	1	1	1	1	1	1	2	4	1

11 Appendix

Medium	Formula	Concentration [%]	High Density Polyethylene HDPE	Polyamide PA	Polycarbonate PC	Polyoxymethylene POM	Polypropylene PP	Polysulfone PSU	Polyvinyl chloride, hard PVC	Polyvinyl chloride, soft PVC	Polytetrafluorethylene PTFE	Acrylonitrile-butadiene- caoutchouc NBR	Aluminium AL
Sodium bisulphite	NaHSO ₃	10	1	1	2	4	1	-	-	-	1	1	1
Sodium carbonate	Na ₂ CO ₃	10	1	1	1	1	1	-	-	-	1	-	3
Sodium chloride	NaCl	30	1	1	1	1	1	1	1	1	1	1	3
Sodium hydroxide	NaOH	30	1	1	4	1	1	1	1	1	1	2	4
Sodium hydroxide	NaOH	50	1	1	4	1	1	1	1	-	1	2	4
Sodium sulfate	Na ₂ SO ₄	10	1	1	1	1	1	1	1	1	1	1	1
Spirits	C ₂ H ₆ O	96	1	1	1	1	1	1	1	3	1	-	1
Styrene	C ₈ H ₈	100	4	1	4	1	3	-	4	4	1	4	1
Sulphuric acid	H ₂ SO ₄	6	1	4	1	4	1	1	1	-	1	2	3
Sulphuric acid	H ₂ SO ₄	fuming	4	4	4	4	4	4	4	4	1	4	3
Tallow	—	100	1	1	1	1	1	-	1	1	1	1	1
Tetrahydrofuran (THF)	C ₄ H ₈ O	100	3	1	4	1	3	4	4	4	1	3	1
Tetrahydronaphthalene	C ₁₀ H ₁₂	100	3	1	4	1	4	4	4	4	1	-	1
Thionyl chloride	Cl ₂ SO	100	4	4	4	2	4	4	4	4	1	-	3
Tin chloride	SnCl ₂	10	1	4	2	2	1	-	-	-	1	1	4
Transformer oil	—	100	1	1	3	3	1	1	1	-	1	1	1
Trichloroethane	C ₂ H ₃ Cl ₃	100	3	3	4	2	4	4	4	4	1	4	4
Urea	CH ₄ N ₂ O	10	1	1	1	1	1	-	-	-	1	1	1
Urine	—	100	1	1	1	1	1	-	1	1	1	-	2
Vinegar	C ₂ H ₄ O ₂	10	1	4	1	1	1	1	1	1	1	2	1
Vinegar	C ₂ H ₄ O ₂	90	1	4	4	4	1	3	1	4	1	-	1
Wax	—	100	-	1	1		1	-	-	-	1	-	1
Wines	—	100	1	1	1	2	1	1	1	1	1	-	4
Xylene	C ₈ H ₁₀	100	3	1	4	1	4	4	4	4	1	4	1

11 Appendix

11.5 EC declaration of conformity



EC Declaration of Conformity

in accordance with the EC Machinery Directive 2006/42/EC, annex II, part 1, section A

Product designation Laboratory centrifuge
Product name Sigma 6-16S, Sigma 6-16HS
Order number 10390, 10391, 10392, 10393
Manufacturer Sigma Laborzentrifugen GmbH
 An der Unteren Söse 50
 37520 Osterode
 Germany

Authorised representative for CE matters Alexander Hörsting

The product named hereinabove was developed, designed, and manufactured in compliance with the relevant, fundamental safety and health requirements of the listed EC directives and norms.

In the event of modifications that were not authorised by the manufacturer or if the product is used in a manner that is not in line with the intended purpose, this declaration will be rendered void.

Directives	Machinery Directive	2006/42/EC
	Low Voltage Directive	2014/35/EU
	EMC Directive	2014/30/EU
	RoHS Directive	2011/65/EU
	Delegated RoHS Directive	(EU) 2015/863

Norms	EN 61010-2-020:2017
	EN IEC 61326-1:2021

Osterode, 2025-11-04



Dr. Michael Sander, Managing Director

CE_Sigma_6-16S_6-16HS_MRL_2025-11-04_en

1 / 1

11 Appendix

11.6 Declaration of conformity – China RoHS 2



Declaration of Conformity

China RoHS 2 (Administrative Measures for the Restriction of the Use of Hazardous Substances in Electrical and Electronic Products)

Laboratory centrifuge models: Sigma 1-14, 1-14K, 1-7, 1-16, 1-16K, 2-7, 2-16KL, 2-16KHL, 3-16L, 3-16KL, 3-16T, 3-17KT, 3-18KS, 3-18KHS, 3-30KS, 3-30KHS, 4-5L, 4-5KL, 4-5KRL, 4-16S, 4-16KS, 4-16KHS, 6-5KRL, 6-16S, 6-16HS, 6-16KS, 6-16KHS, 8KS, 8KBS

Sigma Laborzentrifugen GmbH has made reasonable effort to avoid the use of hazardous substances in the products it manufactures (laboratory centrifuges).

A Product Conformity Assessment (PCA) was performed in order to determine whether the concentration of harmful substances in all homogeneous materials of the component parts is above or below the MCV limit (Maximum Concentration Value limit) as defined in GB/T 26572:

Mercury and its compounds: 0.1 %

Cadmium (Cd) and its compounds: 0.01 %

Lead (Pb) and its compounds: 0.1 %

Hexavalent chromium (Cr (VI)) and its compounds: 0.1 %

Polybrominated biphenyls (PBB): 0.1 %

Polybrominated diphenyl ethers (PBDE): 0.1 %

表1 产品中有害物质的名称及含量

Table 1: Name and content of hazardous substances in the product

部件名称 Component part (PCA)	有害物质 Hazardous substance					
	铅 Lead (Pb)	汞 Mercury (Hg)	镉 Cadmium (Cd)	六价铬 Hexavalent Chromium (Cr (VI))	多溴联苯 Poly-brominated biphenyls (PBB)	多溴二苯醚 Poly-brominated diphenyl ethers (PBDE)
Electronic PCB, cables	X ¹⁾	O	X ¹⁾	O	O	O
Display	O	O	O	O	O	O
Housing	X ¹⁾	O	O	O	O	O
Base, metal, accessories	X ¹⁾	O	O	O	O	O

本表格依据SJ/T 11364的规定编制。

This table is made according to SJ/T 11364.

11 Appendix




O:	表示该有害物质在该部件所有均质材料中的含量均在GB/T 26572规定的限量要求以下。 Indicates that the content of the harmful substance in all homogeneous materials of the component part is below the limit as defined in GB/T 26572.)
X:	表示该有害物质至少在该部件的某一均质材料中的含量超出GB/T 26572规定的限量要求。 (企业可在此处, 根据实际情况对上表打“×”的技术原因进行进一步说明。) Indicates that the content of the harmful substance in at least one homogeneous material of the component part exceeds the limit as defined in GB/T 26572. (Contact the manufacturer for further technical information according to the actual situation.)
1)	Contains parts in compliance with exemptions 6a, 6b, 6c, 7c.I, 7c.II and 8b of 2011/65/EU RoHS.

Apart from the exemptions given in this table, none of the substances listed above have been intentionally added to the product or metallic coatings.

Sigma Laborzentrifugen GmbH
An der Unteren Söse 50
37520 Osterode
Germany

Osterode, 2026-06-10

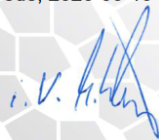

M. Weigoni, Director of Procurement

12 Index

A

Acceleration	46
Acceleration curve	46, 55, 83
Acceleration curve (changing)	55
Acceleration curve (creating)	55
Accessories	
Changes in colour	71
Changes in structure	71
Leaks	71
marking	23
Service life	71
Accessories, cleaning and care	72
Acid	22, 72, 73
Acoustic signal	62
Activating the help function	56
Adapters	39
Adapters, disinfection	70
Alkaline solutions	22, 73
Ambient conditions	80
Ambient temperature	80
Anodised coating	68
Anti-friction coating	73
ASTM standard D91/D96	41
Autoclaving	71
Deformation of vessels	71
Automatic lid opening function	52
Automatic program rotation	60
Automatic rotor identification	46

B

Batch number and serial number of rotors and accessories	23
Blocking a function	51
Blood bag systems	40
Brake	46
Brakeless deceleration	46
Brakeless deceleration from set speed	50
Brief mains power failure	62
Broken glass	75
Broken glass and metal dust, removal from the rotor chamber	75
Bucket	39

Buckets, cleaning and care	72
Buzzer/signal	54

C

Capacity	80
Carrier	39
Carrier systems	39
CE mark in compliance with the directive 2006/42/EC	15
Centrifugation monitoring	49
Centrifugation of blood bag systems	40
Centrifugation of infectious, toxic, radioactive, or pathogenic substances	21
Centrifugation with low capacity	38
Centrifugation with tubes for mineral oil	41
Centrifugation with vessels of various sizes	37
Centrifuge cannot be started	62
Centrifuge decelerates during operation	62
Centrifuge is inclined	62
Centrifuge was moved during run	62
Centrifuges	
definition	9
Changes in colour (accessories)	71
Changes in structure (accessories)	71
Changing an acceleration curve	55
Changing the code	51
Changing the contrast	57
Chemical and biological safety	21
Chemical resistance of plastic	73
China RoHS 2 – Declaration of conformity	93
Cleaning of the centrifuge	67
Cleaning the holes of angle rotors	72
Closing the lid	20, 34
Communication error	64
Condensation	32
Condenser dirty	62
Contamination	21, 72
Continuous run	45
Control system	42
Copyright	10
Corrosion	20, 24, 39, 66, 72, 75
Cost estimate	77

Index

Cracks	68, 73	Exposure to sunlight	66
Creating an acceleration curve.....	55	External signal	54
Curve menu	55	F	
Customer-provided fuses	33	FBRCI-025.....	76
Cycles.....	34, 53	Fire preventions	21
D		Form for the return of defective parts	77
Dangerous goods.....	22	Function.....	52
Date of manufacture	13	Functional and operating elements	12
Deactivating the help function	56	Fuses have tripped	62
Deceleration (brake)	46	G	
Deceleration curve	46, 55, 83	General conditions	10
Declaration of conformity	11, 91	General notes (maintenance).....	66
Declaration of conformity – China RoHS 2.	93	Gravitational field	80
Declaration of decontamination.....	77, 79	Grease for load-bearing bolts.....	72, 73, 74
Decontamination agent	72	H	
Deformation of vessels (autoclaving)	71	Hard running noise during the centrifugation	62
Deleting a program	59	Hazard warnings.....	9, 10, 34, 35, 39
Density.....	13, 22, 49	Hazardous substances	66
Dialog box.....	27	Heavy-duty grease for load-bearing bolts.....	72, 73, 74
Different service life of rotors and accessories	24, 85	Help function.....	56
Dimensions	29	Hexagon socket wrench.....	11
Direct hazard to the life and health.....	16	Highly corrosive substances	21
Directive 2002/96/EC	79	I	
Disinfection of adapters.....	70	Imbalance	37, 39
Disinfection of rotors and buckets	69	Imbalance error.....	64
Display.....	42, 53	Imbalance monitoring system	27
Display deceleration time	52	Importance of the operating manual.....	9
Disposal of the centrifuge.....	79	Important information	16
Disposal of the packaging	79	Improper loading	62
Documentation.....	11	Infectious substances	66
Drive problem	62	Inflammable substances	21
E		Informal safety instructions	19
Earth conductor check	27	Information sheet of the German Berufsgenossenschaft Rohstoffe und chemische Industrie FBRCI-025.....	76
EC declaration of conformity	11, 91	Initial start-up	34
EEPROM error.....	64	Input fuse.....	80
Electrical connection	80	Installation of a rotor	34
Electrical safety.....	20	Installation of accessories	37
Emergency lid release	63	Installation of an angle rotor with a hermetic lid	35
Error correction	62	Installation of rotors and accessories	34
Error message	62		
Executing a program.....	59		
Explanation of the symbols and notes.....	16		
Explosive substances	21		

Index

Installation site	32
Insufficiently greased load-bearing bolts	74
Intended use	9
Interrupting a centrifugation run	43
Interrupting a deceleration process	43
Invert.....	53
IP code	80

K

Kinetic energy	13, 80
----------------------	--------

L

Language.....	53
Layout of the centrifuge.....	12
Leaks (accessories).....	71
Lid cannot be opened	62
Lid error	64
Lid lock device	27
Lid lock has not released	62
Lid seal sticks	62
Linear curves	83
Load-bearing bolts (maintenance).....	73
Loading a process.....	48
Loading a program.....	58
Lock	50
Lye.....	72

M

Mains power switch off.....	62
Mains switch	12
Mains voltage.....	20
Maintenance and service	76
Maintenance of plastic accessories.....	73
Maintenance of the load-bearing bolts	73
Malfunctions and error correction	62
Mandatory testing of centrifuges	76
Manual mode	43
Manufacturer.....	80
Marking of rotors and accessories	23
Marking of the rotors	73
Marking of the unit	15
Maximum speed for tubes	39
Measures in the event of hazards and accidents	28
Mechanical safety	20
Menu Curve	55
Menu Help	56

Menu Parameters	49
Menu Processlibrary	48
Menu Setup	52
Menu Standard	44
Mode of operation	14
Modification mode active.....	43
Motor error	64
Motor shaft.....	34, 35, 72
Multiple carrier	39
Multiple carriers, cleaning and care	72

N

Name plate	12, 13
No indication on the display	62
No power in the mains supply	62
Noise level	80
Nominal voltage	13
Notes on safety and hazards	9
Notes on transport	29

O

Online download of forms	78
Open lid after run	52
Opening the lid.....	34
Operating personnel	19
Operating voltage.....	32
Operational safety.....	72
Option for the input and output of data	54
Options for data input and output.....	61
Overseas shipping	29

P

Packaging	30
Parameter error	64
Parameters menu	49
Pathogenic substances	21, 66
Pick-up order	78
Potential hazard to the life and health	16
Potentially hazardous situation.....	16
Power consumption	13, 80
Power cord is not plugged in.....	62
Power supply	32
Precautions for the care of the accessories	72
Pressure corrosion	
see corrosion	72
Pressure marks.....	73
Prevention of accidents.....	9

Index

Principle of centrifugation	14	Rotors with Sigma "Comfort" rotor coating .	73
Problem description	77	Rotors, cleaning and care	72
Process	49	Runtime	44
Process error	64	Runtime as of set speed	52
Process library	48	S	
Product designation	13	Safety area	22
Program list.....	46, 58, 59	Safety devices	27
Program mode	57	Safety distance	20
Program rotation	52	Safety instructions.....	9, 34, 35, 39
Progress indicator	47	Safety instructions concerning the centrifugation process	22
Protection class	32, 80	Safety of rotors and accessories	23
Q		Safety range	22, 32
Quadratic curves.....	83	Safety, chemical and biological	21
Qualified electrician.....	18	Safety, electrical	20
Quick stop.....	43	Safety, mechanical.....	20
R		Saving a program.....	58
Radioactive substances	21	Scope of supply	11
Radius	14, 49	Screws of the transport safety device are not removed	62
RCF	44	Selection, display, and modification of data	43
Relative centrifugal force (RCF).....	14, 44	Serial number	13, 65
Removal of the transport safety device	30	Serial number of rotors.....	23
Removing a rotor	35, 36	Service contact	65
Removing broken glass and metal dust from the rotor chamber	75	Service life	66
Requirements concerning the personnel....	18	Service life of rotors and accessories... 24, 85	
Responsibility of the operator.....	17	Service life of the accessories.....	71
Return of centrifuges, spare parts, or accessories	77	Set-up and connection	32
Return of defective parts	77	Setup menu	52
Rotor coating Sigma "Comfort"	25	Short run	45
Rotor identification, automatic.....	46	Sigma "Comfort" rotor coating.....	25
Rotor monitoring system	27	Sigma "Comfort" rotor coating.....	73
Rotor radii	82	Signs of corrosion	73
Rotor removal	36	Solvent.....	67, 72
Rotor selection list.....	45	Solvents	22, 73
Rotor tie-down screw	34, 35	Sound signal	27
Rotor wrench	11, 34, 35	Sound signal (warning)	54
Rotors		Spare part enquiries.....	65
marking.....	23	Specialised personnel	19
Rotors and accessories		Speed	13, 14, 44, 80
batch number and serial number	23	Speed/RCF fine	53
Rotors and accessories with a different service life	24, 85	Speedometer error.....	64
Rotors and buckets (disinfection).....	69	Spin-out from set speed	50
Rotors for microtiter plate formats	35	Standard menu	44
		Standards and regulations	11

Index

Standstill monitoring system	27	Thermal stress	66
Start delay.....	52	Time fine.....	53
Starting a centrifugation run	43	Time range.....	80
Starting a process	48	Toxic substances	21
Stopping of the centrifuge due to an imbalance.....	74	Transport safety device.....	30
Storage and transport	29	Tubes.....	39
Storage conditions	29	Type.....	80
Storage locations	80	Type of connection.....	32
Structural changes	21	Type of the centrifuge	65
Supply voltage	32	U	
Supporting rings.....	39	Unblocking a function.....	51
Surface damage	73	Ungreased load- bearing bolts	62
Switching the centrifuge off	61	Units of protection class I	32
Switching the centrifuge on	34	Useful volume	
System.....	53	- volume that is stated for the tube	39
System check	27	User interface	42
System error	62, 64	UV radiation	32
T		UV radiation	66
Table of error codes.....	64	V	
Table of rotors and accessories with a different service life	24, 85	Vents	32
Technical data	80	Vessels	39
Technical documentation	81	W	
Temperature error.....	64	Warning signal	54
Temperature value not reached	62	Warranty and liability.....	10
		Weight	29