

turns more than heads



**Sigma 3-18KT IVD**

**Sigma 3-18KT**



Refrigerated centrifuge  
from serial no. 188383

## Operating manual

Original operating manual

**Read this information carefully before using the centrifuge and retain it for future reference!**

# Revision history

**Document no. 0704802**

**Version 09/2026**

Revision 1.0 - 2026-09-16

Copyright by  
Sigma Laborzentrifugen GmbH  
An der Unteren Söse 50  
37520 Osterode am Harz  
Germany

Telephone: +49 (0) 55 22 / 50 07-0  
E-mail: [info@sigma-zentrifugen.de](mailto:info@sigma-zentrifugen.de)  
[www.sigma-zentrifugen.de](http://www.sigma-zentrifugen.de)

# Contents

|          |                                                                    |           |
|----------|--------------------------------------------------------------------|-----------|
| <b>1</b> | <b>General information .....</b>                                   | <b>11</b> |
| 1.1      | Importance of the operating instructions .....                     | 11        |
| 1.2      | Further applicable documents .....                                 | 11        |
| 1.3      | Representation conventions.....                                    | 11        |
| 1.3.1    | Warning notes .....                                                | 11        |
| 1.3.2    | Typographic conventions .....                                      | 13        |
| 1.3.3    | References to directions.....                                      | 13        |
| 1.4      | Warranty and liability .....                                       | 13        |
| 1.5      | Copyright .....                                                    | 14        |
| 1.6      | Regulations, directives and standards .....                        | 14        |
| 1.7      | Scope of supply .....                                              | 14        |
| <b>2</b> | <b>Safety .....</b>                                                | <b>15</b> |
| 2.1      | Marking of the centrifuge .....                                    | 15        |
| 2.2      | Marking of rotors and accessories .....                            | 16        |
| 2.3      | Intended use.....                                                  | 20        |
| 2.3.1    | Centrifuge .....                                                   | 20        |
| 2.3.2    | Rotors, buckets, carriers and accessories .....                    | 20        |
| 2.3.3    | Starting products .....                                            | 20        |
| 2.4      | Reasonably foreseeable misuse .....                                | 21        |
| 2.5      | Responsibility of the operator.....                                | 21        |
| 2.6      | Safety of products with an MD and IVD mark (medical devices) ..... | 22        |
| 2.7      | Personnel qualification .....                                      | 23        |
| 2.8      | Personal protective equipment.....                                 | 24        |
| 2.9      | Remaining risks.....                                               | 24        |
| 2.9.1    | Hazards due to electrical voltage .....                            | 25        |
| 2.9.2    | Hazards due to mechanical force.....                               | 25        |
| 2.9.3    | Hazards due to hot surfaces .....                                  | 26        |
| 2.9.4    | Hazards due to biological substances.....                          | 26        |
| 2.9.5    | Hazards caused by the operating fluids .....                       | 26        |
| 2.10     | Safety devices and guards .....                                    | 27        |
| <b>3</b> | <b>Equipment description .....</b>                                 | <b>29</b> |
| 3.1      | Principle of centrifugation .....                                  | 29        |
| 3.2      | Centrifuge overview .....                                          | 29        |
| 3.3      | F-gases label .....                                                | 30        |
| 3.4      | Name plate .....                                                   | 31        |
| <b>4</b> | <b>Packaging, transport and storage .....</b>                      | <b>32</b> |
| 4.1      | Safety .....                                                       | 32        |

|          |                                                            |           |
|----------|------------------------------------------------------------|-----------|
| 4.2      | Symbols on the packaging .....                             | 32        |
| 4.3      | Unpacking the centrifuge .....                             | 33        |
| 4.4      | Transporting the centrifuge .....                          | 33        |
| 4.5      | Storing the centrifuge .....                               | 33        |
| <b>5</b> | <b>Setup and connection .....</b>                          | <b>34</b> |
| 5.1      | Setting up the centrifuge .....                            | 34        |
| 5.2      | Establishing the power supply .....                        | 35        |
| 5.3      | Fuses inside the unit .....                                | 35        |
| 5.4      | On-site fuse protection .....                              | 35        |
| <b>6</b> | <b>Operation .....</b>                                     | <b>36</b> |
| 6.1      | Safety .....                                               | 36        |
| 6.2      | Switching the centrifuge on and off .....                  | 36        |
| 6.3      | Opening and closing the lid .....                          | 36        |
| 6.4      | Using rotors, buckets, carriers and accessories .....      | 37        |
| 6.4.1    | Preparing the accessories .....                            | 37        |
| 6.4.2    | Preparing and installing a rotor .....                     | 38        |
| 6.4.3    | Removing the rotor and accessories .....                   | 42        |
| 6.5      | User interface .....                                       | 44        |
| 6.5.1    | Main screen .....                                          | 44        |
| 6.5.2    | Basic functions .....                                      | 45        |
| 6.5.2.1  | Starting a centrifugation .....                            | 45        |
| 6.5.2.2  | Stopping a centrifugation .....                            | 46        |
|          | Performing a quick stop .....                              | 46        |
| 6.5.3    | Process parameters .....                                   | 47        |
| 6.5.3.1  | Entering process parameters .....                          | 47        |
| 6.5.3.2  | Runtime: process parameters .....                          | 48        |
| 6.5.3.3  | RPM and RCF: process parameters .....                      | 50        |
| 6.5.3.4  | Temperature: process parameters .....                      | 50        |
|          | Rapid temperature control function <i>RAPID TEMP</i> ..... | 51        |
| 6.5.4    | Process-related parameters .....                           | 52        |
| 6.5.4.1  | Entering process-related parameters .....                  | 52        |
| 6.5.4.2  | Runtime: process-related parameters .....                  | 54        |
| 6.5.4.3  | RPM: process-related parameters .....                      | 55        |
| 6.5.4.4  | RCF: process-related parameters .....                      | 56        |
| 6.5.4.5  | Temperature: process-related parameters .....              | 57        |
| 6.5.4.6  | Expert mode: Process-related parameters .....              | 58        |
| 6.5.5    | Rotor library .....                                        | 59        |
| 6.5.5.1  | Selecting rotors and buckets .....                         | 59        |
| 6.5.6    | Profiles .....                                             | 61        |

|          |                                                      |           |
|----------|------------------------------------------------------|-----------|
| 6.5.6.1  | Selecting acceleration and deceleration curves ..... | 61        |
| 6.5.7    | Program library .....                                | 62        |
| 6.5.7.1  | Creating or editing a program .....                  | 62        |
| 6.5.7.2  | Creating or editing a custom profile .....           | 64        |
| 6.5.7.3  | Executing a program .....                            | 66        |
| 6.5.7.4  | Deleting a program .....                             | 67        |
| 6.5.8    | Status indicator and status menu .....               | 67        |
| 6.5.8.1  | Opening the status menu .....                        | 67        |
| 6.5.9    | Error log .....                                      | 69        |
| 6.5.9.1  | Deleting error log entries .....                     | 69        |
| 6.5.10   | Settings .....                                       | 70        |
| 6.5.10.1 | Device settings .....                                | 71        |
| 6.5.10.2 | ECO mode .....                                       | 79        |
| 6.5.10.3 | Language selection .....                             | 81        |
| 6.5.10.4 | User management .....                                | 81        |
| 6.5.10.5 | Info .....                                           | 89        |
| 6.5.10.6 | Service .....                                        | 90        |
| 6.5.10.7 | Cycles .....                                         | 92        |
| 6.5.10.8 | Quality Control .....                                | 93        |
| 6.5.11   | Data transfer .....                                  | 95        |
| 6.5.11.1 | Exporting and importing programs .....               | 95        |
| 6.5.11.2 | Exporting QC data .....                              | 97        |
| <b>7</b> | <b>Maintenance .....</b>                             | <b>99</b> |
| 7.1      | Safety .....                                         | 99        |
| 7.2      | General notes .....                                  | 99        |
| 7.3      | Cleaning tasks .....                                 | 101       |
| 7.3.1    | Cleaning the centrifuge .....                        | 101       |
| 7.3.2    | Cleaning the rotors and buckets .....                | 101       |
| 7.3.3    | Cleaning the adapters and carriers .....             | 102       |
| 7.4      | Disinfection .....                                   | 102       |
| 7.4.1    | Disinfecting the centrifuge .....                    | 102       |
| 7.4.2    | Disinfecting the rotors and buckets .....            | 103       |
| 7.4.3    | Disinfecting the adapters and carriers .....         | 103       |
| 7.5      | Autoclaving .....                                    | 104       |
| 7.6      | Maintenance tasks .....                              | 105       |
| 7.6.1    | Maintaining the centrifuge .....                     | 105       |
| 7.6.1.1  | Cleaning the condenser .....                         | 106       |
| 7.6.1.2  | Removing broken glass .....                          | 106       |

---

|           |                                                              |            |
|-----------|--------------------------------------------------------------|------------|
| 7.6.2     | Maintaining the rotors, buckets, adapters and carriers ..... | 107        |
| 7.6.2.1   | Rotors with a Sigma Comfort rotor coating .....              | 108        |
| 7.6.2.2   | Maintaining the load-bearing bolts .....                     | 110        |
| 7.6.2.3   | Maintaining the bayonet lock .....                           | 111        |
| 7.7       | Inspection to be performed by the manufacturer .....         | 112        |
| <b>8</b>  | <b>Troubleshooting .....</b>                                 | <b>113</b> |
| 8.1       | Safety .....                                                 | 113        |
| 8.2       | Identifying and eliminating faults and malfunctions .....    | 113        |
| 8.2.1     | Table of possible faults .....                               | 114        |
| 8.2.2     | Table of error messages .....                                | 115        |
| 8.3       | Opening the lid manually (emergency lid release) .....       | 117        |
| 8.4       | Service contact .....                                        | 118        |
| 8.5       | Returning the centrifuge to the manufacturer .....           | 118        |
| <b>9</b>  | <b>Disposal .....</b>                                        | <b>120</b> |
| 9.1       | Safety .....                                                 | 120        |
| 9.2       | Disposal of the centrifuge .....                             | 120        |
| 9.3       | Disposal of the packaging material .....                     | 120        |
| <b>10</b> | <b>Technical data .....</b>                                  | <b>121</b> |
| 10.1      | Ambient conditions .....                                     | 122        |
| <b>11</b> | <b>Appendix .....</b>                                        | <b>123</b> |
| 11.1      | Rotor program .....                                          | 123        |
| 11.1.1    | Spare parts for the rotors and covers .....                  | 125        |
| 11.2      | Speed, radius and relative centrifugal force (RCF) .....     | 126        |
| 11.3      | Density .....                                                | 126        |
| 11.4      | Acceleration and deceleration curves .....                   | 127        |
| 11.5      | Table of symbols .....                                       | 130        |
| 11.6      | Table of chemical resistance data .....                      | 136        |
| 11.7      | Declarations of conformity .....                             | 142        |
| <b>12</b> | <b>Index .....</b>                                           | <b>146</b> |

## List of figures

|                                                                                                      |    |
|------------------------------------------------------------------------------------------------------|----|
| Fig. 1: Rotor with an engraved batch number and serial number (example) .....                        | 17 |
| Fig. 2: Marking of the set number and the bucket number on the bucket and on the rotor (example).... | 18 |
| Fig. 3: Marking of a bucket with a different service life (example) .....                            | 18 |
| Fig. 4: Marking of a rotor with a different service life (example) .....                             | 19 |
| Fig. 5: Marking of a rotor with a Comfort coating (example).....                                     | 19 |
| Fig. 6: General view, lid open .....                                                                 | 29 |
| Fig. 7: Right side .....                                                                             | 30 |
| Fig. 8: F-gases label in compliance with the implementation regulation (EU) 2024/2729.....           | 30 |
| Fig. 9: Name plate (example).....                                                                    | 31 |
| Fig. 10: Safety range around the centrifuge (example) .....                                          | 34 |
| Fig. 11: Permissible and not permissible loading of angle rotors (example) .....                     | 39 |
| Fig. 12: Biosafe angle rotor and cover with bayonet closure (example).....                           | 39 |
| Fig. 13: Biosafe bucket with a cover (example) .....                                                 | 41 |
| Fig. 14: Permissible and not permissible loading of swing-out rotors (example) .....                 | 41 |
| Fig. 15: Marking on the centrifuge: loading of a swing-out rotor.....                                | 42 |
| Fig. 16: Main screen .....                                                                           | 44 |
| Fig. 17: Message after a quick stop .....                                                            | 46 |
| Fig. 18: Entering process parameters (example).....                                                  | 47 |
| Fig. 19: Runtime area of the main screen (example).....                                              | 48 |
| Fig. 20: Short run – Selecting and holding the Start button .....                                    | 49 |
| Fig. 21: Selecting the $\infty$ button and confirming the input.....                                 | 49 |
| Fig. 22: RPM and RCF areas of the main screen (example).....                                         | 50 |
| Fig. 23: Temperature area of the main screen (example) .....                                         | 50 |
| Fig. 24: Starting the RAPID TEMP function .....                                                      | 51 |
| Fig. 25: Entering process-related parameters (example: Quality Control temperature).....             | 53 |
| Fig. 26: Drop-down menu "Process-related parameters", here: Runtime .....                            | 54 |
| Fig. 27: Main screen with an activated start timer .....                                             | 54 |
| Fig. 28: Drop-down menu "Process-related parameters", here: RPM.....                                 | 55 |
| Fig. 29: Drop-down menu "Process-related parameters", here: RCF .....                                | 56 |
| Fig. 30: Drop-down menu "Process-related parameters", here: Temperature .....                        | 57 |
| Fig. 31: Drop-down menu "Process-related parameters", here: runtime in expert mode .....             | 58 |
| Fig. 32: Drop-down menu "Process-related parameters", here: speed (RMP) in expert mode.....          | 58 |
| Fig. 33: Selecting a rotor or rotor/bucket combination (example).....                                | 60 |
| Fig. 34: Error message of the automatic rotor identification system (example) .....                  | 60 |
| Fig. 35: Curve selection – example: deceleration .....                                               | 61 |

|                                                                                          |     |
|------------------------------------------------------------------------------------------|-----|
| Fig. 36: Creating or editing a program .....                                             | 63  |
| Fig. 37: Creating or editing a profile .....                                             | 65  |
| Fig. 38: Executing a program .....                                                       | 66  |
| Fig. 39: Deleting a program .....                                                        | 67  |
| Fig. 40: Status menu (example) .....                                                     | 68  |
| Fig. 41: Deleting error log entries (example) .....                                      | 69  |
| Fig. 42: "Settings" menu .....                                                           | 70  |
| Fig. 43: "Device settings" menu .....                                                    | 71  |
| Fig. 44: Activating the Automatic lid opening function .....                             | 72  |
| Fig. 45: Activating the Display deceleration time function .....                         | 73  |
| Fig. 46: Setting the time (example) .....                                                | 74  |
| Fig. 47: Entering the signal duration .....                                              | 75  |
| Fig. 48: Selecting a temperature unit .....                                              | 76  |
| Fig. 49: "Display settings" menu for adjusting the brightness .....                      | 77  |
| Fig. 50: "Display settings" menu for activating the screen lock .....                    | 78  |
| Fig. 51: "Display settings" menu with an active "QR code for Sigma Timer" function ..... | 79  |
| Fig. 52: ECO mode function .....                                                         | 80  |
| Fig. 53: Language selection menu (example) .....                                         | 81  |
| Fig. 54: Rights management: creating an administrator account .....                      | 83  |
| Fig. 55: User login .....                                                                | 84  |
| Fig. 56: Creating a user account .....                                                   | 86  |
| Fig. 57: User login via the status menu .....                                            | 88  |
| Fig. 58: User logout .....                                                               | 89  |
| Fig. 59: "Info" menu – entering a centrifuge name .....                                  | 90  |
| Fig. 60: "Service" menu – QR code for downloading an operating manual .....              | 91  |
| Fig. 61: "Cycles current setup" menu .....                                               | 92  |
| Fig. 62: "Quality Control" menu .....                                                    | 93  |
| Fig. 63: "Quality control" dialog box – temperature error message .....                  | 94  |
| Fig. 64: Program export .....                                                            | 96  |
| Fig. 65: "Export programs" dialog box - starting the transfer .....                      | 96  |
| Fig. 66: "Export programs" dialog box - data transfer completed .....                    | 97  |
| Fig. 67: Export of QC data (example) .....                                               | 98  |
| Fig. 68: "QC data transfer" dialog box - data transfer completed .....                   | 98  |
| Fig. 69: Load-bearing bolt with traces of use – no greasing required .....               | 109 |
| Fig. 70: Load-bearing bolt with completely worn-off coating – greasing required .....    | 109 |
| Fig. 71: Applying the grease for load-bearing bolts .....                                | 110 |
| Fig. 72: Cover of an angle rotor with a bayonet lock (example) .....                     | 111 |

---

|                                                                                 |     |
|---------------------------------------------------------------------------------|-----|
| Fig. 73: Measure to be taken when the cover is difficult to open or close. .... | 112 |
| Fig. 74: Position of the opening for the emergency lid release.....             | 117 |
| Fig. 75: Diagram of linear acceleration curves .....                            | 127 |
| Fig. 76: Diagram of quadratic acceleration curves .....                         | 128 |
| Fig. 77: Diagram of the optimised quadratic acceleration curves .....           | 129 |

## List of tables

|                                                                                      |     |
|--------------------------------------------------------------------------------------|-----|
| Table 1: Marking of the centrifuge .....                                             | 16  |
| Table 2: Marking of the accessories .....                                            | 16  |
| Table 3: Symbols on the packaging.....                                               | 32  |
| Table 4: Table of the required activities and their frequency .....                  | 100 |
| Table 5: Autoclaving table.....                                                      | 104 |
| Table 6: Table of possible faults .....                                              | 115 |
| Table 7: Table of error messages .....                                               | 116 |
| Table 8: Product data.....                                                           | 121 |
| Table 9: Connection data.....                                                        | 121 |
| Table 10: Performance parameters .....                                               | 121 |
| Table 11: Other parameters.....                                                      | 121 |
| Table 12: Physical data.....                                                         | 122 |
| Table 13: Sound pressure level .....                                                 | 122 |
| Table 14: Operating fluids and quantities .....                                      | 122 |
| Table 15: Ambient conditions.....                                                    | 122 |
| Table 16: Table of available rotors .....                                            | 124 |
| Table 17: Table indicating the spare parts for the rotors and covers .....           | 125 |
| Table 18: Table of the gradient of linear acceleration curves.....                   | 127 |
| Table 19: Table of the gradient of quadratic acceleration curves.....                | 128 |
| Table 20: Table of the gradient of the optimised quadratic acceleration curves ..... | 129 |
| Table 21: Symbols of the control bar .....                                           | 130 |
| Table 22: Symbols of the configuration bar .....                                     | 131 |
| Table 23: Symbols on the status bar .....                                            | 131 |
| Table 24: Symbols of the controls in the input area.....                             | 132 |
| Table 25: Symbols of the elements of information in the input area.....              | 135 |
| Table 26: Table of chemical resistance data .....                                    | 140 |

# 1 General information

## 1.1 Importance of the operating instructions

Familiarity with the operating instructions ensures the safe and trouble-free operation of the centrifuge. This manual is an integral part of the product.

The operating instructions and, in particular, the safety and hazard notes, must be observed by everyone working with this centrifuge.

- The operating instructions and any other applicable documents must be kept at the location of use of the centrifuge. Ensure that these documents are accessible at all times.
- Any changes, additions or updates received must be added to the operating manual.
- The operating instructions and any other applicable documents must be handed over to every subsequent operator of the centrifuge.

## 1.2 Further applicable documents

The following documents must be observed in addition to this instruction manual:

- Brochure Sigma 3-18KT,
- Cleaning and care leaflet,
- Quality certificate.

## 1.3 Representation conventions

The following symbols and typographic conventions are used throughout this document.

### 1.3.1 Warning notes

#### Warning symbols



This is the hazard symbol. It accompanies all hazard, warning and caution notes. It warns of possible injuries.

- Follow all of the measures marked with a hazard sign to avoid injuries and death.

### Section-related warning notes

Warning of injury hazards that exist permanently or repeatedly during a sequence of actions.

|                                                                                                  |                                                                                                                                                                                           |
|--------------------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
|  <b>DANGER</b>  | <b>Description of a dangerous situation</b><br><br>Will lead to serious or fatal injuries if not avoided. <ul style="list-style-type: none"> <li>• Hazard prevention measures.</li> </ul> |
|  <b>WARNING</b> | <b>Description of a dangerous situation</b><br><br>Can lead to serious or fatal injuries if not avoided. <ul style="list-style-type: none"> <li>• Hazard prevention measures.</li> </ul>  |
|  <b>CAUTION</b> | <b>Description of a dangerous situation</b><br><br>Can lead to minor or moderate injuries if not avoided. <ul style="list-style-type: none"> <li>• Hazard prevention measures.</li> </ul> |
| <b>ATTENTION</b>                                                                                 | <b>Description of damage to property and the environment</b><br><br>Used to indicate damage to property. <ul style="list-style-type: none"> <li>• Hazard prevention measures.</li> </ul>  |

### Integrated warning notes

Warning of injury hazards that exist when performing one or a few related steps within a sequence of actions.

**DANGER** – Type of danger. Instructions to avoid a dangerous situation that will lead to serious or fatal injuries.

**WARNING** – Type of danger. Instructions to avoid a dangerous situation that can lead to serious or fatal injuries.

**CAUTION** – Type of danger. Instructions to avoid a dangerous situation that can lead to minor or moderate injuries.

**ATTENTION** – Type of property damage. Instructions to avoid damage to property and the environment.

---

### 1.3.2 Typographic conventions

#### Instructions for a specific course of action

- ▶ Prerequisite.
- 1. Action (single step).
- 2. Action (single step).
  - Intermediate result.
- 3. Action (single step).
- 4. Action (single step).
  - Final result.

#### Lists

- Bullet point.
  - Subordinate point.

#### Notices

-  This symbol indicates a notice.

#### References

-  This symbol refers to separate applicable documents.

---

### 1.3.3 References to directions

Throughout this document, references to directions refer to the position from the perspective of the equipment user.

## 1.4 Warranty and liability

Our general terms and conditions apply. They have been handed over to the purchaser on the conclusion of the contract.

Warranty and liability claims are excluded if they are due to one or several of the following reasons:

- Improper use
- Non-compliance with the safety and hazard notes in the operating instructions and other applicable documents
- Improper set-up, connection, start-up, operation, transport or maintenance of the centrifuge
- Disasters caused by external influences or force majeure

## 1.5 Copyright

This document is solely intended for the operator and their personnel. It includes instructions and notes that must not be duplicated, distributed or communicated in any other way, neither in full nor in parts.

The copyright concerning this document remains with Sigma Laborzentrifugen GmbH.

Non-compliance may be prosecuted under criminal law.

## 1.6 Regulations, directives and standards

The applicable regulations, directives and standards are listed in the EC declaration of conformity in the appendix (see chapter 11.7 - "Declarations of conformity").

## 1.7 Scope of supply




|                                                         |                  |
|---------------------------------------------------------|------------------|
| 1 laboratory centrifuge Sigma 3-18KT                    |                  |
| 1 power cord with an IEC C13 connection                 |                  |
| 1 socket spanner, size 10 (rotor fastener)              | Ref. no. 930 100 |
| 1 socket spanner, size 5 (emergency lid release system) | Ref. no. 930 053 |
| 1 tube (30 g) heavy-duty grease for load-bearing bolts  | Ref. no. 71 401  |
| 1 operating manual                                      | Ref. no. 0704802 |
| 1 brochure                                              | Ref. no. 0613902 |

### Accessories

According to your order, our order confirmation and our delivery note.

## 2 Safety

### 2.1 Marking of the centrifuge
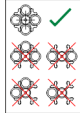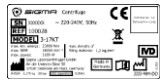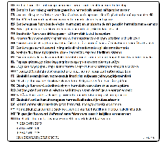
| Marking                                                                             | Description                                                                                                                       |
|-------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------|
|    | ON (mains power)                                                                                                                  |
|    | OFF (mains power)                                                                                                                 |
|    | Arrow indicating the direction of rotation (example)                                                                              |
|    | Information about the loading of the rotor                                                                                        |
|   | Name plate                                                                                                                        |
|  | F-gases label                                                                                                                     |
|  | CE mark in compliance with the regulation (EU) 2017/746 on in vitro diagnostic medical devices or with the directive 2006/42/EC   |
|  | Symbol indicating the need for the separate disposal of waste electrical and electronic equipment as per the directive 2012/19/EU |
|  | Symbol indicating that compliance with the operating instructions is mandatory                                                    |
|  | Symbol indicating that the product is a medical device as per the regulation (EU) 2017/746                                        |
|  | NRTL mark (USA and Canada only)                                                                                                   |

| Marking                                                                           | Description                               |
|-----------------------------------------------------------------------------------|-------------------------------------------|
|  | RCM mark (Australia only)                 |
|  | California Proposition 65 mark (USA only) |
|  | China RoHS 2 mark (China only)            |

*Table 1: Marking of the centrifuge*

-  The markings on the centrifuge must not be covered, removed or altered, and they must be replaced if they become illegible.
-  The marking of the centrifuge varies depending on the version and country of destination.

## 2.2 Marking of rotors and accessories

| Marking                                                                             | Description                                                                                 |
|-------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------|
|  | "Biosafe" symbol (biosafety leak tightness as per DIN EN 61010, part 2-020, appendix AA)    |
|  | "Biosafe" lettering (biosafety leak tightness as per DIN EN 61010, part 2-020, appendix AA) |

*Table 2: Marking of the accessories*

### Batch number and serial number

During production, every rotor and bucket receives a batch number enabling conclusions to be drawn concerning the production process and the subsequent quality inspection.

Some rotors also have an additional serial number providing further detailed information.

The batch number and serial number is engraved on the rotor as follows:



**i** If there are any enquiries concerning the rotor or bucket, please state the batch number and serial number!

### **Set number and bucket number**

Buckets are always produced in sets for a specific rotor. The set number is engraved on the bucket.

In addition, every bucket is assigned a fixed position in the rotor by way of a bucket number. The bucket number is engraved on the bucket and also on the rotor (see the following picture).

- Always install the buckets in their intended positions and align the inscription on the buckets in the same direction for all positions (inwards or outwards).



*Fig. 2: Marking of the set number and the bucket number on the bucket and on the rotor (example)*


- 1 Bucket number
- 2 Set number

### Marking in case of a different service life

Rotors and buckets have a limited service life (see chapter 7.6.2 - "Maintaining the rotors, buckets, adapters and carriers").

If the actual service life differs from the standard service life, the rotors and buckets are marked accordingly.

- A bucket with the engraving "max. cycles = 15,000" has a maximum service life of 15,000 cycles.

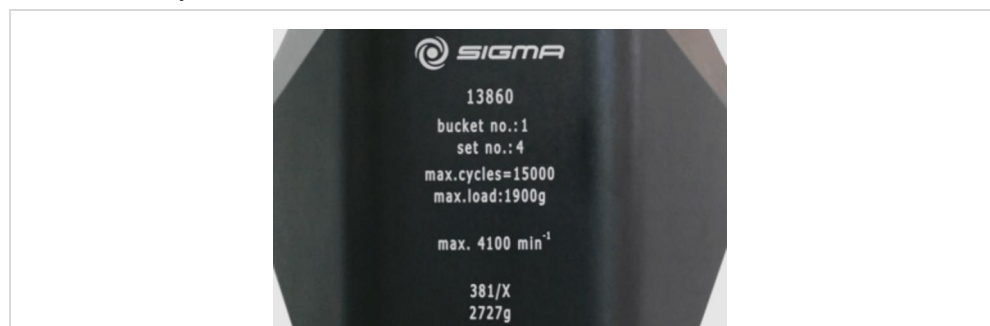

*Fig. 3: Marking of a bucket with a different service life (example)*

## Safety

- A rotor with the engraving "Exp. Date 01/29" (expiry date 01/29) must be scrapped by January 2029 at the latest.



*Fig. 4: Marking of a rotor with a different service life (example)*

### Marking of rotors with Sigma Comfort coating

All swing-out rotors with a Sigma Comfort coating have a "C" shown after their number. These rotors do not need to be greased during the service life of the coating.

-  Rotors with a Sigma Comfort rotor coating must not be used in combination with buckets or carriers with a cleanroom coating.



*Fig. 5: Marking of a rotor with a Comfort coating (example)*

- 1 Marking indicating the Sigma Comfort coating
- 2 Load-bearing bolt

-  A list of the available rotors with a Sigma Comfort coating as well as information about the service life of the coating can be found in the "Sigma Comfort rotor coating" document that comes supplied with every rotor which has a Comfort coating.

---

## 2.3 Intended use

This chapter and its subchapters describe the centrifuge's intended area of application. Any use other than the one described here or any use going beyond this description is regarded as improper use. Sigma Laborzentrifugen GmbH cannot be held liable for any damage resulting from this.

The intended use also includes:

- Observation of all of the notes and instructions that are included in the operating manual,
- Compliance with the care, cleaning and repair instructions.

---

### 2.3.1 Centrifuge

Centrifuges are power-driven machines that use centrifugal force to separate mixtures with a maximum density of 1.2 g/cm<sup>3</sup> into their components. They are intended for this specific purpose only.

Centrifuges with an IVD mark are intended for human biological samples, including donated blood and tissue, in conjunction with diagnostic in-vitro applications. This means that they are medical devices in accordance with the regulation (EU) 2017/746 on in-vitro diagnostic medical devices.

The centrifuge may only be operated by trained and authorised specialist personnel within closed spaces if the following conditions are fulfilled:

- The centrifuge has been properly set up and connected.
- The centrifuge is in a perfect technical state.
- Maintenance and service intervals are adhered to.
- The centrifuge is not located in an area with an explosion hazard.

---

### 2.3.2 Rotors, buckets, carriers and accessories

- The rotors, buckets, carriers and accessories have been produced by Sigma Laborzentrifugen GmbH and are authorised for use in this centrifuge.
- The maintenance intervals are adhered to.
- The service life specifications and, in particular, the specifications concerning the cycle limitations are complied with (see chapter 7.6.2 - "Maintaining the rotors, buckets, adapters and carriers").

---

### 2.3.3 Starting products

The centrifuge is intended for the centrifugation of solid substances as well as mixtures of liquids or solids. The substances to be centrifuged must fulfil the following requirements:

- The substances must not damage the material of the centrifuge, rotors or inserts nor impair its mechanical strength.
- Biological material may only be centrifuged in special, certified containment systems with a biological seal to prevent the release of the material.

## 2.4 Reasonably foreseeable misuse

There are applications where it is foreseeable that malfunctions, accidents or unsafe conditions may occur.

- Do not use the centrifuge in an area with an explosion hazard.
- Do not centrifuge explosive or flammable substances.
- Do not centrifuge substances that are prone to a high-energy chemical reaction.

## 2.5 Responsibility of the operator

The operator is the person who uses the centrifuge for commercial, economic or scientific purposes or lets a third party use it and who bears the legal product responsibility for the protection of personnel or third parties during its operation.

### Personnel

The operator is obliged:

- to comply with the specifications concerning the personnel (see chapter 2.7 - "Personnel qualification"),
- to clearly define the areas of responsibility of the personnel concerning the operation, cleaning and care.
- to ensure that any personnel undergoing training as well as persons with physical, sensory or mental disabilities only work on the centrifuge under the supervision of an experienced person.
- to inform the operating personnel about the procedure concerning serious incidents in connection with products bearing an MD or IVD mark (see chapter 2.6 - "Safety of products with an MD and IVD mark (medical devices)").

### Work environment

The operator is obliged:

- to perform a risk assessment concerning potential accidents in connection with the centrifuge and take design-related countermeasures, if necessary.
- to perform a risk assessment concerning the specific hazards with regard to the processing of the various products in the centrifuge (e.g. leakage of harmful substances, contaminated surfaces, fire hazard) and to take appropriate protective measures, if necessary.
- to perform compatibility tests of the substances to be processed in the centrifuge. Substances that damage the materials of the centrifuge, rotors or inserts or that weaken their mechanical strength must not be used.
- to provide standard operating procedures (SOPs) for the processing of the various products in the centrifuge.
- to check the safety-conscious work of the personnel in compliance with the operating instructions, other applicable documents and SOPs as well as with the relevant EU directives and the national laws concerning health and safety and the prevention of accidents at regular intervals (e.g. every month).
- to ensure that a sufficient quantity of suitable fire-extinguishing equipment is provided and ready for use.

**Maintenance**

The operator is obliged:

- to have the centrifuge checked by authorised specialist personnel at regular intervals to ensure its correct operation.
- to immediately stop using a damaged device and to have defective parts replaced without delay.
- to keep all safety and hazard notes on the centrifuge readable at all times and to replace them, if necessary.
- to check – at their own accord and at regular intervals – the number of cycles performed and the condition of the rotors and accessories.
- to comply with the inspection obligation as per the information FBRCI-025 of the German Berufsgenossenschaft Rohstoffe und chemische Industrie for centrifuges operated in Germany (see chapter 7.7 - “Inspection to be performed by the manufacturer”).

**2.6 Safety of products with an MD and IVD mark (medical devices)**

Centrifuges labelled as medical devices (MD) in accordance with regulation (EU) 2017/745 and those labelled as in-vitro diagnostic medical devices (IVD) in accordance with regulation (EU) 2017/746 are exclusively intended for use by specialist users in their respective field of work.

**Serious incident (definition)**

*Incident* means any malfunction or deterioration in the characteristics or performance of an MD or IVD device which is made available on the market for an intended use as defined by the manufacturer, including user errors due to ergonomic features as well as any inadequacy in the information supplied by the manufacturer or any undesirable side effects.

*Serious incident* means any incident that directly or indirectly led, might have led or might lead to any of the following:

- the death of a patient, user or other person
- the temporary or permanent serious degradation of a patient's, user's or other person's state of health
- a serious public health threat

**Reporting of a serious incident to the manufacturer and the competent authorities within the EU**

A reportable serious incident is an incident that has a causal relationship with an MD or IVD product or where such a causal relationship is reasonably possible. Such an incident must be reported without delay to the manufacturer, Sigma Laborzentrifugen GmbH or the competent national authorities within the EU. If the report is issued to the manufacturer or their authorised representative, they must fulfil their reporting obligation to the authorities. EUDAMED provides forms for the electronic reporting of incidents by manufacturers and authorised representatives via MIR (Manufacturer Incident Report) or by users, operators and distributors (User Report).

## Safety

Authorised distributors, vendors and service providers are required by law to cooperate with the manufacturer in order to ensure the safety of medical devices. Non-compliance may result in substantial fines and penalties.

To report a serious incident to the manufacturer:

1. Please complete the "Vorkommnismeldung an den Hersteller" (incident report for the manufacturer) form in full.
2. Send the form to [support@sigma-zentrifugen.de](mailto:support@sigma-zentrifugen.de)

The form for reporting an IVD incident to the manufacturer is available on our website.

[www.sigma-zentrifugen.de/en.html](http://www.sigma-zentrifugen.de/en.html) → [Service]  
→ [Overhaul and repair]

### Scan QR code

If you have any questions about this topic, kindly contact our service department:



## 2.7 Personnel qualification

Work on the centrifuge may only be performed by specialist personnel that have been authorised by the operator. Specialist personnel can identify and avoid hazards without assistance. Specialist personnel are familiar with the applicable statutory accident prevention regulations as well as the relevant standards and provisions.

The qualifications of the various groups of persons are defined as follows:

### Operating personnel

The centrifuge may only be operated by personnel who, due to their training as clinical or medical laboratory technicians or due to comparable training, are able to perform the assigned tasks in a safe manner.

The tasks of the operating personnel include:

- the operation of the centrifuge
- the cleaning, decontamination and maintenance of the centrifuge
- the notification of the competent specialist in the event of malfunctions
- the notification of the manufacturer or competent authority in the event of serious incidents (see chapter 2.6 - "Safety of products with an MD and IVD mark (medical devices)").

### Logistics personnel

Logistics personnel are in a position to perform their assigned tasks in a safe manner due to their training, experience and familiarity with the relevant regulations.

The tasks of the logistics personnel include:

- the transport and packaging of the centrifuge, rotors and accessories
- the placement and setup of the centrifuge

### Qualified electricians, mechanical specialists and certified refrigeration specialists

Qualified electricians, mechanical specialists and certified refrigeration specialists are in a position to perform their assigned tasks in a safe manner due to their specialised training, experience and familiarity with the relevant regulations.

- The tasks of a **qualified electrician** include the maintenance and repair of the electrical system of the centrifuge.
- The tasks of a **mechanical specialist** include the maintenance and repair of the mechanical components of the centrifuge.
- The tasks of a **certified refrigeration specialist** include the maintenance and repair of the refrigeration system of the centrifuge.

## 2.8 Personal protective equipment

Personal protective equipment protects you from health and safety hazards at work.

The personnel must wear personal protective equipment when performing certain types of work on the centrifuge. This is pointed out in the respective sections of this manual. In addition, personal protective equipment may be required depending on the substances being centrifuged.

## 2.9 Remaining risks

The centrifuge was built in accordance with the state of the art and in compliance with the generally recognised safety rules. However, danger to life and limb of the user, or of third parties, or impairments of the centrifuge or other material assets cannot be completely excluded when the centrifuge is in use.



Read and comply with the following safety information prior to using the centrifuge!

### **2.9.1 Hazards due to electrical voltage**

#### **Danger to life caused by electric shock**

If the centrifuge or mains power cable are damaged, there is a risk of electric shock when touching current-carrying components.

This may lead to ventricular fibrillation, cardiac arrest or respiratory paralysis.

- Ensure that the mains voltage matches the voltage that is stated on the name plate of the centrifuge.
- Only use the centrifuge if the device itself and the mains power cable are undamaged. Damaged or faulty mains power cables must be replaced immediately.
- Ensure that no liquid penetrates the inside of the centrifuge. Spilled liquids or liquids that are used for cleaning can damage the electrical or mechanical components if they penetrate into the centrifuge.
- Only qualified electricians are authorised to perform work on the electrical system of the centrifuge that requires the removal of the housing panels.

#### **Fire hazard due to electrical overload**

Electrical overload due to defective components can cause a fire.

- Have the electrical system of the centrifuge checked by a qualified electrician at regular intervals or at least every 12 months.

---

### **2.9.2 Hazards due to mechanical force**

#### **Hazards caused by rotors, buckets, carriers and accessories**

- Use the centrifuge and the accessories only in a perfect state. Stop using the centrifuge and the accessories immediately if they are damaged.
- In the event of rotor breakage resulting in a crash, there is a risk of injury from ejected parts or movement of the centrifuge. Cutting, impact or crushing injuries may occur.
- Do not move a running centrifuge. If a mobile centrifuge bench is used, lock the locking castors.
- Do not open the lid of a running centrifuge by way of the emergency lid release system.
- Comply with the specified service life of rotors, buckets, carriers and accessories.
- Use only rotors, buckets, carriers and accessories that have been manufactured and authorised for use with this centrifuge by the manufacturer.
- Ensure that the rotor is loaded in a rotationally symmetrical manner and that the weight is evenly distributed. Do not exceed the maximum permissible weight.
- Fasten the rotor correctly in place.
- Operate the rotor at the permissible speed.
- In case of swing-out rotors, ensure that the maximum permissible speed of the buckets is not exceeded.
- Comply with the specified maintenance and replacement intervals for the centrifuge, rotor, carriers and accessories.
- Only use spare parts that have been authorised by the manufacturer.

**Hazards caused by the centrifuge lid**

- In the event of defective lid stays, there is a risk of injury from the centrifuge lid. Crushing injuries may occur.
- Defective lid stays must be replaced immediately by specialist personnel.

---

**2.9.3 Hazards due to hot surfaces****Risk of burns on hot surfaces inside the rotor chamber**

After a centrifugation, the surfaces inside the rotor chamber and, in particular the motor shaft, can be hot.

There is a risk of burns when touching the surfaces.

- Wear suitable protective gloves.
- Do not touch the surfaces without good reason.

---

**2.9.4 Hazards due to biological substances****Health hazard due to biological substances**

The release of biological material can cause damage to health as well as cross-contamination.

- Set the centrifuge up in a room that fulfils the safety requirements for the respective biological substances (as per the directive 2000/54/EC).
- Wear personal protective equipment.
- Use special, certified containment systems with biological seals.
- Fill, close and open the rotors and buckets only at workstations that fulfil the safety requirements for the respective biological substances (as per the directive 2000/54/EC).
- The sealing ring of a rotor cover or bucket is made of NBR. Chemicals may damage the material and cause leaks. It is essential to comply with the information provided in the table of chemical resistance data (see chapter 11.6 - "Table of chemical resistance data").
- Observe any national regulations and the safety data sheets of the substances.

---

**2.9.5 Hazards caused by the operating fluids****Hazards caused by the refrigerant**

If refrigerant leaks from the centrifuge, harmful substances may be released. Contact with, or the aspiration of, the refrigerant can cause severe damage to health.

- Switch the centrifuge off and leave the affected area immediately.
- Ensure adequate ventilation.
- Have the centrifuge checked and repaired by a certified refrigeration specialist.

## 2.10 Safety devices and guards

The safety devices and guards are in compliance with the applicable national and international laws, rules and regulations concerning health and safety and the prevention of accidents.

- Check the safety devices and guards at regular intervals to ensure that they are operating correctly.
- Do not tamper with the safety devices and guards.

### **Mains power plug**

Unplugging the mains power plug disconnects the centrifuge from the power supply.

Position: close to the centrifuge; accessible at all times

### **Mains power switch**

The mains power switch disconnects the centrifuge from the power supply.

Position: right side of the centrifuge; accessible at all times

### **Protective earth conductor**

A qualified electrician can test the protective conductor using an appropriate measuring device.

### **System for stopping the centrifuge due to an imbalance**

If the centrifuge runs with an impermissible imbalance, the drive will be switched off. A visual and audible signal will be issued.

Position: user interface

### **System check**

An internal system check monitors the data traffic and sensor signals for plausibility. The system performs continuous self-monitoring to detect faults. Error messages will be displayed in a dialogue box.

Position: user interface

### **Temperature monitoring system**

If the temperature inside the rotor chamber exceeds +50°C, the drive will be switched off automatically. The centrifuge cannot be restarted until it has cooled down. A corresponding error message will be displayed.

Position: user interface

**Rotor identification system**

The rotor identification system automatically identifies the rotor and adapts the speed and relative centrifugal force (RCF) accordingly. If several buckets are authorised for use with the rotor in question, a dialogue box will show the permissible variants.

Position: user interface

## 3 Equipment description

### 3.1 Principle of centrifugation

Centrifugation is a process for separating heterogeneous mixtures of substances (suspensions, emulsions or gas mixtures) into their components. Rotating on a circular trajectory, the mixture of substances is subjected to centrifugal acceleration many times greater than the acceleration due to gravity.

Centrifuges harness inertia in the rotor chamber to separate substances. Due to their higher inertia, particles or media with a higher density migrate outwards, thereby displacing any components with a lower density which are forced towards the centre.

The centrifugal acceleration of an object in a centrifuge depends on the distance of the object from the axis of rotation and on the angular velocity. The larger the radius of the rotor chamber and the higher the speed are, the greater the centrifugal acceleration will be.


For detailed information about the various areas of application, see [www.sigma-zentrifugen.de/en.html](http://www.sigma-zentrifugen.de/en.html) → [Applications].

### 3.2 Centrifuge overview



*Fig. 6: General view, lid open*

- 1 Lid
- 2 Rotor chamber
- 3 User interface

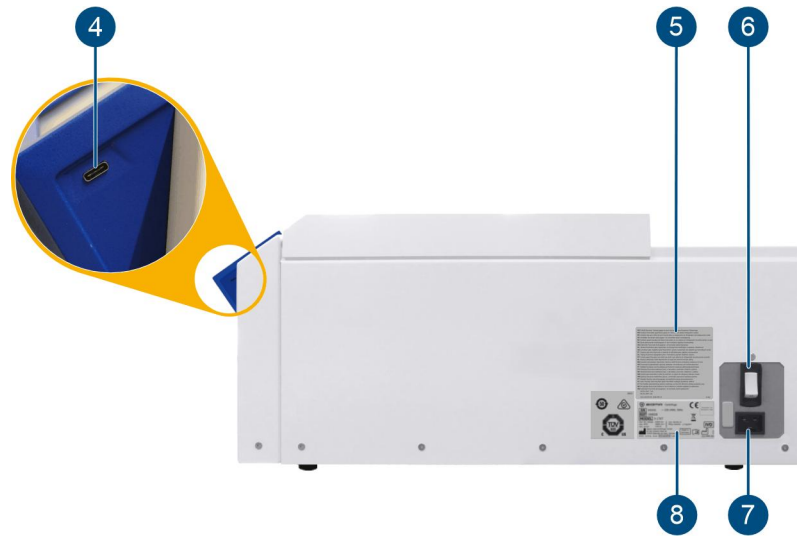

*Fig. 7: Right side*

- |                      |                     |
|----------------------|---------------------|
| 4 USB-C port         | 7 Mains power input |
| 5 F-gases label      | 8 Name plate        |
| 6 Mains power switch |                     |

### 3.3 F-gases label

**DE** Enthält fluorierte Treibhausgase in einer hermetisch geschlossenen Kälteanlage  
**EN** Contains fluorinated greenhouse gases in a hermetically sealed refrigeration system  
**FR** Contient des gaz à effet de serre fluorés dans un équipement de réfrigération hermétiquement scellé  
**SV** Innehåller fluorerade växthusgaser i en hermetiskt sluten kylanläggning  
**ES** Contiene gases fluorados de efecto invernadero en un sistema de refrigeración herméticamente cerrado  
**NL** Bevat gefluoreerde broeikasgassen in een hermetisch gesloten koelinstallatie  
**DA** Indeholder fluorerede drivhusgasser i et hermetisk lukket kølesystem  
**PL** Zawiera fluorowane gazy cieplarniane w hermeticznie zamkniętym urządzeniu chłodniczym  
**GA** Cuimsítear gáis cheaptha teasa fhluaírinthe i gcóras cuisniúcháin atá séalaithe go heirméiteach anseo  
**IT** Contiene gas serra fluorurati in impianto di raffreddamento sigillato ermeticamente  
**SL** Vsebuje fluorirane toplogredne pline v hermetično zaprtem hladilnem sistemu  
**PT** Contém gases fluorados com efeito de estufa num sistema de refrigeração hermeticamente fechado  
**EL** Περιέχει φθοριούχα αέρια θερμοκηπίου σε ερμητικά κλειστό σύστημα ψύξης  
**BG** Съдържа флуорирани парникови газове в херметически затворена хладилна система  
**HU** Fluorozott üvegházhatású gázokat tartalmaz, hermetikusan zárt hűtőrendszerben  
**FI** Sisältää fluorattuja kasvihuonekaasuja ilmatiiviisti suljetussa jäähdytysjärjestelmässä  
**CS** Obsahuje fluorované skleníkové plyny v hermeticky uzavřeném chladičím systému  
**SK** Obsahuje fluórované skleníkové plyny v hermeticky uzavretom chladiacom systéme  
**RO** Conține gaze fluorurate cu efect de seră într-un sistem de refrigerare etanșat ermetic  
**HR** Sadržava fluorirane stakleničke plinove u hermetički zatvorenoj rashladnoj opremi  
**ET** Sisaldab fluuritüd kasvuhoonegaase hermeetiliselt suletud jahutussüsteemis  
**LV** Satur fluorētas siltumnīcefekta gāzes hermētiski noslēgtā dzesēšanas sistēmā  
**LT** Hermetiškai užsandarintoje šaldymo sistemoje yra fluorintų šiltnamio efektą sukeliančių dujų  
**MT** Fih gassijiet fluworurati b'effett ta' serra f'sistema ta' tkessi h'issigillata ermetikament  
**NO** Inneholder fluorinerte drivhusgasser i et hermetisk lukket kjølesystem  
R452A GWP: 2140  
R513A GWP: 630  
R455A GWP: 145  
(EU) 2024/2729; 2028-DEC-31

31402 a

*Fig. 8: F-gases label in compliance with the implementation regulation (EU) 2024/2729*

### 3.4 Name plate

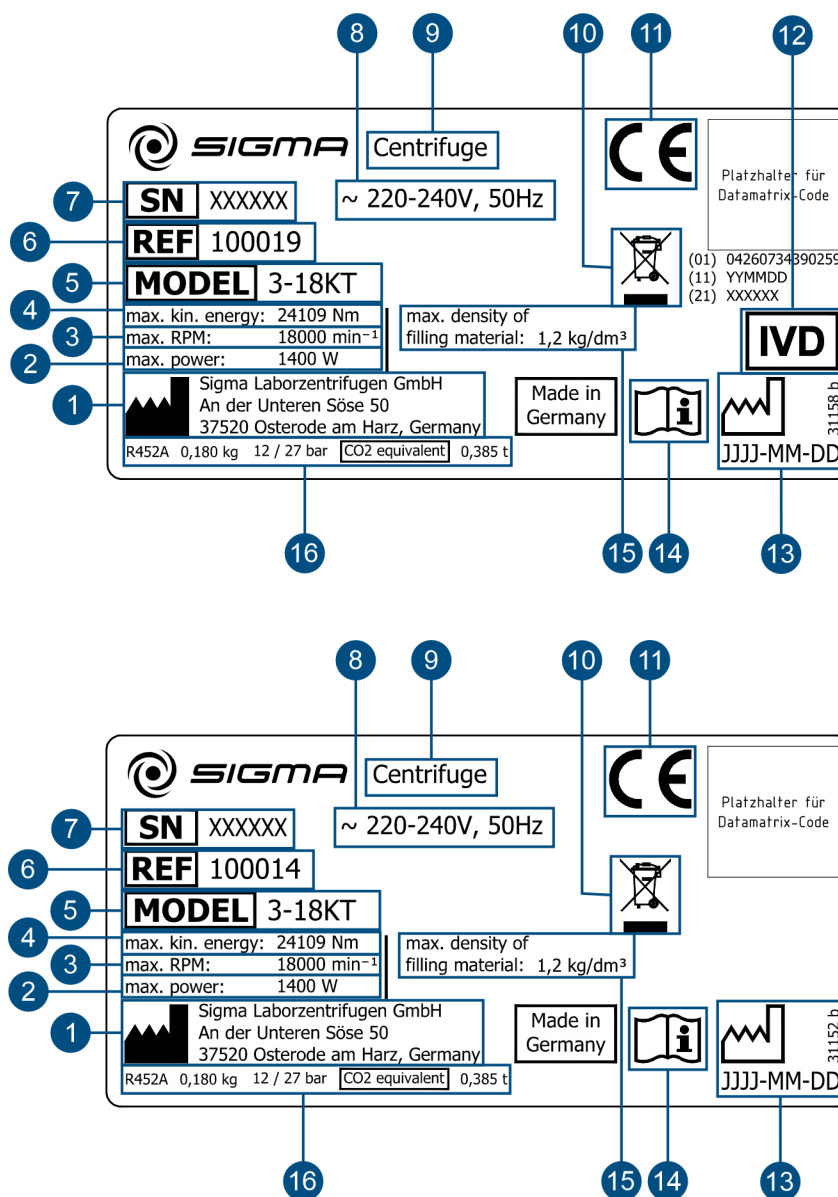


Fig. 9: Name plate (example)

|   |                     |    |                                                                                |
|---|---------------------|----|--------------------------------------------------------------------------------|
| 1 | Manufacturer        | 10 | Symbol for separate disposal                                                   |
| 2 | Power consumption   | 11 | CE mark                                                                        |
| 3 | Max. speed          | 12 | IVD mark as per the regulation (EU) 2017/746 (if applicable)                   |
| 4 | Max. kinetic energy | 13 | Date of manufacture                                                            |
| 5 | Product name        | 14 | Symbol indicating that compliance with the operating instructions is mandatory |
| 6 | Part number         | 15 | Max. permissible density                                                       |
| 7 | Serial number       |    |                                                                                |
| 8 | Nominal voltage     |    |                                                                                |
| 9 | Product designation |    |                                                                                |

## 4 Packaging, transport and storage

### 4.1 Safety



#### **WARNING**

#### **Risk of injury when lifting or carrying the centrifuge**

Serious injuries can occur when the centrifuge falls down or tips over.

- Use a suitable lifting device to lift the centrifuge.
- Comply with the fundamental health and safety rules and regulations as well as the rules and regulations for the prevention of accidents.
- Wear personal protective equipment (e.g. safety shoes and work gloves).

### 4.2 Symbols on the packaging

The following symbols, which are described hereinbelow, are used on the packaging of the centrifuge:

| Symbol                                                                              | Description                                                                                |
|-------------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------|
|   | Fragile goods                                                                              |
|  | Keep away from moisture                                                                    |
|  | Top                                                                                        |
|  | Temperature limit                                                                          |
|  | Recycling symbol for corrugated cardboard                                                  |
|  | RECY symbol (recycling)                                                                    |
|  | Symbol indicating that the product is a medical device as per the regulation (EU) 2017/746 |

Table 3: Symbols on the packaging

### 4.3 Unpacking the centrifuge

The centrifuge is packaged in a corrugated cardboard container (consisting of a cover, frame and base on a pallet). The rotors and accessories come supplied in a separate cardboard box that is fastened to the corrugated cardboard container.

1. Take the separate cardboard box off the corrugated cardboard container.
2. Remove the cover.
3. Remove the packaging material.
4. Lift the frame upwards to remove it.
5. Lift the centrifuge off the base by way of a suitable lifting device.

 Retain the original packaging for any possible future transport of the centrifuge.

### 4.4 Transporting the centrifuge

 Remove the rotor whenever the centrifuge needs to be transported.

#### **In-house transport**

The centrifuge can be transported without packaging over short distances within a building.

#### **Onward transport or return transport**

- The centrifuge must be packaged in its original or equivalent packaging.
- The symbols on the original packaging, indicating the necessary transport and environmental conditions, must be observed at all times.

 Please note/use the information and forms for returns (see chapter 8.5 - "Returning the centrifuge to the manufacturer").

### 4.5 Storing the centrifuge

The centrifuge can be stored for up to one year under the following conditions:

- The centrifuge is stored in its original packaging.
- It is stored in dry conditions and free from dust, mechanical load and strong temperature fluctuations.
- Comply with the permissible storage temperature of -20°C to +60°C.
- The requirements concerning the location of use as well as the ambient conditions must be fulfilled (see chapter 5.1 - "Setting up the centrifuge", chapter 10.1 - "Ambient conditions").

 If the centrifuge needs to be stored for more than one year, consult with the manufacturer.

## 5 Setup and connection

### 5.1 Setting up the centrifuge

#### Location requirements

To ensure safe operation, the location of use must fulfil the following requirements:

- The workbench is stable and resonance-free.
- If a mobile centrifuge bench is used, the locking castors are locked.
- The worktop is level and horizontal.
- The centrifuge stands with all four feet firmly on the worktop.
- The centrifuge is not subjected to thermal stress (e.g. caused by heat generators) or UV radiation (e.g. exposure to direct sunlight).


**ATTENTION** – When transporting the centrifuge from a cold to a warmer environment, condensation forms. Let the centrifuge dry completely (24 hours minimum) before using it.

#### Ensuring sufficient ventilation

The entire energy supplied to the centrifuge is converted into heat and released into the ambient air.

**ATTENTION** – The samples may overheat if the air inlet is blocked. Observe the following instructions to ensure adequate ventilation.

- Ensure the fulfilment of the specified climatic conditions (see chapter 10.1 - “Ambient conditions”).
- Set the centrifuge up so that it is freely accessible. Keep a safety range of at least 30 cm free around the centrifuge as well as with regard to walls or other devices so that the vents of the centrifuge remain unobstructed and fully effective.

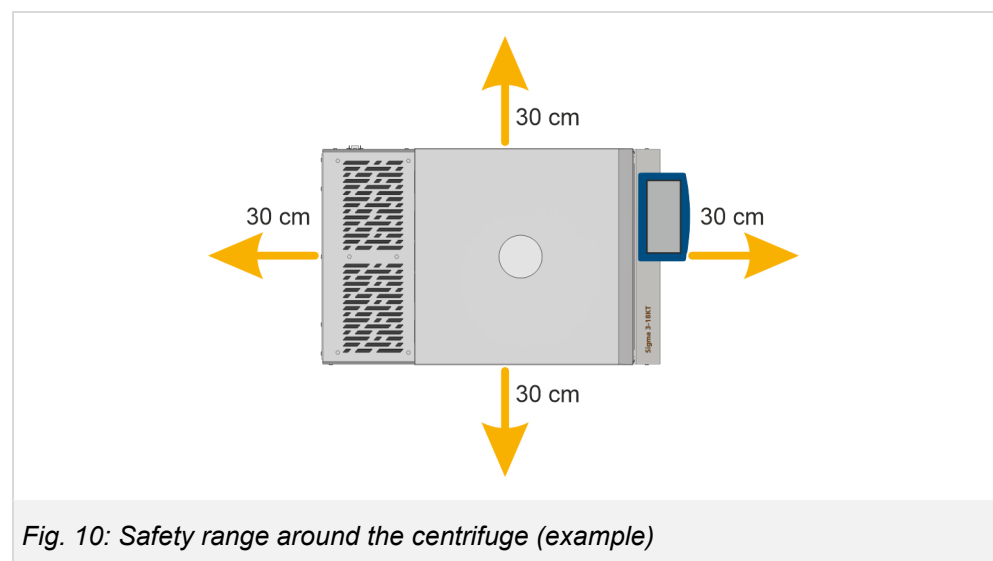

## 5.2 Establishing the power supply



The mains power plug serves as a disconnecter.  
It must be accessible at all times!

**ATTENTION** – The centrifuge may be damaged if the operating voltage is incorrect. Ensure that the operating voltage on the name plate corresponds to the local supply voltage.

Sigma laboratory centrifuges are units of protection class I. The units of this model series have a three-wire power cord with an IEC C13 connection.

- ▶ The removable power cord must not be longer than 3 m.
- ▶ The original power cord must not be replaced with an insufficiently rated one.
  1. Ensure that the power cord and on-site socket are not visibly damaged.
  2. Plug the power plug into the socket.

## 5.3 Fuses inside the unit

The centrifuge is equipped with a mains power switch with an integrated thermal circuit breaker.

After the fuses have tripped:

1. Switch the centrifuge off by actuating the mains power switch and let it cool down for approx. 2 minutes.
2. Switch the centrifuge back on.
  - The centrifuge is now ready for operation.

## 5.4 On-site fuse protection

- The centrifuge must be typically protected with an on-site type B 16 A fuse.
- An AC/DC-sensitive RCD (residual current device) is required in the building's electrical installation to ensure a safe shutdown in the event of a fault.

## 6 Operation

### 6.1 Safety



In an emergency, switch the centrifuge off immediately and disconnect the mains power plug!

### 6.2 Switching the centrifuge on and off

#### Power on

► The centrifuge has been properly set up and installed.

1. Press the mains power switch.
  - The touch panel lights up. The centrifuge is ready for operation.

#### Power off

► The centrifuge is at a standstill.

1. Open the lid so that any potential moisture can escape.
2. Press the mains power switch.
  - The centrifuge is off.



After the centrifuge has been switched off, wait at least 30 seconds before switching it back on.




If the centrifuge remains switched off for a longer period of time (> 3 months), its date and time must be checked and reset when the centrifuge is switched on again.

### 6.3 Opening and closing the lid

#### Opening the lid

► The centrifuge is at a standstill.

1. Select the "Lid" button.
  - The lid opens approximately halfway.
2. Push the lid manually further into its opening direction.
  - The lid opens completely.

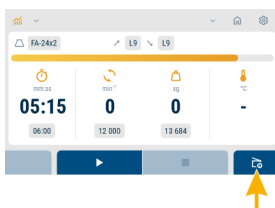

**WARNING** – Risk of injury caused by defective gas struts. The lid can fall down on its own, potentially crushing fingers or limbs. Have the gas struts checked regularly and replaced, if necessary.

### Closing the lid

**CAUTION** – Risk of injury when closing the lid. Fingers or limbs may be crushed. Never reach between the lid and the housing.

1. Push the lid down until you hear the lid lock engage.
  - The lid is now closed correctly.

## 6.4 Using rotors, buckets, carriers and accessories

---

### 6.4.1 Preparing the accessories

#### Vessels

- The vessels have been authorised by Sigma Laborzentrifugen GmbH for use in this centrifuge.
1. Fill the vessels outside of the centrifuge.
  2. Fill the vessels carefully and only up to the specified maximum fill level.
  3. Seal the vessels.

#### Maximum speed for tubes

Some tubes, such as centrifuge glass tubes, microtubes, culture tubes, PTFE tubes and especially high-volume tubes can be used in the rotors, buckets, and adapters of Sigma Laborzentrifugen GmbH at higher speed than their breaking limit.

**ATTENTION** – When using glass vessels, do not exceed 4,000 xg. The only exceptions are heavy-duty centrifuge vials, in which case the manufacturer's instructions apply.

#### Adapters

Adapters enable the use of vessels of different sizes in a rotor.

- The adapters have been produced by Sigma Laborzentrifugen GmbH and are authorised for use in this centrifuge.
- Load the adapters with the same number of vessels and the same weights.

**ATTENTION** – When using rotors for microtiter plates, never place the removal tool in the bucket without the plate.

## 6.4.2 Preparing and installing a rotor



### WARNING

#### Hazard caused by insufficient fastening of the rotor

If a rotor is not properly fastened, it may detach itself from the shaft, thereby causing a crash.

- Prior to starting a centrifugation run, ensure that the rotor tie-down screw has been properly tightened.
- Loosen and tighten the rotor tie-down screw at regular intervals (depending on the frequency of use) to ensure a proper connection between the rotor holder and the motor shaft.



### WARNING

#### Risk of toxic, pathogenic or radioactive contamination

Toxic, pathogenic or radioactive substances can cause contamination if they are not handled properly.

It may be necessary to clean and/or decontaminate the laboratory centrifuge, rotors and accessories for reasons of safety before they can be maintained, repaired or handed over.

- Disinfect and clean the centrifuge, rotors, buckets, carriers and accessories after a contamination.
- Comply with the relevant safety instructions concerning the handling of toxic, pathogenic or radioactive substances.
- Wear personal protective equipment.
- Document any measures taken (see chapter 8.5 - "Returning the centrifuge to the manufacturer" / Declaration of decontamination)



### CAUTION

#### Hazard due to the improper loading of the rotor

Rotors with insufficiently balanced loads can lead to overload as a result of imbalance during operation.

- Always load the rotor in a rotationally symmetrical manner with even weight distribution.



The design of Sigma accessories with the "Biosafe" marking ensures that no liquids or aerosols can leave the rotor or bucket if a vessel breaks during the centrifugation.


#### Angle rotor (Biosafe type)

##### Loading an angle rotor


- ▶ The rotor and covers, including all of the rubber seals, are clean and undamaged.
- ▶ All of the seals are properly in place.
- ▶ The areas of the bioseals of the cover are absolutely dry.
- ▶ The sample vessels are filled and sealed.

## Operation

1. Insert the sample vessels into the holes of the angle rotor.
  - Load the rotationally symmetric compartments of the rotors with the same accessories and content.
  - When performing a centrifugation run at low capacity, distribute the sample vessels so that the load is rotationally symmetrical. Loading only one axis of the rotor is not permissible.

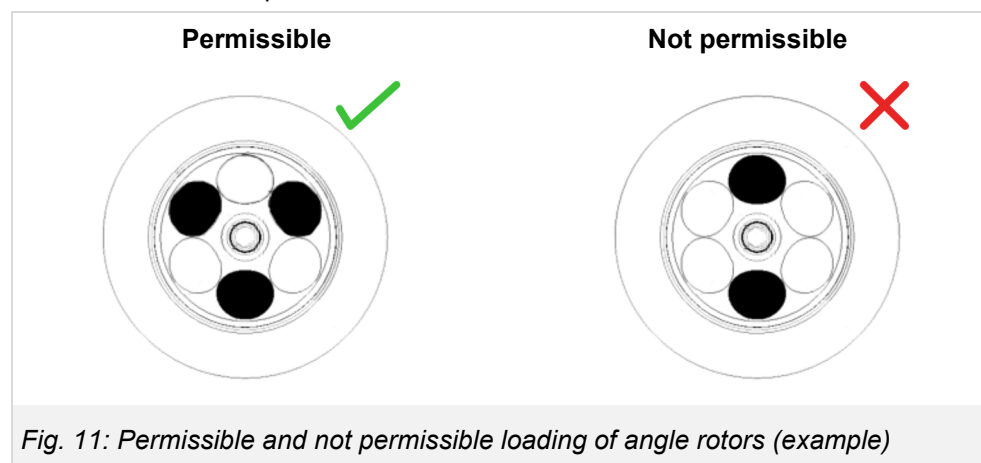

2. Position the cover horizontally on top of the rotor.
3. For sealing the rotor, turn the cover screw of the bayonet closure clockwise until the mark on the cover screw and the mark on the rotor are aligned (see the following picture).
  - The angle rotor is loaded.

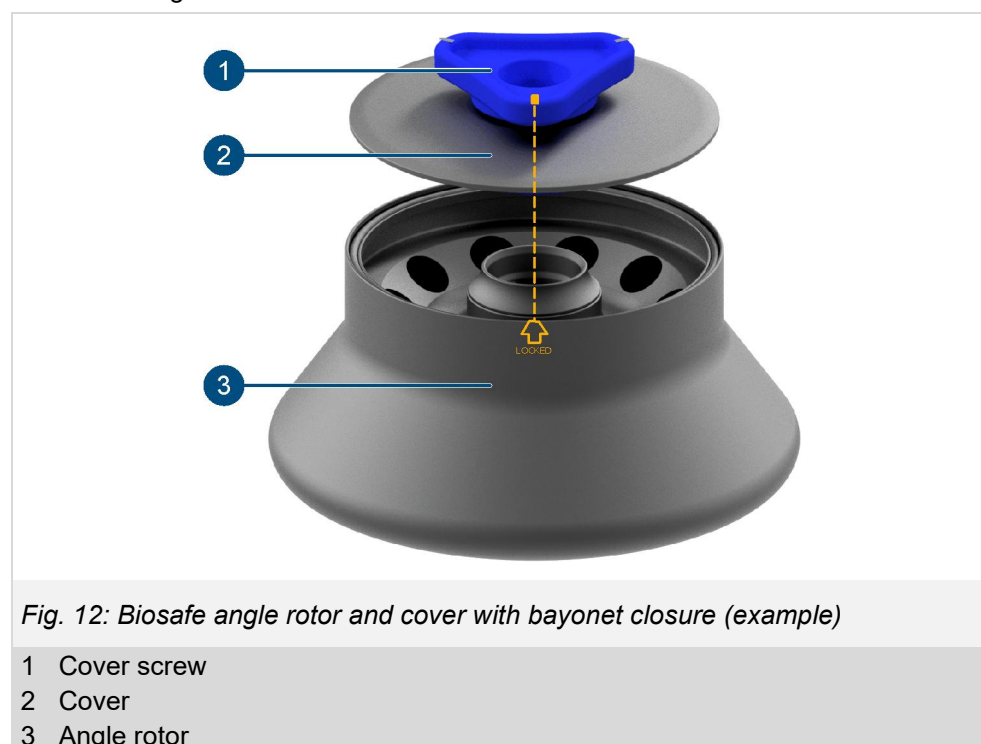

**ATTENTION** – The cover screw is solely intended for securing the cover on the rotor and not for fastening the rotor to the motor shaft.

#### Installing an angle rotor

- ▶ The centrifuge lid is open.
- ▶ The rotor chamber and the lid seal are clean and dry.
  1. Unscrew the rotor tie-down screw from the motor shaft (anticlockwise).
  2. Lower the rotor with its central bore onto the motor shaft straight from above.
  3. Hold the rotor on its outer rim and tighten the rotor tie-down screw clockwise by way of the rotor wrench. Tightening torque: 7.5 Nm.
  4. Close the centrifuge lid.
    - The angle rotor is installed.
    - The centrifuge is ready for operation.

 Using an angle rotor without a cover will result in increased noise levels.

 Angle rotors can also be used without a rotor cover. However, in this case, they do not fulfil the biosafety leak tightness in accordance with DIN EN 61010, part 2-020, appendix AA.

#### **Swing-out rotor (Biosafe type)**

##### Installing a rotor cross

- ▶ The centrifuge lid is open.
- ▶ The rotor chamber and the lid seal are clean and dry.
- ▶ The rotor cross is clean and undamaged.
  1. Unscrew the rotor tie-down screw from the motor shaft (anticlockwise).
  2. Lower the rotor cross without any buckets with its central bore onto the motor shaft straight from above.
  3. Hold the rotor on its outer rim and tighten the rotor tie-down screw clockwise by way of the rotor wrench. Tightening torque: 7.5 Nm.

##### Loading and installing the buckets


- ▶ The buckets as well as the rubber seals, bucket caps and adapters are clean and undamaged.
- ▶ All of the seals are properly in place.
- ▶ The areas of the bioseals of the buckets are absolutely dry.
- ▶ Only buckets of the same set are used (see the set number on the buckets).
- ▶ The load-bearing bolts of the rotor cross are greased or they have a Sigma Comfort rotor coating.
  1. Place the sample vessels in the bucket or in the adapter in the bucket.
    - As a general rule, all of the compartments must be filled with buckets.
    - Load the rotationally symmetric compartments of the rotors with the same accessories and content.
    - When performing a centrifugation run at low capacity, distribute the sample vessels so that the load is rotationally symmetrical.
  2. Place the caps on the buckets and screw them on in a clockwise direction.

 The buckets can also be used without a cover. However, in this case, they do not fulfil the biosafety leak tightness in accordance with DIN EN 61010, part 2-020, appendix AA.

## Operation




3. Place all of the buckets in the rotor cross. If the buckets have numbers, place them in the corresponding compartments of the rotor (the compartments are numbered as well).

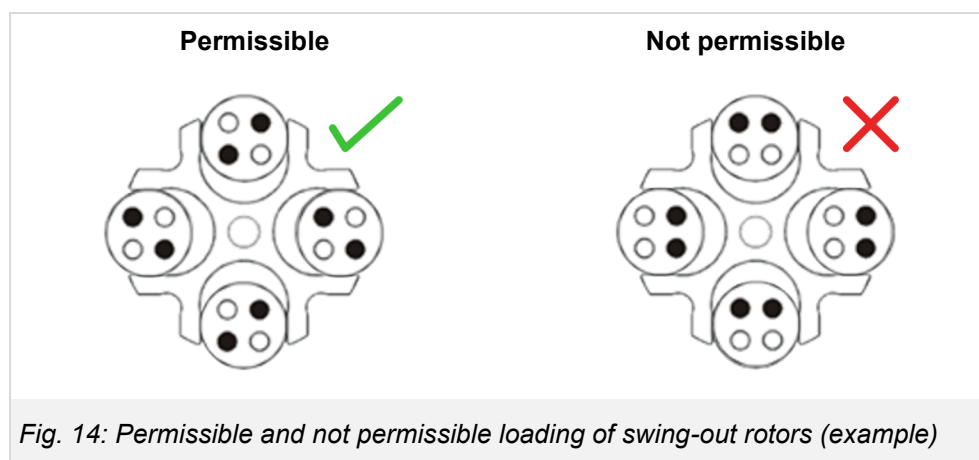

4. Close the centrifuge lid.



Comply with the instructions on the label in the centrifuge lid (see the following illustration)! The adhesive labels must be kept readable at all times. If necessary, they must be replaced.

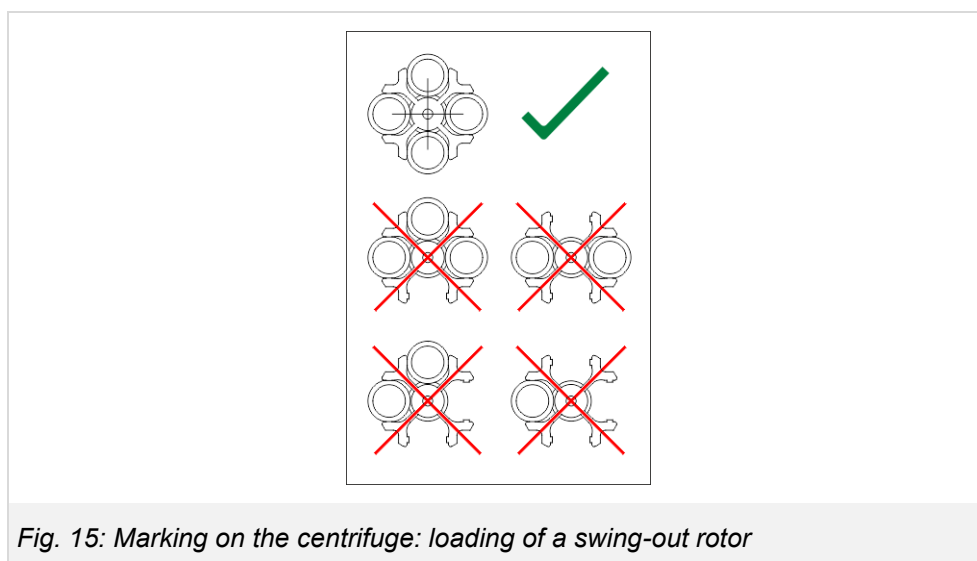

### 6.4.3 Removing the rotor and accessories



#### **WARNING**

#### **Risk of toxic, pathogenic or radioactive contamination when opening the vessels**

During centrifugation, overpressure may build up in the vessels due to heat generated by air friction.

- This is why it must be ensured that there is no liquid in the cap/sealing area prior to opening the vessels.



#### **CAUTION**

#### **Risk of burns on hot surfaces inside the rotor chamber**

After a centrifugation, the surfaces inside the rotor chamber and, in particular the motor shaft, can be hot.

- Wear suitable protective gloves.
- Do not touch the surfaces without good reason.

#### **Angle rotor (Biosafe type) removal**

► The centrifuge is at a standstill and the lid is open.

1. Unscrew the rotor tie-down screw by turning it anticlockwise by way of the rotor wrench.
2. Lift the rotor with its cover perpendicularly off the motor shaft.
3. Loosen the cover by turning the bayonet lock anticlockwise and then remove it.
4. Remove the sample vessels carefully in order to prevent their contents from being mixed.

## Operation

---

### Swing-out rotor (Biosafe type) removal

- The centrifuge is at a standstill and the lid is open.
  - 1. Remove the buckets with their caps from the rotor cross.
  - 2. Remove the caps from the buckets.
  - 3. Remove the sample vessels carefully in order to prevent their contents from being mixed.
  - 4. Unscrew the rotor tie-down screw (anticlockwise) by way of the rotor wrench and lift the rotor cross perpendicularly off the motor shaft.
- Or:
- Fill the buckets with sample vessels again, close them and reinstall them.


## 6.5 User interface

The centrifuge is operated by way of a touch screen. Buttons are activated by touch. Depending on the specific function, an additional window with options or input fields will open. Numerous symbols are provided for the navigation within the menus.

### 6.5.1 Main screen

The main screen is the default display of the control system. This is where all of the transactions are started.

The main screen includes several areas:

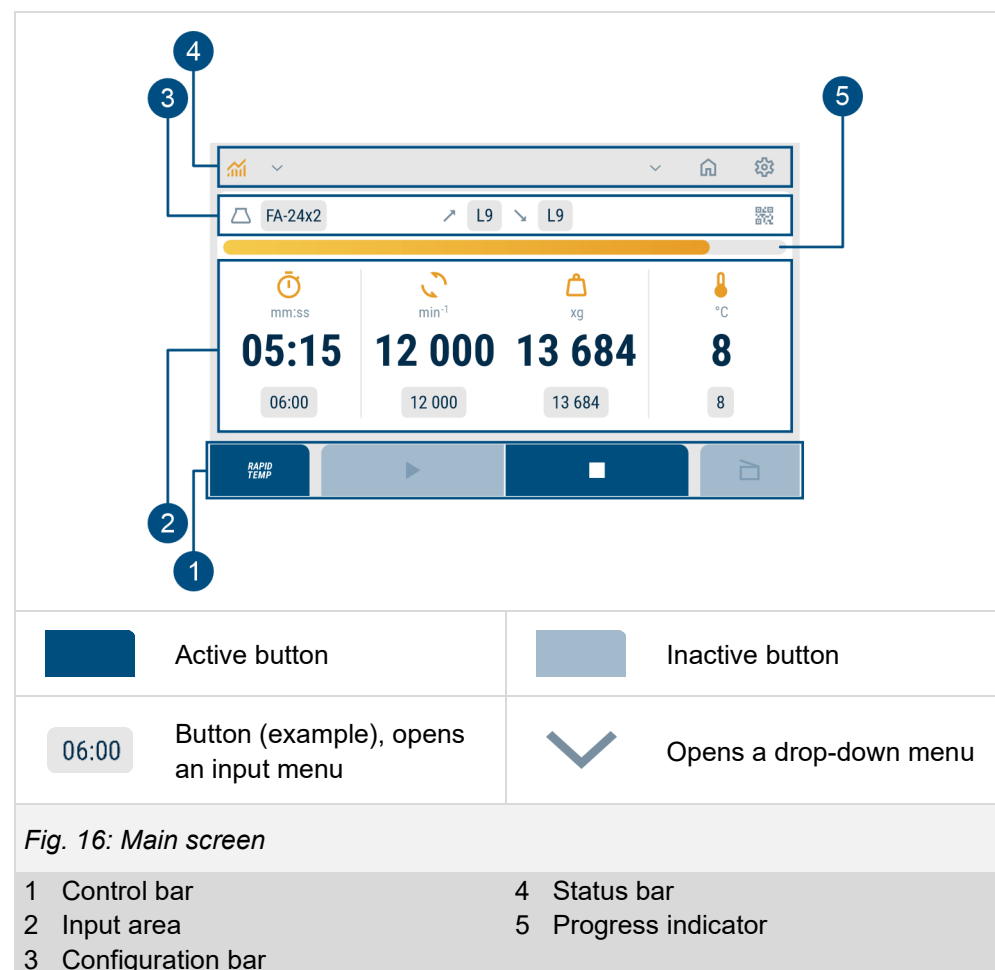

#### Control bar (1)

The buttons of the control bar are used to control the basic functions of the centrifuge as well as the rapid temperature control function:

- Starting a centrifugation (see chapter 6.5.2.1)
- Stopping a centrifugation (see chapter 6.5.2.2)
- Opening the lid (see chapter 6.3)
- *RAPID TEMP* (see chapter 6.5.3.4)

## Operation

### Input area (2)

The input area of the main screen is used to configure the following process parameters:

- Runtime (see chapter 6.5.3.2)
- RPM (see chapter 6.5.3.3)
- Relative centrifugal force RCF (see chapter 6.5.3.3)
- Temperature (see chapter 6.5.3.4)

The individual fields show the current value of the parameter and a button for configuring the set value of the corresponding process parameter.

Inputs are also performed in this area in every subsequent menu.

### Configuration bar (3)

The configuration bar includes the following buttons:

- Rotor and bucket selection (see chapter 6.5.5)
- Profiles (see chapter 6.5.6)
- QR code generation for the Sigma Timer app (see chapter 6.5.10.1)


### Status bar (4)

In addition to the status indicator, the status bar includes the following buttons:

- Program (see chapter 6.5.7)
- Status menu (see chapter 6.5.8)
- Settings (see chapter 6.5.10)
- Home (this button is used to return to the main screen from any other screen).

## 6.5.2 Basic functions


- The centrifuge is on.

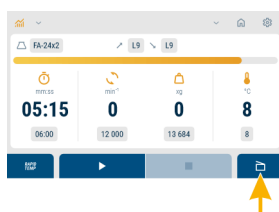

### Opening and closing the lid

See chapter 6.3 - "Opening and closing the lid"


### 6.5.2.1 Starting a centrifugation

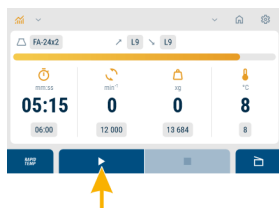

- The rotor has been prepared for centrifugation.
  1. Open the lid.
  2. Install the rotor and briefly turn it manually in the direction of rotation (anticlockwise).
    - The system identifies the rotor and displays it on the main screen.
  3. Close the lid.

4. Configuring the centrifugation parameters:
  - Runtime (see chapter 6.5.3.2 - “Runtime: process parameters”)
  - Speed or RCF (see chapter 6.5.3.3 - “RPM and RCF: process parameters”)
  - Temperature (see chapter 6.5.3.4 - “Temperature: process parameters”).
5. Select the *Start* button.
  - The centrifugation run starts.

### 6.5.2.2 Stopping a centrifugation

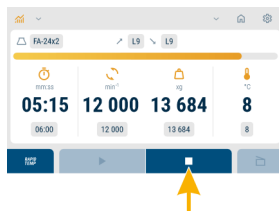

1. Wait for the set time to elapse.
- Or:
2. Select the *Stop* button.
    - The centrifugation will be stopped prematurely.

### Performing a quick stop

During a quick stop, the centrifuge decelerates with maximum deceleration power (curve L9).

- The centrifuge is on.
  - A centrifugation run is performed.
1. Select the *Stop* button and hold it for a minimum of three seconds.
    - The centrifugation will be stopped prematurely.
    - A message will be displayed:

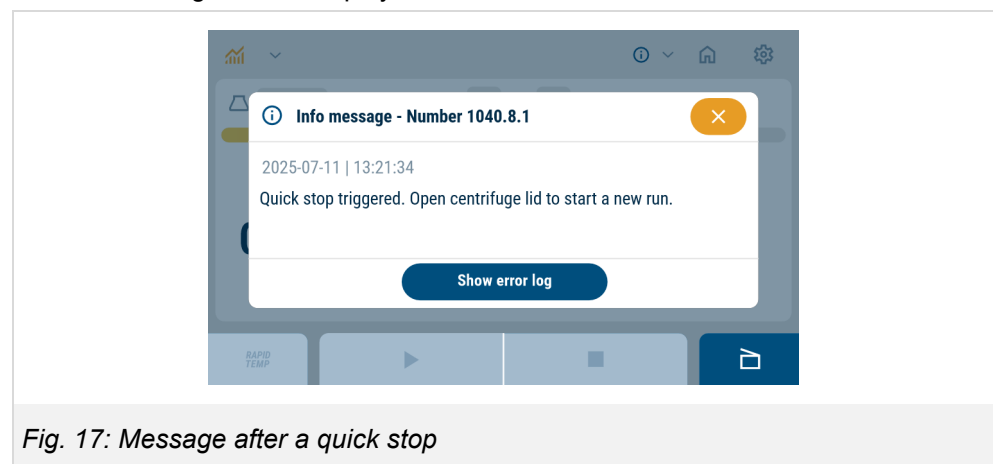


Fig. 17: Message after a quick stop

### 6.5.3 Process parameters

The process parameters area shows the set and actual values of the runtime, speed, RCF and temperature as well as information about the activated functions.

#### 6.5.3.1 Entering process parameters

Process for entering a process parameter (example: runtime)

► The user interface shows the main screen.


1. Select the *Set value* button (A).

– The input menu opens.

2. Enter the set value by way of the keyboard (B).

3. Confirm the input (C).

– The main screen will be displayed.

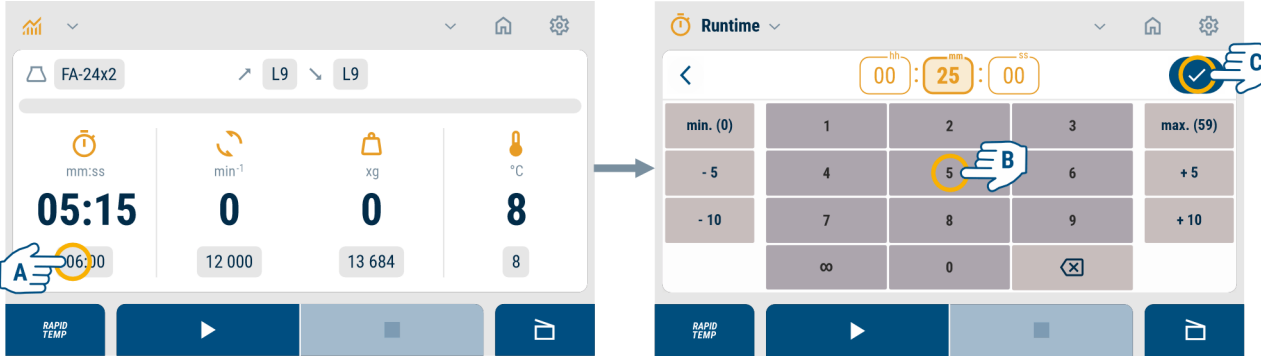

|                                                                                     |                                                     |                                                                                      |                                                |
|-------------------------------------------------------------------------------------|-----------------------------------------------------|--------------------------------------------------------------------------------------|------------------------------------------------|
|  | Set value button                                    |   | Active input field (example)                   |
|  | Inactive input field (example)                      |   | Input error (example)                          |
|  | Minimum (min.) and maximum (max.) permissible value |  | Adds 5 or 10 units to the selected value       |
|  | Subtracts 5 or 10 units from the selected value     |   | Confirms the input and closes the current view |
|  | Opens the previous view without saving              |   | Opens the main screen                          |

*Fig. 18: Entering process parameters (example)*

### 6.5.3.2 Runtime: process parameters

-  In the runtime area, the upper value indicates the remaining runtime (current value) with the preselected centrifugation time (set value) below. Values up to 99 h:59 min:59 s maximum can be entered. The runtime is counted down from the set value. From 59 min:59 s downwards, the unit switches from "h:min" to "min:s".

In the factory settings, the function *Runtime at set speed* is activated so that the runtime is not counted until the set speed is reached. Then, the countdown continues until the deceleration phase starts. The corresponding symbol is displayed in the runtime area on the main screen. The function can be deactivated (see chapter 6.5.4.2 - "Runtime: process-related parameters").

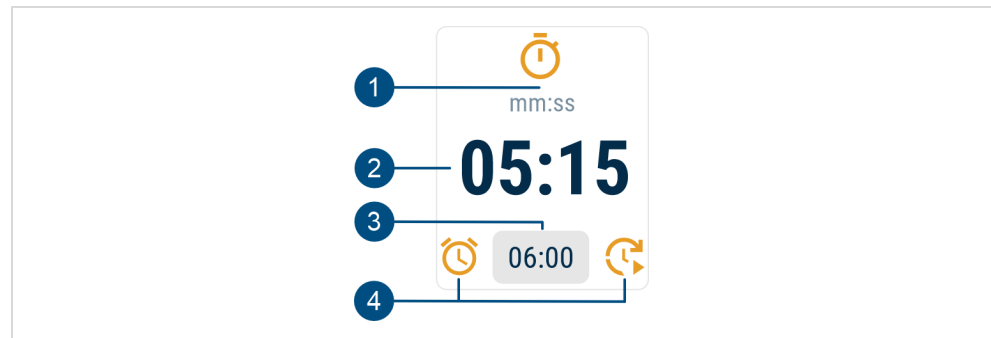



Fig. 19: Runtime area of the main screen (example)

|                  |                                                 |
|------------------|-------------------------------------------------|
| 1 Runtime symbol | 3 <i>Set value</i> button for the runtime       |
| 2 Current value  | 4 Symbols of the activated functions (examples) |

#### Performing a short run

- The centrifuge is at a standstill.
  - The user interface shows the main screen.
1. Select the *Start* button and hold it for the desired duration of the short run (A).
    - The short run starts when the colour of the button has completely changed from blue to yellow (B).
    - The centrifuge accelerates with maximum possible power (curve L9) up to the set speed.
    - The runtime is counted upwards.

## Operation

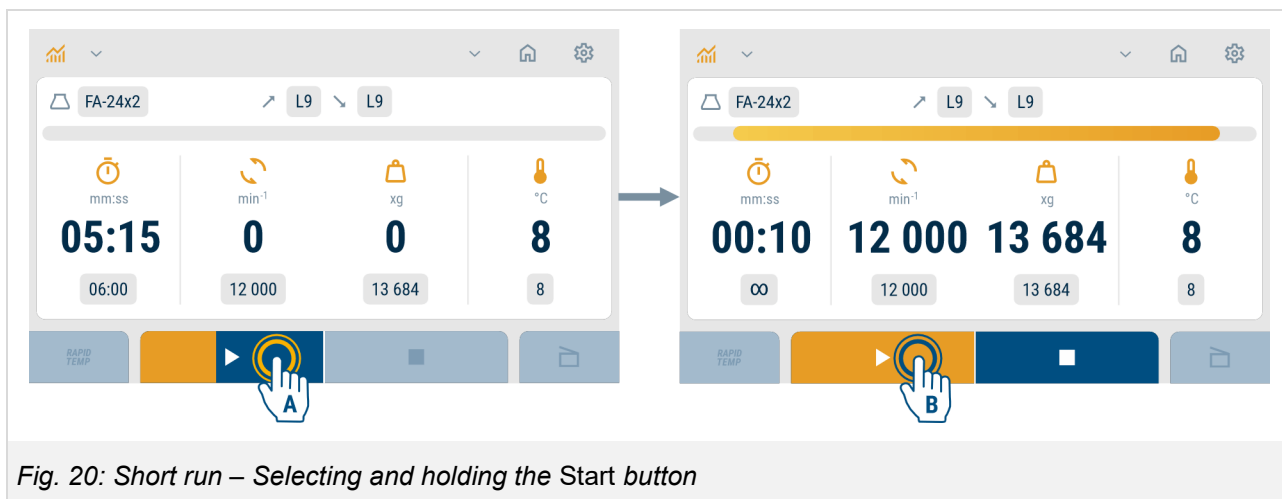


### 2. Release the *Start* button.

- The centrifuge decelerates with maximum deceleration power (curve L9) down to a standstill.

## Performing a continuous run

During the continuous run, the runtime of the centrifuge is unlimited and must be stopped manually. The centrifuge accelerates during the continuous run until the set speed is reached.

- The centrifuge is at a standstill.
  - The user interface shows the main screen.
1. Select the *Set value* button for the runtime.
    - The input menu opens.
  2. Select the ∞ button on the keyboard (A) and confirm the input (B).

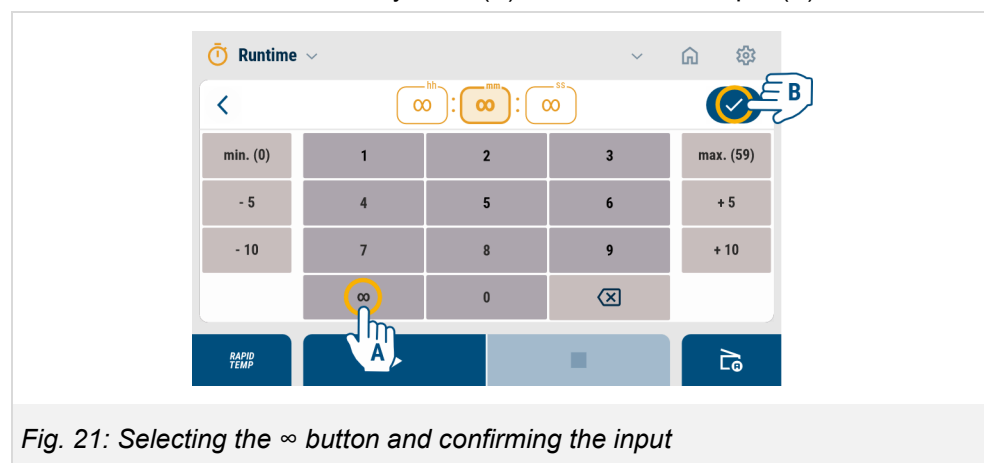

- The main screen will be displayed.
  - The ∞ symbol will be displayed for the set runtime.
3. Select the *Start* button to start a centrifugation run.
    - The continuous run is activated and the runtime will be counted upwards.

4. Select the Stop button to stop the continuous run.
  - The centrifuge decelerates to a standstill.

### 6.5.3.3 RPM and RCF: process parameters



In the RPM and RCF fields, the upper value indicates the current RPM or RCF value. The values below are the preselected (set) values.

The RPM values are indicated in revolutions per minute ( $\text{min}^{-1}$ ) and the RCF values as multiples of the acceleration due to gravity (xg). The values are interdependent. The maximum RPM/RCF values depend on the rotor that is used.



Fig. 22: RPM and RCF areas of the main screen (example)

- |                         |                                         |
|-------------------------|-----------------------------------------|
| 1 RPM/RCF symbol        | 3 Set value button of the RPM/RCF value |
| 2 Current RPM/RCF value |                                         |

### 6.5.3.4 Temperature: process parameters



In the temperature area, the upper value indicates the current temperature (current value) with the preselected temperature (set value) below. Temperatures between  $-10^{\circ}\text{C}$  and  $+40^{\circ}\text{C}$  can be preselected.



The centrifuge is designed for spinning samples in the available rotors in a temperature range from  $+4^{\circ}\text{C}$  to  $23^{\circ}\text{C}$ . It does not have a heater, which is why temperatures above the rotor temperature depend on the air friction of the spinning rotor.

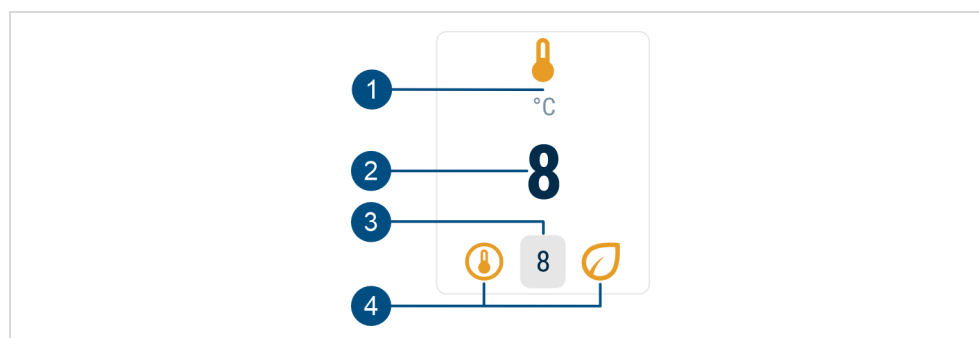


Fig. 23: Temperature area of the main screen (example)

- |                             |                                                 |
|-----------------------------|-------------------------------------------------|
| 1 Temperature symbol        | 3 Set value button for the temperature          |
| 2 Current temperature value | 4 Symbols of the activated functions (examples) |

**RAPID  
TEMP**

**Rapid temperature control function *RAPID TEMP***

The *RAPID TEMP* function is used to precool or preheat the system consisting of the rotor chamber and rotor or rotor/bucket combination. The cooling process is controlled by the refrigeration unit while the waste heat of the motor as well as the air friction of the spinning rotor are used for heating.

The system is designed for centrifuging samples at a defined temperature and holding the original sample temperature. The centrifuge is not intended for cooling or heating the samples.

**Starting the *RAPID TEMP* function**

- ▶ The centrifuge is on.
  - ▶ For rapid cooling, the current temperature of the rotor chamber is above the set value.
  - ▶ For heating, the current temperature of the rotor chamber is below the set value.
1. Select the *RAPID TEMP* button (A).
  2. Select the *Start* button (B).

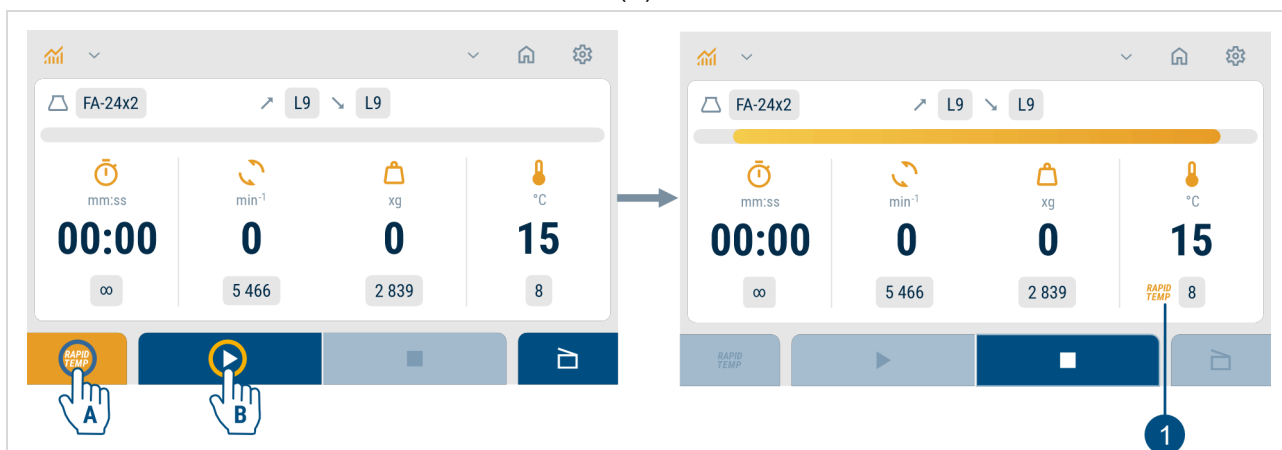


Fig. 24: Starting the *RAPID TEMP* function

1 *RAPID TEMP* symbol

- The *RAPID TEMP* function starts with the following parameters:
  - The set speed for quick cooling corresponds to 1/3 of the maximum rotor speed. The RCF value will be adapted accordingly.
  - The set value for heating corresponds to the maximum permissible rotor speed.
  - The acceleration and deceleration curves correspond to the curve L9.
  - The symbol *RAPID TEMP* is displayed in the temperature area.
  - The symbol for a continuous run is displayed in the runtime area. The runtime is counted upwards.
- The *RAPID TEMP* function stops when the set value is reached. The *RAPID TEMP* symbol disappears.
- After quick cooling, the standstill cooling function will be activated (see chapter 6.5.4.5 - “Temperature: process-related parameters”).
- The automatic lid opening function is suppressed in order to prevent repeated heating or cooling.

**Aborting the *RAPID TEMP* function**

1. Select the *Stop* button.

Or:

2. Make any type of entry.
  - The *RAPID TEMP* function will be aborted. No message will be issued when the set temperature is reached.

---

**6.5.4 Process-related parameters**

The drop-down menus of the process parameters can be used to configure additional process-related parameters.


**6.5.4.1 Entering process-related parameters**

Input based on the Quality Control temperature example


► The user interface shows the "Temperature" input menu.

1. Open the "Temperature" drop-down menu (A).
2. Activate the *Quality Control* function (B).
  - The colour of the *Value input* button changes from grey to yellow.
3. Select the *Value input* button (C) to enter the set values.
  - The input menu opens.
4. Enter the set values (D) and confirm the input (E).
  - The main screen will be displayed.
  - The QC temperature symbol will be displayed in the temperature area.

Operation

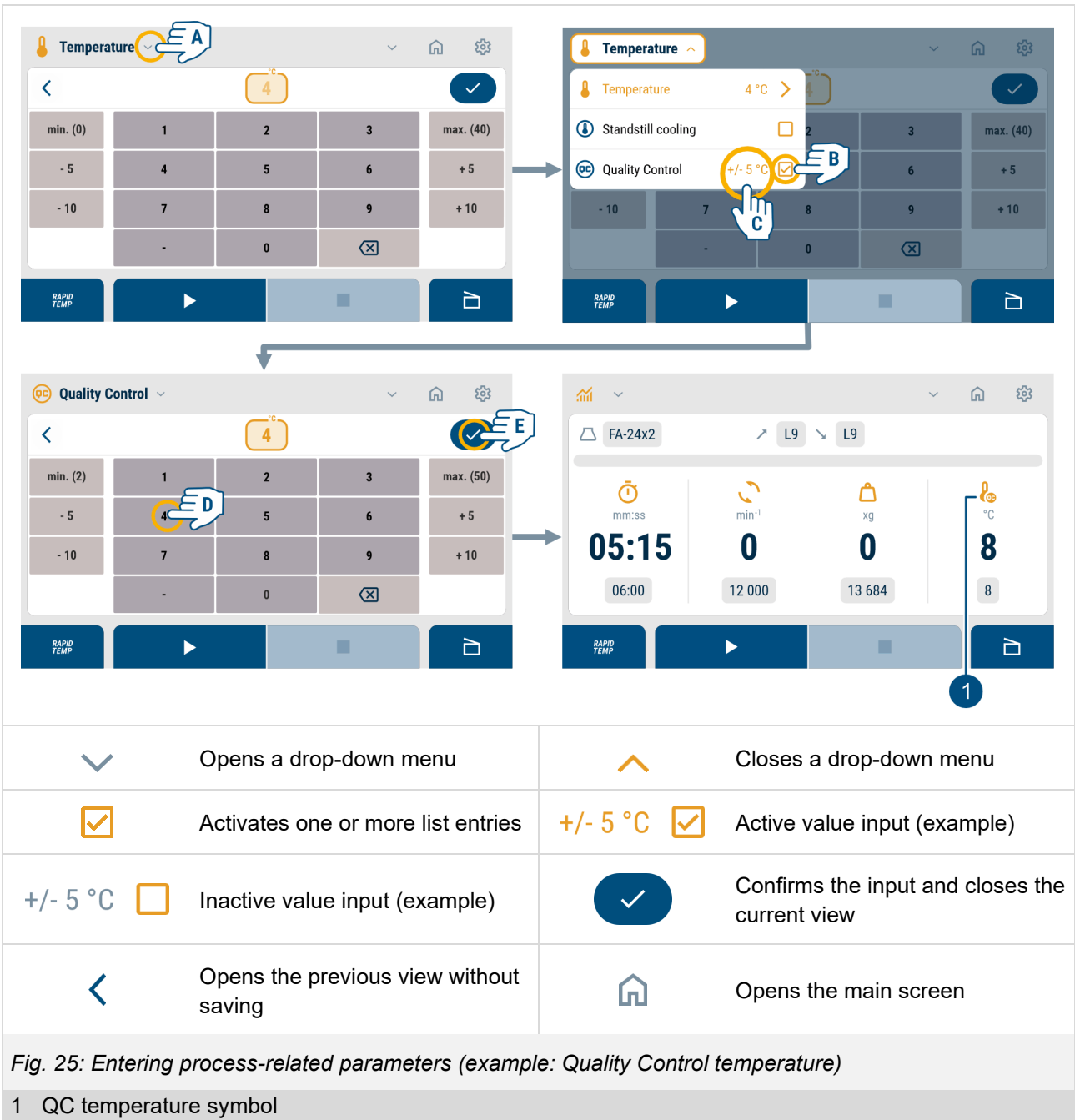


#### 6.5.4.2 Runtime: process-related parameters

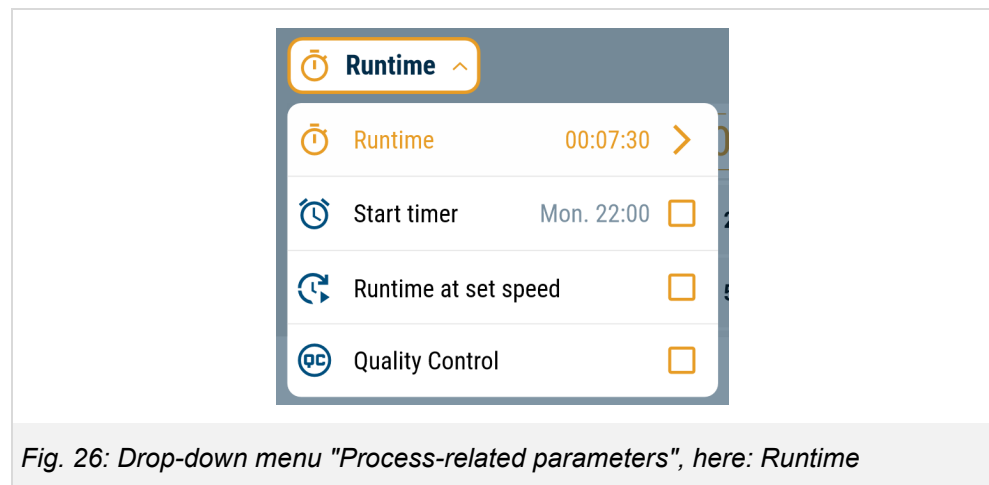

#### Start timer

When this function is active, the centrifugation run will start a specified time. The start time can be programmed up to 6 days, 23 hours and 59 minutes in advance.

- The centrifuge is ready for operation.
  1. Activate the start timer.
    - The *Value input* button will be enabled.
  2. Select the *Value input* button.
  3. Enter the start time and confirm the entry.
    - The main screen will be displayed.
    - In the runtime area, the start timer symbol will be displayed.
    - The *Stop* button flashes to indicate that the start timer is active.
    - The progress bar indicates the remaining time up to the start.
    - It is not possible to change parameters or perform any type of input.
    - The centrifuge will start at the programmed time.

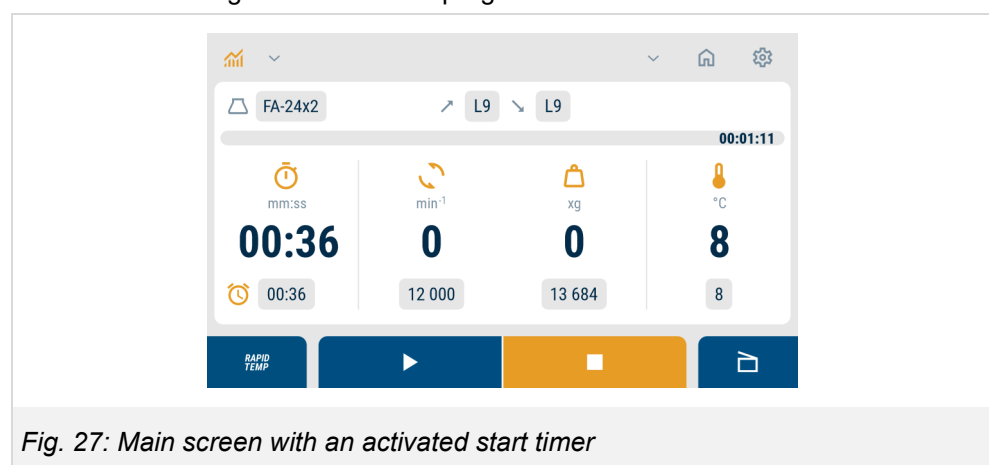

## Operation

The start timer mode will be stopped:

- when the centrifugation run starts as scheduled,
- when the *Stop* button is selected,
- if the start conditions are not fulfilled at the start time (e.g. power failure or error message).

In this case, the start timer must be reprogrammed.



### Runtime at set speed

If this function is active (factory setting), the runtime will not be counted until the set speed is reached. The "Runtime at set speed" symbol is displayed in the runtime area of the main screen.

If the function is deactivated, the runtime will be counted when the centrifuge starts. The symbol will disappear.



### Quality Control

If the function is active, the runtime will be continuously checked (in combination with the speed/RCF parameters) during the centrifugation run. QC is added to the runtime symbol on the main screen.

For detailed information about the Quality Control function, see chapter 6.5.10.8 - "Quality Control".



### Display deceleration time

If this function is active, the deceleration time will be displayed instead of the runtime during the deceleration process. The "Display deceleration time" symbol will be displayed instead of the runtime symbol in the runtime area of the main screen. In addition, the time will be displayed in grey and italics.

## 6.5.4.3 RPM: process-related parameters

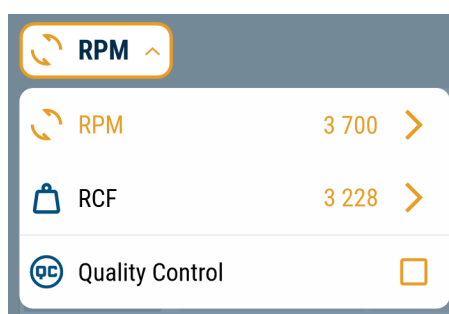


Fig. 28: Drop-down menu "Process-related parameters", here: RPM



### RCF

The RPM and RCF values are interdependent. The maximum values depend on the rotor that is used. This menu can be used to directly access the RCF input menu.

### **Quality Control**

If the function is active, the RPM will be continuously checked (in combination with the runtime and RCF parameters) during the centrifugation run. QC is added to the RPM symbol on the main screen.

For detailed information about the Quality Control function, see chapter 6.5.10.8 - "Quality Control".

#### 6.5.4.4 RCF: process-related parameters

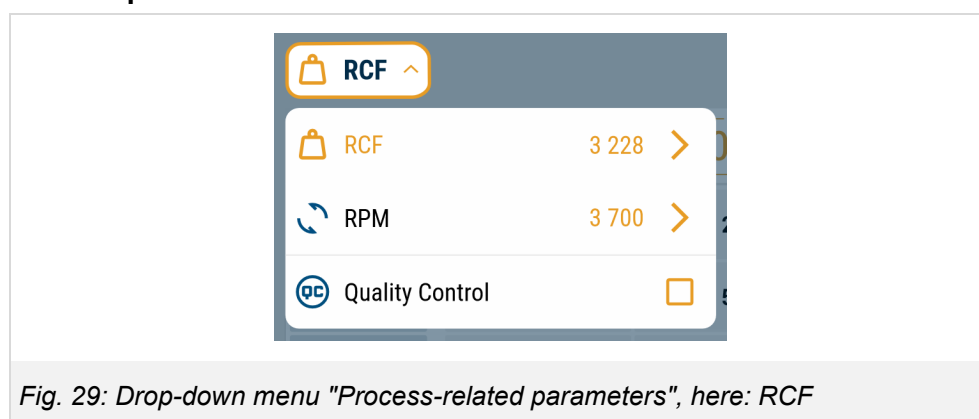


Fig. 29: Drop-down menu "Process-related parameters", here: RCF

### **RPM**

The RPM and RCF values are interdependent. The maximum values depend on the rotor that is used. This menu can be used to directly access the RPM input menu.

### **Quality Control**

If this function is active, the RPM will be continuously checked (in combination with the runtime parameter) during the centrifugation run and the corresponding RCF value will be calculated. QC is added to the RCF symbol on the main screen.

For detailed information about the Quality Control function, see chapter 6.5.10.8 - "Quality Control".

#### 6.5.4.5 Temperature: process-related parameters

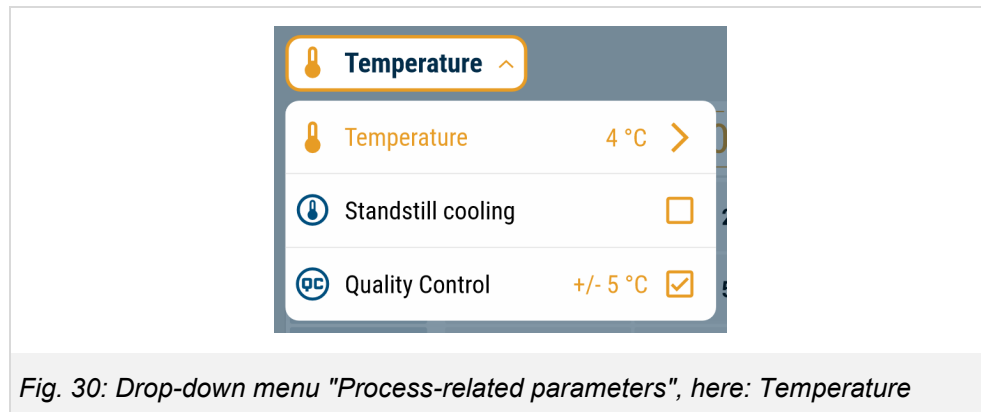


Fig. 30: Drop-down menu "Process-related parameters", here: Temperature

##### **Standstill cooling**

The standstill cooling function is activated in the factory settings so that the cooling process will be started when the centrifuge is switched on and the lid is closed and if the set temperature is below the current value. After the centrifugation run, the cooling function remains active in order to hold the samples on the set temperature.

The "Standstill cooling" symbol will be displayed on the main screen.

**ATTENTION** – Damage to property caused by freezing of the compressor, inaccurate measurement and control behaviour and freezing of the product. Do not cool the rotor below 0°C at a standstill.

##### **Quality Control**


If the function is active, the temperature will be continuously checked during the centrifugation run. QC is added to the temperature symbol on the main screen.

For detailed information about the Quality Control function, see chapter 6.5.10.8 - "Quality Control".

#### 6.5.4.6 Expert Mode: Process-related parameters


After Expert Mode has been activated (see chapter 6.5.10.1 - "Device settings"), additional process-related parameters will be enabled.

##### Runtime: Integral

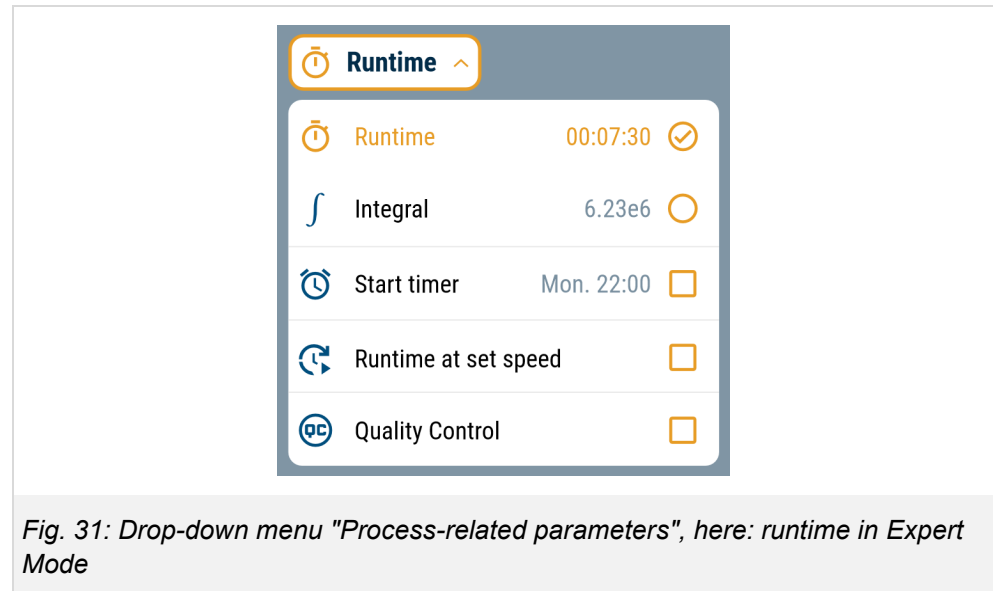

∫ The "Integral" function ensures that the centrifugal force acting on the samples over time (xg\*s) can be compared across different centrifuges, rotors and rotor loads. The function is used instead of a time value and the "Integral" area is shown instead of the "Runtime" area.

##### RPM and RCF: Radius reduction

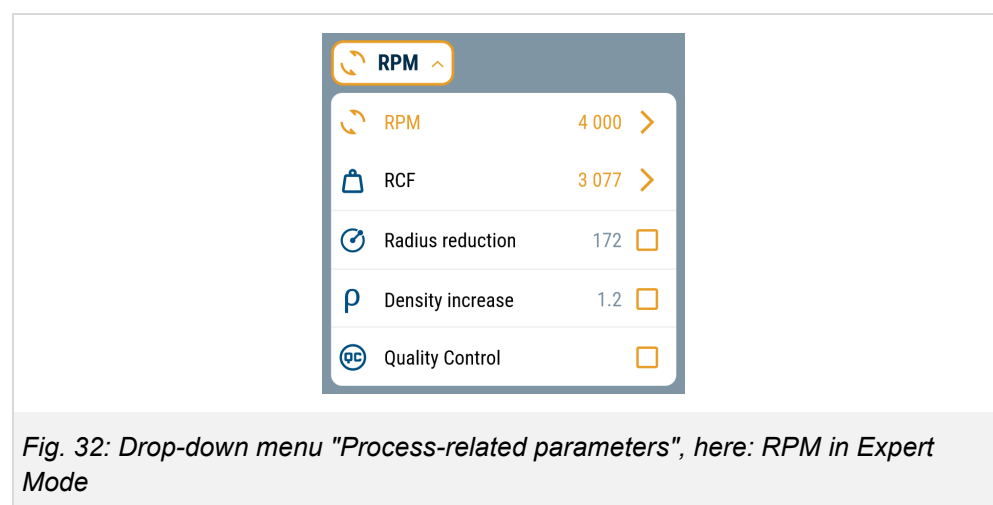

⌚ The value of the maximum possible radius is preset for every rotor and bucket. This function is used to reduce the radius of a rotor or bucket. This is relevant if the radius changes because a radius-reducing adapter is used. After the adaptation, the RCF acting on the sample will be displayed correctly in the RCF

area. In the "RPM" area, the radius reduction is indicated by the radius reduction symbol.

### **RPM and RCF: Density increase**

 The laboratory centrifuge is preset for the separation of substances of different densities in mixtures with a maximum density of 1.2 g/cm<sup>3</sup>. If the density of the liquid to be centrifuged is higher, the value must be adapted manually. When the density is increased, the RPM and RCF values will be automatically reduced to the maximum permissible value. In the "RCF" area, the density increase is indicated by the density increase symbol.

---

## **6.5.5 Rotor library**




The rotor library is a list of all the available rotors and buckets. The rotor or rotor/bucket combination that is used can be selected manually or it will be identified by the automatic rotor identification system.

For detailed information about all of the available rotors and rotor/bucket combinations, see chapter 11.1 - "Rotor program".

### **6.5.5.1 Selecting rotors and buckets**

Selection of a rotor based on the example of a rotor/bucket combination


- ▶ The user interface shows the main screen.
- ▶ A rotor/bucket combination is installed.
  1. Select the *Rotor & buckets* button (A).
    - The selection menu opens.
  2. Select the current rotor/bucket combination (B) and confirm the input (C).
    - The main screen will be displayed.

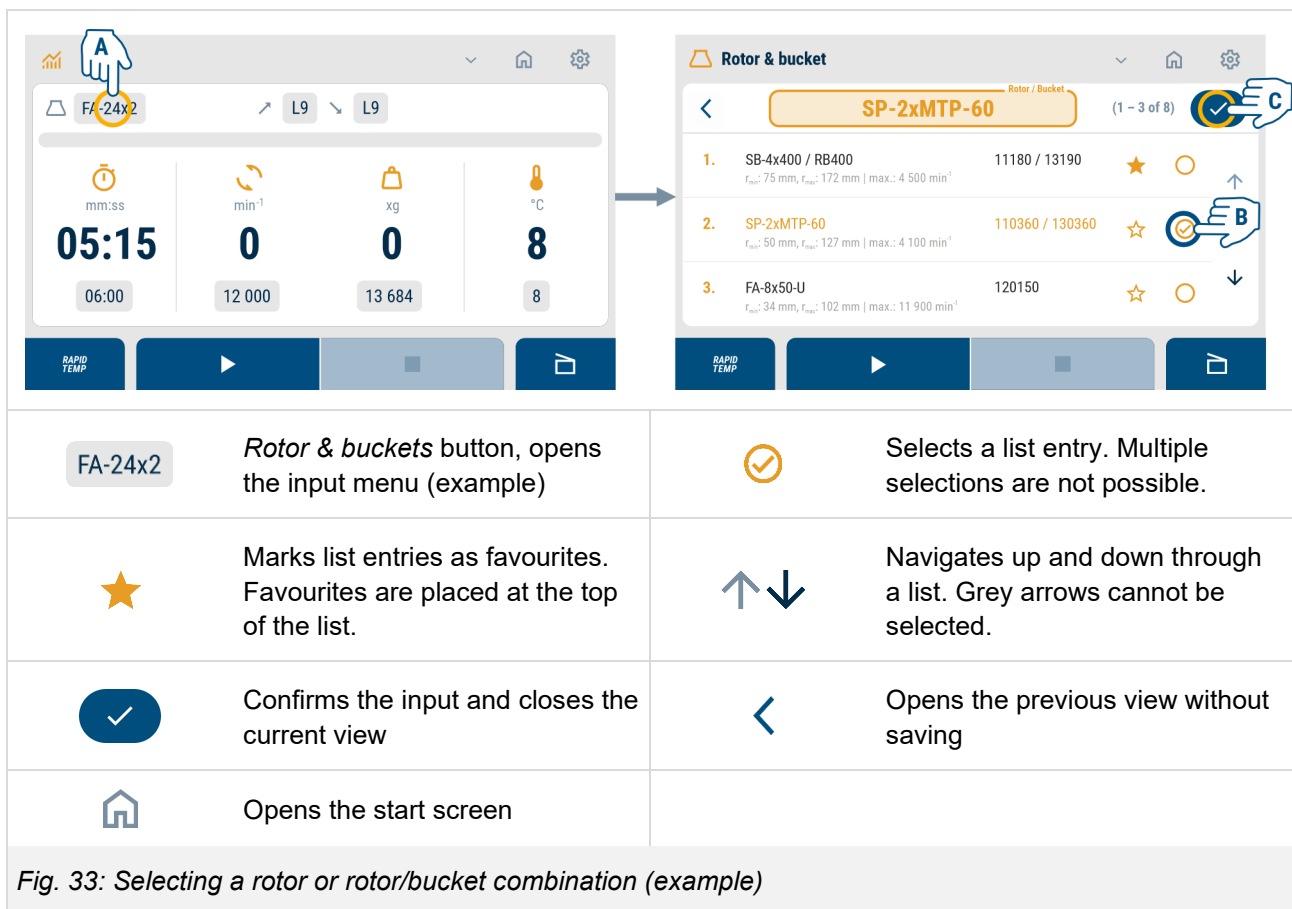

### Automatic rotor identification

If the rotor is turned manually in the direction of rotation (anticlockwise) or if a centrifugation run is started, the rotor code of the current rotor will be checked against the selected rotor.

- If the system identifies a rotor other than the selected one, the rotor input will be automatically corrected and an error message will be displayed.
- An error message will also be displayed if the system cannot identify the rotor. The rotor cannot be used in the centrifuge.

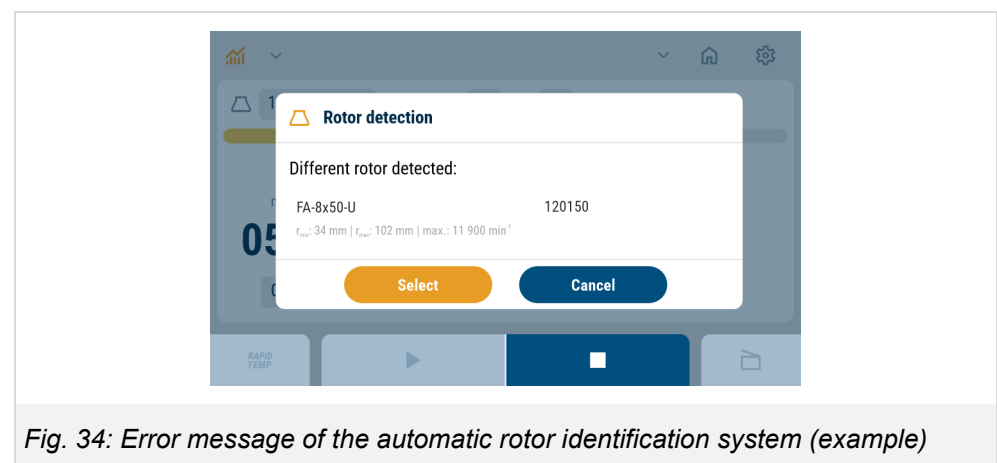

6.5.6 Profiles


↗ The *Curve selection (Profiles)* buttons open the menu for selecting the acceleration and deceleration curves. The system offers nine linear, nine quadratic and nine optimised quadratic acceleration and deceleration curves as well as a spin-out option.

For detailed information about the programmed curves, see chapter 11.4 - “Acceleration and deceleration curves”.

6.5.6.1 Selecting acceleration and deceleration curves

Selection of curves based on the example of a deceleration

- ▶ The user interface shows the main screen.
  1. Select the *Deceleration curve* button (A).
    - The input menu opens.
  2. Select a deceleration curve (B) and confirm the input (C).
    - The main screen will be displayed.

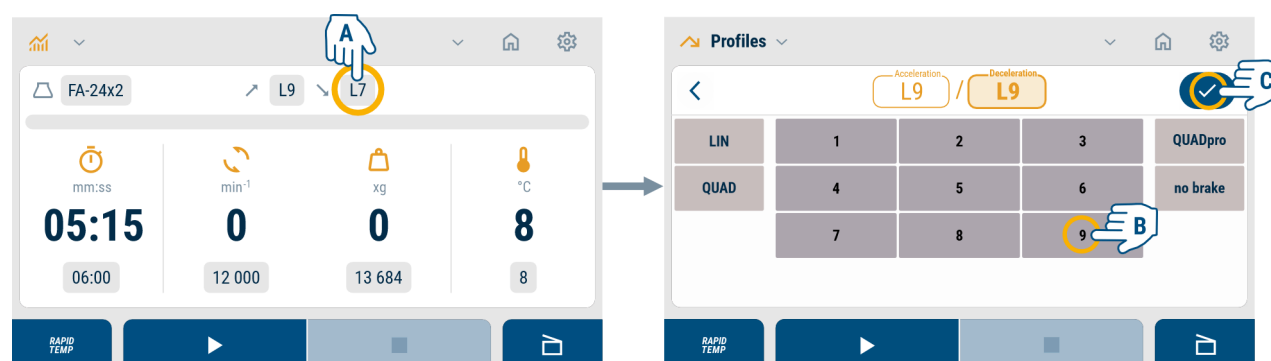

|                                                                                     |                                                 |                                                                                     |                                                |
|-------------------------------------------------------------------------------------|-------------------------------------------------|-------------------------------------------------------------------------------------|------------------------------------------------|
|  | Deceleration curve button, opens the input menu |  | Active input field (example)                   |
|  | Inactive input field (example)                  |  | Selects a linear curve (L)                     |
|  | Selects a quadratic curve (Q)                   |  | Selects an optimised quadratic curve (P)       |
|  | Selects the brake-free spin-out option          |  | Confirms the input and closes the current view |
|  | Opens the previous view without saving          |  | Opens the main screen                          |

Fig. 35: Curve selection – example: deceleration

## 6.5.7 Program library



In the program library, up to 99 programs can be managed. Programs can be created, selected, edited or deleted. The list shows every program with its name and program-specific parameters (runtime, speed and temperature).


Programs can be marked as favourites so that they will be placed at the top of the program library.

The three most recently used programs can be selected via a quick selection feature.

The functions *Program import* and *Program export* (see chapter 6.5.11.1 - "Exporting and importing programs") can be used to exchange programs between centrifuges with a Spincontrol T control system.

### 6.5.7.1 Creating or editing a program

- ▶ The user interface shows the main screen.
  1. Select the *Program* button.
    - The "Programs" drop-down menu opens.
  2. Select the *Show all programs* button.
    - The user interface shows the program library.
  3. Select the *New* button (A1) to create a new program.
- Or:
  4. Select the *Next* button (A2) to edit an existing program.
    - The input menu of the program will be displayed.
  5. Select the *Edit* button on the status bar to edit the program name (B).
    - The input menu of the program name will be displayed.
  6. Enter a program name (C) and confirm the input (D).
    - The input menu of the program will be displayed.
  7. Enter the process parameters, process-related parameters and profiles (if applicable) (see chapter 6.5.7.2 - "Creating or editing a custom profile") for the program and select Enter (E).
    - The program will be saved or edited. The program library will be displayed.
  8. Select the *Confirm*, *Back* or *Home* button.
    - The main screen will be displayed.

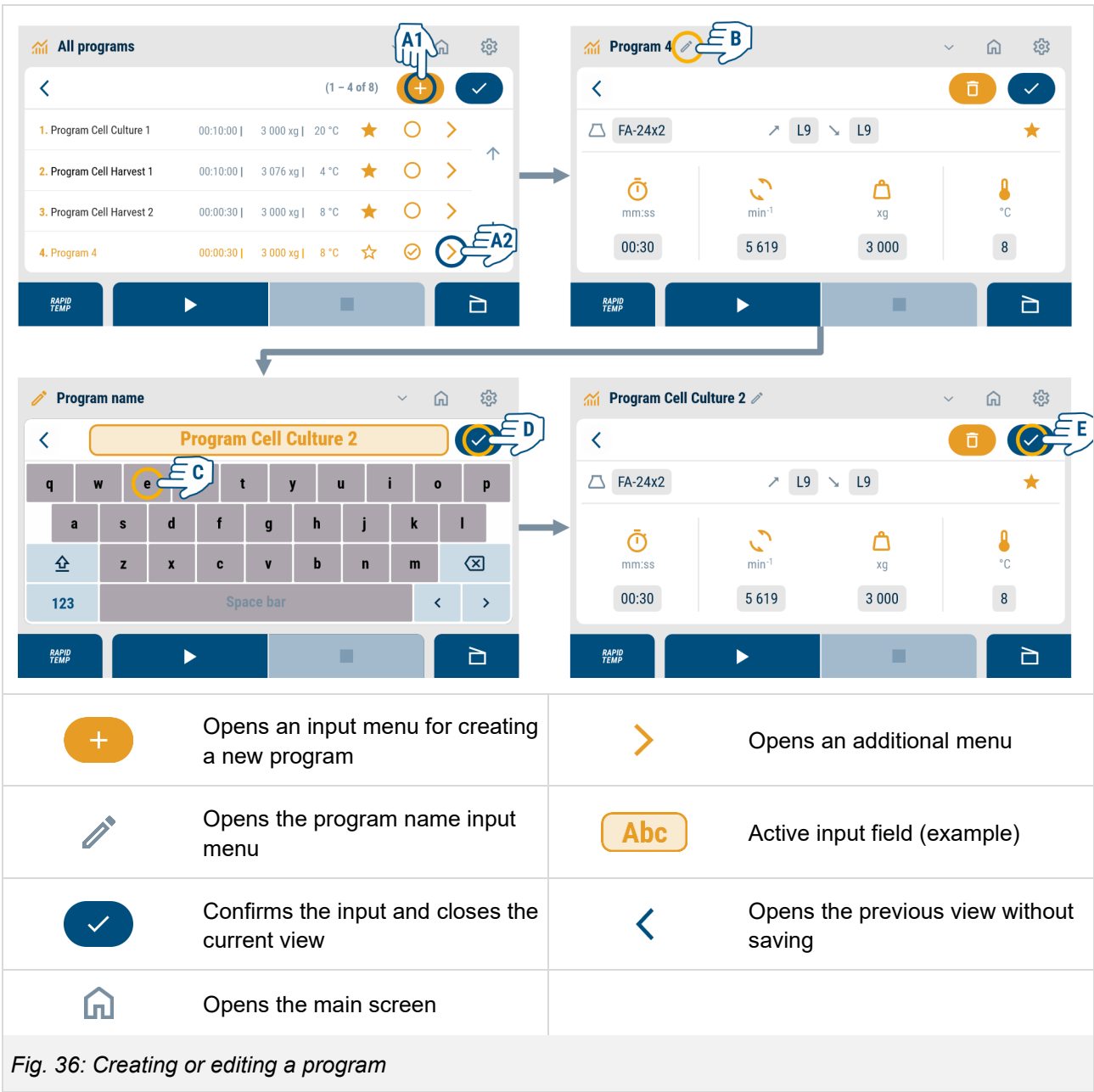

### 6.5.7.2 Creating or editing a custom profile



Profiles can only be created within a program.



Profiles (custom profiles) are freely programmable centrifugation runs composed of various intervals. A profile always starts with an acceleration interval (start interval) and ends with the deceleration interval (stop interval). Up to 19 additional intervals can be configured between the acceleration and deceleration intervals.

#### Creating a profile


- ▶ The input menu of the program is displayed.
- 1. Select the *Acceleration curve* or *Deceleration curve* button (A).
  - The "Curve selection" (Profiles) menu will be opened.
- 2. Open the "Curve selection" drop-down menu.
- 3. Select the "Profiles" function (Custom profiles) (B).
  - The "Profiles" menu opens.
- 4. Use the toggle feature to select the desired representation (RPM or RCF).
- 5. Define the intervals for the profile (see hereinbelow).

#### Editing a profile

- ▶ The input menu of the program is displayed.
- 1. Select the *Profiles* button.
  - The "Profiles" menu opens.
- 2. Edit the intervals of the profile (see hereinbelow).

#### Defining intervals

- ▶ A profile based on the current parameters is displayed.
- 1. Select the *New* button (C) to add a new interval. The new interval will be added after the selected intervals (continue with point 4).
  - Or:
- 2. Select the *Next* or *Previous* button to select an existing interval.
  - The selected interval will be highlighted.
- 3. Select the "Edit" button (D).
  - The input menu of the interval will be displayed.
- 4. Use the toggle feature to select the desired representation (RPM or RCF).
- 5. Enter the target speed or target RCF value that should be reached within a user-defined time.
  - Or:
- 6. Enter the desired speed or RCF value plus the time (E) during which this value should be reached. Then, confirm your input (F).
  - The "Profiles" menu will be displayed.
- 7. Configure the desired number of intervals in this way and then select Enter.
  - The profile will be stored within the selected program.
  - The main screen will be displayed.

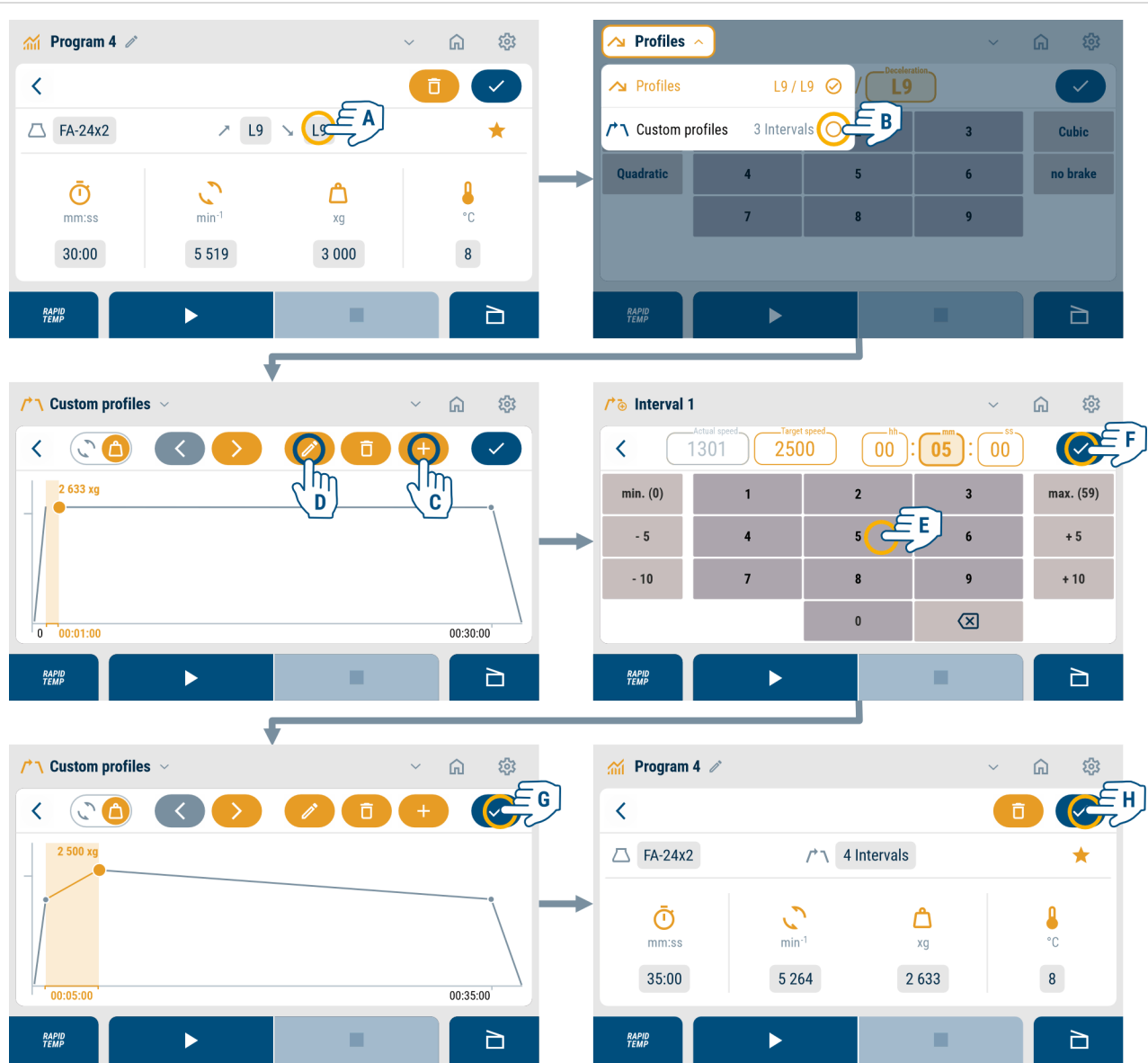


|  |                                                 |  |                                                            |
|--|-------------------------------------------------|--|------------------------------------------------------------|
|  | Deceleration curve button, opens the input menu |  | Selects a list item. Multiple selections are not possible. |
|  | Active input field (example)                    |  | Selects the desired representation (RPM or RCF)            |
|  | Selects the previous interval (inactive)        |  | Selects the next interval (active)                         |
|  | Opens the input menu of an existing interval    |  | Deletes the selected interval                              |
|  | Opens an input menu for creating a new interval |  | Confirms the input and closes the current view             |

Fig. 37: Creating or editing a profile

### 6.5.7.3 Executing a program

- The user interface shows the main screen.
- 1. Select the *Program* button (A).
- The "Programs" drop-down menu opens.
- 2. Select one of the three most recently used programs via the quick selection feature (B1).
- Or:  
Select the *Show all programs* button (B2).
- The program library opens.
- 3. Select a program (C) and confirm the input (D).
- The program will be loaded and the main screen will be displayed.
- 4. Select the *Start* button to execute the program.

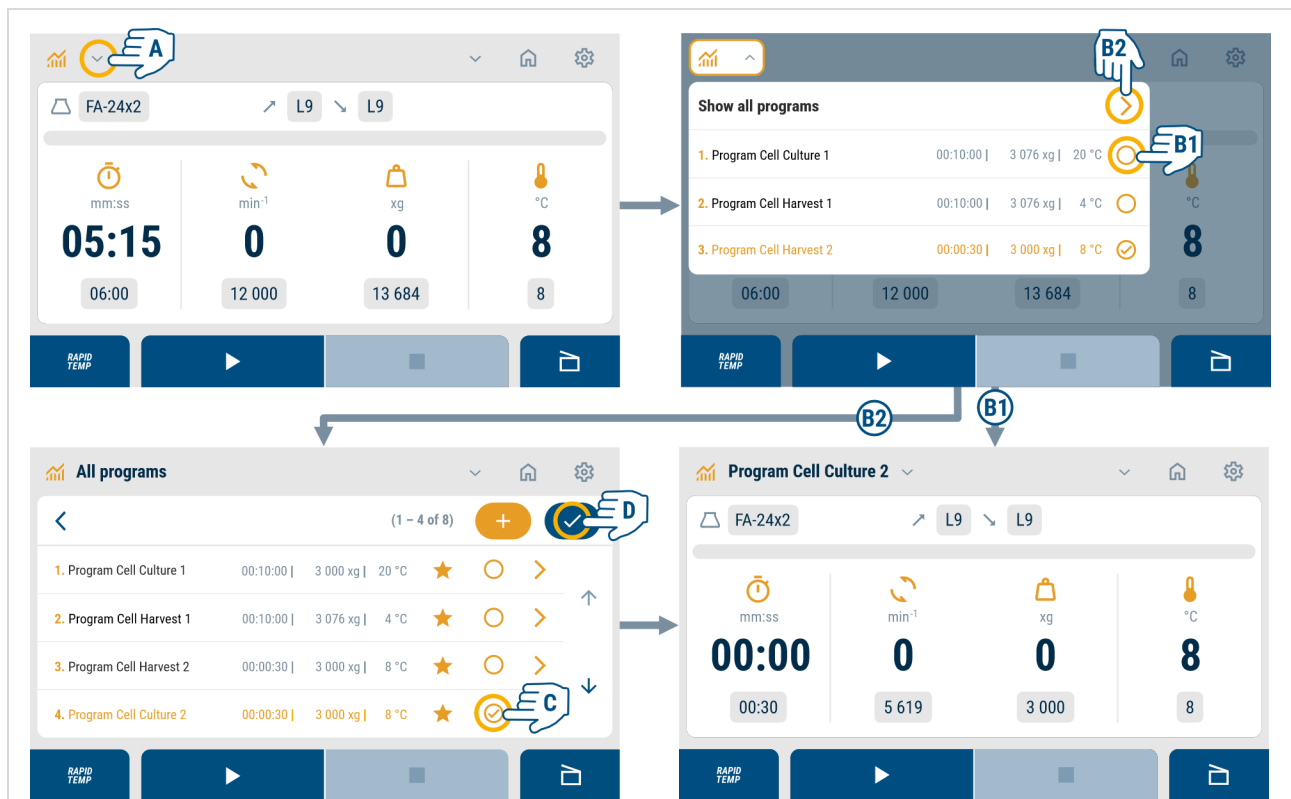


|                                                                                     |                                                             |                                                                                     |                                                                                 |
|-------------------------------------------------------------------------------------|-------------------------------------------------------------|-------------------------------------------------------------------------------------|---------------------------------------------------------------------------------|
|  | Opens a drop-down menu                                      |  | Opens an additional menu                                                        |
|  | Selects a list entry. Multiple selections are not possible. |  | Marks list entries as favourites. Favourites are placed at the top of the list. |
|  | Confirms the input and closes the current view              |  | Opens the previous view without saving                                          |
|  | Opens the main screen                                       |                                                                                     |                                                                                 |

**Fig. 38: Executing a program**

6.5.7.4 Deleting a program

- ▶ The user interface shows the main screen.
- 1. Select the *Program* button.
  - The "Programs" drop-down menu opens.
- 2. Select the *Show all programs* button.
  - The user interface shows the program library.
- 3. Select a program and then the *Next* button (A).
- 4. Select the *Delete* button (B).
  - The program will be deleted and the program library will be displayed.
- 5. Select the *Confirm*, *Back* or *Home* button to return to the main screen.

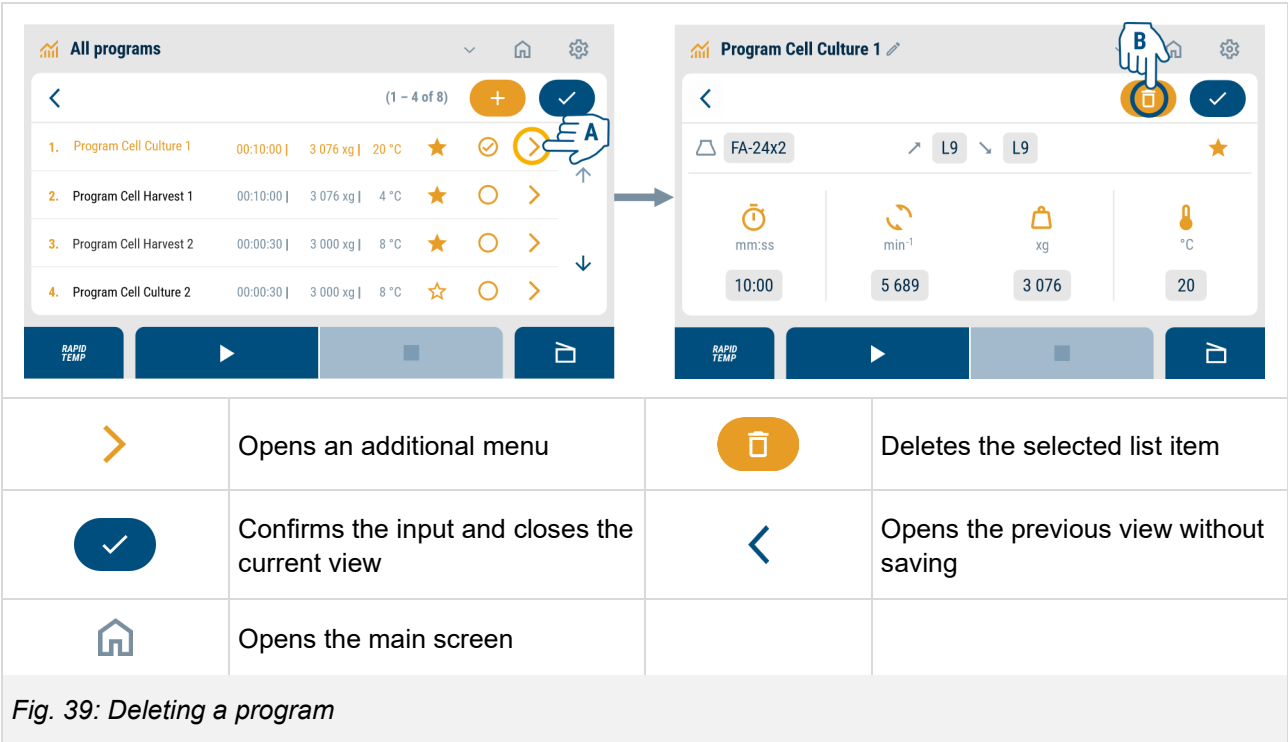


6.5.8 Status indicator and status menu

- The status indicator includes a variable number of symbols informing the user about the status of events and functions.
- The status menu provides additional information about the status.

6.5.8.1 Opening the status menu

- ▶ The user interface shows the main screen.
- 1. Select the *Status menu* button (A).
  - The status menu opens.
- 2. Select the desired function.

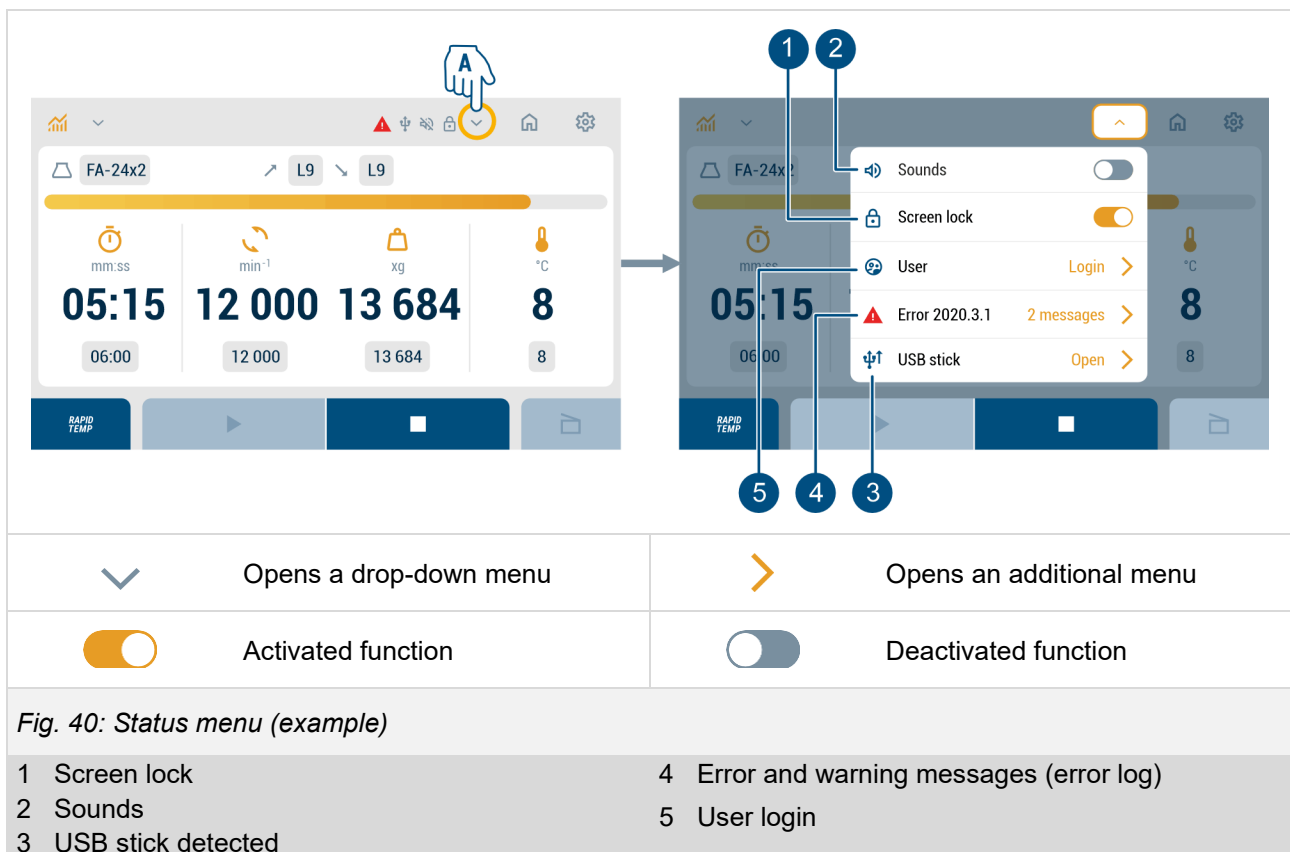

### Activating or deactivating the screen lock (1)

This element is used to immediately activate or deactivate the screen lock. For more information about the configuration of the screen lock, see chapter 6.5.10.1 - “Device settings”/ Display settings / Screen lock.

### Activating or deactivating sounds (2)

This element is used to activate or deactivate all of the sound signals in their entirety. For information about the configuration of individual sounds as well as their duration and volume, see chapter 6.5.10.1 - “Device settings”/ Signals.

### USB stick (3)

If a USB stick has been plugged into the port of the control panel and if it has been detected by the system, a corresponding line will be displayed. The *Next* button can be used to display the data stored on the USB stick. Programs can be imported or exported and QC data can be transferred (see chapter 6.5.9 - “Error log”).

### Error and warning messages (4)

The drop-down menu displays the most recent message as well as the total number of active messages. The *Next* button opens the error log (see chapter 6.5.11 - “Data transfer”).

### User login (5)

If the user login function is active, users can log in or out. For a detailed description of the function, see chapter 6.5.10.4 - "User management"/ "Account management by the user".

## 6.5.9 Error log




Any faults that occur are listed in the error log together with an error number, time stamp and description. The error log displays the last 100 entries.

### 6.5.9.1 Deleting error log entries

- The status menu is open.
  1. Select the *Next* button in the message line (A).
    - The error log opens.
  2. Select the *Delete* button (B).
    - The messages in the list will be permanently deleted from the error log without further notice.



The messages remain accessible to authorised specialist personnel in the error history of the service menu.

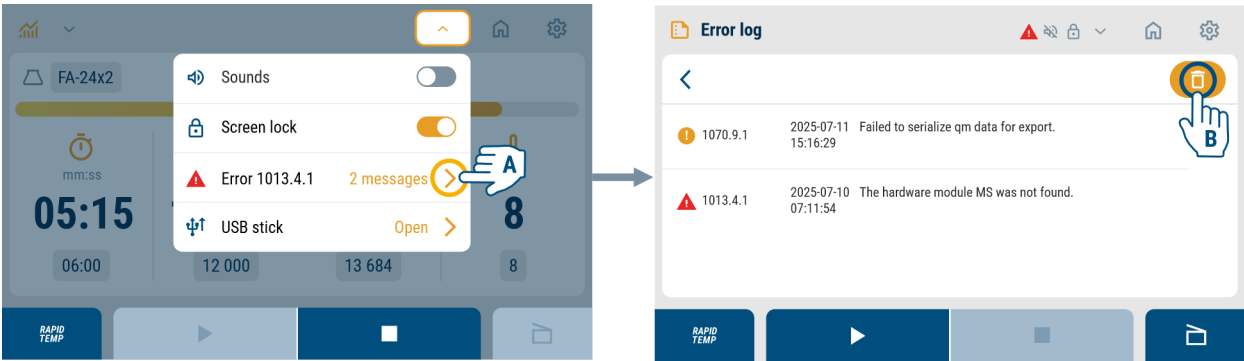

|                                                                                     |                                        |                                                                                     |                               |
|-------------------------------------------------------------------------------------|----------------------------------------|-------------------------------------------------------------------------------------|-------------------------------|
|  | Opens an additional menu               |  | Deletes all of the list items |
|  | Error message                          |  | Warning message               |
|  | Opens the previous view without saving |  | Opens the main screen         |


*Fig. 41: Deleting error log entries (example)*

## 6.5.10 Settings

 The "Settings" menu can be used to configure device settings and to view information about the centrifuge.

Opening the "Settings" menu:


1. Select the *Settings* button (A).
- The "Settings" menu opens.

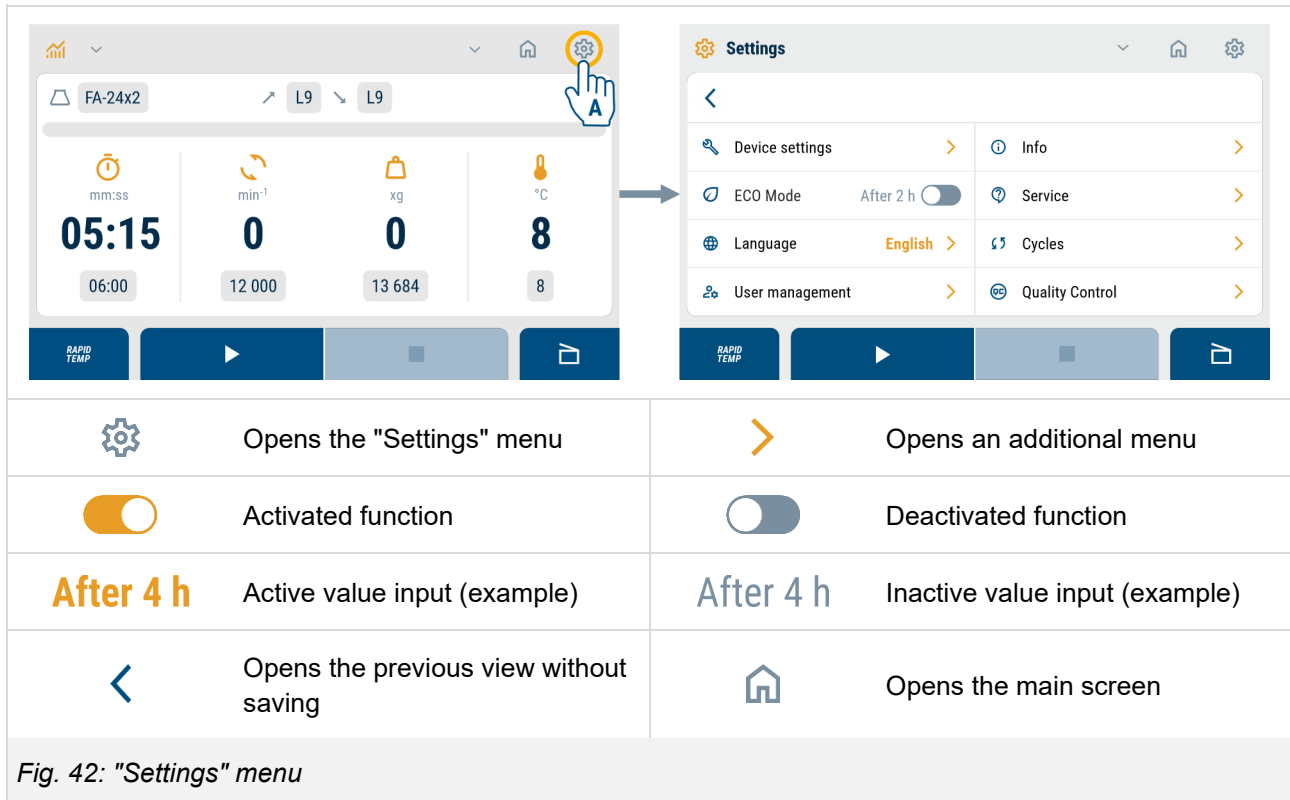

### 6.5.10.1 Device settings



Device settings are universal. They cannot be adapted in a program-specific manner.

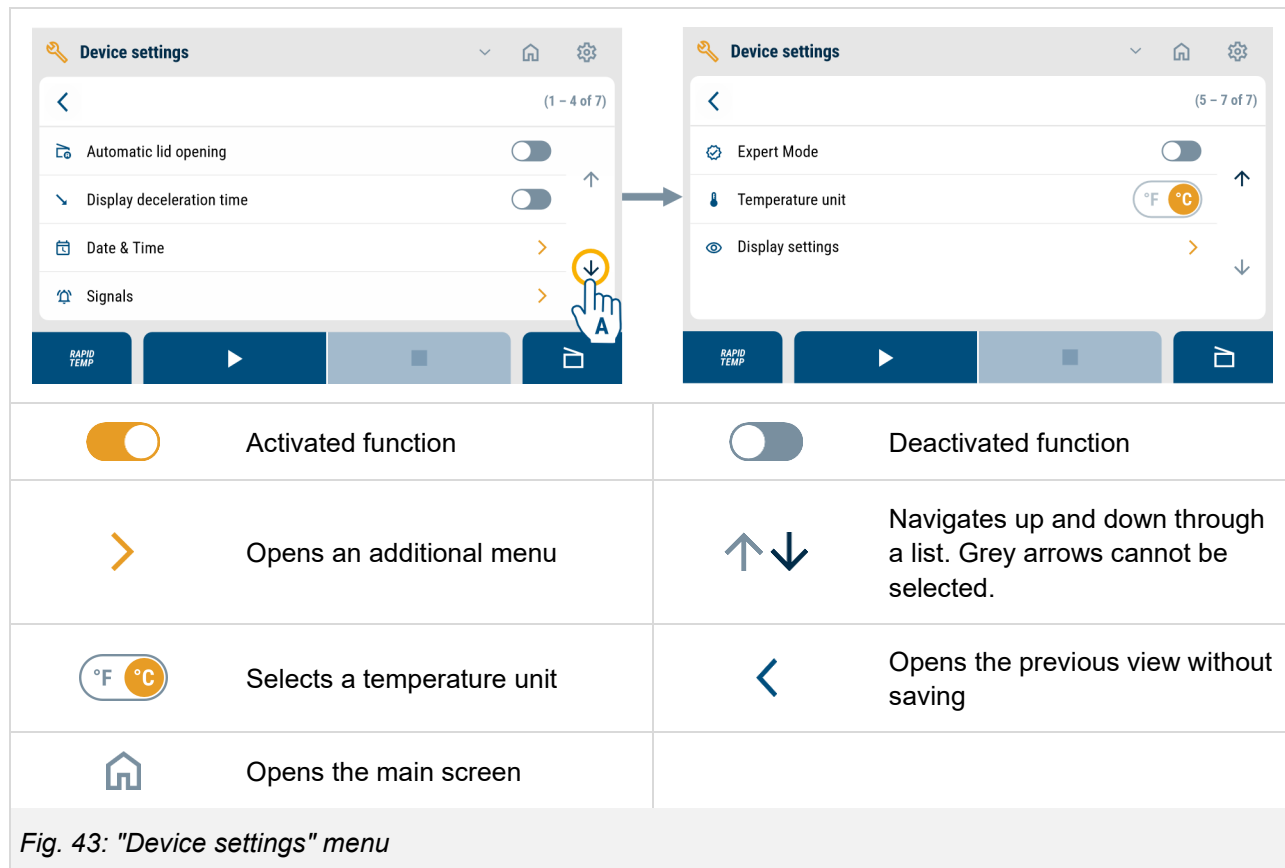

#### Automatic lid opening function

This function activates or deactivates the automatic opening of the centrifuge lid after the completion of a centrifugation run.

In the factory settings, this function is deactivated.


► The user interface shows the "Device settings" menu.

1. Activate the *Automatic lid opening* function (A).
2. Select the *Home* button (B).

– The main screen will be displayed.

– The *Lid* button on the control bar is marked with an A.

After the centrifugation run, the lid lock will be released and the lid opens automatically.

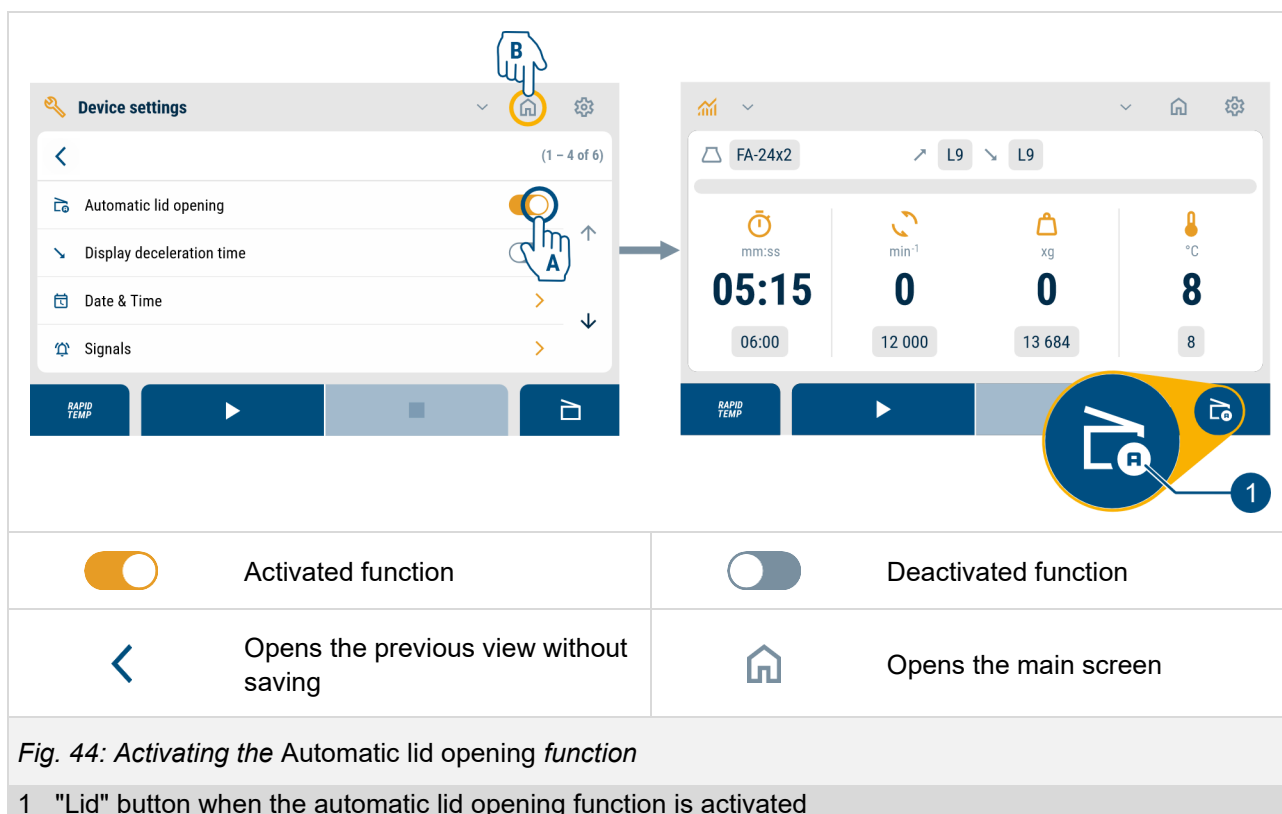

**CAUTION** – When the lid is open, the cooling system is inactive. The samples may heat up.

In the following cases, the lid will remain closed even though the automatic lid opening function has been activated:

- The standstill cooling function has been activated.
- A RAPID-TEMP run has been completed.

This prevents the samples from heating up.

### **Display deceleration time**

This function activates or deactivates the indication of the deceleration time during and after the deceleration. The deceleration time will be displayed in the "Runtime" area.


In the factory settings, this function is deactivated.

► The user interface shows the "Device settings" menu.

1. Activate the *Display deceleration time* function (A).
2. Select the *Back* or *Home* button (B).
  - The main screen will be displayed.
  - During the deceleration, the symbol for the deceleration time indication will be displayed in the "Runtime" area.
  - The deceleration time will be displayed in grey and in italics.

Once the deceleration is complete, the system will once again display the runtime.

## Operation

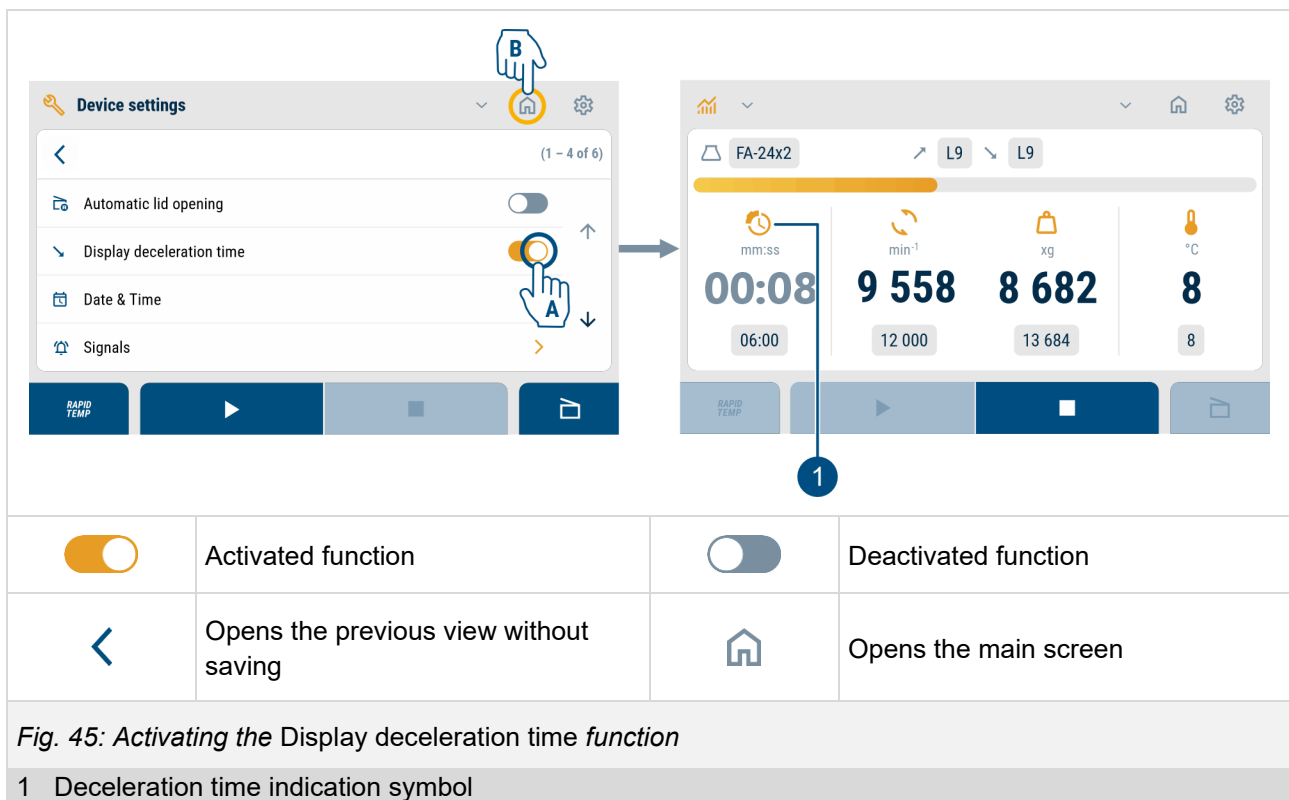

### Date and time

This function can be used to set the date and time. In addition, the time zone of the location of use of the centrifuge can be specified.


In the factory settings, the time zone for the centrifuge is Central European Time (CET).



If the centrifuge remains switched off for a longer period of time (> 3 months), its date and time must be checked and reset when the centrifuge is switched on again.

Process for setting the time (example):

- The user interface shows the "Device settings" menu.
  1. Select the *Next* button to open the "Date & time" menu.
  2. Select the *Next* button to open the "Time" input menu (A).
  3. Enter the time (B) and confirm the input (C).
  4. Select the *Home* button.
    - The main screen will be displayed.

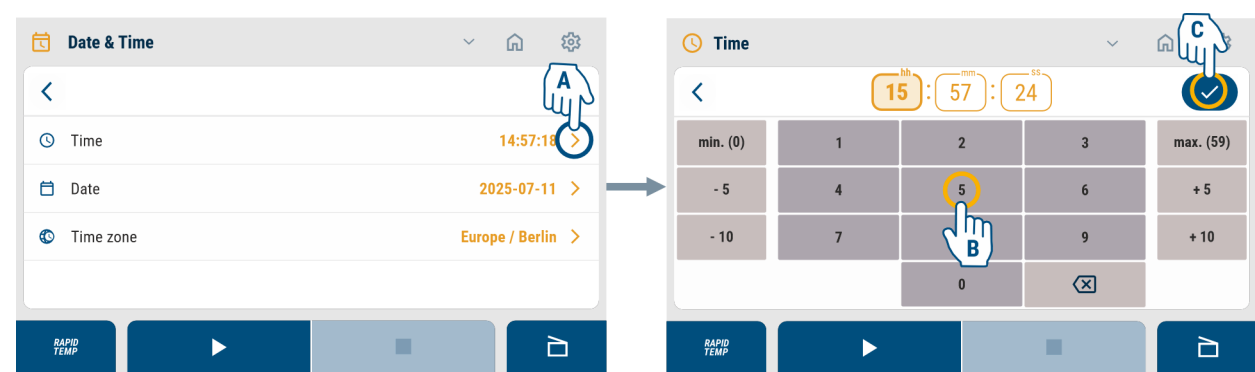

|                                                                                     |                                                |                                                                                    |                                                     |
|-------------------------------------------------------------------------------------|------------------------------------------------|------------------------------------------------------------------------------------|-----------------------------------------------------|
|    | Opens an additional menu                       |   | Active input field (example)                        |
|    | Inactive input field (example)                 |   | Minimum (min.) and maximum (max.) permissible value |
|    | Adds 5 or 10 units to the selected value       |   | Subtracts 5 or 10 units from the selected value     |
|   | Confirms the input and closes the current view |  | Opens the previous view without saving              |
|  | Opens the main screen                          |                                                                                    |                                                     |

**Fig. 46: Setting the time (example)**

## Signals

This function is used to activate or deactivate a sound signal at the end of a centrifugation run or in the event of imbalance or errors.

A signal duration between 1 and 360 seconds and a volume level between 1 and 5 can be set.

-  The function *Sound off* in the status menu must be deactivated so that the signal sounds are audible. The corresponding symbol must be displayed.

Entering the signal duration (example):

- The user interface shows the "Device settings" menu.
  1. Deactivate the *Sound off* function if it is still active.
  2. Select the *Next* button (A).
    - The "Signals" menu opens.
  3. Activate the desired signals (B).
  4. Scroll down through the list (C) and open the "Signal duration" menu (D).
  5. Enter the signal duration and confirm the input (E).
  6. Select the *Home* button.
    - The main screen will be displayed.



|                                                                                     |                                                                         |                                                                                      |                                                |
|-------------------------------------------------------------------------------------|-------------------------------------------------------------------------|--------------------------------------------------------------------------------------|------------------------------------------------|
|  | Sound off function activated                                            |   | Opens an additional menu                       |
|  | Activated function                                                      |   | Deactivated function                           |
|  | Navigates up and down through the list. Grey arrows cannot be selected. |   | Active input field (example)                   |
|  | Minimum (min.) and maximum (max.) permissible value                     |  | Adds 10 or 50 units to the selected value      |
|  | Subtracts 10 or 50 units from the selected value                        |   | Confirms the input and closes the current view |
|  | Opens the previous view without saving                                  |   | Opens the main screen                          |

Fig. 47: Entering the signal duration


### Expert mode

The function activates or deactivates the extended process-related parameters (see chapter 6.5.4.6 - "Expert Mode: Process-related parameters"):

- Integral
- Radius reduction
- Density increase

### Temperature unit

This function can be used to switch between °C (degree Celsius) and °F (degree Fahrenheit) for the temperature unit by way of a toggle button.

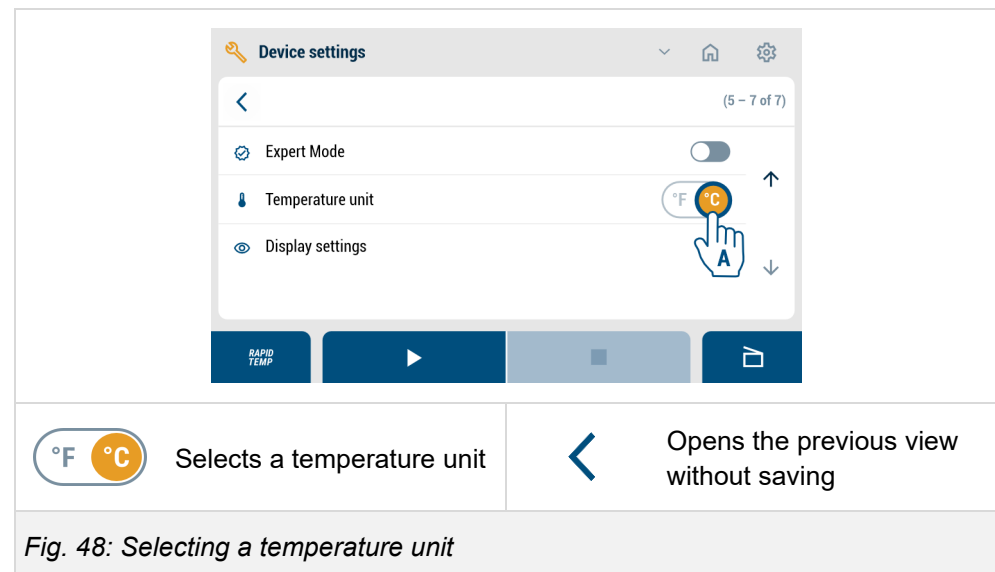


### Display settings


#### Brightness

Configuration of display settings (here: brightness):

- The user interface shows the "Device settings" menu.
  1. Select the *Next* button to open the "Display settings" menu.
  2. Select the *Next* button (A) to open the "Brightness" menu.
    - The input menu opens.
  3. Select a value between level 1 (dark) and level 10 (bright) and confirm your selection (B).
  4. Select the *Home* button (C).
    - The main screen will be displayed.

## Operation

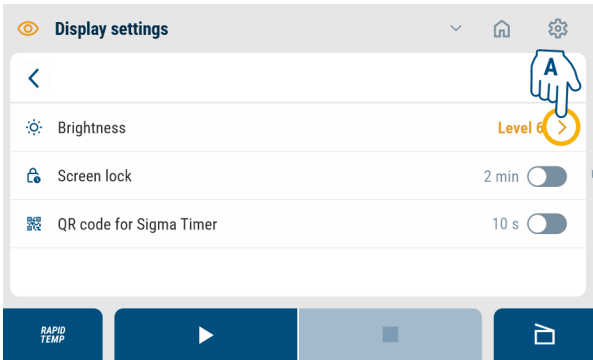

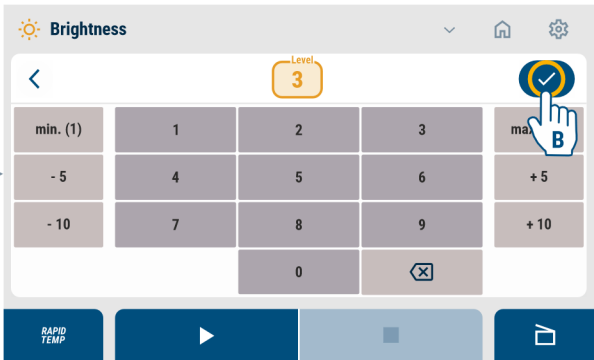


|                                                                                    |                                                     |                                                                                    |                                                |
|------------------------------------------------------------------------------------|-----------------------------------------------------|------------------------------------------------------------------------------------|------------------------------------------------|
|   | Opens an additional menu                            |   | Active input field (example)                   |
|   | Minimum (min.) and maximum (max.) permissible value |   | Adds 5 or 10 units to the selected value       |
|   | Subtracts 5 or 10 units from the selected value     |   | Confirms the input and closes the current view |
|  | Opens the previous view without saving              |  | Opens the main screen                          |

*Fig. 49: "Display settings" menu for adjusting the brightness*

### Screen lock

The automatic screen lock can be activated for a duration between 1 and 60 minutes. If no input is made during the set period:

- the user interface will be locked automatically,
- the system switches to the main screen,
- all of the submenus will be closed, any unconfirmed inputs will be rejected and
- the corresponding symbol will be displayed on the status bar (see the illustration).

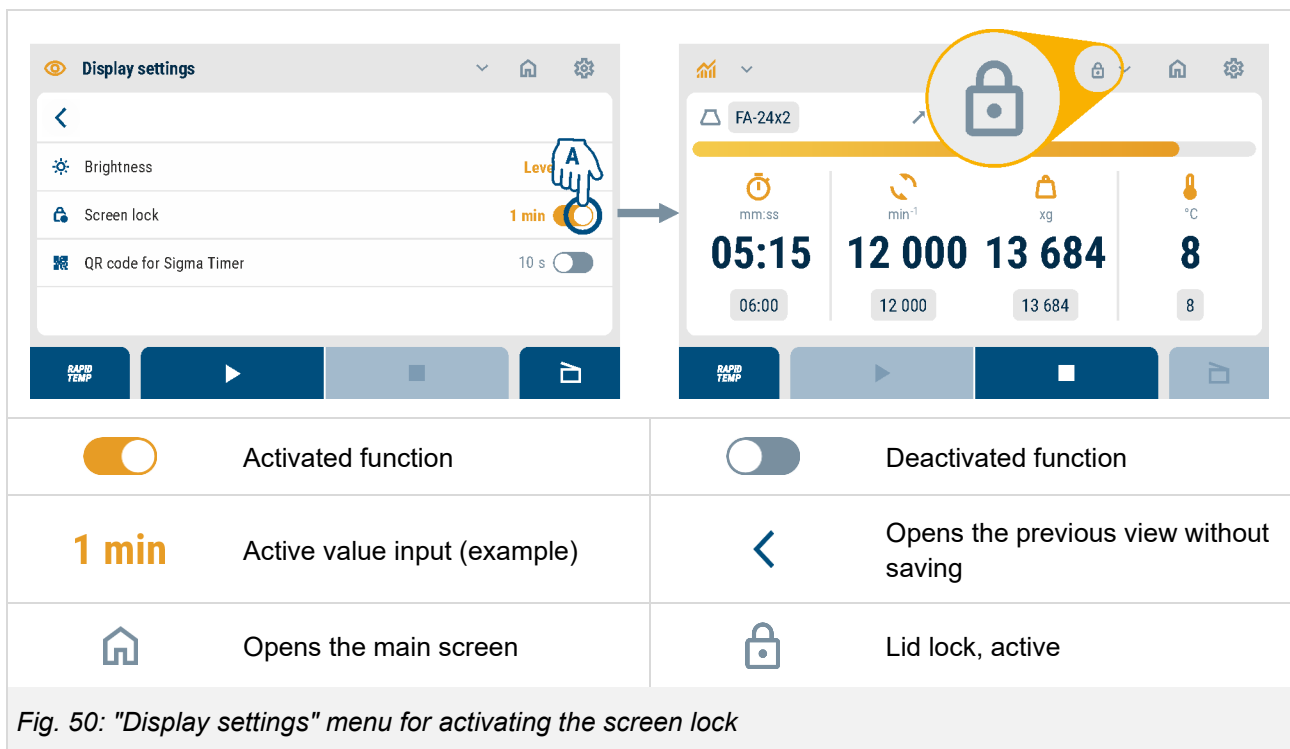

When the screen is locked, the following transactions are still possible:

- Touching of buttons. In this case, a prompt will be displayed asking whether the screen should be unlocked.
- Deactivation of the screen lock in the status menu. The screen will be unlocked, but the general setting in the "Display settings" menu will remain unchanged.

#### QR code for Sigma Timer

When this function is active, the centrifuge generates a QR code for the Sigma Timer app. With this app, the remaining runtime of the centrifuge can be monitored by way of a mobile device.

The following happens when the function is active:


- A QR code will be generated at the start of a centrifugation run. This QR code can then be scanned with the Sigma Timer app.
- The *QR code* button will be displayed on the configuration bar in the main screen during an active centrifugation run. When this button is selected, the QR code will be displayed again for a preset time.

A display time between 10 and 60 seconds can be configured. In the factory settings, this value is set to 10 seconds.

Generating a QR code:

- The user interface shows the main screen.
- 1. Select the *Generate QR code* button (A).
  - The "QR code" dialog box opens.

## Operation

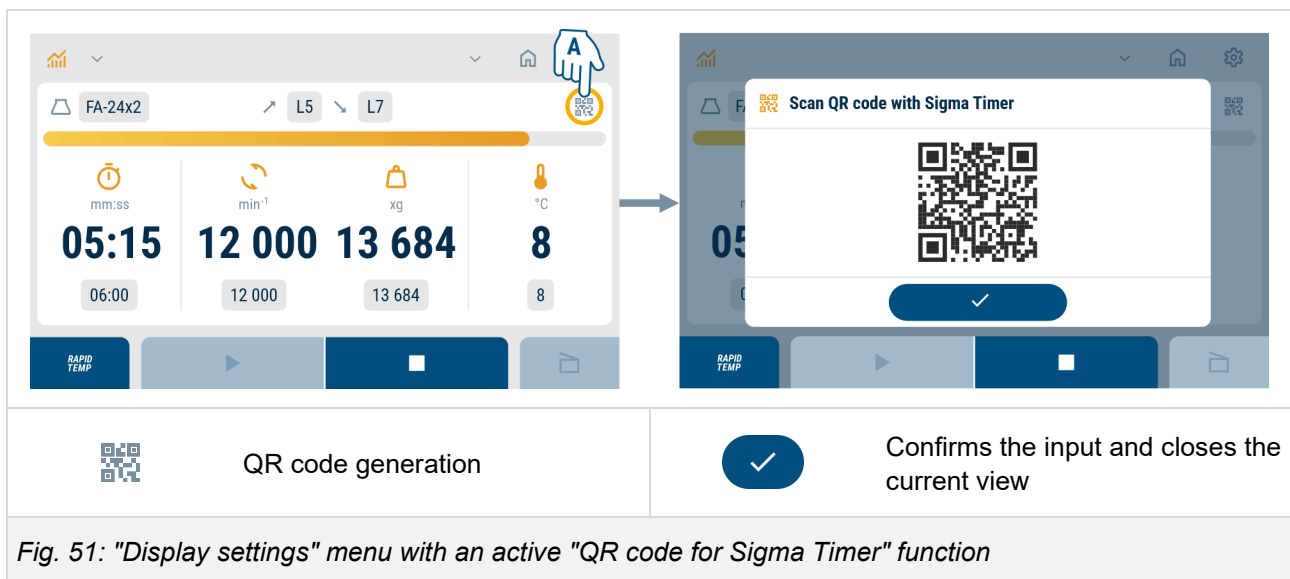

### 6.5.10.2 ECO mode



If the centrifuge is not used for a defined period of time, the *ECO mode* function will switch off the refrigeration unit. In addition, the brightness of the user interface will be set to level 1 (dark) and the centrifuge lid will be opened so that any moisture can escape.

The setting for this function is universal, i.e. it cannot be activated or deactivated for individual programs.

Activating the ECO mode:

1. Activate the *ECO mode* function.
2. Select the *Value input* button (A).
  - The input menu opens.
3. Enter a set value between 1 and 8 hours in steps of 60 minutes.
4. Confirm the input (B) and select the *Home* button (C).
  - The main screen will be displayed.



|                                                                                     |                                                     |                                                                                      |                                                |
|-------------------------------------------------------------------------------------|-----------------------------------------------------|--------------------------------------------------------------------------------------|------------------------------------------------|
|  | Activated function                                  |   | Deactivated function                           |
| <b>After 4 h</b>                                                                    | Active value input                                  |   | Active input field (example)                   |
|  | Minimum (min.) and maximum (max.) permissible value |  | Adds 1 or 4 units to the selected value        |
|  | Subtracts 1 or 4 units from the selected value      |   | Confirms the input and closes the current view |
|  | Opens the previous view without saving              |   | Opens the main screen                          |

**Fig. 52: ECO mode function**

- 1 The ECO mode is activated.
- 2 The ECO mode will start in 30 minutes.

The following will happen when the ECO mode is activated:

- The turn-off countdown will start in the background after the last keystroke while the centrifuge is at a standstill.
- The ECO mode symbol will be displayed in the temperature area of the main screen. In the last 30 minutes before the end of the countdown, the symbol will flash green.

Touching the user interface will reset the countdown to its original value so that it will restart.

In the factory settings, the countdown of the ECO mode is set to 4 hours.


### 6.5.10.3 Language selection

 This menu can be used to select the language of the user interface.

Selecting another language:

► The user interface shows the "Settings" menu.

1. Select the *Next* button (A).
  - The selection menu opens.
2. Select the desired language (B) and confirm the input (C).
3. Select the *Home* button (D).
  - The main screen will be displayed.

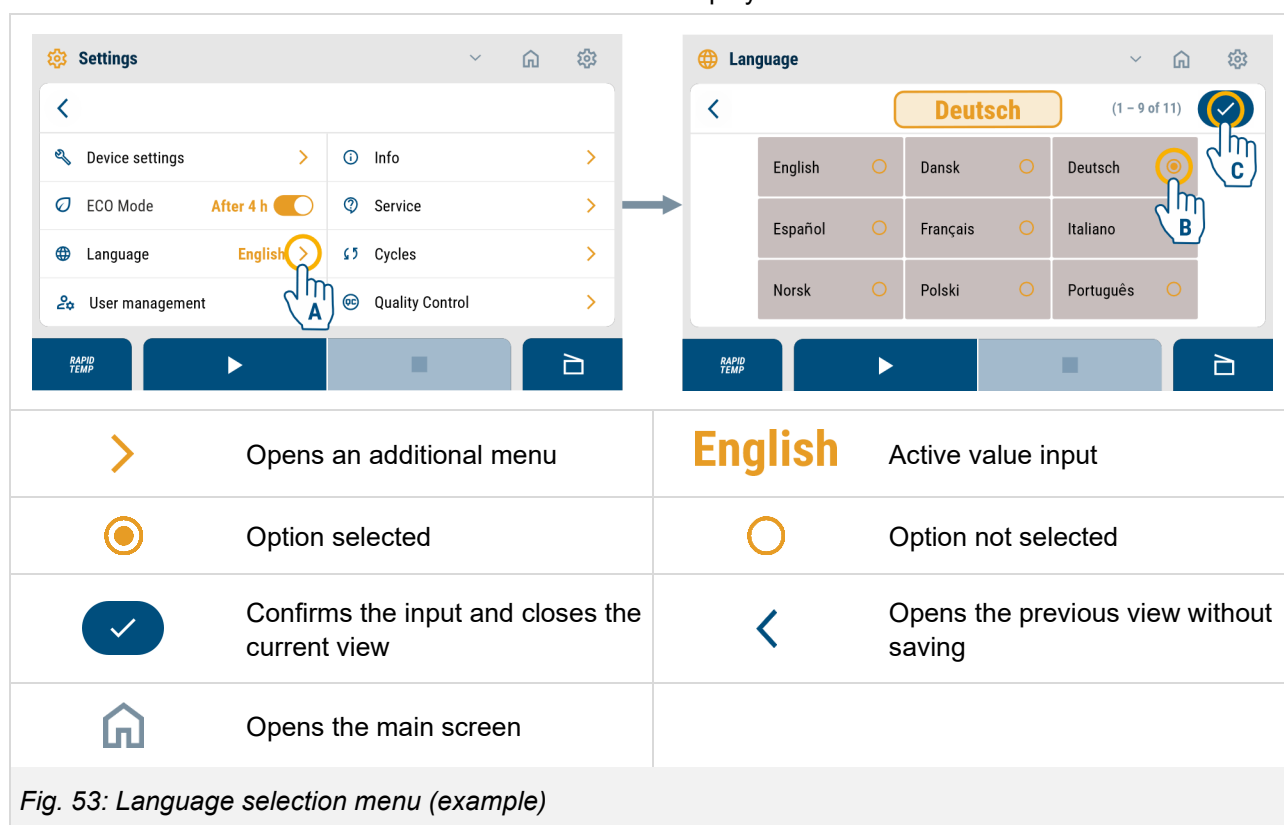

### 6.5.10.4 User management

 The user management function enables an administrator to globally assign permissions to operate the centrifuge. Optionally, user login with or without a password can be configured to enable the operation of the centrifuge. Up to 99 users can be set up. User passwords are initially assigned by the administrator. However, the users can change their password afterwards.

The administrator can assign the following rights:

- Loading and starting of already existing programs
- Creation, storage and execution of programs (the loading and starting of programs will be automatically activated as well)
- Editing of process parameters
- Editing of global device settings, e.g. ECO mode
- Export and import of data via a USB stick, either individually or combined

If the rights management function is active, but no rights are enabled and the user login is deactivated (default setting after rights management activation), only a centrifugation run with the current parameters can be started.

### **Rights management**

The rights management function comprises several options for restricting access to the centrifuge:

- Unrestricted accessed – the rights management function is deactivated and the centrifuge can be operated freely without any restrictions. An administrator is not necessary.
- Global assignment of rights – the rights management function is active and the functions of the centrifuge are enabled or disabled for all users.
- User login with a defined scope of functions – the rights management function is active and an administrator has created a user account for every user.

#### Activating the rights management function


- The user interface shows the "Settings" menu.
  1. Select the *Next* button to open the "User management" menu.
  2. Select the *Next* button to open the "Rights management" menu.
  3. Activate the rights management function (A). If the rights management function is being created for the first time, an administrator account must be created first.

#### Creating an administrator account

The administrator has unrestricted access to all the functions. An administrator account can be created only once. Deleting or renaming the account is not possible.

- The user interface shows the "Rights management" menu.
  1. Activate the rights management function (A).
    - A "Create admin" dialog box will be opened.
  2. Select the *Create admin* button (B).
    - The "New password" menu for entering an administrator password will be opened.
  3. Enter an administrator password (C), repeat it (D) and confirm the input.
    - The password will be saved and the administrator account will be created.
    - The "User management" menu will be displayed.


## Operation

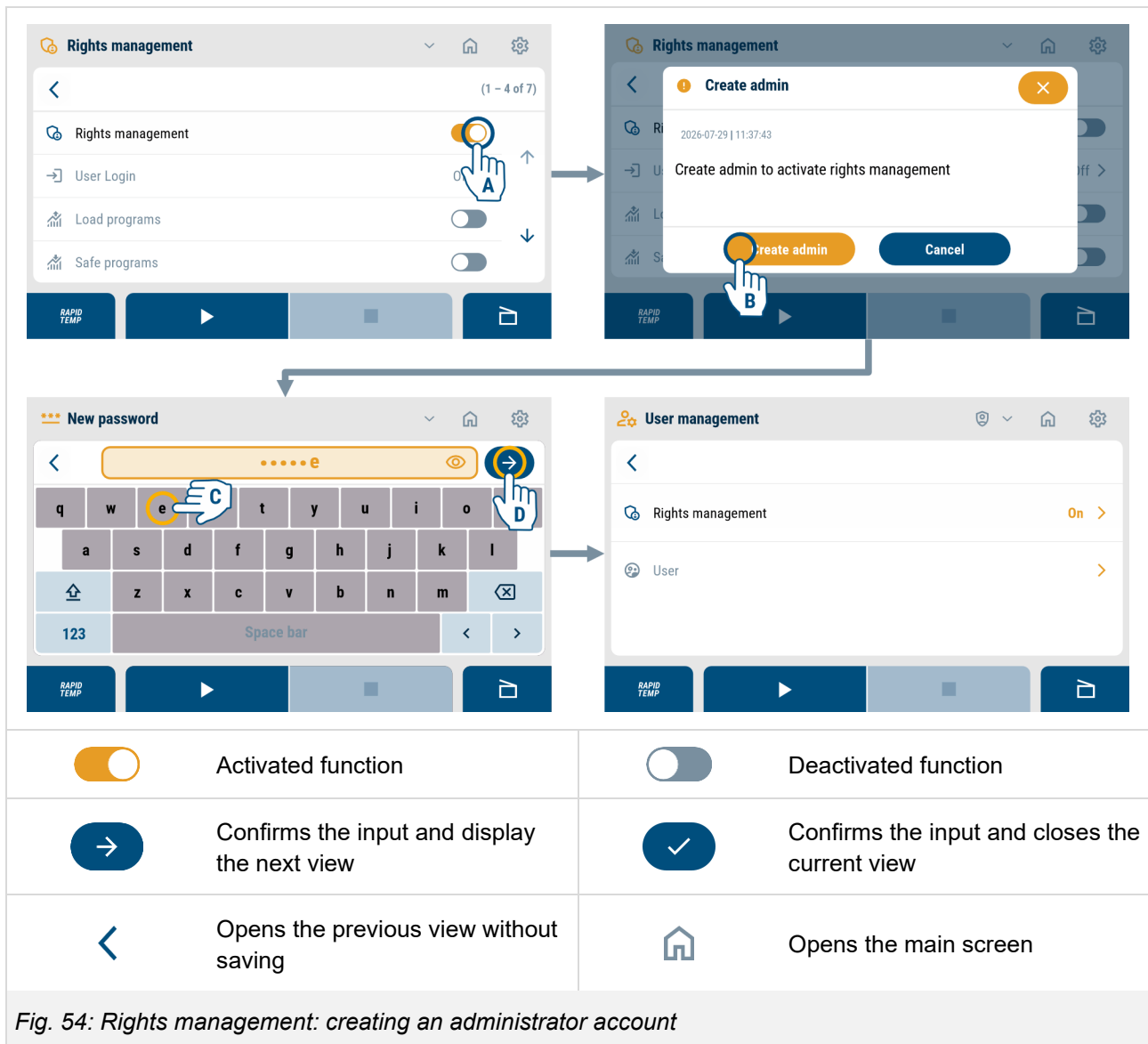

### User login

The administrator uses the "User login" menu to define whether user login is required, whether a password is required for logging in and/or whether users are automatically logged out after a specified period of time.

- ▶ The user interface shows the "Rights management" menu.
  - ▶ The administrator is logged in.
  - ▶ The rights management function is active.
1. Select the *Next* button (A) to open the "User login" menu.
  2. Activate the functions "User login" (B), "Login with password" and/or "Automatic logout".
  3. If the "Automatic logout" function is active, select the *Value input* button to change the logout time (C).
    - The input menu opens.
  4. Enter a value between 1 minute and 99 hours and confirm your input.
    - The settings will be adopted and the "User login" menu will be displayed.

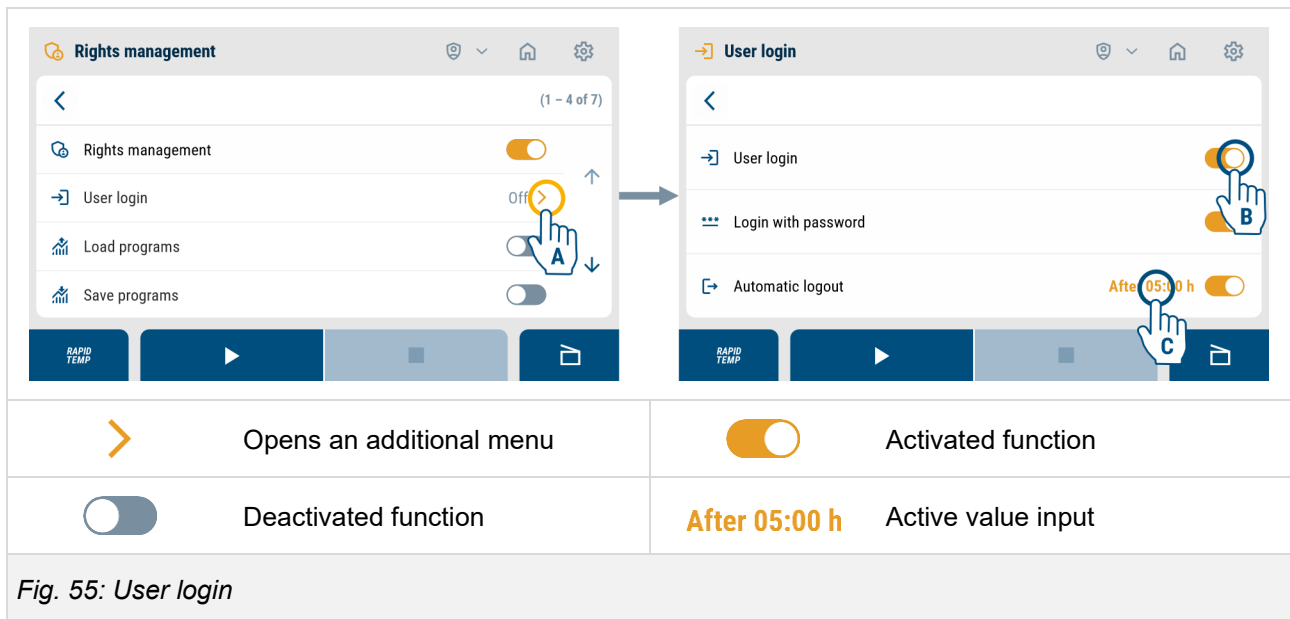

#### Assigning global rights

The administrator can assign rights for the following functions. Functions that have been activated by the administrator can be used by all users.

- ▶ The user interface shows the "Rights management" menu.
- ▶ The administrator is logged in.
- ▶ The rights management function is active.

##### 1. Activating or deactivating the desired functions:

- Load programs
- Save programs: When this function is activated, the "Load programs" functions will automatically be activated as well.
- Change process parameters
- Change settings
- USB access



If all global rights are deactivated and no user has been created, the centrifuge can only be operated by the administrator.

#### Deactivating the rights management function

- ▶ The user interface shows the "Rights management" menu.
- ▶ The rights management function is active.

##### 1. Deactivate the rights management function.

- The centrifuge can be operated without any restrictions.
- The user accounts are deactivated. However, they will be available once again when the rights management function is reactivated.

## Operation

---

### User

The "User" menu enables the administrator to create, edit and delete individual user accounts.

#### Creating a user account

- ▶ The user interface shows the "User management" menu.
- ▶ The administrator is logged in.
- ▶ The rights management function is active.
- ▶ The "User login" function is active.
- 1. Select the *Next* button (A) to open the "User" menu.
- 2. Select the *New* button (B) to create a new user account.
  - The input menu for the user name opens.
- 3. Enter a user name (B) and select the *Next* button (C) to confirm the input and to open the next menu.
  - The input menu for the password opens.
- 4. Enter a password and confirm your input.
  - The menu for the created user will be displayed.

5. Select the *Previous* button.

- The "User management" menu will be displayed.
- The creation of the user account is complete.

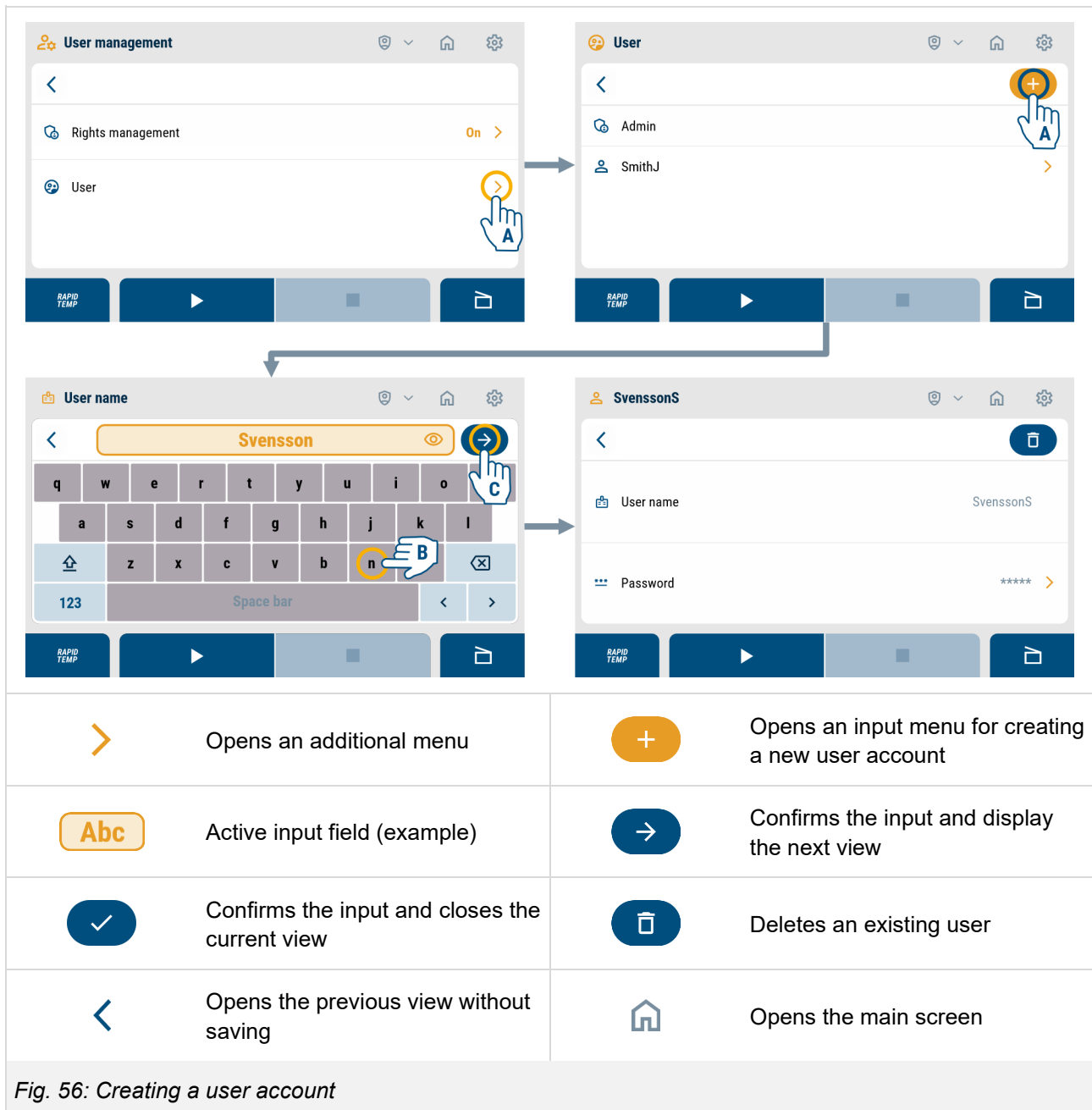


Fig. 56: Creating a user account

#### Changing a password (by the administrator)

- The user interface shows the "User management" menu.
  - The administrator is logged in.
  - The rights management function is active.
1. Select the *Next* button to open the "User" menu.
  2. Select the *Next* button to open an existing user account.
    - The menu for the created user will be displayed.
  3. Select the *Next* button to access the "Change password" function.
    - The input menu for a new password opens.
  4. Enter the new password and confirm your input.

## Operation

---

5. Select the *Previous* button.
  - The "User management" menu will be displayed.
  - The password has been changed.

### Deleting a user account


- ▶ The user interface shows the "User management" menu.
  - ▶ The administrator is logged in.
  - ▶ The rights management function is active.
1. Select the *Next* button to open the "User" menu.
  2. Select the *Next* button to open an existing user account.
    - The input menu for the password opens.
  3. Enter the current password and confirm your input.
    - The menu for the created user will be displayed.
  4. Select the *Delete* button.
    - A dialog box for removing a user will be opened.
  5. Confirm.
  6. Select the *Previous* button.
    - The "User management" menu will be displayed.
    - The user account has been deleted.

### **Account management by the user**

If the user login with a password has been activated by the administrator, every user must enter their user name and password to log in and be able to operate the centrifuge.

### User login

- ▶ The centrifuge is at a standstill.
1. Select the *Status menu* button.
    - The status menu opens.
  2. Select the *Next* button (A) to log in.
  3. Select the *Next* button (B) to select the user account.
    - The input menu for the password opens.
  4. Enter the password (C) and confirm (D).
    - The user is logged in.
    - The main screen will be displayed.
    - The user symbol will be displayed in the status line.

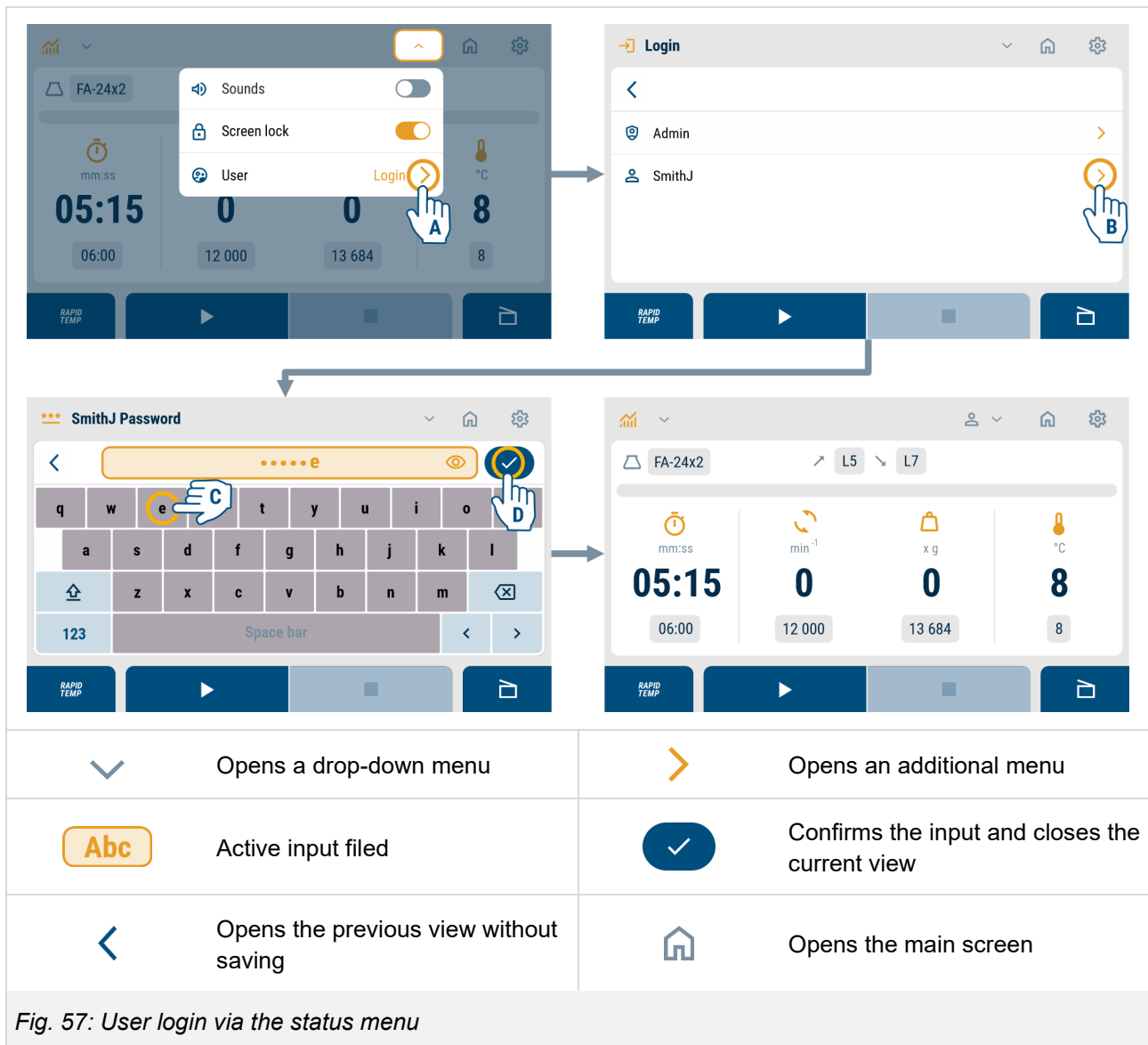

#### Changing a password (by the user)


- ▶ The centrifuge is at a standstill.
- ▶ The user is logged in.
- 1. Select the *Status menu* button.
  - The status menu opens.
- 2. Select the *Next* button to open the "User" menu.
- 3. Select the *Next* button to access the "Change password" function.
- 4. Enter the current password and confirm your input.
- 5. Enter new password, repeat it and confirm your input.
  - The "User" menu will be displayed.
- 6. Select the *Previous* button.
  - The main screen will be displayed.
  - The password has been changed.

#### User logout

- ▶ The centrifuge is at a standstill.
- ▶ The user is logged in.

## Operation

1. Select the *Status menu* button.
  - The status menu opens.
2. Select the *Next* button (A) to open the "User" menu.
3. Select the *Logout* button (B).
  - The main screen will be displayed.
  - The user is logged out.

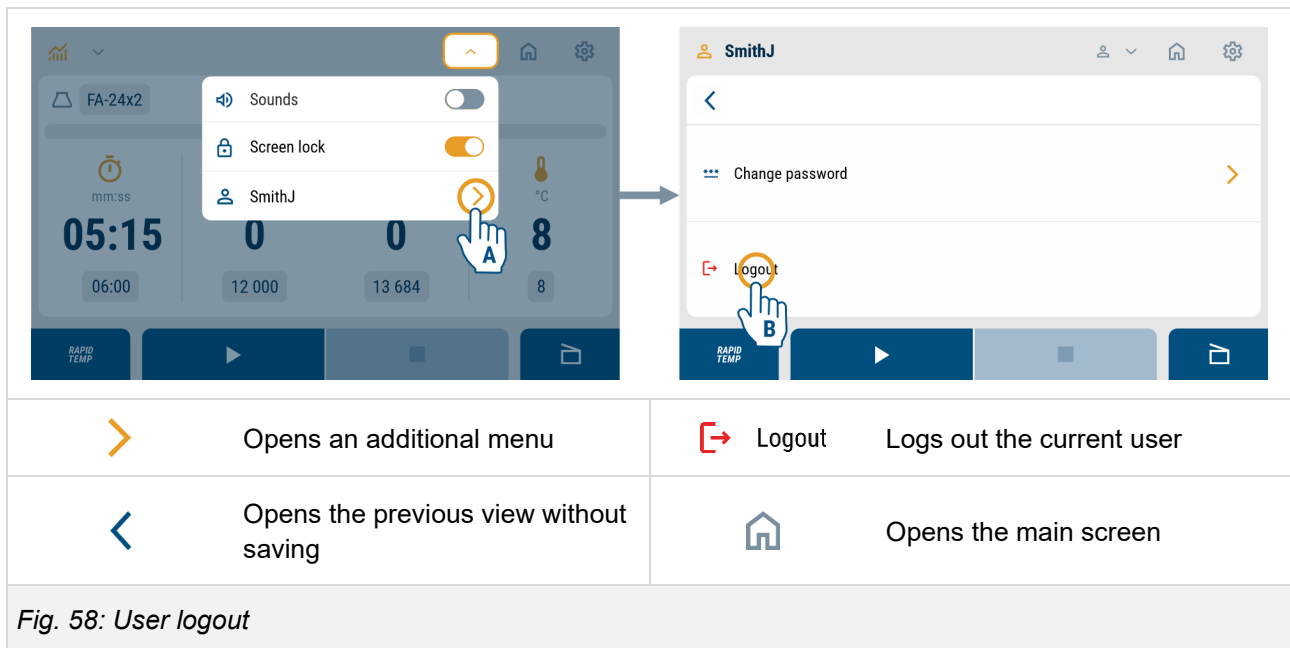


### 6.5.10.5 Info

-  This menu shows all of the system and software information, including the serial number.

#### Centrifuge name

This menu can be used to specify a name for the centrifuge to enable its identification.

- The user interface shows the "Settings" menu.
1. Select the *Next* button to open the "Info" menu.
  2. Select the *Next* button to open the "Centrifuge name" menu (A).
    - The input menu opens.
  3. Enter the centrifuge name. 20 alphanumeric characters are available.
  4. Confirm the input (B).
  5. Select the *Home* button (C).
    - The main screen will be displayed.

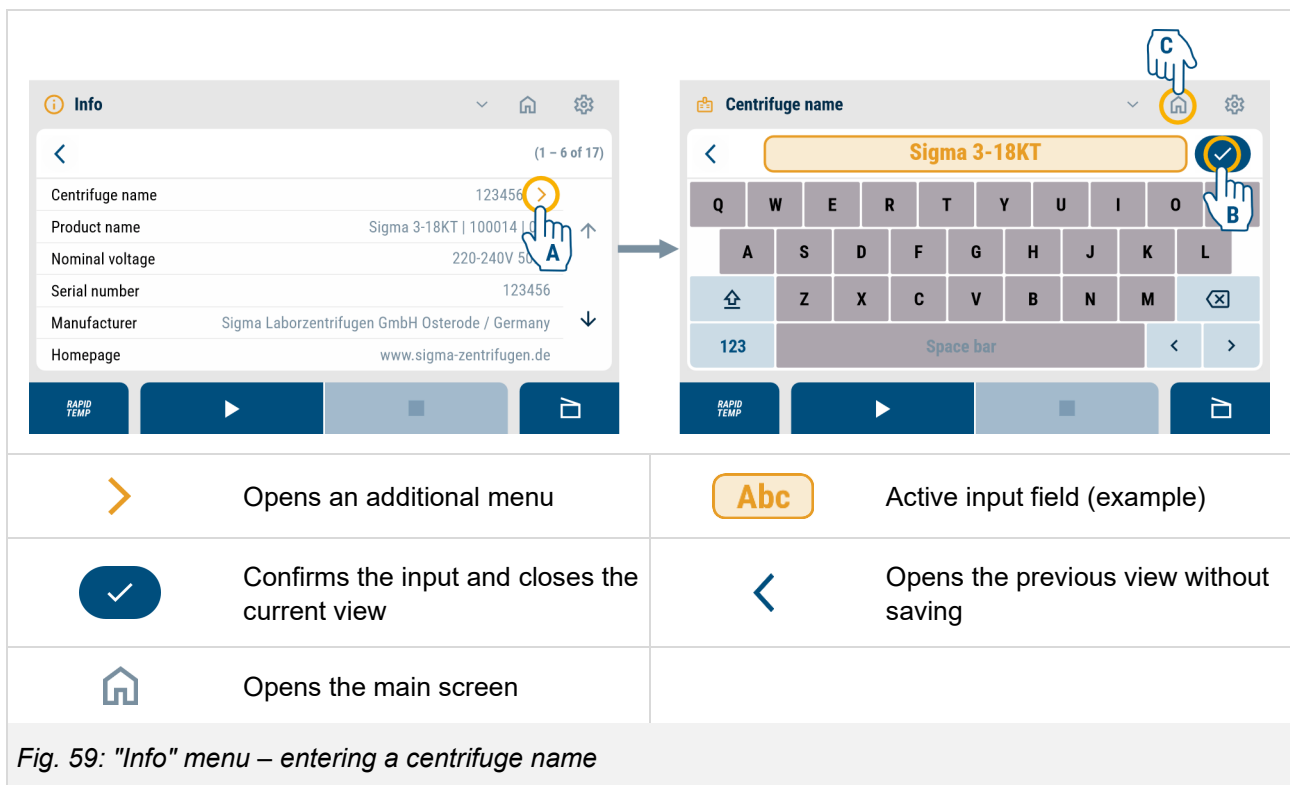

#### 6.5.10.6 Service


 This menu can be used to retrieve various types of service information.

Retrieving information (here: operating manual):

► The user interface shows the "Settings" menu.

1. Select the *Next* button to open the "Service" menu.
2. Select the *Next* button (A) to open the "Operating manual" menu.
  - A dialog box with a QR code for downloading the operating manual will be displayed.
3. Scan the QR code with a mobile device.
4. Confirm (B) to close the dialog box.
5. Select the *Home* button (C).
  - The main screen will be displayed.

## Operation

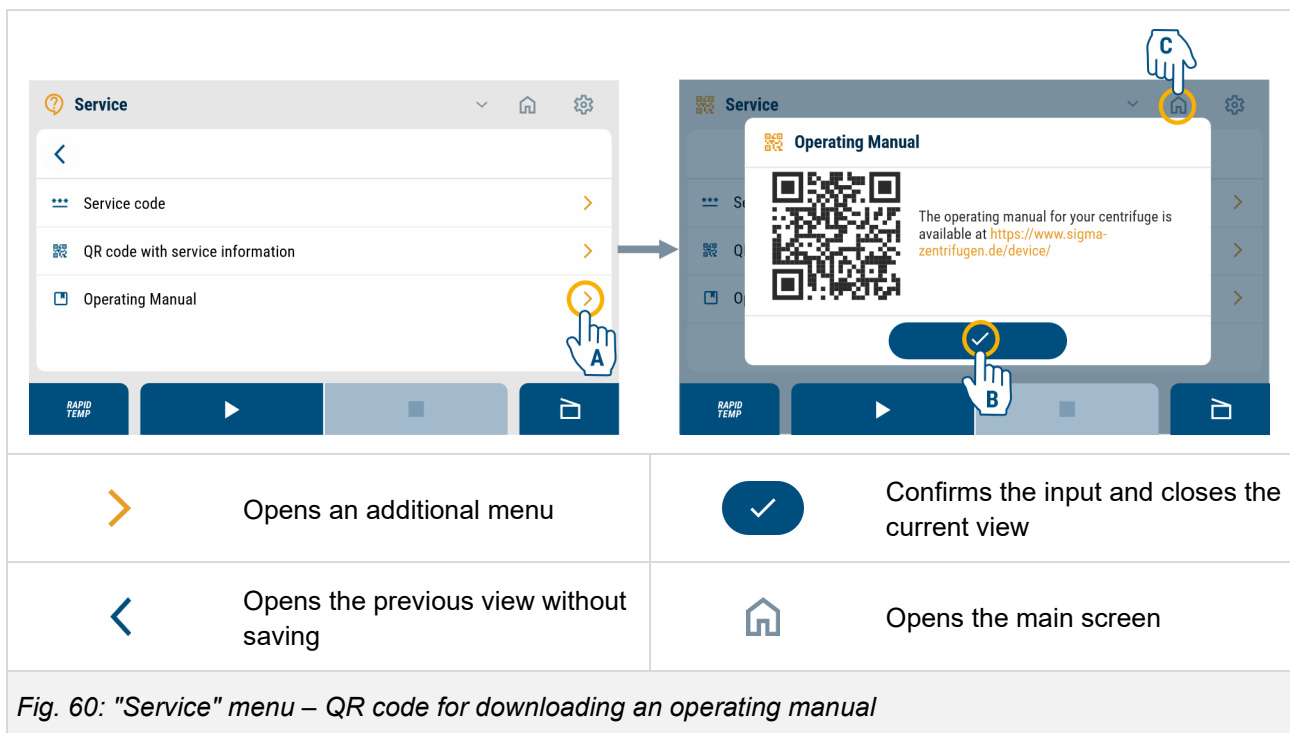

### \*\*\* Service code

This menu opens a keypad for entering a service code for maintenance and repairs



Only authorised specialist personnel are authorised to enter service codes.



### QR code with service information

This menu generates a QR code to access all of the relevant service information in the event of errors or problems. When the QR code is scanned by way of the camera of a mobile phone, the centrifuge data and a list of errors will be copied into an e-mail so that they can be sent to a certified service partner ([www.sigma-zentrifugen.de/en/sales-partners](http://www.sigma-zentrifugen.de/en/sales-partners)).



### Operating manual

This menu generates a QR code for downloading an operating manual ([www.sigma-zentrifugen.de/device](http://www.sigma-zentrifugen.de/device)).


### 6.5.10.7 Cycles

 This menu provides an overview of the cycle data.

Retrieving an overview (here: cycles current setup)

► The user interface shows the "Settings" menu.

1. Select the *Next* button to open the "Cycles" menu.
2. Select the *Next* button (A) to open the "Cycles current setup" menu.
  - An overview of the cycle data of the current setup will be displayed.
3. Select the *Home* button (B).

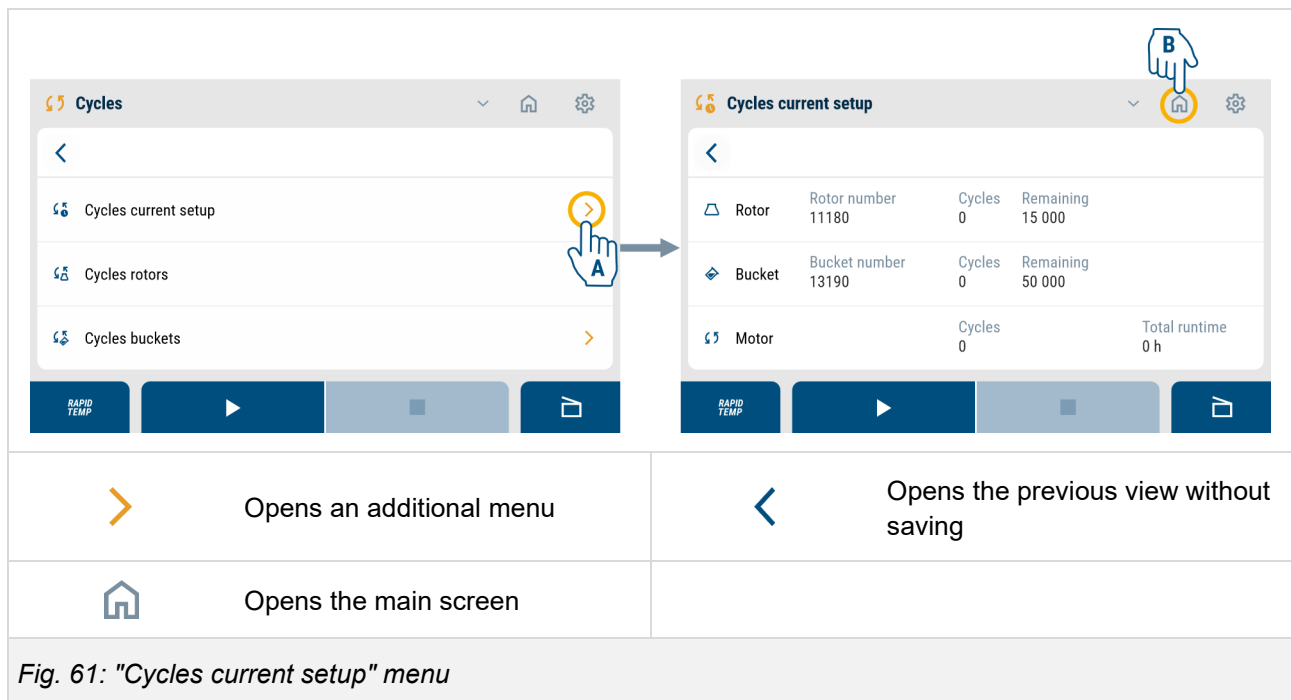

 The information engraved on the rotor or bucket regarding the service life is binding. The information provided by the control system is for informational purposes only and, thereby, non-binding.

#### **Cycles current setup**

This menu shows the rotor and bucket number as well the completed and remaining cycles of the current rotor and buckets. In addition, the number of cycles and the total runtime of the motor are displayed.

#### **Rotor cycles**

This menu shows the completed and remaining cycles for all the rotors that are used.

## Operation



### Bucket cycles

This menu shows the completed and remaining cycles for all the buckets that are used.

### 6.5.10.8 Quality Control




When the *Quality control* function is active, the runtime, speed and temperature will be continuously checked during the centrifugation.

The QC data can be exported into a USB stick for external evaluation (see chapter 6.5.11.2 - "Exporting QC data").

Opening the "Quality control" overview

- The user interface shows the "Settings" menu.
  1. Select the *Next* button to open the "Quality Control" menu.
    - An overview of the overall status of the last centrifugation run will be displayed.
  2. Select the *Home* button (B).
    - The main screen will be displayed.

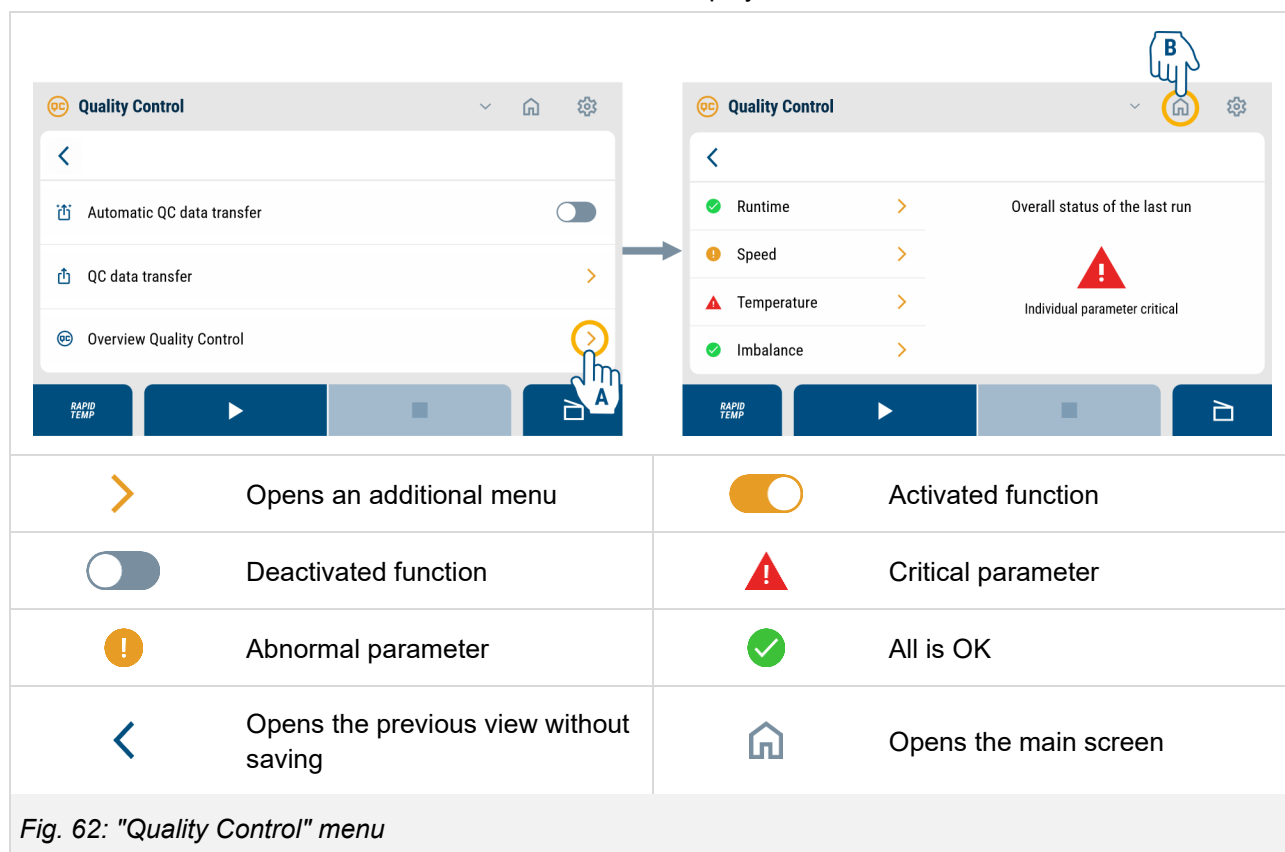



### Automatic QC data transfer

For the automatic transfer of the QC data, a USB stick must have been plugged in and detected by the system. The *Toggle* button can be used to activate the automatic transfer so that the QC data will be automatically transferred to the USB stick after every completed centrifugation run. In the factory settings, this function is deactivated.



### QC data transfer

This menu can be used to transfer the QC data to a USB stick. When the submenu is opened, a list of all of the available QC data will be displayed. The user can then select specific data for the transfer.



### Quality control overview

#### Runtime and speed (RPM)

The *Quality control* function is activated jointly for the runtime and speed. The RCF value will be calculated based on the speed value.

- The runtime will be considered faulty if the centrifugation run is stopped prematurely, e.g. by way of the *Stop* button.
- The speed will be considered faulty if the set speed is not reached or if the speed deviates from the set value by  $\pm 50$  rpm.

#### Temperature

The *Quality control* function must be activated separately for the temperature parameter. Values between  $\pm 2^{\circ}\text{C}$  and  $\pm 8^{\circ}\text{C}$  can be set for the permissible temperature deviation.

If the current value is outside of the permissible temperature deviation range, an error message will be displayed.

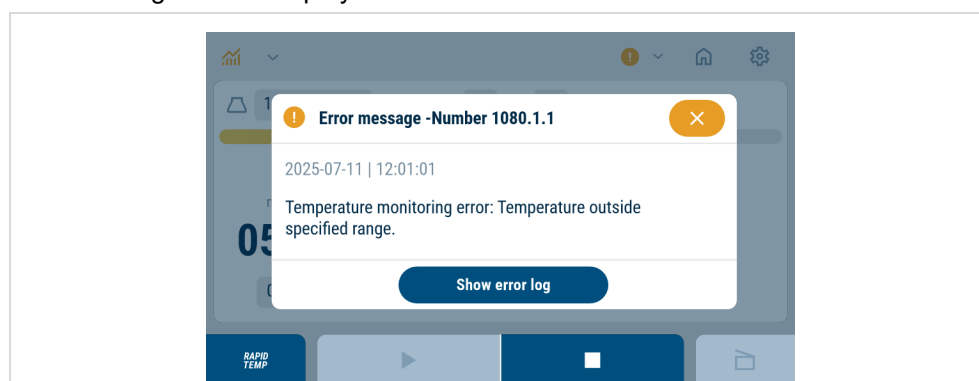


Fig. 63: "Quality control" dialog box – temperature error message

### 6.5.11 Data transfer

If a USB stick is plugged into the USB port of the control panel and the system detects the USB stick, the USB symbol will be displayed on the status bar.

The following transactions are possible:

- Export and import of programs (see chapter 6.5.11.1 - "Exporting and importing programs")
- Export of QC data for external evaluation (see chapter 6.5.11.2 - "Exporting QC data")

#### 6.5.11.1 Exporting and importing programs



Saved programs can be exchanged between centrifuges with a Spincontrol T control system:

- During the transfer from a refrigerated centrifuge (KT) to a non-refrigerated centrifuge (T), the temperature values will not be transferred.
- During the transfer from a non-refrigerated centrifuge (T) to a refrigerated centrifuge (KT), a temperature of 20°C will be automatically added to the programs. This value must be adapted manually.

For the program import, the centrifuge must be at a standstill.

Procedure (here: program export):

- ▶ A USB stick has been plugged in and detected by the unit.
- ▶ The status menu is open.
- 1. Select the *Next* button (A) to access the USB stick.
- 2. Select the *Drop-down* button (B) to open the drop-down menu for selecting a program export or import.
- 3. Select "Export programs" (C).
  - A list of the available programs will be displayed.
- 4. Select the programs for the export (D).
- 5. Select the *Data transfer* button (E).

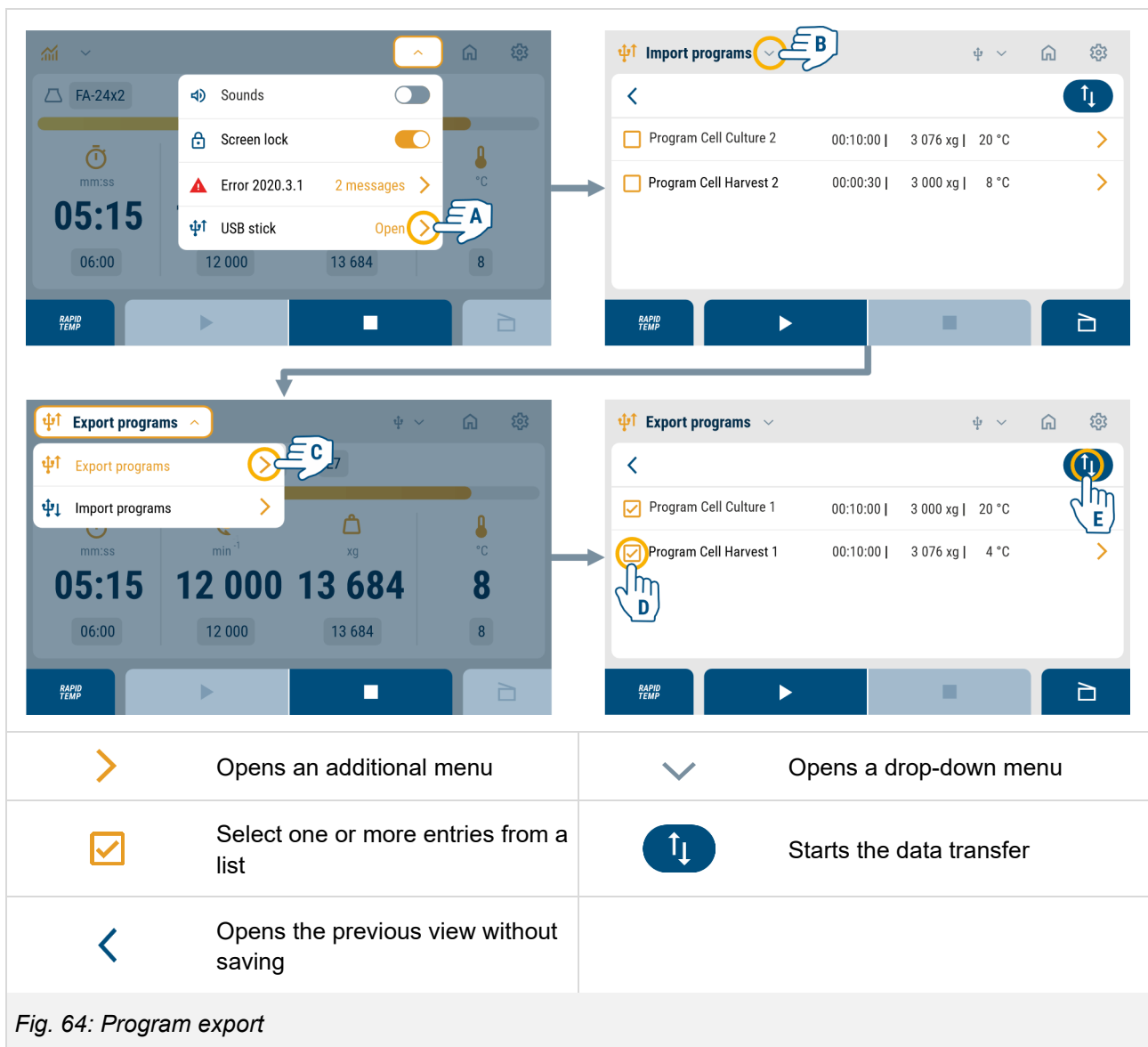


Fig. 64: Program export

- A dialog box will be displayed.
- 6. Select the *Start transfer* button.
- The data will be transferred.

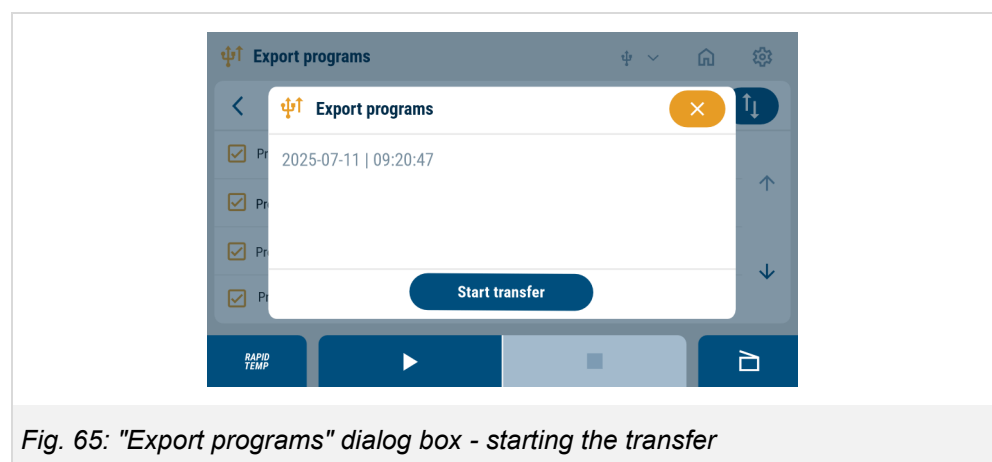


Fig. 65: "Export programs" dialog box - starting the transfer

## Operation

- Following the completion of the data transfer, an additional dialog box will be displayed.

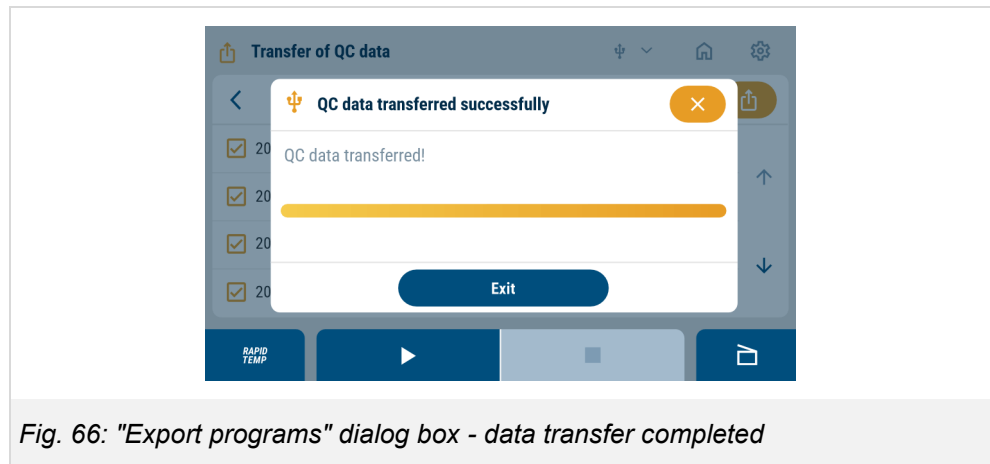


Fig. 66: "Export programs" dialog box - data transfer completed

7. Select the *Exit* button to exit the data transfer.
8. Unplug the USB stick from the unit.
  - The USB symbol disappears from the status bar.

### 6.5.11.2 Exporting QC data

- ▶ Saved QC data can be exported for external evaluation.
  - ▶ A USB stick has been plugged in and detected by the unit.
  - ▶ The user interface shows the main screen.
    1. Select the *Settings* button (A) to open the menu.
    2. Select the *Next* button (B) to open the "Quality control" menu.
    3. Select the *Next* button (C) to open the "QC data transfer" menu.
      - A list of the available data sets will be displayed.
    4. Select the QC data for the export (D).
    5. Select the *Data export* button (E).
      - The data will be transferred.

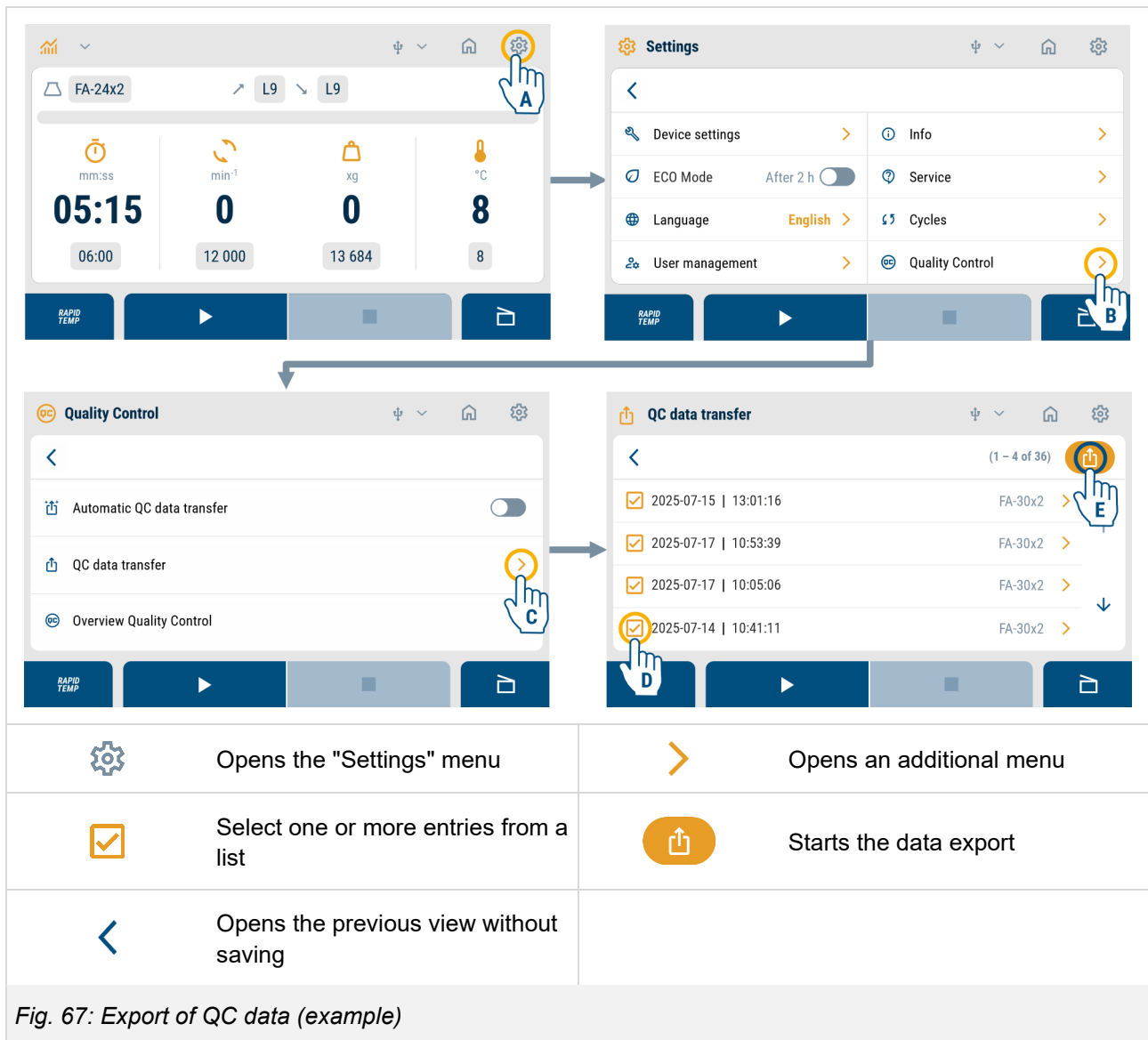


Fig. 67: Export of QC data (example)

- Following the completion of the data export, a dialog box will be displayed.

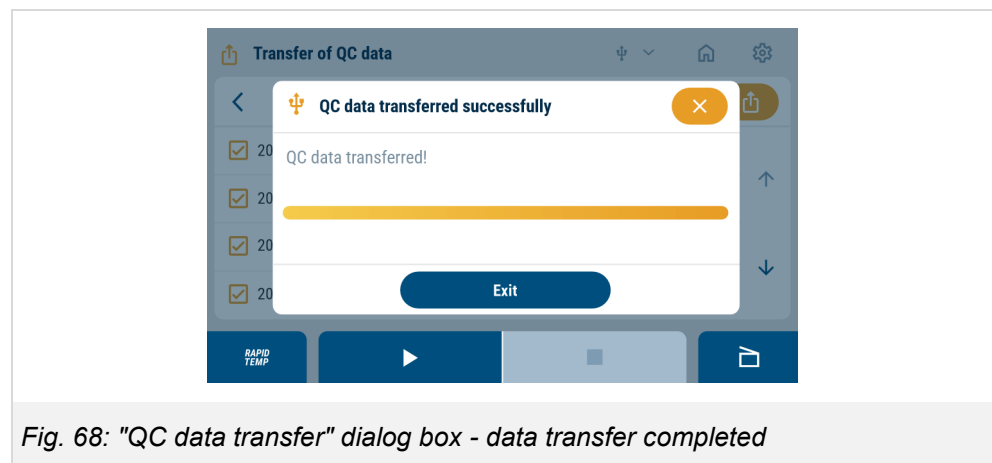


Fig. 68: "QC data transfer" dialog box - data transfer completed

6. Select the *Exit* button to exit the data export.
7. Unplug the USB stick from the unit.
  - The USB symbol disappears from the status bar.

## 7 Maintenance

### 7.1 Safety



#### **DANGER**

#### **Risk of injury from electric shock or sharp-edged components when removing the housing panels**

In the case of maintenance tasks requiring the removal of the housing panels, there is a risk of electric shock or mechanical injury.

- Maintenance tasks requiring the removal of the housing panels must only be carried out by authorised specialists.
- After the completion of maintenance tasks that involved the removal of the housing panels, a final standard-compliant check must be carried out by a qualified electrician.



#### **WARNING**

#### **Risk of toxic, pathogenic or radioactive contamination**

Toxic, pathogenic or radioactive substances can cause contamination if they are not handled properly.

It may be necessary to clean and/or decontaminate the laboratory centrifuge, rotors and accessories for reasons of safety before they can be maintained, repaired or handed over.

- Disinfect and clean the centrifuge, rotors, buckets, carriers and accessories after a contamination.
- Comply with the relevant safety instructions concerning the handling of toxic, pathogenic or radioactive substances.
- Wear personal protective equipment.
- Document any measures taken (see chapter 8.5 - "Returning the centrifuge to the manufacturer" / Declaration of decontamination)

### 7.2 General notes

The centrifuge and its accessories are subject to high mechanical stress. In order to withstand the continuous, extreme levels of stress within their area of application with high centrifugal force, the products are manufactured with the highest levels of precision. The service life of the centrifuge depends on its cleaning, maintenance, correct storage, frequency of use and other specific conditions of use. Diligent care and maintenance by the user will extend the service life of the products and prevent premature wear.

The centrifuge and accessories are made from different materials. Please check for any possible incompatibilities.

Prior to using cleaning agents and disinfectants other than those recommended, the user must consult with the manufacturer to ensure that they will not damage the centrifuge or its accessories.

For autoclaving, the long-term temperature resistance of the individual materials must be observed (see chapter 7.5 - "Autoclaving").

Protect the centrifuge and its accessories from intensive UV radiation (e.g. direct exposure to sunlight) and thermal stress (e.g. caused by heat generators).

**Frequency of the required activities**

| Chapter         | Activity                                                                | Frequency                      |
|-----------------|-------------------------------------------------------------------------|--------------------------------|
|                 | <b>Cleaning</b>                                                         |                                |
| chapter 7.3.1   | Cleaning of the centrifuge                                              | As required                    |
| chapter 7.3.2   | Cleaning of the rotors and buckets                                      | As required                    |
| chapter 7.3.3   | Cleaning of the adapters and carriers                                   | As required                    |
|                 |                                                                         |                                |
|                 | <b>Disinfection</b>                                                     |                                |
| chapter 7.4.1   | Disinfection of the centrifuge                                          | As required                    |
| chapter 7.4.2   | Disinfection of the rotors and buckets                                  | As required                    |
| chapter 7.4.3   | Disinfection of the adapters and carriers                               | As required                    |
|                 |                                                                         |                                |
|                 | <b>Maintenance</b>                                                      |                                |
| chapter 7.6.1   | Inspection of the centrifuge for signs of damage                        | Monthly                        |
| chapter 7.6.1   | Greasing of the motor shaft                                             | After cleaning                 |
| chapter 7.6.1.1 | Check of the condenser for soiling                                      | Monthly                        |
| chapter 7.6.1.1 | Cleaning of the condenser                                               | As required                    |
| chapter 7.6.2   | Check of the rotors, buckets, adapters and carriers for signs of damage | Monthly                        |
| chapter 7.6.2   | Greasing of the rotor tie-down screw                                    | After cleaning                 |
| chapter 7.6.2   | Greasing of the load-bearing bolts                                      | As required/<br>after cleaning |

*Table 4: Table of the required activities and their frequency*



If corrosion or other damage occurs due to lack of care, the manufacturer cannot be held liable and no warranty claims can be made.

## 7.3 Cleaning tasks

For cleaning the centrifuge and its accessories:

- Use soapy water or other water-soluble, mild cleaning agents with a pH value between 6 and 8.
- Clean the centrifuge by hand only.
- The water temperature must be < 40°C.
- Avoid corrosive and aggressive cleaning agents.
- Do not use any solvents.
- Do not use any agents with abrasive substances.

---

### 7.3.1 Cleaning the centrifuge

- ▶ The centrifuge is switched off and disconnected from the power supply.
  - ▶ The lid is open.
1. Remove the rotor and buckets (see chapter 6.4.3 - "Removing the rotor and accessories").
  2. Clean the housing and the rotor chamber with the specified cleaning agents and a soft, lint-free cloth.
  3. Remove any cleaning agent residues completely with water and a soft, lint-free cloth.
  4. After cleaning, dry the surfaces with a soft, lint-free and absorbent cloth.
    - The centrifuge has been cleaned.

---

### 7.3.2 Cleaning the rotors and buckets

- ▶ The rotors and buckets are located outside of the centrifuge.
  - ▶ The covers have been removed from the rotors.
  - ▶ Any existing carriers or adapters have been removed.
1. Remove the seals of the caps/covers, rotors or buckets from the groove with the aid of the seal removal tool (ref. no. 170360).
  2. Remove any old grease and deposits from the load-bearing bolts of the swing-out rotors and from the grooves of the buckets.

**ATTENTION** – Rotors and buckets with an anodised coating must never be cleaned in a dishwasher. The dishwasher will remove the anodised coating of the rotors and buckets, thereby leading to cracks in areas under stress.

3. Clean the rotor, buckets and the associated seals with the specified cleaning agents and a soft, lint-free cloth.
4. To clean narrow holes, use a test tube brush with a wool tip and then rinse with water.
5. Remove any cleaning agent residues completely with water and a soft, lint-free cloth.
6. Dry the rotor, buckets and seals with a soft, lint-free and absorbent cloth or in a drying cabinet at 50°C maximum.
7. To dry the holes, let the rotor dry lying upside down on a surface.
  - Rotors and buckets have been cleaned.

**ATTENTION** – The chemical resistance of rotors, buckets and carriers made of plastic decreases at higher temperatures (e.g. during drying). See the table in chapter 11.6 - "Table of chemical resistance data".

### 7.3.3 Cleaning the adapters and carriers

- ▶ The adapters and carriers have been removed from the rotor or bucket.
- 1. Clean the adapters and carriers with the specified cleaning agents and a soft, lint-free cloth.
- 2. To clean narrow holes, use a test tube brush with a wool tip and then rinse the holes complete with water.
- 3. Remove any cleaning agent residues completely with water and a soft, lint-free cloth.
- 4. Dry the adapters with a soft, lint-free and absorbent cloth.
- 5. To dry the holes, let the adapter dry lying upside down on a surface.
  - The adapters and carriers have been cleaned.

**ATTENTION** – Adapters and carriers with an anodised coating must never be cleaned in a dishwasher. The dishwasher will remove the anodised coating of the adapters and carriers, thereby leading to cracks in areas under stress.

## 7.4 Disinfection

-  Prior to a disinfection, the cleaning steps described in chapter 7.3 - "Cleaning tasks" must be performed.

### 7.4.1 Disinfecting the centrifuge

- To disinfect the centrifuge (control panel, contact surfaces and rotor chamber), use an alcohol-based surface disinfectant with a concentration between 30% and 75% (e.g. Bacillo® AF).
  - Disinfect the lid window and the display (acrylic glass), motor cover (silicone) and (screen-printed) labelling with a disinfectant that is suitable for the material (e.g. Bacillo® 30 sensitive tissues).  
Never use a surface disinfectant with a high alcohol concentration!
  - Follow the manufacturer's instructions regarding the application time and concentration of the disinfectant.
- 
- ▶ The centrifuge is off and disconnected from the power supply.
  - ▶ The lid is open.
  - ▶ There is no rotor inside the rotor chamber.
  - ▶ The centrifuge has been cleaned as described in chapter 7.3.1 - "Cleaning the centrifuge".
1. Wipe down the housing and rotor chamber with the above-mentioned disinfectant and a soft, lint-free cloth.
  2. Comply with the specified application time.
  3. Let the disinfected surfaces dry completely.
    - The centrifuge has been disinfected.

#### **7.4.2 Disinfecting the rotors and buckets**

- To disinfect the rotors, buckets and their seals, use an alcohol-based surface disinfectant with a concentration between 30% and 75% (e.g. Bacillol® AF).
  - Follow the manufacturer's instructions regarding the application time and concentration of the disinfectant.
- The rotors and buckets have been cleaned as described in chapter 7.3.2 - "Cleaning the rotors and buckets".
1. Wipe down the rotor, buckets and their seals with the above-mentioned disinfectant and a soft, lint-free cloth.
  2. To disinfect narrow holes, fill them with the disinfectant without entrapping any air.
  3. Ensure complete wetting.
  4. Comply with the specified application time.
  5. Let the disinfected rotors, buckets and seals dry completely.
    - The rotors and buckets have been disinfected.

---

#### **7.4.3 Disinfecting the adapters and carriers**

- To disinfect the adapters and carriers, use an alcohol-based surface disinfectant with a concentration between 30% and 75% (e.g. Bacillol® AF).
  - Follow the manufacturer's instructions regarding the application time and concentration of the disinfectant.
- The adapters and carriers have been cleaned as described in chapter 7.3.3 - "Cleaning the adapters and carriers".
1. Wipe down the adapters and carriers with the above-mentioned disinfectant and a soft, lint-free cloth.
  2. To disinfect narrow holes, fill them with the disinfectant without entrapping any air.
  3. Ensure complete wetting.
  4. Comply with the specified application time.
  5. Let the disinfected adapters and carriers dry completely.
    - The adapters and carriers have been disinfected.

## 7.5 Autoclaving

The service life of the accessories depends on their frequency of use and on the autoclaving frequency.

- The caps have been removed from the vessels in order to prevent the vessels from being deformed.
  - The seals have been removed from the caps/covers and buckets.
1. Position the accessories in the autoclave and run an autoclave cycle.
  2. After each sterilisation, check the accessories for colour changes, structural changes or leaks and replace the affected parts immediately if any signs are detected.
  3. After the autoclave run, reinstall the seals of the caps/covers and buckets (see chapter 7.3.2 - "Cleaning the rotors and buckets").

**i** The specified autoclaving parameters exclusively refer to the material resistance and not to the degree of sterility.

| Category                    | Name                                                 | Material abbreviation | 121°C<br>20 minutes | Notes               |
|-----------------------------|------------------------------------------------------|-----------------------|---------------------|---------------------|
| <b>Rotors and covers</b>    | Aluminium rotors                                     | AL                    | yes                 |                     |
|                             | Polypropylene rotors                                 | PP                    | no                  |                     |
|                             | Aluminium covers (without seals) for angle rotors    | AL                    | yes                 | 50 cycles max.      |
|                             | Polycarbonate covers for angle rotors                | PC                    | no                  |                     |
|                             | Polyallomer covers for angle rotors                  | PA                    | no                  |                     |
|                             | Polysulfone covers for angle rotors                  | PSU                   | yes                 | 100 cycles max.     |
| <b>Buckets and caps</b>     | Aluminium buckets (without seals)                    | AL                    | yes                 |                     |
|                             | Polyamide buckets                                    | PA                    | no                  | 13035, 13296, 13299 |
|                             | Polyphenylsulfone caps                               | PPSU                  | yes                 | 100 cycles max.     |
|                             | Polysulfone caps                                     | PSU                   | yes                 | 100 cycles max.     |
| <b>Adapters</b>             | Polyallomer racks                                    | PA                    | no                  |                     |
|                             | Polyallomer racks                                    | PC                    | no                  |                     |
|                             | Polyallomer racks                                    | PP                    | no                  |                     |
| <b>Vessels</b>              | Stainless steel vessels and bottles                  | --                    | yes                 |                     |
|                             | Glass vessels                                        | --                    | yes                 |                     |
|                             | Polyethylene vessels                                 | PE                    | no                  |                     |
|                             | Polyflor vessels                                     | PF                    | yes                 | 100 cycles max.     |
|                             | Polycarbonate vessels                                | PC                    | no                  |                     |
|                             | Polypropylene copolymer vessels                      | PPCO                  | yes                 | 20 cycles max.      |
|                             | Polystyrene vessels                                  | PS                    | no                  |                     |
| <b>Additional equipment</b> | Seals of the aluminium covers and buckets            | NBR                   | yes                 | 50 cycles max.      |
|                             | Stainless steel taring weights for blood bag systems | --                    | yes                 |                     |

Table 5: Autoclaving table

## **7.6 Maintenance tasks**

---

### **7.6.1 Maintaining the centrifuge**

The service life of the centrifuge depends on several factors, e.g. the type and frequency of use, the area of application and the performance of care, service and maintenance tasks.

- ▶ The centrifuge is off and disconnected from the power supply.
- ▶ There is no rotor inside the rotor chamber.

#### **Post-centrifugation measures**

- In the event of a contamination with toxic, pathogenic or radioactive substances, clean all of the internal surfaces immediately with a suitable decontamination agent depending on the type of contamination (see chapter 7.3 - "Cleaning tasks" and chapter 7.4 - "Disinfection").
- Liquids, such as water, solvents, acids and lyes, must be removed from the rotor chamber immediately using a cloth in order to prevent damage to the motor bearings.

#### **Post-cleaning measures**

- Grease the motor shaft and rotor tie-down screw slightly with a small amount of heavy-duty grease for load-bearing bolts and distribute the grease with a cloth so that it forms a thin layer.
- Rub talcum powder into the dry lid seal, if necessary.

#### **Visual inspection**

- Perform a visual inspection of the centrifuge at least monthly to detect the following signs of damage:
  - Cracks
  - Visible surface damage
  - Pressure marks
  - Signs of corrosion
  - Other changes
- Open the centrifuge lid completely and ensure that it remains in this position.

#### **Service life**

The centrifuge has an expected service life of 10 years if all of the specified maintenance intervals are observed and all of the necessary maintenance tasks are performed without delay. Non-compliance will shorten the service life of the centrifuge accordingly.

The availability of spare parts can no longer be guaranteed after 10 years from the date of manufacture of the centrifuge.

### 7.6.1.1 Cleaning the condenser

Refrigerated centrifuges with an air-cooled refrigeration system use a finned condenser to cool the refrigerant that is compressed by the refrigeration unit.

Dust and dirt on the condenser's ventilation opening impede cooling by the air flow and reduce the heat transfer, thereby compromising the performance of the refrigeration system.

- Check the condenser at least monthly for soiling and clean it whenever necessary.
- In case of dirt build-up on the condenser fins and tubes, contact the Sigma Service (see chapter 8.4 - "Service contact").

### 7.6.1.2 Removing broken glass



#### **CAUTION**

#### **Risk of injury from broken glass**

Broken glass in the rotor chamber can cause cuts.

- Broken glass must be removed immediately and completely.
  - Wear personal protective equipment when removing the broken glass.
- 

#### **Potential problems caused by broken glass in the rotor chamber**

- Broken glass damages the surface coating of the buckets, thereby causing corrosion.
- Broken glass in the adapters of the rotors and buckets causes more glass breakage.
- Broken glass in the swivel bearings of the load-bearing bolts prevents the buckets from swinging out evenly, thereby causing imbalance.
- Broken glass in the rotor chamber causes metallic abrasion due to the strong air circulation. The fine metal dust contaminates the rotor chamber, rotor and samples. At the same time, it also damages the surfaces of the rotor chamber, rotor, buckets, adapters, carriers and accessories.

#### **Removing broken glass and metal dust from the rotor chamber**

1. Apply a thick layer of petroleum jelly or similar to the upper third of the rotor chamber.
2. Switch the centrifuge on and let the rotor rotate for a few minutes at medium speed (approx. 2,000 min<sup>-1</sup>). The broken glass and metal dust will be bound in the grease.
3. Remove the grease with a cloth.
4. Repeat this process until all the broken glass and metal dust has been removed.

## **7.6.2 Maintaining the rotors, buckets, adapters and carriers**

Chemical reactions and pressure corrosion (a combination of alternating pressure and chemical reactions) can attack or destroy the structure of metals. Barely detectable cracks on the surface expand and weaken the material without leaving any clearly visible signs.

In the case of plastic materials, the resistance against chemicals decreases when the temperature increases (see the table in chapter 7.3.2 - "Cleaning the rotors and buckets").

### **Post-centrifugation measures**

- Liquids, such as water, solvents, acids and lyes, must be removed immediately and completely from the rotor, buckets, adapters, carriers and plastic accessories with a cloth.

### **Post-cleaning measures**

- Grease the rotor tie-down screw slightly with a small amount of heavy-duty grease for load-bearing bolts and distribute the grease with a cloth so that it forms a thin layer.
- Grease the bayonet lock in the cover of an angle rotor (see chapter 7.6.2.3 - "Maintaining the bayonet lock").
- When using aluminium buckets, slightly grease the load-bearing bolts of the rotors with a thin layer of grease for load-bearing bolts (see chapter 7.6.2.2 - "Maintaining the load-bearing bolts").

**CAUTION** – The load-bearing bolts must not be greased if plastic buckets are used.

-  For rotors with a Sigma Comfort coating, it is not necessary to grease the load-bearing bolts until the coating wears off.

### **Visual inspection**

- The products must be checked regularly, at least monthly, for the following:
  - Cracks
  - Visible surface damage
  - Pressure marks
  - Signs of corrosion
  - Other changes
- Check the holes of the rotors, buckets, adapters and carriers.
- Damaged parts must be replaced immediately for your own safety.

### **Service life**

Rotors and buckets have a limited service life.

- The rotors and buckets must be put out of service after 10 years. Any use after this period may be permissible in individual cases after an inspection is performed by the manufacturer.
- For safety reasons, rotors must be scrapped after a maximum of 50,000 cycles.

- Aluminium covers of Biosafe angle-rotors must be replaced after 3 years or 50 autoclaving cycles to guarantee the required tightness.
- If different information regarding the service life is engraved on a rotor or bucket, this information applies accordingly (see chapter 2.2 - "Marking of rotors and accessories").
- If a specification concerning the maximum number of cycles **and** a specification concerning the service life (i.e. a specific date) are provided, the specification that occurs first shall apply.

 The information engraved on the rotor or bucket regarding the service life is binding. The information provided by the control system is for informational purposes only and, thereby, non-binding.

### **Disposal of wear parts**

When replacing wear parts, e.g. a cover or bucket seal, contamination of the material must be assumed.

1. Dispose of all wear parts in compliance with the special provisions of the waste disposal regulations.

## **7.6.2.1 Rotors with a Sigma Comfort rotor coating**

### **Service life**

The service life of the coating varies and depends on the degree of utilisation of the rotor. Tests have shown a service life of 7,000 to 40,000 cycles. The service life of the Sigma Comfort rotor coating can be increased by way of these specific measures:

The following factors have an influence on the service life:

- **Speed/load:**  
The service life of the coating can be increased by a factor of 3.5 if the speed or load is decreased by approximately 15%. In case of a load decrease by 20%, the service life increases approximately by a factor of 5.5.
- **Temperature:**  
The service life of the coating decreases by approximately 30% at lower temperatures (approximately +4°C) compared to room temperature (approximately +20°C).
- **Use of buckets:**  
The buckets and rotor form a joint unit. For the coating to reach the maximum possible service life, ensure to always use the same set of buckets and to always install the individual buckets in their specific position within the rotor.
- **Cleaning frequency:**  
Regular cleaning of the rotors is indispensable. The more frequently a rotor needs to be cleaned, the shorter the service life of the coating will be.

 Traces of use on the black coating do not affect the anti-friction properties.




*Fig. 69: Load-bearing bolt with traces of use – no greasing required*

End of the service life of the Sigma Comfort rotor coating

When the coating has worn off nearly completely and the metallic surface of the load-bearing bolt becomes visible (see the following illustration), the anti-friction effect decreases. As a result, the buckets will swing out irregularly, thereby potentially causing the centrifuge to stop due to an imbalance.

From this moment on, the rotor must be used with greased load-bearing bolts until it reaches the end of its service life.

Mark the rotor accordingly to prevent it from being used without grease on the load-bearing bolts.

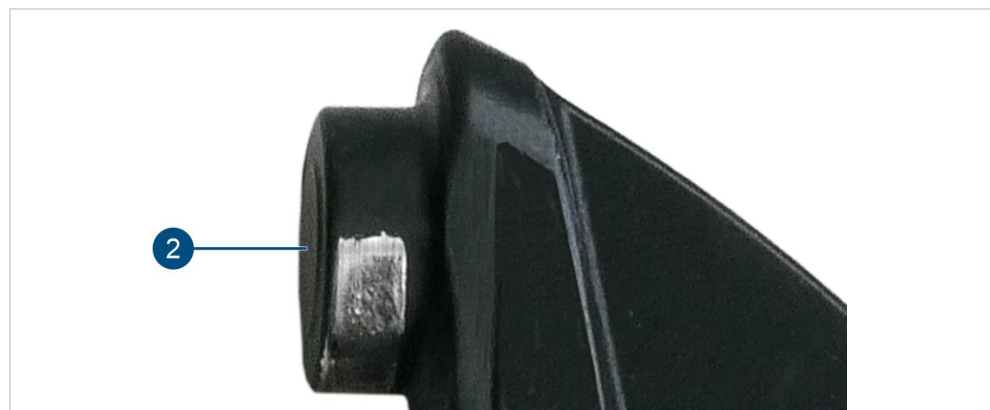

*Fig. 70: Load-bearing bolt with completely worn-off coating – greasing required*

### 7.6.2.2 Maintaining the load-bearing bolts

#### Rotors with Sigma Comfort coating

The load-bearing bolts of some types of swing-out rotors have an anti-friction coating. This coating prevents friction between the buckets and bolts.

- i** The bolts do not require any greasing during the service life of the Sigma Comfort rotor coating.

All swing-out rotors with a Sigma Comfort coating have a “C” shown after their number.

- 📄** A list of the available rotors with a Sigma Comfort coating as well as information about the service life of the coating can be found in the “Sigma Comfort rotor coating” document that comes supplied with every rotor which has a Comfort coating.

#### Rotors without Sigma Comfort coating

The load-bearing bolts of rotors that do not (or no longer) have anti-friction coating must be greased. This is the only way to ensure a uniform swing-out motion of the buckets and, therefore, the smooth operation of the centrifuge. Load-bearing bolts that are insufficiently greased may cause the centrifuge to stop due to an imbalance.

1. Clean the load-bearing bolts and bucket groove in order to remove the old grease.
2. Apply a small amount of grease for load-bearing bolts to each of the load-bearing bolts of the rotor cross (see the following illustration).



*Fig. 71: Applying the grease for load-bearing bolts*

1 Sufficient quantity of grease for one bolt

3. Install a bucket and swing it manually back to its swing-out position in order to distribute the grease.
4. Repeat this process with all the other buckets.


### 7.6.2.3 Maintaining the bayonet lock

#### Greasing the groove of the bayonet lock

The bayonet lock of the cover of an angle rotor must be greased after every cleaning process to ensure the smooth operation of the lock.


- ▶ The cover has been cleaned as described in chapter 7.3.2 - “Cleaning the rotors and buckets”.
  - ▶ The cover has been disinfected as described in chapter 7.4.3 - “Disinfecting the adapters and carriers”.
1. Place the cover upside down on a worktop.
  2. Apply a small amount of heavy-duty grease for load-bearing bolts (included in the scope of supply) to the groove of the cover screw (see the following illustration), e.g. by way of a cotton bud.

**CAUTION** – When greasing the groove, ensure that no grease reaches the seal. The area of the bioseal must be absolutely dry and free from grease.

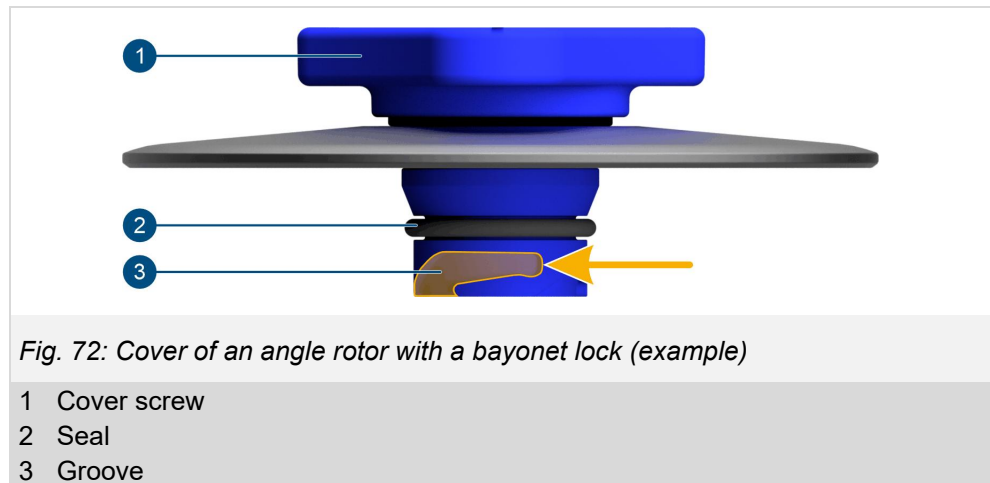

#### Adjusting the cover of an angle rotor

- ▶ The cover of an angle rotor is difficult to open and close.
  - ▶ The rotor with cover is located outside of the centrifuge.
1. Lay the rotor with its cover down on a firm surface.
  2. Place both hands (one above the other) on the cover screw and press down from above with your own bodyweight.
    - As a result, the cover can then be opened and closed easily.

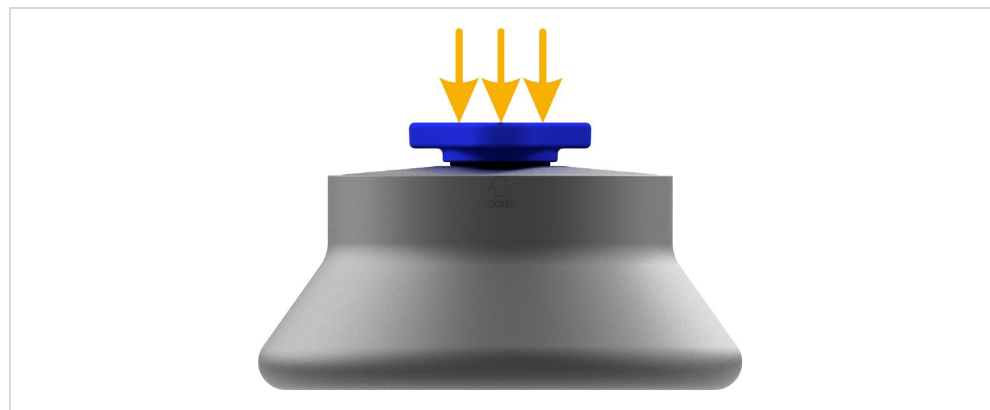

*Fig. 73: Measure to be taken when the cover is difficult to open or close.*

## 7.7 Inspection to be performed by the manufacturer

Centrifuges used in Germany are subject to mandatory testing as per the information sheet of the German Berufsgenossenschaft Rohstoffe und chemische Industrie FBRCI-025 (November 2023). The test intervals must ensure that the centrifuge can be used safely between tests. The operator is responsible for specifying the test intervals.

When specifying the test intervals, the following time frames should not be exceeded:

- Every 1-2 years in operating condition
- Every 3-4 years in dismantled condition

We recommend the following:

- The motor damping elements (rubber parts of the motor suspension) should be replaced every 3 to 4 years or every 10,000 cycles.

Sigma Laborzentrifugen GmbH performs the inspection and maintenance of centrifuges as per FBRCI-025 and in accordance with the German DGUV regulation 3 (DGUV = German statutory accident insurance) (electrical systems and equipment) .

### Maintenance contract

These services can be agreed as part of a maintenance contract.

[www.sigma-zentrifugen.de/en.html](http://www.sigma-zentrifugen.de/en.html) →  
 [Service] → [Maintenance, calibration and  
 maintenance service contracts]

### Scan QR code

If you have any questions  
 about this topic,  
 kindly contact  
 our service  
 department:



## 8 Troubleshooting

### 8.1 Safety



#### **DANGER**

#### **Risk of injury from electric shock or sharp-edged components when removing the housing panels**

In the case of maintenance tasks requiring the removal of the housing panels, there is a risk of electric shock or mechanical injury.

- Maintenance tasks requiring the removal of the housing panels must only be carried out by authorised specialists.
- After the completion of maintenance tasks that involved the removal of the housing panels, a final standard-compliant check must be carried out by a qualified electrician.



#### **WARNING**

#### **Risk of toxic, pathogenic or radioactive contamination**

Toxic, pathogenic or radioactive substances can cause contamination if they are not handled properly.

It may be necessary to clean and/or decontaminate the laboratory centrifuge, rotors and accessories for reasons of safety before they can be maintained, repaired or handed over.

- Disinfect and clean the centrifuge, rotors, buckets, carriers and accessories after a contamination.
- Comply with the relevant safety instructions concerning the handling of toxic, pathogenic or radioactive substances.
- Wear personal protective equipment.
- Document any measures taken (see chapter 8.5 - "Returning the centrifuge to the manufacturer" / Declaration of decontamination)

### 8.2 Identifying and eliminating faults and malfunctions

When a fault occurs, a dialogue box opens. If audible signals are activated, they sound at the same time as when the error or warning message is displayed.

Every fault is listed in the error log together with the error number and a time stamp. The error log shows the last 200 entries.

#### **Eliminating a fault**

► A dialogue box showing an error or warning message is displayed.

1. Eliminate the source of the error (see chapter 8.2.1 - "Table of possible faults" and chapter 8.2.2 - "Table of error messages").
2. Close the dialogue box or open the error log.



If an error message is acknowledged although the source of the error has not been eliminated, the dialogue box will be closed. The error message will not be cleared, but the centrifuge can be used.

3. Contact the service department if the recommended measures are not successful.

### Opening the error log

The error log can be opened from the status bar via a drop-down menu.

- The status bar displays an error or warning.
- 1. Open the drop-down menu from the status bar.
- 2. Select the "Error" menu.
  - The error log will be displayed.

### 8.2.1 Table of possible faults

| Error type                                                                                    | Possible cause                                             | Measures                                                                                                                                                                                                                               |
|-----------------------------------------------------------------------------------------------|------------------------------------------------------------|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| The display remains dark, the user interface is not displayed                                 | Mains power failure                                        | <ul style="list-style-type: none"> <li>Check the on-site fuse.</li> </ul>                                                                                                                                                              |
|                                                                                               | The mains power connector is not connected.                | <ul style="list-style-type: none"> <li>Connect the mains power connector.</li> </ul>                                                                                                                                                   |
|                                                                                               | The fuses have tripped.                                    | <ul style="list-style-type: none"> <li>Switch the centrifuge off and let it cool down for approx. 2 minutes.</li> <li>Switch the centrifuge back on.</li> </ul>                                                                        |
|                                                                                               | The mains power switch is switched off.                    | <ul style="list-style-type: none"> <li>Switch the centrifuge on.</li> </ul>                                                                                                                                                            |
| The centrifuge cannot be started: The "Start" button is inactive.                             | Miscellaneous                                              | <ul style="list-style-type: none"> <li>Switch the centrifuge off and back on. If the error occurs again, contact the Sigma service department.</li> </ul>                                                                              |
|                                                                                               | The lid lock has not closed properly.                      | <ul style="list-style-type: none"> <li>Open and close the lid. If the error occurs again, contact the Sigma service department.</li> </ul>                                                                                             |
| The centrifuge decelerates during a centrifugation run.                                       | Brief mains power failure                                  | <ul style="list-style-type: none"> <li>Select the <i>Start</i> button to restart the centrifugation run.</li> </ul>                                                                                                                    |
|                                                                                               | System error                                               | <ul style="list-style-type: none"> <li>Switch the centrifuge off and back on. If the error occurs again, contact the Sigma service department.</li> </ul>                                                                              |
| The centrifuge decelerates during a centrifugation run and an imbalance message is displayed. | Uneven distribution of the load                            | <ul style="list-style-type: none"> <li>Eliminate the imbalance and restart the centrifuge. If the error occurs again, contact the Sigma service department.</li> </ul>                                                                 |
|                                                                                               | The centrifuge is tilted.                                  |                                                                                                                                                                                                                                        |
|                                                                                               | Malfunction of the drive                                   |                                                                                                                                                                                                                                        |
|                                                                                               | The centrifuge has been moved during a centrifugation run. |                                                                                                                                                                                                                                        |
| The centrifuge decelerates during a centrifugation run and a temperature error is displayed.  | Unsuitable location                                        | <ul style="list-style-type: none"> <li>Adapt the conditions at the location of use or change the location (see chapter 5.1 - "Setting up the centrifuge"). If the error occurs again, contact the Sigma service department.</li> </ul> |

## Troubleshooting

| Error type                                                                      | Possible cause                            | Measures                                                                                                                                |
|---------------------------------------------------------------------------------|-------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------|
|                                                                                 | Soiled condenser                          | <ul style="list-style-type: none"> <li>Clean the condenser. If the error occurs again, contact the Sigma service department.</li> </ul> |
|                                                                                 | Insufficient refrigerant supply           | <ul style="list-style-type: none"> <li>Contact the Sigma service department.</li> </ul>                                                 |
| The lid cannot be opened.                                                       | The lid locks have not unlocked properly. | <ul style="list-style-type: none"> <li>Unlock the lid manually and contact the Sigma service department.</li> </ul>                     |
|                                                                                 | The lid seal is sticky.                   | <ul style="list-style-type: none"> <li>Clean the lid seal and apply some talcum powder.</li> </ul>                                      |
| Insufficient sedimentation of the product                                       | Insufficient fastening of the rotor       | <ul style="list-style-type: none"> <li>Fasten the rotor correctly and restart the centrifugation.</li> </ul>                            |
|                                                                                 | Defective sensor                          | <ul style="list-style-type: none"> <li>Contact the Sigma service department.</li> </ul>                                                 |
| The error message "Inspection by the manufacturer is recommended" is displayed. | Maintenance interval reached or exceeded  | <ul style="list-style-type: none"> <li>Contact the Sigma service department.</li> </ul>                                                 |

Table 6: Table of possible faults

### 8.2.2 Table of error messages

The error code displayed consists of the error group code and two additional digits to further categorise the error for authorised service personnel. An imbalance error, for example, is indicated by way of the error code 48.1.4.

| Error type          | Error group | Measures                                                                                                                                                                                              |
|---------------------|-------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Software problems   | 1<br>2      | <ul style="list-style-type: none"> <li>Let the rotor coast down.</li> <li>Switch the centrifuge off and back on.</li> <li>If the error occurs again: Contact the Sigma service department.</li> </ul> |
| Hardware problems   | 3           | <ul style="list-style-type: none"> <li>Switch the centrifuge off and back on.</li> <li>If the error occurs again: Contact the Sigma service department.</li> </ul>                                    |
| Speed sensor error  | 10          | <ul style="list-style-type: none"> <li>Let the rotor coast down.</li> <li>Switch the centrifuge off and back on.</li> <li>If the error occurs again: Contact the Sigma service department.</li> </ul> |
| Power supply error  | 18<br>19    | <ul style="list-style-type: none"> <li>Check the mains power supply and mains power connector.</li> <li>Check the fuses.</li> </ul>                                                                   |
| Drive control error | 20          | <ul style="list-style-type: none"> <li>Let the rotor coast down.</li> <li>Switch the centrifuge off.</li> <li>Ensure sufficient ventilation.</li> </ul>                                               |

| Error type                        | Error group                            | Measures                                                                                                                                                                                                                                                                                   |
|-----------------------------------|----------------------------------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Inverter error                    | 22                                     | <ul style="list-style-type: none"> <li>Switch the centrifuge off and back on.</li> <li>If the error occurs again: Contact the Sigma service department.</li> </ul>                                                                                                                         |
| Faulty input signal               | 30                                     | <ul style="list-style-type: none"> <li>Switch the centrifuge off and back on.</li> <li>If the error occurs again: Contact the Sigma service department.</li> </ul>                                                                                                                         |
| Data storage error                | 33                                     | <ul style="list-style-type: none"> <li>Let the rotor coast down.</li> <li>Switch the centrifuge off and back on.</li> </ul>                                                                                                                                                                |
| Temperature error                 | 40<br>41<br>43<br>44<br>46<br>47<br>48 | <ul style="list-style-type: none"> <li>Let the rotor coast down.</li> <li>Switch the centrifuge off.</li> <li>Let the centrifuge cool down.</li> <li>Air-cooled centrifuges: Ensure better ventilation.</li> <li>Water-cooled centrifuges: Ensure a sufficient water flow rate.</li> </ul> |
| Imbalance error                   | 48                                     | <ul style="list-style-type: none"> <li>Let the rotor coast down.</li> <li>Switch the centrifuge off.</li> <li>Eliminate the imbalance.</li> </ul>                                                                                                                                          |
| Lid error                         | 50<br>51<br>52                         | <ul style="list-style-type: none"> <li>Select the "Lid" button.</li> <li>Check whether there are any foreign objects in the opening for the lid hook.</li> <li>Close the lid.</li> </ul>                                                                                                   |
| Error of the refrigeration system | 53<br>54                               | <ul style="list-style-type: none"> <li>Contact the Sigma service department.</li> </ul>                                                                                                                                                                                                    |
| Communication error               | 71<br>72<br>73                         | <ul style="list-style-type: none"> <li>Let the rotor coast down.</li> <li>Switch the centrifuge off and back on.</li> </ul>                                                                                                                                                                |
| Rotor error                       | 80                                     | <ul style="list-style-type: none"> <li>Check the rotor.</li> <li>Check the rotor cycles.</li> </ul>                                                                                                                                                                                        |

*Table 7: Table of error messages*

### 8.3 Opening the lid manually (emergency lid release)



#### **WARNING**

#### **Risk of injury due to a spinning rotor**

If the emergency lid release is activated during a centrifugation run, the centrifuge will be switched off and then the rotor freely coasts to a standstill. There is a risk of injury due to the spinning rotor when reaching into the rotor chamber.

- Do not activate the emergency lid release until the rotor has come to a standstill.

In the event of a power supply failure, the centrifuge lid can be opened manually in order to withdraw the samples.

► The centrifuge is off and disconnected from the mains voltage.

1. Look through the window in the lid to verify that the rotor has come to a complete standstill.
2. Prise the plug out of the front (see the following illustration, item 1), e.g. by way of a screwdriver.



*Fig. 74: Position of the opening for the emergency lid release*

1 Plug

3. Push the supplied hexagon key horizontally into the opening as far as it will go.
4. Unlock the lid lock by turning the hexagon key anticlockwise.
5. Remove the key and seal the opening with the plug.
6. Reconnect the mains power connector.
7. Switch the centrifuge back on and acknowledge the error message.

## 8.4 Service contact

### In Germany

Sigma Laborzentrifugen GmbH  
An der Unteren Söse 50  
37520 Osterode (Germany)  
Telephone: +49 (0) 55 22 / 50 07-44 44  
E-mail: [support@sigma-zentrifugen.de](mailto:support@sigma-zentrifugen.de)

### Outside of Germany

The contact data of our international sales and service partners can be found at  
[www.sigma-zentrifugen.de/en.html](http://www.sigma-zentrifugen.de/en.html) → [Sales Partners]

#### Scan QR code

If you have any questions about this topic, kindly contact our service department:



## 8.5 Returning the centrifuge to the manufacturer

Complete and comprehensive information concerning the matter is required in order to ensure that all returns of centrifuges, rotors, buckets, adapters, carriers, accessories or spare parts are processed quickly and efficiently.

The necessary forms can be downloaded at  
[www.sigma-zentrifugen.de/en.html](http://www.sigma-zentrifugen.de/en.html) → [Service] → [Overhaul and repair].

### Completion of the declaration of decontamination

Due to the legal regulations for the protection of our employees and the environment, the manufacturer is obliged to document the harmlessness of all incoming goods.

- The form must be filled in completely as well as signed by authorised and specialised personnel.
- Affix the original document in a clearly visible manner to the outside of the packaging.



We will return the part/unit at your expense if no declaration of decontamination is provided.

### Completion of the return declaration

Product-specific data are required to ensure that the return can be processed without delay.

- Enclose the completed return declaration, including a detailed error description, with the defective part/unit.
- Fill in a separate form for every defective part even if several parts are returned in one package.



We will return the part/unit at your expense if no error description is provided.

### Cost estimate

Upon request, a cost estimate will be prepared prior to performing the repair.

-  If the cost estimate has still not been confirmed after 4 weeks, the part/unit will be returned at the customer's expense.

### **Pick-up order**

Upon request, we will commission a forwarding agent to collect the unit

- Complete the pick-up order and send it to [support@sigma-zentrifugen.de](mailto:support@sigma-zentrifugen.de)

### **Returning a product**

- Package the part/unit safely for transport, and in its original packaging, if possible. The original packaging can be supplied by the manufacturer for a fee.
- Return the product with the completed forms to  
Sigma Laborzentrifugen GmbH  
Service laboratory centrifuges  
An der Unteren Söse 50  
37520 Osterode  
Germany

-  If the part/unit is sent to the manufacturer in unsuitable packaging, it will be repackaged for return at your expense.

## 9 Disposal

### 9.1 Safety

Disposal regulations vary among the member states of the European Union.

- Comply with the national rules and regulations.

### 9.2 Disposal of the centrifuge



In accordance with the directive 2012/19/EU about waste electrical and electronic equipment (WEEE), Sigma centrifuges are marked with a symbol showing a crossed-out waste bin. It means:

- that the centrifuge was delivered after 13 August 2005.
- that the centrifuge must not be disposed of via the municipal waste system.

Products whose disposal is subject to the WEEE directive must be properly decontaminated and disposed of via a return scheme that ensures their proper collection, treatment, reuse or recycling.

Sigma Laborzentrifugen GmbH will take back appropriately labelled discarded units free of charge and ensure their proper disposal.

- Prior to returning the discarded unit, the operator must ensure that the unit has been decontaminated. A duly completed decontamination declaration must be enclosed (see chapter 8.5 - "Returning the centrifuge to the manufacturer").

### 9.3 Disposal of the packaging material

- Use the packaging for returning the centrifuge to the manufacturer.
- Or:
- Dispose of the packaging separately by material type.

## 10 Technical data

| Product data        |                                                                                   |                           |
|---------------------|-----------------------------------------------------------------------------------|---------------------------|
| Manufacturer        | Sigma Laborzentrifugen GmbH<br>An der Unteren Söse 50<br>37520 Osterode (Germany) |                           |
| Product designation | Laboratory centrifuge                                                             |                           |
| Product name        | Sigma 3-18KT                                                                      | Sigma 3-18KT IVD          |
| Part numbers:       | 100014; 100016;<br>100017                                                         | 100019; 100021;<br>100022 |

Table 8: Product data

| Connection data       |                    |
|-----------------------|--------------------|
| Power supply          | See the name plate |
| Protection class      | I                  |
| IP code               | 20                 |
| Power consumption (W) | 1 400              |

Table 9: Connection data

| Performance parameters                    |                  |
|-------------------------------------------|------------------|
| Speed                                     |                  |
| – Set speed (min <sup>-1</sup> )          | 100 ... 18 000   |
| – Display range (min <sup>-1</sup> )      | 0 ... 99 999     |
| – Resolution (min <sup>-1</sup> )         | 1                |
| – Measuring accuracy (min <sup>-1</sup> ) | ± 30 (at 18 000) |
| Max. capacity (ml)                        | 1 600            |
| Max. RCF (x g)                            | 31 655           |
| Max. kinetic energy (Nm)                  | 22 218           |

Table 10: Performance parameters

| Other parameters               |                                                    |
|--------------------------------|----------------------------------------------------|
| Time                           |                                                    |
| – Set time (m:s) or (h:m)      | 10 s ... 99 h 59 min,<br>Short run, continuous run |
| – Display range (m:s) or (h:m) | 00:00 ... 99:59                                    |
| – Resolution (s)               | 1                                                  |
| – Measuring accuracy (‰)       | ± 1                                                |
| Program storage locations      | 99                                                 |

Table 11: Other parameters

| Physical data             |     |
|---------------------------|-----|
| Height (mm)               | 337 |
| Height with open lid (mm) | 755 |
| Width (mm)                | 460 |
| Depth (mm)                | 716 |
| Weight without rotor (kg) | 80  |

Table 12: Physical data

| Sound pressure level |                  |
|----------------------|------------------|
| Noise level (dB(A))  | 61 at max. speed |

Table 13: Sound pressure level

| Operating fluids and quantities |       |
|---------------------------------|-------|
| Refrigerant                     | R452A |
| Global warming potential (GWP)  | 2 140 |
| Quantity (kg)                   | 0.180 |
| Max. permissible pressure (bar) | 27    |
| CO2 equivalent (t)              | 0.385 |

Table 14: Operating fluids and quantities

## 10.1 Ambient conditions

The information applies to an ambient temperature of +23°C and the rated voltage  $\pm 10\%$ . For a nominal voltage of 100 V or 200 V, tolerances of +10%/-5% apply.








The minimum temperatures are  $\leq +4^\circ\text{C}$  and they depend on the rotor type, speed and ambient temperature.









| Ambient conditions                       |                                                              |
|------------------------------------------|--------------------------------------------------------------|
| Use                                      | Indoor                                                       |
| Altitude above sea level (m)             | $\leq 2\,000$                                                |
| Pollution degree                         | 2                                                            |
| Ambient temperature ( $^\circ\text{C}$ ) | +5 ... +32                                                   |
| Humidity (%)                             | 80 at +5°C to 31°C,<br>linearly decreasing to 76.75 at +32°C |

Table 15: Ambient conditions

## 11 Appendix

### 11.1 Rotor program
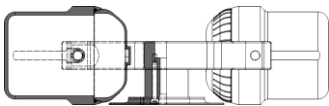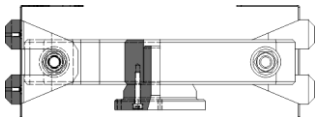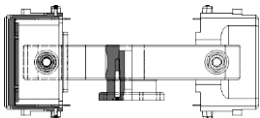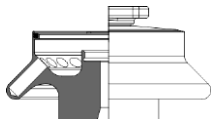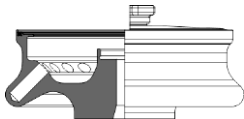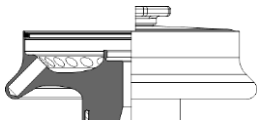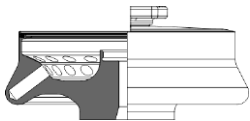
| Product                                                                             | Description                                                                            | Ref. no.        | Max. speed (min <sup>-1</sup> ) | Max. RCF (x g) | Max. capacity                     |
|-------------------------------------------------------------------------------------|----------------------------------------------------------------------------------------|-----------------|---------------------------------|----------------|-----------------------------------|
|    | SB-4x400 RB400<br>Swing-out rotor for round buckets (optionally Biosafe with a cap)    | 11180 / 13190   | 4 700                           | 4 248          | 4 x 400 ml                        |
|    | SP-2x MTP-60<br>Swing-out rotor, including carriers for microtiter plates              | 110360 / 130360 | 4 300                           | 2 625          | 2 x 60 mm                         |
|   | SPC-2x MTP-66C55<br>Swing-out rotor, including carriers for microtiter plates, Biosafe | 11241 / 13241   | 4 700                           | 3 285          | 2 x 66 mm, incl. cap<br>2 x 55 mm |
|  | FA-24x2<br>Angle rotor, Biosafe                                                        | 120154          | 18 000                          | 30 790         | 24 x 2 ml                         |
|  | FA-24x2 spin column<br>Angle rotor, Biosafe                                            | 120324          | 15 500                          | 26 054         | 24 x 2 ml spin column             |
|  | FA-30x2<br>Angle rotor, Biosafe                                                        | 120331          | 16 500                          | 31 655         | 30 x 2 ml                         |
|  | FA-48x2<br>Angle rotor, Biosafe                                                        | 120349          | 15 500                          | 26 860         | 48 x 2 ml                         |
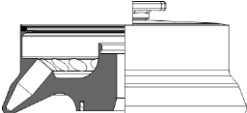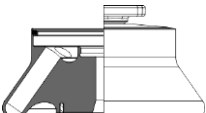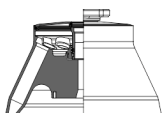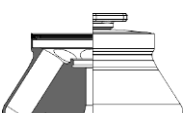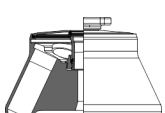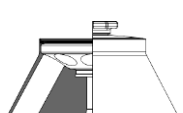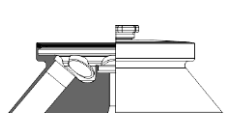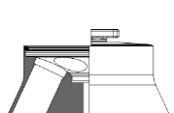
| Product                                                                             | Description                        | Ref. no. | Max. speed (min <sup>-1</sup> ) | Max. RCF (x g) | Max. capacity |
|-------------------------------------------------------------------------------------|------------------------------------|----------|---------------------------------|----------------|---------------|
|    | FA-20x5-V<br>Angle rotor, Biosafe  | 120322   | 15 400                          | 27 575         | 20 x 5 ml     |
|    | FA-10x10-U<br>Angle rotor, Biosafe | 120111   | 17 600                          | 30 475         | 10 x 10 ml    |
|    | FA-12x15-V<br>Angle rotor, Biosafe | 120215   | 13 500                          | 20 172         | 12 x 15 ml    |
|   | FA-8x30-U<br>Angle rotor, Biosafe  | 120158   | 17 500                          | 31 500         | 8 x 30 ml     |
|  | FA-6x50-V<br>Angle rotor, Biosafe  | 120250   | 14 500                          | 22 096         | 6 x 50 ml     |
|  | FA-8x50-U<br>Angle rotor, Biosafe  | 120150   | 13 000                          | 19 272         | 8 x 50 ml     |
|  | FA-8x50-V<br>Angle rotor, Biosafe  | 120309   | 9 200                           | 11 828         | 8 x 50 ml     |
|  | FA-6x94-U<br>Angle rotor, Biosafe  | 120395   | 12 500                          | 17 119         | 6 x 94 ml     |

V = Conical

U = Round bottom

**Table 16: Table of available rotors**

## Appendix

### 11.1.1 Spare parts for the rotors and covers

| Angle rotor (Biosafe) |                        | Cover    | Seal set | Rotor seal | Cover seal | Tool for O-ring removal |
|-----------------------|------------------------|----------|----------|------------|------------|-------------------------|
| Ref. no.              | Name                   | Ref. no. | Ref. no. | Ref. no.   | Ref. no.   | Ref. no.                |
| 120154                | FA-24x2                | 170154   | 160154   | 480400     | 481099     | 170360                  |
| 120324                | FA-24x2<br>spin column | 170324   | 160324   | 480401     | 481099     | 170360                  |
| 120331                | FA-30x2                | 170331   | 160331   | 480401     | 481099     | 170360                  |
| 120349                | FA-48x2                | 170349   | 160349   | 480401     | 481099     | 170360                  |
| 120322                | FA-20x5-V              | 170322   | 160322   | 480401     | 481099     | 170360                  |
| 120111                | FA-10x10-U             | 170111   | 160111   | 480400     | 481099     | 170360                  |
| 120215                | FA-12x15-V             | 170215   | 160215   | 480400     | 481099     | 170360                  |
| 120158                | FA-8x30-U              | 170158   | 160158   | 480400     | 481099     | 170360                  |
| 120250                | FA-6x50-V              | 170250   | 160250   | 480400     | 481099     | 170360                  |
| 120150                | FA-8x50-U              | 170150   | 160150   | 480400     | 481099     | 170360                  |
| 120309 <sup>1</sup>   | FA-8x50-V              | 170309   | 160309   | 480401     | 481099     | 170360                  |
| 120395                | FA-6x94-U              | 170395   | 160395   | 480402     | 481099     | 170360                  |

<sup>1</sup> approved for 3-18KT only

*Table 17: Table indicating the spare parts for the rotors and covers*

## 11.2 Speed, radius and relative centrifugal force (RCF)

The acceleration  $g$  that the samples are subjected to can be increased by expanding the radius in the rotor chamber and by increasing the rotational speed. The three parameters are interdependent and linked with one another by the following formula:

$$\text{Relative centrifugal force (RCF)} = 11.18 \times 10^{-6} \times r \times n^2$$

$r$  = radius in cm

$n$  = speed in  $\text{min}^{-1}$

RCF is dimensionless

When two values are entered, the third is defined by the specified formula. If the speed or radius is changed thereafter, the resulting relative centrifugal force will be automatically recalculated by the centrifuge control system. If the RCF value is changed, the speed will be adapted accordingly using the specified radius.

## 11.3 Density

The laboratory centrifuge is suitable for the separation of substances of different densities in mixtures with a maximum density of  $1.2 \text{ g/cm}^3$ . All the information on the speed of rotors, buckets and accessories refers to liquids with a density corresponding to this specification. If the density of the liquid is higher than this value, the maximum permissible speed of the centrifuge must be reduced by way of the following formula:

$$n = n_{\text{max}} \times \sqrt{(1.2/\rho)}$$

$\rho$  = density in  $\text{g/cm}^3$

## 11.4 Acceleration and deceleration curves

The curves are numbered in relation to increasing acceleration (from right to left). With curve 9, the centrifuge accelerates with maximum power. The time required to reach the set speed solely depends on the inertia of the rotor.

The deceleration curves are inverted images of the acceleration curves and are labelled with identical numbers. Curve 0 is an exception. It ensures a free coast-down.

### Linear curves (LIN)

| Curve | Type   | Gradient                     |
|-------|--------|------------------------------|
| L1    | Linear | 6 [min <sup>-1</sup> /s]     |
| L2    | Linear | 8 [min <sup>-1</sup> /s]     |
| L3    | Linear | 17 [min <sup>-1</sup> /s]    |
| L4    | Linear | 25 [min <sup>-1</sup> /s]    |
| L5    | Linear | 33 [min <sup>-1</sup> /s]    |
| L6    | Linear | 50 [min <sup>-1</sup> /s]    |
| L7    | Linear | 100 [min <sup>-1</sup> /s]   |
| L8    | Linear | 200 [min <sup>-1</sup> /s]   |
| L9    | Linear | 1,000 [min <sup>-1</sup> /s] |

Table 18: Table of the gradient of linear acceleration curves

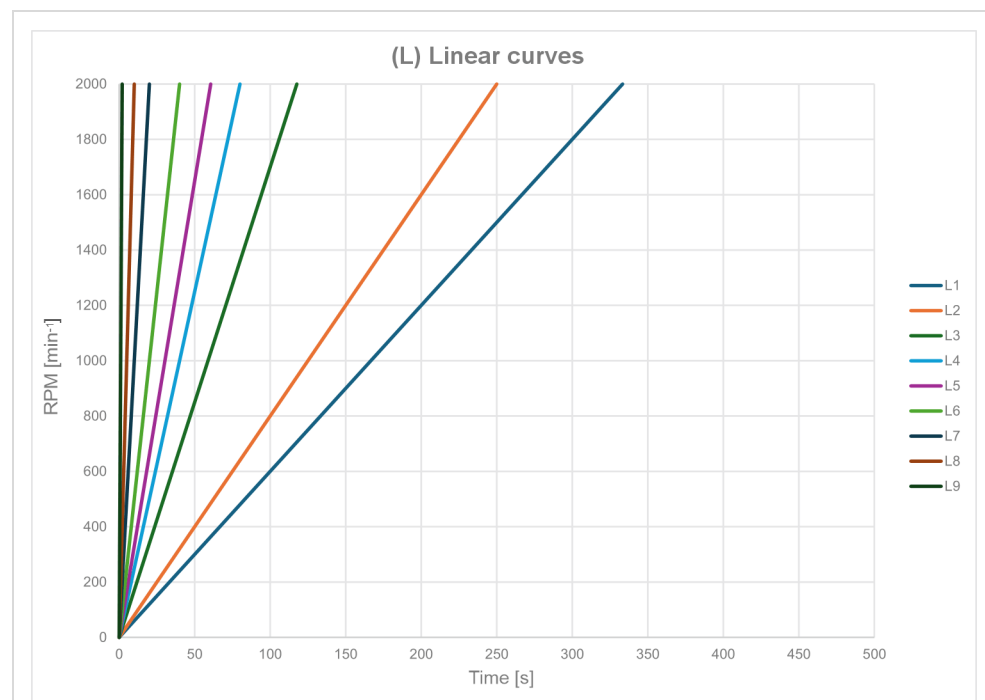


Fig. 75: Diagram of linear acceleration curves

### Quadratic curves (QUAD)

**i** Above 1,000 min<sup>-1</sup>, square curves will continue linearly.

| Curve | Type      | Time to 1,000 min <sup>-1</sup> | Gradient                     |
|-------|-----------|---------------------------------|------------------------------|
| Q1    | Quadratic | 333 [s]                         | 6 [min <sup>-1</sup> /s]     |
| Q2    | Quadratic | 250 [s]                         | 8 [min <sup>-1</sup> /s]     |
| Q3    | Quadratic | 118 [s]                         | 17 [min <sup>-1</sup> /s]    |
| Q4    | Quadratic | 80 [s]                          | 25 [min <sup>-1</sup> /s]    |
| Q5    | Quadratic | 61 [s]                          | 33 [min <sup>-1</sup> /s]    |
| Q6    | Quadratic | 40 [s]                          | 50 [min <sup>-1</sup> /s]    |
| Q7    | Quadratic | 20 [s]                          | 100 [min <sup>-1</sup> /s]   |
| Q8    | Quadratic | 10 [s]                          | 200 [min <sup>-1</sup> /s]   |
| Q9    | Quadratic | 2 [s]                           | 1,000 [min <sup>-1</sup> /s] |

Table 19: Table of the gradient of quadratic acceleration curves


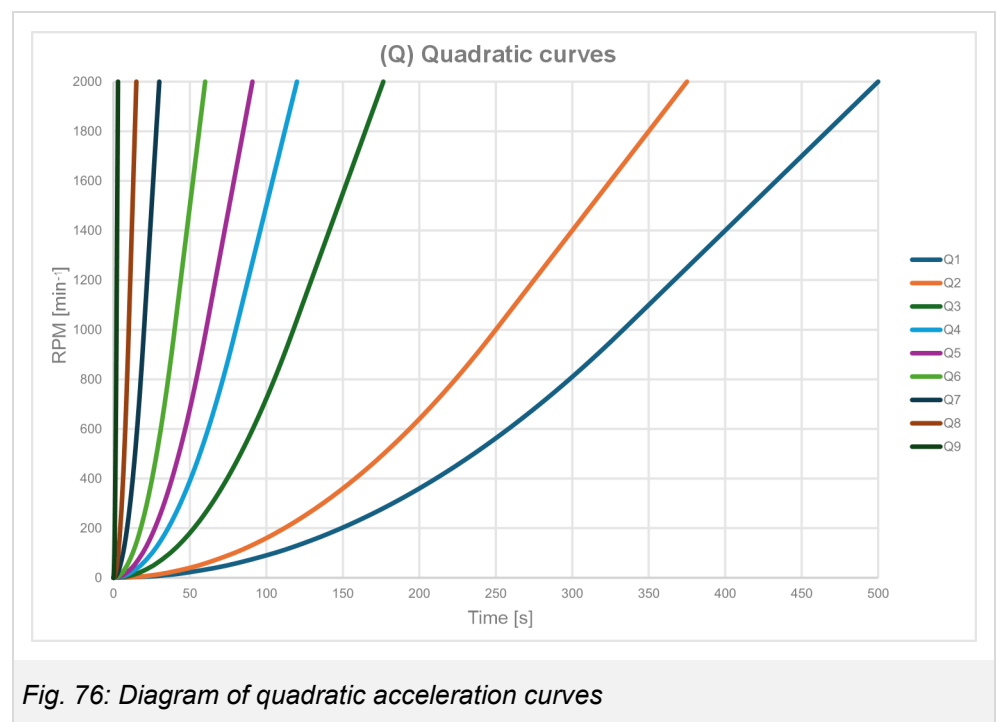


Fig. 76: Diagram of quadratic acceleration curves

### Optimised quadratic curves (QUADpro)

| Curve | Type    | Time to 1,000 rpm |
|-------|---------|-------------------|
| P1    | QUADpro | 333 [s]           |
| P2    | QUADpro | 250 [s]           |
| P3    | QUADpro | 118 [s]           |
| P4    | QUADpro | 80 [s]            |
| P5    | QUADpro | 61 [s]            |
| P6    | QUADpro | 40 [s]            |
| P7    | QUADpro | 20 [s]            |
| P8    | QUADpro | 10 [s]            |
| P9    | QUADpro | 2 [s]             |

Table 20: Table of the gradient of the optimised quadratic acceleration curves

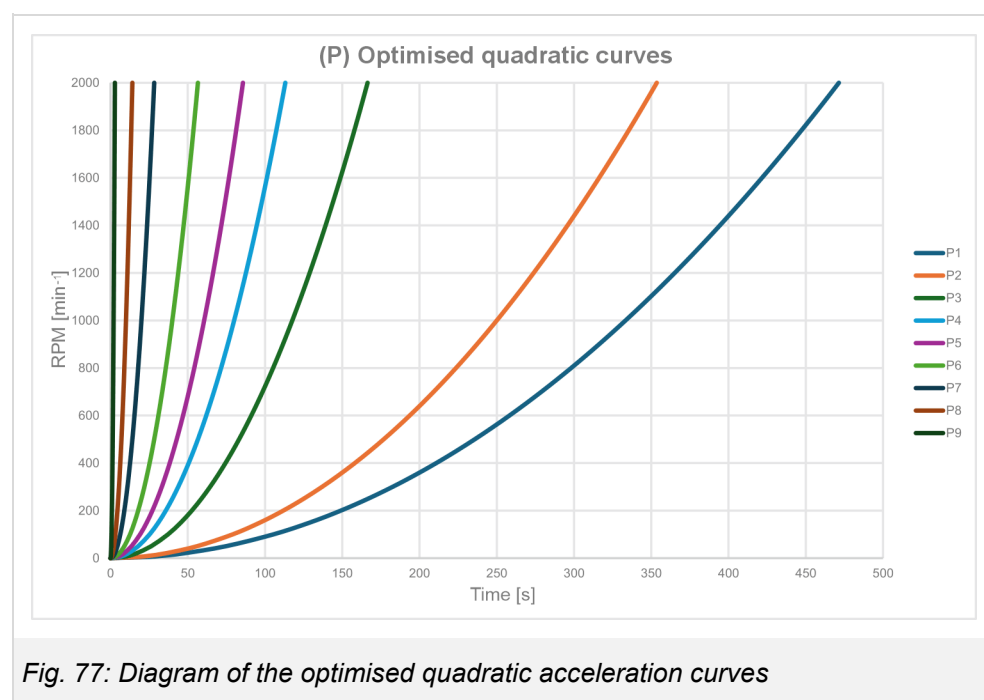

## 11.5 Table of symbols

<sup>1</sup> Symbols with this mark are used only in the Spincontrol T Expert control system.

### Control bar

| Symbol                                                                            | Description                              | Symbol                                                                             | Description                         |
|-----------------------------------------------------------------------------------|------------------------------------------|------------------------------------------------------------------------------------|-------------------------------------|
|  | "Start" button, active                   |  | "Start" button, inactive            |
|  | "Stop" button, active                    |  | "Stop" button, inactive             |
|  | "Lid" button, active                     |   | "Lid" button, inactive              |
|  | Automatic lid opening function activated |   | <i>RAPID TEMP</i> button, active    |
|  | <i>RAPID TEMP</i> button, inactive       |   | <i>RAPID TEMP</i> button, activated |

Table 21: Symbols of the control bar

### Configuration bar

| Symbol                                                                              | Description                                         | Symbol                                                                              | Description                                        |
|-------------------------------------------------------------------------------------|-----------------------------------------------------|-------------------------------------------------------------------------------------|----------------------------------------------------|
|  | "Rotor & buckets" button, opens an input menu       |  | "Curve" button, opens an input menu                |
|  | Generates a QR code for the Sigma Timer app         |  | Opens the previous view without saving             |
|  | Confirms the input and closes the current view      |  | Opens a window for creating a new list item        |
|  | Selects the previous interval (active) <sup>1</sup> |  | Selects the next intervals (inactive) <sup>1</sup> |
|  | Opens the input menu of the interval <sup>1</sup>   |  | Selects a profile view <sup>1</sup>                |
|  | Deletes the selected files                          |  | Closes a dialog box                                |
|  | Starts the program export                           |  | Starts the QC data transfer                        |

## Appendix

| Symbol                                                                            | Description                          | Symbol                                                                            | Description                               |
|-----------------------------------------------------------------------------------|--------------------------------------|-----------------------------------------------------------------------------------|-------------------------------------------|
|  | Numeric input field active (example) |  | Alphanumeric input field active (example) |
|  | Inactive input field (example)       |  | Input error (example)                     |

Table 22: Symbols of the configuration bar

### Status bar

| Symbol                                                                              | Description                          | Symbol                                                                              | Description                       |
|-------------------------------------------------------------------------------------|--------------------------------------|-------------------------------------------------------------------------------------|-----------------------------------|
|    | Opens a drop-down menu               |    | Closes a drop-down menu           |
|    | Opens the main screen                |    | Opens the "Settings" menu         |
|  | Lid lock, inactive                   |  | Lid lock, active                  |
|  | Error message                        |  | Warning message                   |
|  | Information                          |  | Sound off function activated      |
|  | USB stick plugged in and detected    |  | Opens the program name input menu |
|  | Administrator logged in <sup>1</sup> |  | User logged in <sup>1</sup>       |
|  | User logout <sup>1</sup>             |                                                                                     |                                   |

Table 23: Symbols on the status bar

### Input area, controls

| Symbol                                                                              | Description                                        | Symbol                                                                               | Description                                                   |
|-------------------------------------------------------------------------------------|----------------------------------------------------|--------------------------------------------------------------------------------------|---------------------------------------------------------------|
|  | "Set value" button (example), opens an input menu  |  | Minimum (min.) and maximum (max.) permissible value (example) |
|  | Adds 5 or 10 units to the selected value (example) |  | Subtracts 5 or 10 units from the selected value (example)     |

| Symbol                                                                              | Description                                                                   | Symbol                                                                              | Description                                                |
|-------------------------------------------------------------------------------------|-------------------------------------------------------------------------------|-------------------------------------------------------------------------------------|------------------------------------------------------------|
|    | Selects a linear curve (L)                                                    |    | Selects a quadratic curve (Q)                              |
|    | Selects an optimised quadratic curve (P) <sup>1</sup>                         |    | Selects the brake-free spin-out option (F)                 |
|    | Navigates up and down through a list. Grey arrows cannot be selected.         |    | Activates one or more list items                           |
|    | List item not selected                                                        |    | Selects a list item. Multiple selections are not possible. |
|    | List item not selected                                                        |    | Option selected                                            |
|    | Option not selected                                                           |    | Opens an additional menu                                   |
|  | Marks list items as favourites. Favourites are placed at the top of the list. |  | Activated function                                         |
|  | Deactivated function                                                          | <b>After 4 h</b>                                                                    | Active value input (example)                               |
| <b>After 4 h</b>                                                                    | Inactive value input (example)                                                |  | Selects a temperature unit                                 |

Table 24: Symbols of the controls in the input area

### Input area, elements of information

| Symbol                                                                              | Description                               | Symbol                                                                              | Description                                 |
|-------------------------------------------------------------------------------------|-------------------------------------------|-------------------------------------------------------------------------------------|---------------------------------------------|
|  | Runtime                                   |  | Runtime, Quality Control function activated |
|  | Start timer                               |  | Runtime at set speed                        |
|  | Deceleration time indication              |  | Speed                                       |
|  | Speed, Quality Control function activated |  | Radius reduction <sup>1</sup>               |
|  | Density increase <sup>1</sup>             |  | RCF                                         |

## Appendix

| Symbol                                                                              | Description                                     | Symbol                                                                              | Description                                               |
|-------------------------------------------------------------------------------------|-------------------------------------------------|-------------------------------------------------------------------------------------|-----------------------------------------------------------|
|    | RCF, Quality Control function activated         |    | Temperature                                               |
|    | Temperature, Quality Control function activated |    | <i>RAPID TEMP</i> rapid temperature control function      |
|    | Standstill cooling                              |    | Rotor library                                             |
|    | Integral <sup>1</sup>                           |    | Integral, Quality Control function activated <sup>1</sup> |
|    | Curve selection                                 |    | Profiles <sup>1</sup>                                     |
|    | Program library                                 |    | Error log                                                 |
|  | Error                                           |  | Warning                                                   |
|  | Settings                                        |  | Device settings                                           |
|  | Automatic lid opening function                  |  | Display deceleration time                                 |
|  | Date and time                                   |  | Time                                                      |
|  | Date                                            |  | Time zone                                                 |
|  | Signals                                         |  | Signal at end of run                                      |
|  | Signal in case of imbalance                     |  | Signal in case of error                                   |
|  | Signal duration                                 |  | Volume                                                    |
|  | Expert mode <sup>1</sup>                        |  | Temperature unit                                          |
|  | Display settings                                |  | Brightness                                                |
|  | Screen lock                                     |  | QR code for Sigma Timer                                   |
|  | ECO mode                                        |  | Language selection                                        |

| Symbol                                                                              | Description                                           | Symbol                                                                              | Description                                                     |
|-------------------------------------------------------------------------------------|-------------------------------------------------------|-------------------------------------------------------------------------------------|-----------------------------------------------------------------|
|    | User management <sup>1</sup>                          |    | Rights management <sup>1</sup>                                  |
|    | User login <sup>1</sup>                               |    | Loading of programs (rights management) <sup>1</sup>            |
|    | Saving of programs (rights management) <sup>1</sup>   |    | Changing of process parameters (rights management) <sup>1</sup> |
|    | Changing of settings (rights management) <sup>1</sup> |    | USB access (rights management) <sup>1</sup>                     |
|    | Admin <sup>1</sup>                                    |    | User <sup>1</sup>                                               |
|    | User logged in <sup>1</sup>                           |    | User name <sup>1</sup>                                          |
|  | Password <sup>1</sup>                                 |  | Info                                                            |
|  | Centrifuge name                                       |  | Service                                                         |
|  | Service code                                          |  | QR code with service information                                |
|  | Operating manual                                      |  | Cycles                                                          |
|  | Cycles current setup                                  |  | Rotor cycles                                                    |
|  | Bucket cycles                                         |  | Rotor                                                           |
|  | Bucket                                                |  | Motor                                                           |
|  | Quality control                                       |  | Critical parameter                                              |
|  | Abnormal parameter                                    |  | System OK                                                       |
|  | Automatic QC data transfer                            |  | QC data transfer                                                |
|  | Quality control overview                              |  | Program export                                                  |

## Appendix

| Symbol                                                                            | Description    | Symbol                                                                            | Description    |
|-----------------------------------------------------------------------------------|----------------|-----------------------------------------------------------------------------------|----------------|
|  | Program import |  | QC data export |

Table 25: Symbols of the elements of information in the input area

## 11.6 Table of chemical resistance data

**i** The information refers to the resistance at 20°C.

The information in the following tables applies to accessories with an anodised coating.

| – no data<br>1 resistant<br>2 practically resistant<br>3 partially resistant<br>4 not resistant |                                                 | Concentration | High-density polyethylene | Polyamide | Polycarbonate | Polyoxymethylene | Polypropylene | Polysulfone | Polyvinyl chloride, hard | Polyvinyl chloride, soft | Polytetrafluorethylene | Acrylonitrile butadiene caoutchouc | Aluminium |
|-------------------------------------------------------------------------------------------------|-------------------------------------------------|---------------|---------------------------|-----------|---------------|------------------|---------------|-------------|--------------------------|--------------------------|------------------------|------------------------------------|-----------|
| Medium                                                                                          | Formula                                         | [%]           | HDPE                      | PA        | PC            | POM              | PP            | PSU         | PVC                      | PVC                      | PTFE                   | NBR                                | AL        |
| Acetaldehyde                                                                                    | C <sub>2</sub> H <sub>4</sub> O                 | 40            | 3                         | 2         | 4             | 2                | 3             | 4           | 4                        | -                        | 1                      | 4                                  | 1         |
| Acetamide                                                                                       | C <sub>2</sub> H <sub>5</sub> NO                | saturated     | 1                         | 1         | 4             | 1                | 1             | 4           | 4                        | -                        | 1                      | -                                  | 1         |
| Acetic acid                                                                                     | C <sub>2</sub> H <sub>4</sub> O <sub>2</sub>    | 10            | 1                         | 4         | 1             | 1                | 1             | 1           | 1                        | 1                        | 1                      | 2                                  | 1         |
| Acetic acid                                                                                     | C <sub>2</sub> H <sub>4</sub> O <sub>2</sub>    | 90            | 1                         | 4         | 4             | 4                | 1             | 3           | 1                        | 4                        | 1                      | -                                  | 1         |
| Acetone                                                                                         | C <sub>3</sub> H <sub>6</sub> O                 | 100           | 1                         | 1         | 4             | 1                | 1             | 4           | 4                        | -                        | 1                      | 4                                  | 1         |
| Acrylonitrile                                                                                   | C <sub>3</sub> H <sub>3</sub> N                 | 100           | 1                         | 1         | 4             | 3                | 3             | 4           | 4                        | 4                        | 1                      | 4                                  | 1         |
| Allyl alcohol                                                                                   | C <sub>3</sub> H <sub>6</sub> O                 | 96            | 1                         | 3         | 3             | 2                | 2             | 2           | 2                        | 4                        | 1                      | 1                                  | 1         |
| Aluminium chloride                                                                              | AlCl <sub>3</sub>                               | saturated     | 1                         | 3         | 2             | 4                | 1             | -           | 1                        | -                        | 1                      | 1                                  | 4         |
| Aluminium sulphate                                                                              | Al <sub>2</sub> (SO <sub>4</sub> ) <sub>3</sub> | 10            | 1                         | 1         | 1             | 3                | 1             | 1           | 1                        | 1                        | 1                      | 1                                  | 1         |
| Ammonium chloride                                                                               | (NH <sub>4</sub> )Cl                            | aqueous       | 1                         | 1         | 1             | 2                | 1             | 1           | 1                        | 1                        | 1                      | 1                                  | 3         |
| Ammonium hydroxide                                                                              | NH <sub>3</sub> +H <sub>2</sub> O               | 30            | 1                         | 3         | 4             | 1                | 1             | 2           | 1                        | -                        | 1                      | -                                  | 1         |
| Aniline                                                                                         | C <sub>6</sub> H <sub>7</sub> N                 | 100           | 1                         | 3         | 4             | 1                | 2             | 4           | 4                        | 4                        | 1                      | 4                                  | 1         |
| Antimony trichloride                                                                            | SbCl <sub>3</sub>                               | 90            | 1                         | 4         | 1             | 4                | 1             | -           | 1                        | -                        | 1                      | -                                  | 4         |
| Benzaldehyde                                                                                    | C <sub>7</sub> H <sub>6</sub> O                 | 100           | 1                         | 3         | 4             | 1                | 1             | 3           | 4                        | 4                        | 1                      | 4                                  | 1         |
| Benzene                                                                                         | C <sub>6</sub> H <sub>6</sub>                   | 100           | 3                         | 2         | 4             | 1                | 3             | 4           | 4                        | -                        | 1                      | 4                                  | 1         |
| Boric acid                                                                                      | H <sub>3</sub> BO <sub>3</sub>                  | aqueous       | 1                         | 3         | 1             | 2                | 1             | -           | -                        | -                        | 1                      | 1                                  | 1         |
| Butyl acrylate                                                                                  | C <sub>7</sub> H <sub>12</sub> O <sub>2</sub>   | 100           | 1                         | 2         | 4             | 2                | 3             | 4           | 4                        | 4                        | 1                      | -                                  | 1         |
| Butyl alcohol                                                                                   | C <sub>4</sub> H <sub>10</sub> O                | 100           | 1                         | 1         | 2             | 1                | 1             | 2           | 2                        | 4                        | 1                      | 1                                  | 1         |
| Calcium chloride                                                                                | CaCl <sub>2</sub>                               | alcoholic     | 1                         | 4         | 2             | 3                | 1             | -           | -                        | 4                        | 1                      | 1                                  | 3         |
| Carbon disulfide                                                                                | CS <sub>2</sub>                                 | 100           | 4                         | 3         | 4             | 2                | 4             | 4           | 4                        | 4                        | 1                      | 3                                  | 1         |

**Appendix**

| – no data<br>1 resistant<br>2 practically resistant<br>3 partially resistant<br>4 not resistant |                                                          |               |                           |           |               |                  |               |             |                          |                          |                        |                                    |           |
|-------------------------------------------------------------------------------------------------|----------------------------------------------------------|---------------|---------------------------|-----------|---------------|------------------|---------------|-------------|--------------------------|--------------------------|------------------------|------------------------------------|-----------|
|                                                                                                 |                                                          | Concentration | High-density polyethylene | Polyamide | Polycarbonate | Polyoxymethylene | Polypropylene | Polysulfone | Polyvinyl chloride, hard | Polyvinyl chloride, soft | Polytetrafluorethylene | Acrylonitrile butadiene caoutchouc | Aluminium |
| Medium                                                                                          | Formula                                                  | [%]           | HDPE                      | PA        | PC            | POM              | PP            | PSU         | PVC                      | PVC                      | PTFE                   | NBR                                | AL        |
| Carbon tetrachloride (TETRA)                                                                    | CCl <sub>4</sub>                                         | 100           | 4                         | 4         | 4             | 2                | 4             | 4           | 4                        | 4                        | 1                      | 3                                  | 1         |
| Chlorine                                                                                        | Cl <sub>2</sub>                                          | 100           | 4                         | 4         | 4             | 4                | 4             | 4           | 4                        | 4                        | 1                      | -                                  | 3         |
| Chlorine water                                                                                  | Cl <sub>2</sub> x H <sub>2</sub> O                       |               | 3                         | 4         | 4             | 4                | 3             | -           | 3                        | 3                        | 1                      | -                                  | 4         |
| Chlorobenzene                                                                                   | C <sub>6</sub> H <sub>5</sub> Cl                         | 100           | 3                         | 4         | 4             | 1                | 3             | 4           | 4                        | 4                        | 1                      | 4                                  | 1         |
| Chloroform                                                                                      | CHCl <sub>3</sub>                                        | 100           | 3                         | 3         | 4             | 4                | 3             | 4           | 4                        | 4                        | 1                      | 4                                  | 3         |
| Chromic acid                                                                                    | CrO <sub>3</sub>                                         | 10            | 1                         | 4         | 2             | 4                | 1             | 4           | 1                        | -                        | 1                      | 4                                  | 1         |
| Chromic potassium sulphate                                                                      | KCr(SO <sub>4</sub> ) <sub>2</sub> x 12 H <sub>2</sub> O | saturated     | 1                         | 2         | 1             | 3                | 1             | -           | 1                        | -                        | 1                      | -                                  | 3         |
| Citric acid                                                                                     | C <sub>6</sub> H <sub>8</sub> O <sub>7</sub>             | 10            | 1                         | 1         | 1             | 2                | 1             | 1           | 1                        | 1                        | 1                      | 1                                  | 1         |
| Citric acid                                                                                     | C <sub>6</sub> H <sub>8</sub> O <sub>7</sub>             | 50            | 1                         | 3         | 1             | 2                | 1             | -           | -                        | -                        | 1                      | 1                                  | 1         |
| Copper sulphate                                                                                 | CuSO <sub>4</sub> x 5H <sub>2</sub> O                    | 10            | 1                         | 1         | 1             | 1                | 1             | 1           | 1                        | 1                        | 1                      | 1                                  | 4         |
| Cyclohexanol                                                                                    | C <sub>6</sub> H <sub>12</sub> O                         | 100           | 1                         | 1         | 3             | 1                | 1             | 1           | 1                        | 4                        | 1                      | 2                                  | 1         |
| Decane                                                                                          | C <sub>10</sub> H <sub>22</sub>                          | 100           | -                         | 1         | 2             | 1                | 3             | -           | -                        | -                        | 1                      | 2                                  | 1         |
| Diaminoethane                                                                                   | C <sub>2</sub> H <sub>8</sub> N <sub>2</sub>             | 100           | 1                         | 1         | 3             | 1                | 1             | -           | 3                        | 4                        | 1                      | 1                                  | 1         |
| Diesel fuel                                                                                     | –                                                        | 100           | 1                         | 1         | 3             | 1                | 1             | -           | 1                        | 3                        | 1                      | 1                                  | 1         |
| Dimethyl formamide (DMF)                                                                        | C <sub>3</sub> D <sub>7</sub> NO                         | 100           | 1                         | 1         | 4             | 1                | 1             | 4           | 3                        | -                        | 1                      | 3                                  | 1         |
| Dimethyl sulfoxide (DMSO)                                                                       | C <sub>2</sub> H <sub>6</sub> SO                         | 100           | 1                         | 2         | 4             | 1                | 1             | 4           | 4                        | -                        | 1                      | -                                  | 1         |
| Dimethylaniline                                                                                 | C <sub>8</sub> H <sub>11</sub> N                         | 100           | -                         | 3         | 4             | 2                | 4             | -           | -                        | -                        | 1                      | -                                  | 1         |
| Dioxane                                                                                         | C <sub>4</sub> H <sub>8</sub> O <sub>2</sub>             | 100           | 2                         | 1         | 4             | 1                | 3             | 2           | 3                        | 4                        | 1                      | 3                                  | 1         |
| Dipropylene glycol monomethyl ether                                                             | C <sub>4</sub> H <sub>10</sub> O                         | 100           | 3                         | 1         | 4             | 1                | 4             | 4           | 4                        | 4                        | 1                      | -                                  | 1         |
| Ethanol                                                                                         | C <sub>2</sub> H <sub>6</sub> O                          | 96            | 1                         | 1         | 1             | 1                | 1             | 1           | 1                        | 3                        | 1                      | -                                  | 1         |
| Ethyl acetate                                                                                   | C <sub>4</sub> H <sub>8</sub> O <sub>2</sub>             | 100           | 1                         | 1         | 4             | 1                | 1             | 4           | 4                        | 4                        | 1                      | 4                                  | 1         |

| – no data<br>1 resistant<br>2 practically resistant<br>3 partially resistant<br>4 not resistant |                                                                  |               |                           |           |               |                  |               |             |                          |                          |                        |                                    |           |
|-------------------------------------------------------------------------------------------------|------------------------------------------------------------------|---------------|---------------------------|-----------|---------------|------------------|---------------|-------------|--------------------------|--------------------------|------------------------|------------------------------------|-----------|
|                                                                                                 |                                                                  | Concentration | High-density polyethylene | Polyamide | Polycarbonate | Polyoxymethylene | Polypropylene | Polysulfone | Polyvinyl chloride, hard | Polyvinyl chloride, soft | Polytetrafluorethylene | Acrylonitrile butadiene caoutchouc | Aluminium |
| Medium                                                                                          | Formula                                                          | [%]           | HDPE                      | PA        | PC            | POM              | PP            | PSU         | PVC                      | PVC                      | PTFE                   | NBR                                | AL        |
| Ethylene chloride                                                                               | C <sub>2</sub> H <sub>4</sub> Cl <sub>2</sub>                    | 100           | 3                         | 3         | 4             | 1                | 3             | 4           | 4                        | 4                        | 1                      | -                                  | 1         |
| Ferrous chloride                                                                                | FeCl <sub>2</sub>                                                | saturated     | 1                         | 3         | 1             | 3                | 1             | 1           | 1                        | 1                        | 1                      | -                                  | 4         |
| Formaldehyde solution                                                                           | CH <sub>2</sub> O                                                | 30            | 1                         | 3         | 1             | 1                | 1             | -           | -                        | -                        | 1                      | 2                                  | 1         |
| Formic acid                                                                                     | CH <sub>2</sub> O <sub>2</sub>                                   | 100           | 1                         | 4         | 3             | 4                | 1             | 3           | 3                        | 1                        | 1                      | 2                                  | 1         |
| Furfural                                                                                        | C <sub>5</sub> H <sub>4</sub> O <sub>2</sub>                     | 100           | 1                         | 3         | 3             | 2                | 4             | -           | -                        | -                        | 1                      | 4                                  | 1         |
| Gasoline                                                                                        | C <sub>5</sub> H <sub>12</sub> - C <sub>12</sub> H <sub>26</sub> | 100           | 2                         | 1         | 3             | 1                | 3             | 3           | 2                        | -                        | 1                      | 1                                  | 1         |
| Glycerol                                                                                        | C <sub>3</sub> H <sub>8</sub> O <sub>3</sub>                     | 100           | 1                         | 1         | 3             | 1                | 1             | 1           | 1                        | 2                        | 1                      | 1                                  | 1         |
| Heptane, n-                                                                                     | C <sub>7</sub> H <sub>16</sub>                                   | 100           | 2                         | 1         | 1             | 1                | 2             | 1           | 2                        | 4                        | 1                      | 1                                  | 1         |
| Hexane, n-                                                                                      | C <sub>6</sub> H <sub>14</sub>                                   | 100           | 2                         | 1         | 2             | 1                | 2             | 1           | 2                        | 4                        | 1                      | 1                                  | 1         |
| Hydrochloric acid                                                                               | HCl                                                              | 5             | 1                         | 4         | 1             | 4                | 1             | 1           | 1                        | -                        | 1                      | 2                                  | 4         |
| Hydrochloric acid                                                                               | HCl                                                              | concentrated  | 1                         | 4         | 4             | 4                | 1             | 1           | 2                        | 3                        | 1                      | 4                                  | 4         |
| Hydrogen peroxide                                                                               | H <sub>2</sub> O <sub>2</sub>                                    | 3             | 1                         | 3         | 1             | 1                | 1             | 1           | 1                        | -                        | 1                      | 3                                  | 3         |
| Hydrogen peroxide                                                                               | H <sub>2</sub> O <sub>2</sub>                                    | 30            | 1                         | 4         | 1             | 4                | 1             | 1           | 1                        | -                        | 1                      | 3                                  | 3         |
| Hydrogen sulphide                                                                               | H <sub>2</sub> S                                                 | 10            | 1                         | 1         | 1             | 1                | 1             | 1           | 1                        | 3                        | 1                      | 3                                  | 1         |
| Iodine, tincture of                                                                             | I <sub>2</sub>                                                   |               | 1                         | 4         | 3             | 1                | 1             | -           | 4                        | 4                        | 1                      | 1                                  | 1         |
| Isopropyl alcohol                                                                               | C <sub>3</sub> H <sub>8</sub> O                                  | 100           | 1                         | 1         | 1             | 1                | 1             | 1           | 1                        | 4                        | 1                      | -                                  | 2         |
| Lactic acid                                                                                     | C <sub>3</sub> H <sub>6</sub> O <sub>3</sub>                     | 3             | 1                         | 3         | 1             | 2                | 1             | 1           | 2                        | -                        | 1                      | 1                                  | 1         |
| Magnesium chloride                                                                              | MgCl <sub>2</sub>                                                | 10            | 1                         | 1         | 1             | 1                | 1             | 1           | 1                        | 1                        | 1                      | 1                                  | 1         |
| Mercuric chloride                                                                               | HgCl <sub>2</sub>                                                | 10            | 1                         | 4         | 1             | 3                | 1             | 1           | 1                        | 1                        | 1                      | 1                                  | 4         |
| Mercury                                                                                         | Hg                                                               | 100           | 1                         | 1         | 1             | 1                | 1             | 1           | 1                        | 3                        | 1                      | 1                                  | 3         |
| Methanol                                                                                        | CH <sub>4</sub> O                                                | 100           | 1                         | 2         | 4             | 1                | 1             | 3           | 1                        | 3                        | 1                      | 2                                  | 1         |
| Methyl acetate                                                                                  | C <sub>3</sub> H <sub>6</sub> O <sub>2</sub>                     | 100           | 1                         | 1         | 4             | 2                | 1             | -           | 4                        | 4                        | 1                      | -                                  | 1         |
| Methyl benzene                                                                                  | C <sub>7</sub> H <sub>8</sub> O                                  | 100           | 3                         | 4         | 4             | 1                | 4             | 4           | 2                        | -                        | 1                      | 4                                  | 1         |

## Appendix

| – no data<br>1 resistant<br>2 practically resistant<br>3 partially resistant<br>4 not resistant |                                                                  |               |                           |           |               |                  |               |             |                          |                          |                        |                                    |           |
|-------------------------------------------------------------------------------------------------|------------------------------------------------------------------|---------------|---------------------------|-----------|---------------|------------------|---------------|-------------|--------------------------|--------------------------|------------------------|------------------------------------|-----------|
|                                                                                                 |                                                                  | Concentration | High-density polyethylene | Polyamide | Polycarbonate | Polyoxymethylene | Polypropylene | Polysulfone | Polyvinyl chloride, hard | Polyvinyl chloride, soft | Polytetrafluorethylene | Acrylonitrile butadiene caoutchouc | Aluminium |
| Medium                                                                                          | Formula                                                          | [%]           | HDPE                      | PA        | PC            | POM              | PP            | PSU         | PVC                      | PVC                      | PTFE                   | NBR                                | AL        |
| Methyl ethyl ketone (MEK)                                                                       | C <sub>4</sub> H <sub>8</sub> O                                  | 100           | 1                         | 1         | 4             | 1                | 1             | 4           | 4                        | 4                        | 1                      | 4                                  | 1         |
| Methylene chloride                                                                              | CH <sub>2</sub> Cl <sub>2</sub>                                  | 100           | 4                         | 3         | 4             | 3                | 3             | 4           | 4                        | 4                        | 1                      | -                                  | 1         |
| Mineral oil                                                                                     | –                                                                | 100           | 1                         | 1         | 1             | 1                | 1             | 1           | 1                        | -                        | 1                      | 1                                  | 1         |
| Nitric acid                                                                                     | HNO <sub>3</sub>                                                 | 10            | 1                         | 4         | 1             | 4                | 1             | 1           | 1                        | -                        | 1                      | 4                                  | 3         |
| Nitric acid                                                                                     | HNO <sub>3</sub>                                                 | 100           | 4                         | 4         | 4             | 4                | 4             | -           | 4                        | -                        | 1                      | 4                                  | 1         |
| Nitrobenzene                                                                                    | C <sub>6</sub> H <sub>5</sub> NO <sub>2</sub>                    | 100           | 3                         | 4         | 4             | 3                | 2             | 4           | 4                        | 4                        | 1                      | 4                                  | 1         |
| Oleic acid                                                                                      | C <sub>18</sub> H <sub>34</sub> O <sub>2</sub>                   | 100           | 1                         | 1         | 1             | 2                | 1             | -           | 1                        | -                        | 1                      | 3                                  | 1         |
| Oxalic acid                                                                                     | C <sub>2</sub> H <sub>2</sub> O <sub>4</sub> x 2H <sub>2</sub> O | 100           | 1                         | 3         | 1             | 4                | 1             | 1           | 1                        | 1                        | 1                      | 2                                  | 1         |
| Ozone                                                                                           | O <sub>3</sub>                                                   | 100           | 3                         | 4         | 1             | 4                | 3             | 1           | 1                        | -                        | 1                      | 4                                  | 2         |
| Petroleum                                                                                       | –                                                                | 100           | 1                         | 1         | 3             | 1                | 1             | 1           | 1                        | 3                        | 1                      | 1                                  | 1         |
| Phenol                                                                                          | C <sub>6</sub> H <sub>6</sub> O                                  | 10            | 1                         | 4         | 4             | 4                | 1             | 4           | 1                        | 3                        | 1                      | 3                                  | 1         |
| Phenol                                                                                          | C <sub>6</sub> H <sub>6</sub> O                                  | 100           | 2                         | 4         | 4             | 4                | 1             | 3           | 4                        | 3                        | 1                      | 3                                  | 1         |
| Phosphoric acid                                                                                 | H <sub>3</sub> PO <sub>4</sub>                                   | 20            | 1                         | 4         | 2             | 4                | 1             | -           | -                        | -                        | 1                      | 2                                  | 4         |
| Phosphorus pentachloride                                                                        | PCl <sub>5</sub>                                                 | 100           | -                         | 4         | 4             | 4                | 1             | -           | 4                        | 4                        | 1                      | -                                  | 1         |
| Potassium hydrogen carbonate                                                                    | CHKO <sub>3</sub>                                                | saturated     | 1                         | 1         | 2             | 1                | 1             | -           | -                        | -                        | 1                      | -                                  | 4         |
| Potassium hydroxide                                                                             | KOH                                                              | 30            | 1                         | 1         | 4             | 3                | 1             | 1           | 1                        | 1                        | 1                      | -                                  | 4         |
| Potassium hydroxide                                                                             | KOH                                                              | 50            | 1                         | 1         | 4             | 3                | 1             | 1           | 1                        | 1                        | 1                      | -                                  | 4         |
| Potassium nitrate                                                                               | KNO <sub>3</sub>                                                 | 10            | 1                         | 1         | 1             | 1                | 1             |             | -                        | -                        | 1                      | 1                                  | 1         |
| Potassium permanganate                                                                          | KMnO <sub>4</sub>                                                | 100           | 1                         | 4         | 1             | 1                | 1             | -           | 1                        | -                        | 1                      | 3                                  | 1         |
| Pyridine                                                                                        | C <sub>5</sub> H <sub>5</sub> N                                  | 100           | 1                         | 1         | 4             | 1                | 3             | 4           | 4                        | 4                        | 1                      | 4                                  | 1         |
| Resorcinol                                                                                      | C <sub>6</sub> H <sub>6</sub> O <sub>2</sub>                     | 5             | 1                         | 4         | 2             | 3                | 1             | 4           | 2                        | -                        | 1                      | -                                  | 2         |
| Silver nitrate                                                                                  | AgNO <sub>3</sub>                                                | 100           | 1                         | 1         | 1             | 1                | 1             | 1           | 1                        | 1                        | 1                      | 2                                  | 4         |

| – no data<br>1 resistant<br>2 practically resistant<br>3 partially resistant<br>4 not resistant |                                               |               |                           |           |               |                  |               |             |                          |                          |                        |                                    |           |
|-------------------------------------------------------------------------------------------------|-----------------------------------------------|---------------|---------------------------|-----------|---------------|------------------|---------------|-------------|--------------------------|--------------------------|------------------------|------------------------------------|-----------|
|                                                                                                 |                                               | Concentration | High-density polyethylene | Polyamide | Polycarbonate | Polyoxymethylene | Polypropylene | Polysulfone | Polyvinyl chloride, hard | Polyvinyl chloride, soft | Polytetrafluorethylene | Acrylonitrile butadiene caoutchouc | Aluminium |
| Medium                                                                                          | Formula                                       | [%]           | HDPE                      | PA        | PC            | POM              | PP            | PSU         | PVC                      | PVC                      | PTFE                   | NBR                                | AL        |
| Sodium bisulphite                                                                               | NaHSO <sub>3</sub>                            | 10            | 1                         | 1         | 2             | 4                | 1             | -           | -                        | -                        | 1                      | 1                                  | 1         |
| Sodium carbonate                                                                                | Na <sub>2</sub> CO <sub>3</sub>               | 10            | 1                         | 1         | 1             | 1                | 1             | -           | -                        | -                        | 1                      | -                                  | 3         |
| Sodium chloride                                                                                 | NaCl                                          | 30            | 1                         | 1         | 1             | 1                | 1             | 1           | 1                        | 1                        | 1                      | 1                                  | 3         |
| Sodium hydroxide                                                                                | NaOH                                          | 30            | 1                         | 1         | 4             | 1                | 1             | 1           | 1                        | 1                        | 1                      | 2                                  | 4         |
| Sodium hydroxide                                                                                | NaOH                                          | 50            | 1                         | 1         | 4             | 1                | 1             | 1           | 1                        | -                        | 1                      | 2                                  | 4         |
| Sodium sulphate                                                                                 | Na <sub>2</sub> SO <sub>4</sub>               | 10            | 1                         | 1         | 1             | 1                | 1             | 1           | 1                        | 1                        | 1                      | 1                                  | 1         |
| Styrene                                                                                         | C <sub>8</sub> H <sub>8</sub>                 | 100           | 4                         | 1         | 4             | 1                | 3             | -           | 4                        | 4                        | 1                      | 4                                  | 1         |
| Sulphuric acid                                                                                  | H <sub>2</sub> SO <sub>4</sub>                | 6             | 1                         | 4         | 1             | 4                | 1             | 1           | 1                        | -                        | 1                      | 2                                  | 3         |
| Sulphuric acid                                                                                  | H <sub>2</sub> SO <sub>4</sub>                | fuming        | 4                         | 4         | 4             | 4                | 4             | 4           | 4                        | 4                        | 1                      | 4                                  | 3         |
| Tallow                                                                                          | –                                             | 100           | 1                         | 1         | 1             | 1                | 1             | 1           | 1                        | 1                        | 1                      | 1                                  | 1         |
| Tetrahydrofuran (THF)                                                                           | C <sub>4</sub> H <sub>8</sub> O               | 100           | 3                         | 1         | 4             | 1                | 3             | 4           | 4                        | 4                        | 1                      | 3                                  | 1         |
| Tetrahydronaphthalene                                                                           | C <sub>10</sub> H <sub>12</sub>               | 100           | 3                         | 1         | 4             | 1                | 4             | 4           | 4                        | 4                        | 1                      | -                                  | 1         |
| Thionyl chloride                                                                                | Cl <sub>2</sub> SO                            | 100           | 4                         | 4         | 4             | 2                | 4             | 4           | 4                        | 4                        | 1                      | -                                  | 3         |
| Tin chloride                                                                                    | SnCl <sub>2</sub>                             | 10            | 1                         | 4         | 2             | 2                | 1             | -           | -                        | -                        | 1                      | 1                                  | 4         |
| Toluene                                                                                         | C <sub>7</sub> H <sub>8</sub>                 | 100           | 3                         | 1         | 4             | 1                | 3             | 4           | 4                        | 4                        | 1                      | 4                                  | 1         |
| Transformer oil                                                                                 | –                                             | 100           | 1                         | 1         | 3             | 3                | 1             | 1           | 1                        | -                        | 1                      | 1                                  | 1         |
| Trichloroethane                                                                                 | C <sub>2</sub> H <sub>3</sub> Cl <sub>3</sub> | 100           | 3                         | 3         | 4             | 2                | 4             | 4           | 4                        | 4                        | 1                      | 4                                  | 4         |
| Urea                                                                                            | CH <sub>4</sub> N <sub>2</sub> O              | 10            | 1                         | 1         | 1             | 1                | 1             | -           | -                        | -                        | 1                      | 1                                  | 1         |
| Urine                                                                                           | –                                             | 100           | 1                         | 1         | 1             | 1                | 1             | -           | 1                        | 1                        | 1                      | -                                  | 2         |
| Wax                                                                                             | –                                             | 100           | -                         | 1         | 1             | -                | 1             | -           | -                        | -                        | 1                      | -                                  | 1         |
| Wines                                                                                           | –                                             | 100           | 1                         | 1         | 1             | 2                | 1             | 1           | 1                        | 1                        | 1                      | -                                  | 4         |
| Xylene                                                                                          | C <sub>8</sub> H <sub>10</sub>                | 100           | 3                         | 1         | 4             | 1                | 4             | 4           | 4                        | 4                        | 1                      | 4                                  | 1         |

Table 26: Table of chemical resistance data



## 11.7 Declarations of conformity



### EC Declaration of Conformity

**Product designation** Laboratory centrifuge  
**Product name** Sigma 3-18KT IVD  
**Part numbers** 100019; 100021; 100022  
**Manufacturer** Sigma Laborzentrifugen GmbH  
An der Unteren Söse 50  
37520 Osterode  
Germany

**Basic UDI as referred to in Part C of Annex VI** 426073439IVD01001JQCJ4  
**Single Registration Number (SRN)** DE-MF-000009414

As the manufacturer of the unit(s), we assume full responsibility and hereby declare that the product(s) mentioned hereinabove comply with the requirements as set out in the following regulation(s)/directive(s).

**Regulations** Regulation on in vitro diagnostica (EU) 2017/746  
Implementation Regulation for Regulation (EU) 2024/573 (F-gases Regulation) (EU) 2024/2729

**Directives** Machinery Directive 2006/42/EC  
Low Voltage Directive 2014/35/EU  
EMC Directive 2014/30/EU  
RoHS Directive 2011/65/EU  
Delegated RoHS Directive (EU) 2015/863

**Norms** EN 61010-2-020:2017  
EN IEC 61010-2-011:2021  
EN IEC 61326-1:2021  
EN IEC 61010-2-101:2023

**Risk class in accordance with Annex VIII** class A

Osterode, 2026-06-10

*Michael Sander*

Dr. Michael Sander, Managing Director

CE\_Sigma\_3-18KT\_IVD\_2026-06-10\_en

1 / 1



## ***EC Declaration of Conformity***

in accordance with the EC Machinery Directive 2006/42/EC, annex II, part 1, section A

|                            |                                                                                    |
|----------------------------|------------------------------------------------------------------------------------|
| <b>Product designation</b> | Laboratory centrifuge                                                              |
| <b>Product name</b>        | Sigma 3-18KT                                                                       |
| <b>Order number</b>        | 100014; 100016                                                                     |
| <b>Manufacturer</b>        | Sigma Laborzentrifugen GmbH<br>An der Unteren Söse 50<br>37520 Osterode<br>Germany |

|                                                 |                    |
|-------------------------------------------------|--------------------|
| <b>Authorised representative for CE matters</b> | Alexander Hörsting |
|-------------------------------------------------|--------------------|

The product named hereinabove was developed, designed, and manufactured in compliance with the relevant, fundamental safety and health requirements of the listed EC regulations, directives and norms. In the event of modifications that were not authorised by the manufacturer or if the product is used in a manner that is not in line with the intended purpose, this declaration will be rendered void.

|                    |                                                                             |                                                                       |
|--------------------|-----------------------------------------------------------------------------|-----------------------------------------------------------------------|
| <b>Regulations</b> | Implementation Regulation for Regulation (EU) 2024/573 (F-gases Regulation) | (EU) 2024/2729                                                        |
| <b>Directives</b>  | Machinery Directive                                                         | 2006/42/EC                                                            |
|                    | Low Voltage Directive                                                       | 2014/35/EU                                                            |
|                    | EMC Directive                                                               | 2014/30/EU                                                            |
|                    | RoHS Directive                                                              | 2011/65/EU                                                            |
|                    | Delegated RoHS Directive                                                    | (EU) 2015/863                                                         |
| <b>Norms</b>       |                                                                             | EN 61010-2-020:2017<br>EN IEC 61010-2-011:2021<br>EN IEC 61326-1:2021 |

Osterode, 2025-11-04



Dr. Michael Sander, Managing Director

CE\_Sigma\_3-18KT\_MRL\_2025-11-04\_en

1 / 1



## Declaration of Conformity

### China RoHS 2 (Administrative Measures for the Restriction of the Use of Hazardous Substances in Electrical and Electronic Products)

Laboratory centrifuge models: Sigma 1-14, 1-14K, 1-7, 1-16, 1-16K, 2-7, 2-16KL, 2-16KHL, 3-16L, 3-16KL, 3-16T, 3-17KT, 3-18KT, 3-18KS, 3-18KHS, 3-30KS, 3-30KHS, 4-5L, 4-5KL, 4-5KRL, 4-16S, 4-16KS, 4-16KHS, 6-5KRL, 6-16S, 6-16HS, 6-16KS, 6-16KHS, 8KS, 8KBS

Sigma Laborzentrifugen GmbH has made reasonable effort to avoid the use of hazardous substances in the products it manufactures (laboratory centrifuges).

A Product Conformity Assessment (PCA) was performed in order to determine whether the concentration of harmful substances in all homogeneous materials of the component parts is above or below the MCV limit (Maximum Concentration Value limit) as defined in GB/T 26572:

Mercury and its compounds: 0.1 %

Cadmium (Cd) and its compounds: 0.01 %

Lead (Pb) and its compounds: 0.1 %

Hexavalent chromium (Cr (VI)) and its compounds: 0.1 %

Polybrominated biphenyls (PBB): 0.1 %

Polybrominated diphenyl ethers (PBDE): 0.1 %

表1 产品中有害物质的名称及含量

Table 1: Name and content of hazardous substances in the product

| 部件名称<br>Component part (PCA) | 有害物质<br>Hazardous substance |                   |                   |                                      |                                         |                                                 |
|------------------------------|-----------------------------|-------------------|-------------------|--------------------------------------|-----------------------------------------|-------------------------------------------------|
|                              | 铅<br>Lead (Pb)              | 汞<br>Mercury (Hg) | 镉<br>Cadmium (Cd) | 六价铬<br>Hexavalent Chromium (Cr (VI)) | 多溴联苯<br>Poly-brominated biphenyls (PBB) | 多溴二苯醚<br>Poly-brominated diphenyl ethers (PBDE) |
| Electronic PCB, cables       | X <sup>1)</sup>             | O                 | X <sup>1)</sup>   | O                                    | O                                       | O                                               |
| Display                      | O                           | O                 | O                 | O                                    | O                                       | O                                               |
| Housing                      | X <sup>1)</sup>             | O                 | O                 | O                                    | O                                       | O                                               |
| Base, metal, accessories     | X <sup>1)</sup>             | O                 | O                 | O                                    | O                                       | O                                               |

本表格依据SJ/T 11364的规定编制。

This table is made according to SJ/T 11364.

Appendix



|    |                                                                                                                                                                                                                                                                                                                                                    |
|----|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| O: | 表示该有害物质在该部件所有均质材料中的含量均在GB/T 26572规定的限量要求以下。<br>Indicates that the content of the harmful substance in all homogeneous materials of the component part is below the limit as defined in GB/T 26572.)                                                                                                                                                |
| X: | 表示该有害物质至少在该部件的某一均质材料中的含量超出GB/T 26572规定的限量要求。<br>(企业可在此处, 根据实际情况对上表打“X”的技术原因进行进一步说明。)<br>Indicates that the content of the harmful substance in at least one homogeneous material of the component part exceeds the limit as defined in GB/T 26572. (Contact the manufacturer for further technical information according to the actual situation.) |
| 1) | Contains parts in compliance with exemptions 6a, 6b, 6c, 7c.I, 7c.II and 8b of 2011/65/EU RoHS.                                                                                                                                                                                                                                                    |

Apart from the exemptions given in this table, none of the substances listed above have been intentionally added to the product or metallic coatings.

Sigma Laborzentrifugen GmbH  
An der Unteren Söse 50  
37520 Osterode  
Germany

Osterode, 2026-09-16



M. Weigoni, Director of Procurement

## 12 Index

|                                                           |     |                                                                   |     |
|-----------------------------------------------------------|-----|-------------------------------------------------------------------|-----|
| <b>A</b>                                                  |     |                                                                   |     |
| Acceleration curve selection .....                        | 61  | Technical data .....                                              | 122 |
| Acceleration curves .....                                 | 128 | Continuous run .....                                              | 49  |
| Activating the rights management<br>function .....        | 82  | Copyright .....                                                   | 14  |
| Ambient conditions                                        |     | Cost estimate .....                                               | 120 |
| Technical data .....                                      | 123 | Creating a profile .....                                          | 64  |
| Angle rotor removal .....                                 | 42  | Creating an administrator account                                 | 82  |
| Annual inspection and maintenance<br>.....                | 113 | Curve selection.....                                              | 61  |
| Area of application .....                                 | 20  | Cycles                                                            |     |
| Autoclaving .....                                         | 105 | Control system .....                                              | 92  |
| Automatic QC data transfer .....                          | 94  | <b>D</b>                                                          |     |
| Automatic rotor identification .....                      | 60  | Danger to life caused by electric<br>shock .....                  | 25  |
| <b>B</b>                                                  |     | Data transfer .....                                               | 95  |
| Basic functions                                           |     | Deceleration curve selection .....                                | 61  |
| Spincontrol T .....                                       | 45  | Deceleration curves.....                                          | 128 |
| Batch number and serial number ..                         | 16  | Deceleration time indication:Process-<br>related parameters ..... | 55  |
| Broken glass in the rotor chamber<br>.....                | 107 | Declaration of conformity.....                                    | 143 |
| <b>C</b>                                                  |     | Declaration of conformity – China<br>RoHS2 .....                  | 145 |
| Centrifugation                                            |     | Defining intervals .....                                          | 64  |
| Principle .....                                           | 29  | Deleting a program .....                                          | 67  |
| Centrifugation start .....                                | 45  | Density .....                                                     | 127 |
| Centrifugation stop .....                                 | 46  | Density increase .....                                            | 59  |
| Centrifuge overview .....                                 | 29  | Device settings .....                                             | 71  |
| Certified refrigeration specialists ...                   | 24  | Different service life (marking).....                             | 18  |
| Chemical resistance table .....                           | 137 | Directives .....                                                  | 14  |
| China RoHS2 declaration of<br>conformity .....            | 145 | Disinfectant for the centrifuge ....                              | 103 |
| Cleaning agent .....                                      | 102 | Disinfecting the adapters .....                                   | 104 |
| Cleaning tasks .....                                      | 102 | Disinfecting the buckets .....                                    | 104 |
| Cleaning the adapters .....                               | 103 | Disinfecting the carriers.....                                    | 104 |
| Cleaning the buckets .....                                | 102 | Disinfecting the centrifuge .....                                 | 103 |
| Cleaning the carriers .....                               | 103 | Disinfecting the rotors.....                                      | 104 |
| Cleaning the centrifuge .....                             | 102 | Disinfection                                                      |     |
| Cleaning the condenser                                    |     | Adapters and carriers .....                                       | 104 |
| with an air-cooled refrigeration<br>system .....          | 107 | Centrifuge .....                                                  | 103 |
| Cleaning the rotors .....                                 | 102 | Rotors and buckets .....                                          | 104 |
| Closing the lid .....                                     | 37  | Disposal of the centrifuge .....                                  | 121 |
| Completion of the declaration of<br>decontamination ..... | 119 | Disposal of the packaging material<br>.....                       | 121 |
| Completion of the return declaration<br>.....             | 119 | Dotors with Sigma Comfort coating<br>Marking.....                 | 19  |
| Connection data                                           |     | <b>E</b>                                                          |     |
|                                                           |     | EC declaration of conformity .....                                | 143 |
|                                                           |     | ECO mode .....                                                    | 79  |

## Index

|                                                                   |     |                                                                     |     |
|-------------------------------------------------------------------|-----|---------------------------------------------------------------------|-----|
| Editing a profile .....                                           | 64  | und chemische Industrie FBRCI-025 .....                             | 113 |
| Emergency lid release .....                                       | 118 | In-house transport .....                                            | 33  |
| Error log .....                                                   | 69  | Inspection to be performed by the manufacturer .....                | 113 |
| Error log entries, deletion .....                                 | 69  | Installing a rotor cross .....                                      | 40  |
| Error log, deletion .....                                         | 69  | Installing a swing-out rotor .....                                  | 40  |
| Error messages                                                    |     | Installing an angle rotor .....                                     | 40  |
| Identifying and eliminating faults and malfunctions .....         | 114 | Installing the buckets .....                                        | 40  |
| Executing a program .....                                         | 66  | Instructions for a specific course of action (representation) ..... | 13  |
| Expert mode .....                                                 | 76  | Integral .....                                                      | 58  |
| Expert Mode: Process-related parameters .....                     | 58  | Integrated warning notes .....                                      | 12  |
| Expert Mode: Density increase .....                               | 59  | Intended use                                                        |     |
| Expert mode: Integral .....                                       | 58  | Centrifuge .....                                                    | 20  |
| Expert Mode: Radius reduction .....                               | 58  | Rotors, buckets, carriers and accessories .....                     | 20  |
| Export of programs .....                                          | 95  | Starting products .....                                             | 20  |
| Export of QC data .....                                           | 97  | <b>L</b>                                                            |     |
| Exposure to direct sunlight .....                                 | 100 | Label (F-gases) .....                                               | 30  |
| <b>F</b>                                                          |     | Language selection .....                                            | 81  |
| FBRCI-025 .....                                                   | 113 | Liability claims .....                                              | 13  |
| F-gases label .....                                               | 30  | Linear curves (LIN) .....                                           | 128 |
| Filling the vessels .....                                         | 37  | Lists (representation) .....                                        | 13  |
| Fire hazard due to electrical overload .....                      | 25  | Loading an angle rotor .....                                        | 38  |
| Further applicable documents .....                                | 11  | Loading the adapters .....                                          | 37  |
| Fuses inside the unit .....                                       | 35  | Loading the buckets .....                                           | 40  |
| <b>H</b>                                                          |     | Location of use                                                     |     |
| Hazards cause by the centrifuge lid .....                         | 26  | Climatic conditions .....                                           | 34  |
| Hazards caused by rotors, buckets, carriers and accessories ..... | 25  | Location requirements .....                                         | 34  |
| Hazards caused by the operating fluids .....                      | 26  | Logistics personnel .....                                           | 24  |
| Hazards caused by the refrigerant .....                           | 26  | <b>M</b>                                                            |     |
| Hazards due to biological substances .....                        | 26  | Main screen .....                                                   | 44  |
| Hazards due to electrical voltage .....                           | 25  | Mains power plug                                                    |     |
| Hazards due to hot surfaces .....                                 | 26  | Safety device .....                                                 | 27  |
| Hazards due to mechanical force .....                             | 25  | Mains power switch                                                  |     |
| <b>I</b>                                                          |     | Safety device .....                                                 | 27  |
| Identifying and eliminating faults and malfunctions .....         | 114 | Maintaining the centrifuge                                          |     |
| Import of programs .....                                          | 95  | Post-centrifugation measures .....                                  | 106 |
| Importance of the operating instructions .....                    | 11  | Post-cleaning measures .....                                        | 106 |
| Info                                                              |     | Visual inspection .....                                             | 106 |
| Control system .....                                              | 89  | Maintaining the load-bearing bolts .....                            | 111 |
| Information sheet of the German Berufsgenossenschaft Rohstoffe    |     | Maintaining the rotors, buckets, adapters and carriers              |     |
|                                                                   |     | Post-centrifugation measures .....                                  | 108 |
|                                                                   |     | Post-cleaning measures .....                                        | 108 |

|                                                           |     |                                                         |            |
|-----------------------------------------------------------|-----|---------------------------------------------------------|------------|
| Service life .....                                        | 108 | Power cord with an IEC C13 connection .....             | 35         |
| Visual inspection.....                                    | 108 | Power off .....                                         | 36         |
| Maintenance .....                                         | 22  | Power on .....                                          | 36         |
| Maintenance contract .....                                | 113 | Power supply .....                                      | 35         |
| Mandatory testing of centrifuges .                        | 113 | Principle of centrifugation .....                       | 29         |
| Marking in case of a different service life .....         | 18  | Process parameters .....                                | 47         |
| Marking of rotors with Sigma Comfort coating .....        | 19  | Process parameters:Runtime .....                        | 48         |
| Marking of the accessories.....                           | 16  | Process parameters:Temperature                          | 50         |
| Marking of the centrifuge .....                           | 15  | Process-related parameters .....                        | 52         |
| Maximum speed for tubes .....                             | 37  | RCF .....                                               | 56         |
| Mechanical specialists .....                              | 24  | Process-related parameters:RPM .                        | 55         |
| Misuse .....                                              | 21  | Process-related parameters:Runtime .....                | 54         |
| <b>N</b>                                                  |     | Process-related parameters:Temperature .....            | 57         |
| Name plate .....                                          | 31  | Product data                                            |            |
| Notices (representation) .....                            | 13  | Technical data .....                                    | 122        |
| <b>O</b>                                                  |     | Program creation .....                                  | 62         |
| On-site fuse protection .....                             | 35  | Program editing .....                                   | 62         |
| Onward transport.....                                     | 33  | Program export.....                                     | 95         |
| Opening the error log                                     |     | Program import.....                                     | 95         |
| Identifying and eliminating faults and malfunctions ..... | 115 | Program library .....                                   | 62         |
| Opening the lid .....                                     | 36  | Protective earth conductor                              |            |
| Opening the lid manually                                  |     | Safety device .....                                     | 27         |
| Emergency lid release .....                               | 118 | Protective equipment.....                               | 24         |
| Operating fluids                                          |     | <b>Q</b>                                                |            |
| Technical data .....                                      | 123 | QC data export .....                                    | 97         |
| Operating personnel .....                                 | 23  | QC data transfer .....                                  | 94         |
| Operator                                                  |     | Quadratic curves (QUAD) .....                           | 129        |
| Definition .....                                          | 21  | Qualification of the personnel.....                     | 23         |
| Operator responsibility.....                              | 21  | Qualified electricians .....                            | 24         |
| Optimised quadratic curves (QUADpro).....                 | 130 | Quality control .....                                   | 93         |
| Other parameters                                          |     | Quality control overview .....                          | 94         |
| Technical data .....                                      | 122 | Quality Control:Process-related parameters .....        | 55, 56, 57 |
| Overview of the centrifuge.....                           | 29  | Quantities                                              |            |
| <b>P</b>                                                  |     | Technical data .....                                    | 123        |
| Packaging symbols.....                                    | 32  | Quick stop .....                                        | 46         |
| Performance parameters                                    |     | <b>R</b>                                                |            |
| Technical data .....                                      | 122 | Radius reduction .....                                  | 58         |
| Personal protective equipment.....                        | 24  | Radius, speed and relative centrifugal force (RCF)..... | 127        |
| Personnel qualification .....                             | 23  | RCF: Process-related parameters                         | 56         |
| Personnell.....                                           | 21  | RCF:Process-related parameters .                        | 55         |
| Physical data                                             |     | Reasonably foreseeable misuse ...                       | 21         |
| Technical data .....                                      | 123 | References (representation) .....                       | 13         |
| Pick-up order .....                                       | 120 | Refrigeration specialists .....                         | 24         |
|                                                           |     | Regulations .....                                       | 14         |

## Index

|                                                                                        |     |                                                                                  |     |
|----------------------------------------------------------------------------------------|-----|----------------------------------------------------------------------------------|-----|
| Relative centrifugal force (RCF),<br>speed and radius .....                            | 127 | Safety of in-vitro diagnostic devices<br>.....                                   | 22  |
| Relative centrifugal force: Process<br>parameters .....                                | 50  | Safety range around the centrifuge                                               | 34  |
| Remaining risks .....                                                                  | 24  | Scope of supply .....                                                            | 14  |
| Removing broken glass .....                                                            | 107 | Section-related warning notes .....                                              | 12  |
| Removing broken glass and metal<br>dust from the rotor chamber ....                    | 107 | Serious incident (MD or IVD)<br>Definition .....                                 | 22  |
| Repair .....                                                                           | 113 | Service<br>Control system .....                                                  | 90  |
| Representation conventions .....                                                       | 11  | Service contact .....                                                            | 119 |
| Requalification (repeat test) .....                                                    | 113 | Setting up the centrifuge .....                                                  | 34  |
| Responsibility of the operator .....                                                   | 21  | Settings .....                                                                   | 70  |
| Return transport .....                                                                 | 33  | Short run .....                                                                  | 48  |
| Returning a product .....                                                              | 120 | Sigma Comfort rotor coating<br>Service life .....                                | 109 |
| Returning the centrifuge to the<br>manufacturer .....                                  | 119 | Sound pressure level<br>Technical data .....                                     | 123 |
| Rights management .....                                                                | 82  | Speed, radius and relative centrifugal<br>force (RCF) .....                      | 127 |
| Risk of burns on hot surfaces inside<br>the rotor chamber .....                        | 26  | Standards .....                                                                  | 14  |
| Rotor and bucket selection .....                                                       | 59  | Standstill cooling: Process-related<br>parameters .....                          | 57  |
| Rotor identification system<br>Safety device .....                                     | 28  | Start timer: Process-related<br>parameters .....                                 | 54  |
| Rotor library .....                                                                    | 59  | Starting a centrifugation .....                                                  | 45  |
| Rotor program .....                                                                    | 124 | Status indicator .....                                                           | 67  |
| Rotors without Sigma Comfort<br>coating<br>Maintaining the load-bearing bolts<br>..... | 111 | Status menu .....                                                                | 67  |
| RPM: Process parameters .....                                                          | 50  | Stopping a centrifugation .....                                                  | 46  |
| RPM: Process-related parameters                                                        | 55  | Storing the centrifuge .....                                                     | 33  |
| RPM: Process parameters .....                                                          | 50  | Swing-out rotor removal .....                                                    | 43  |
| RPM: Process-related parameters                                                        | 56  | Switching the centrifuge on and off                                              | 36  |
| Runtime at set speed: Process-<br>related parameters .....                             | 55  | Symbols on the packaging .....                                                   | 32  |
| Runtime: Process parameters .....                                                      | 48  | System check<br>Safety device .....                                              | 27  |
| Runtime: Process-related parameters<br>.....                                           | 54  | System for stopping the centrifuge<br>due to an imbalance<br>Safety device ..... | 27  |
| <b>S</b>                                                                               |     | <b>T</b>                                                                         |     |
| Safety device<br>Mains power plug .....                                                | 27  | Table of chemical resistance data                                                | 137 |
| Mains power switch .....                                                               | 27  | Table of error messages .....                                                    | 116 |
| Protective earth conductor .....                                                       | 27  | Table of possible faults .....                                                   | 115 |
| Rotor identification system .....                                                      | 28  | Table of symbols .....                                                           | 131 |
| System check .....                                                                     | 27  | Temperature monitoring system<br>Safety device .....                             | 27  |
| System for stopping the centrifuge<br>due to an imbalance .....                        | 27  | Temperature: Process parameters                                                  | 50  |
| Temperature monitoring system                                                          | 27  | Temperature: Process-related<br>parameters .....                                 | 57  |
| Safety devices and guards .....                                                        | 27  | Thermal stress .....                                                             | 100 |

Tripped fuses ..... 35

Typographic conventions..... 13

**U**

Unit fuses..... 35

Unpacking the centrifuge ..... 33

User interface ..... 44

User login..... 83

User management ..... 81

UV radiation ..... 100

**V**

Ventilation..... 34

**W**

Warning messages

Identifying and eliminating faults

and malfunctions ..... 114

Warning notes ..... 11

Warranty claims..... 13

Work environment ..... 21